

Public Health Environmental Surveillance Open Data Model (PHES-ODM) Documentation

Version 2.0.0-rc.3.18

2023-03-09

Table of contents

Introduction	1
How to use the documentation	2
Future documentation	2
Quick start	3
1. Get oriented to the ODM tables where you record results	3
2. View mandatory and optional data fields in the tables reference guide	4
3. Try entering example data in ODM Excel templates . . .	4
4. Next steps	5
Getting help	5
Background	9
Getting to know the ODM	11
What are “Parts” and “Part Types”?	11
“Sets” and the rationale behind sets	12
Why category sets are treated differently	13
The “what” and “why” of Specimen IDs	13
The “what” and “why” of Groups and Classes	14
Tables and table types within the ODM	14
Program description tables	15
Results tables	15
Look-up tables	15
Measures, Methods, and Attributes: Key parts	15
Measures	16
Methods	16

Table of contents

Attributes	16
Data quality and reportability	17
Time periods for samples and measures	17
Dates within the ODM	18
Translation and language capabilities	18
Questions and online community	19
How to use the ODM	21
How-To guides	21
1) How to use templates	21
2) How to report Measures and Measure Sets:	23
Quick Start	24
Measures Template	24
Measure Sets Template	25
Detailed description	25
Measures template	25
Measure Sets Template	30
3) How to report Samples and Sample Relationships:	31
Quick Start	31
Samples Template	32
Sample Relationships Template	33
Detailed Description	33
Samples Template	33
Sample Relationships Template	35
4) How to report Protocols, Protocol Steps and Protocol Relationships:	37
Quick Start	37
Protocol Steps Template	37
Protocols Template	38
Protocol Relationships Template	39
Detailed Description	40
Protocol Steps Template	40
Protocols Template	41
Protocol Relationships Template	43

Table of contents

5) How to report Metadata:	44
Quick Start	45
Addresses Template	45
Contacts Template	46
Datasets Template	46
Instruments Template	47
Organizations Template	48
Polygons Template	48
Quality Reports Template	49
Sites Template	49
Detailed Description	50
Addresses Template	50
Contacts Template	51
Datasets Template	53
Instruments Template	54
Organizations Template	56
Polygons Template	57
Quality Reports Template	58
Sites Template	59
Example data for each of the 11 report tables.	61
 I. Reference guide	 63
 Parts	 67
a1306s delta-variant gene target	68
a1918v delta-variant gene target	68
a2710t omicron-variant gene target	68
a63t omicron-variant gene target	69
a67v omicron-variant gene target	69
Atline Analyzer	69
Absolute humidity	70
Absolute humidity unit set	70
Academic institution	70

Table of contents

Number of active cases	71
Active	71
Analysis date	71
Analysis date end	72
Analysis date start	72
Address Line 1	72
Address Line 2	73
Address table	73
Addresses table column order	74
Address header req	74
Address ID	74
Administrative regions	75
Air filter	75
Aggregation scale	75
Aggregation	76
Aggregations	76
Aggregation scales	77
Aggregation set	77
Aggregation sets	77
Air compartment	78
Air compartment set	78
Airplane	78
Airport	79
Airport lists worksheet	79
Air and surface compartment set	79
Environmental temperature	80
Air and water compartment set	80
Alias depreciated	80
Alias ID depreciated	81
All domains	81
Alleles	82
Allele unit set	82
Access to all org	82
Aloh3 precipitation	83

Table of contents

Assay method id default depreciated	83
Amicon ultrafiltration	84
Amicon filter, extract with MP96	84
Amplicon sequencing	84
AND Boolean aggregation	85
Any compartment set	85
Any specimen set	85
Area proportional sample	86
Arrival temperature	86
Use of the articV3 primer	87
Use of the articV4 primer	87
Assay ID	87
Attributes	88
Bacteria unit set	88
Miscellaneous bacteria group	88
Bead beating, extract with MP96	89
bcov	89
bcov culture spike target	90
bcov vaccine spike target	90
Beef extract flocculation	90
Is Before	91
Beta	91
Biologic	91
Digital droplet emulsification PCR	92
Binary Large Object (BLOB) data type	92
Breadth of coverage ($\geq 5\times$ depth)	92
5-Day carbonaceous biochemical oxygen demand	93
5-day total biochemical oxygen demand	93
Boolean data type	94
Boolean aggregation set	94
Boolean value set	94
Bovine respiratory syncytial virus group	95
brsv culture spike target	95
BRSV-N	95

Table of contents

brsv vaccine spike target	96
Branch sewer pipeline	96
Building cleanout	96
Capacity unit set	97
Case report date	97
Categorical data type	98
Categories	98
Long-term acute care hospital	98
Child day care	99
Degrees Celcius	99
Centrifugation	100
Centricon ultrafiltration	100
Ceres nanotrap	100
CFU per 100 ml	101
Changes column	101
Chemical	101
Chemagic viral dna/rna 300 kit	102
Child relationship	102
City	102
Clarified sample	103
Class	103
Classes	104
Cloudy	104
Centimetres	104
Carbon Dioxide	105
Chemical Oxygen Demand	105
Collection group	105
Collection date time	106
Collection date time end	106
Collection date time start	106
Sample collection set	107
Collection number	107
Collection number and period	108
Collection period	108

Table of contents

Sample collection type	109
Co-located sample	109
Column name depreciated	110
Composite sample - archival	110
Composite grab sample of 3	111
Composite grab sample of 3 depreciate	111
Composite 3hr grab sample	111
Composite 3hr grab sample depreciate	112
Composite 8hr grab sample depreciated	112
Composite 8hr grab sample deprecated	112
Compartment	113
Compartments	113
Compartment sets	114
Sample concentrate	114
Confirmed Cases Case report date	114
Is concurrent to	115
Water conductivity	115
Conductivity	116
Conductivity unit set	116
Number of confirmed cases	116
Confirmed cases episode date	117
Confirmed Cases Onset date	117
Contact ID	117
Contact name	118
Contact table	118
Contacts table column order	118
Contact header req	119
Contact table set	119
Confirmed cases test date	120
Cphd ID	120
COSCa ball	120
Country	121
Countries or sovereign states	121
Districts, counties, regions	121

Table of contents

Covid-19	122
SARS-CoV-2 measure	122
SARS-CoV-2-B.1.1.7	123
SARS-CoV-2-B.1.351	123
SARS-CoV-2-E	123
SARS-CoV-2-N1	124
SARS-CoV-2-N2	124
SARS-CoV-2-N3	124
SARS-CoV-2-P.1	125
SARS-CoV-2-RdRp	125
Covid-19 population measurement date	125
Cphd ID	126
crAssphage-N	126
Crosswal table set	126
Cycle threshold or quantification cycle (Ct or Cq)	127
Cumulative count	127
Custody Contact ID	128
Data custodian ID	128
d377y delta-variant gene target	128
d3g omicron-variant gene target { #d3g }	129
d63g delta-variant gene target	129
d796y omicron-variant gene target { #d796y }	129
d950n delta-variant gene target	130
Daily count	130
Dataset creation date	131
Dataset ID	131
Dataset table	131
Datasets table column order	132
Dataset header req	132
Data types	132
Date	133
Datetime data type	133
ddcov_e sars-cov-2 gene target	134
ddcov_n sars-cov-2 gene target { #ddcovN }	134

Table of contents

Deaths	134
del143/145	135
del 157/158	135
del2084/2084	136
del212/212	136
del3674/3676	136
del69/70	137
Delta	137
Depreciated	137
Derived sample	138
Description	138
Description of change	139
Desk or counter	139
Dewatered solids	139
Detected	140
Access to details deprecated	140
Development	140
Dictionary reference table set	141
Additional dictionary sheets	141
Dictionary set type	142
Dilution factor	142
Point dillutions	142
Diseases (human)	143
Dissolved gas concentration unit set	143
Depth of coverage	143
Domain	144
Domains	144
Higher education domitory or residential building	145
e156g delta-variant gene target	145
e484a omicron-variant gene target { #e484a }	145
Escherichia coli	146
Education	146
Efficiency	146
Contact email	147

Table of contents

European Nucleotide Association	147
ENA - notes	147
Rainfall	148
Ground snow depth	148
Snowfall	148
Episode date	149
ERD table set	149
Sarbecovirus-specific E sars-cov-2 gene target	150
Estimated frequency of reads	150
Estuary, natural water body	150
Experiment Failed	151
4s method	151
Nucleic acid extraction method	151
Nucleic Acid Extraction set	152
Extraction volume of sample	152
Fecal matter	153
FALSE	153
Field sample	153
Field sample replicate	154
File location of polygon	154
Filtration	154
First Nation	155
First name of contact	155
First released version	156
Foreign key	156
AI - Inhibition present but addressed	156
B - Trace levels of contamination	157
FI - Inhibition present and unaddressed	157
J - Weak signal extrapolation	157
ND - Non-detect	158
UJ - Trace signal extrapolation	158
UQ - Unquantifiable	159
Float data type	159
Flow-normalized mean	160

Table of contents

Floor	160
Flow rate	160
Flow	161
Flow proportional 24hr sample depreciated	161
Flow proportional sample	161
Volume flow rate unit set	162
Flow volume	162
Influenza virus measure	163
Influenza virus A1	163
influenza virus A2	163
Influenza virus B	164
Fluidigm digital PCR	164
Foggy	164
Fraction analyzed	165
Sample fraction set	165
Freyja Script	165
F-Specific RNA bacteriophages	166
F-Specific RNA bacteriophages, G2	166
Sample frozen	167
Field sample temperature	167
Full dictionary sheets set	167
Funder Contact ID	168
Funding agency ID	168
g215c delta-variant gene target	168
g339d omicron-variant gene target { #g339d }	169
g496s omicron-variant gene target { #g496s }	169
g662s delta-variant gene target	170
Gamma	170
Gene copies per copy of crAssphage	170
gene copies per day per 100,000	171
Gene copies per gram solids	171
Gene copies per L	171
Gene copies per mL	172
Gene copies per PMMoV copy	172

Table of contents

Genetics unit set	172
General missingness set	173
Generic quality flag set	173
European Petroleum Survey Group Coordinates	173
Latitude	174
Longitude	174
Type of geography	175
Geographic set	175
Well-known text	175
Gram per cubic metre	176
Geometric mean	176
Government agency	176
Grams	177
Gravity settling	177
Grab sample	177
Gromstole 1.0 Script	178
Group	178
Groups	179
Header	179
Heat inactivated sars-cov-2 virus spike target	179
hep g armored rna	180
hep g armored rna spike target	180
High severity	180
Sewer Network Health Region	181
Hand Measurement	181
Hollow fiber dead end ultrafiltration	182
General Hospital Admissions	182
Hospital Census	182
Hospital	183
Holding tank wastewater depreciated	183
Holding tank wastewater	183
Holding tank	184
Human compartment	184
Humidity	184

Table of contents

Human compartment set	185
i1566v omicron-variant gene target { #i1566v }	185
i3758v omicron-variant gene target { #i3758v }	185
i82t delta-variant gene target	186
International Classification of Diseases	186
Intensive care unit patients	186
Index	187
Inhibition measure	187
Inhibition method	188
Innovaprep ultrafiltration	188
Input	188
Instrument ID	189
Instrument table	189
Instruments table column order	190
Instrutment header req	190
Other instrument	190
Instrument Type	191
Describe other instrument type, if applicable	191
Instrument set	191
Integer data type	192
Intercept	192
Combined ip2 and ip4 sars-cov-2 gene target	192
ISO639-1	193
ISO639-2B	193
ISO639-2T	194
ISO639-3	194
ISO639-6	195
k856r omicron-variant gene target { #k856r }	195
Kilogram per second	195
Kilolitres	196
Litres	196
l452r delta-variant gene target	196
l981f omicron-variant gene target { #l981f }	197
Lab ID default depricated	197

Table of contents

Laboratory duplicate	198
Lab ID	198
Lagoon system	198
Lake, natural water body	199
Lambda	199
Language ID	199
Language family	200
Language name	200
Language Look-up table	201
Language table column order	201
Language header req	201
Last edited	202
Last name of contact	202
Last updated version	203
LaTeX expression	203
Low breadth of coverage	203
Laboratory control sample duplicate	204
Low depth of coverage	204
Leaked sample	204
Link	205
License	205
Linear scale aggregation set	206
Liquid fraction	206
List worksheet	206
List set	207
Access to local ha	207
Limit of detection (LOD)	207
Limit of detection (LOD) - sewage	208
Logarithmic scale aggregation set	208
Dictionary tables	209
Limit of quantification (LOQ)	209
Low severity	209
Lower limit of a 95% confidence interval	210
Low-volume sample	210

Table of contents

Long-term care facility	210
Long-term care - assisted living or retirement home	211
Other long-term care	211
Life technologies digital PCR	212
Lysis buffer	212
Cubic metres per day	212
Cubic metres per hour	213
Cubic metres per second	213
Manufacturer	213
Mate pair layout	214
Maximum length	214
Maximum value part support	214
Maximum Value	215
Control-plate Flag - NTCs { #Maybe not something we're even reporting, so.. no include? }	215
Arithmetic mean	215
Measurement group	216
Other Measure	216
Measures quality set	216
Measure	217
Measures	217
Report ID	218
Measure report table	218
Measure set report table	218
Measure sets table column order	219
Measure Set header req	219
Measures table column order	219
Measures header req	220
Report set ID	220
Median	220
Normalized arithmetic mean	221
Other Measure Description	221
Measure license	222
Method	222

Table of contents

Concentration method	222
Concentration method set	223
Methods	223
Membrane filtration with acidification and mgcl2	223
Membrane filtration with acidification and mgcl2, mem- brane recombined with separated solids	224
Membrane filtration with addition of mgcl2	224
Membrane filtration with mgcl2 and separated solids	225
Membrane filtration with no amendment	225
Membrane filtration with no amendment, membrane recom- bined with separated solids	225
Membrane filtration with sample acidification	226
Membrane filtration with sample acidification, membrane recombined with separated solids	226
Millions of gallons per day (MG/D)	226
Milligrams per litre	227
Murine Hepatitis Virus	227
Multiple issues	228
Mid-level severity	228
Minimum length	228
Minimum value part support	229
Minimum Value	229
Miscellaneous attribute group	229
Miscellaneous measure group	230
Missingness	230
Missingness sets	230
Mixed/homogenized sample	231
Millilitres	231
Megalitres per day (ML/d)	232
Millimetres	232
Measure, method, or attribute	232
Measure, method, or attribute sets	233
Model	233
Moore swab passive sample	233

Table of contents

Metres per second	234
ms2 coliphage	234
ms2 coliphage spike target	234
Matrix spike duplicate	235
Moore Swab Passive sample depreciated	235
Major sewer pipeline	235
Mu	236
murine coronavirus	236
murine coronavirus spike target	236
Multiple fraction	237
Multiple Purpose	237
Municipalities or communes	238
Municipality	238
Mutations	238
N sars-cov-2 gene target	239
Combined N1 and N2 sars-cov-2 gene target	239
n211i omicron-variant gene target { #n211i }	239
n679k omicron-variant gene target { #n679k }	240
n856k omicron-variant gene target { #n856k }	240
n969k omicron-variant gene target { #n969k }	240
Not applicable	241
Aggregation not applicable	241
Aggregation scale not applicable	242
Aggregation set not applicable	242
Class not applicable	242
Compartment not applicable	243
Compartment set not applicable	243
Domain not applicable	243
Group not applicable	244
Name	244
Not a number	245
Nomenclature not applicable	245
Specimen not applicable	245
Specimen set not applicable	246

Table of contents

Native Name	246
Unit not applicable	246
Unit set not applicable	247
National Center for Biotechnology Information	247
NCBI - notes	248
Neighbourhood	248
Local administrative units or neighborhoods	248
Population of newly vaccinated persons	249
NextClade nomenclature	249
Normalized geometric mean	249
Ammonium Nitrogen	250
niid_2019-ncov_n sars-cov-2 gene target	250
No quality concerns	250
Sample not labelled	251
No liquid concentration, liquid recombined with separated solids	251
Nomenclature	251
Nomenclatures	252
NORMAN	252
NORMAN - notes	253
Norovirus G1	253
Norovirus G2	253
Notes	254
Sample missing time stamp	254
Not reported	254
Not reported, Not applicable missing set	255
Sarbecovirus-specific N sars-cov-2 gene target	255
Nephelometric turbidity unit	255
Nuclisens automated magnetic bead extraction kit	256
Nuclisens manual magnetic bead extraction kit	256
Null	257
National Wastewater Surveillance System	257
NWSS - notes	257
NWSS - process	258

Table of contents

Water	258
coronavirus OC43	258
oc43 spike target	259
Ocean, natural water body	259
Other Collection depreciated	259
Lab Analysis	260
Omicron BA.1	260
Omicron BA.2	261
Omicron BA.2.75	261
Omicron BA.4	261
Omicron BA.5	262
Onset date	262
Online Sensor	262
Ontology reference	263
Orthophosphates	263
Other residential building	264
OR Boolean aggregation	264
ORF1a sars-cov-2 gene target	264
ORF1ab sars-cov-2 gene target	265
ORF1b sars-cov-2 gene target	265
Organization ID	265
Organization table	266
Organizations table column order	266
Organziation header req	266
Organization level	267
Organization level set	267
Organization sector	268
Organization sector set	268
Organization Type	268
Organization type set	269
Sample origin	269
Sample Origin set	269
Other Collection depreciated	270
Other aggregation scale	270

Table of contents

Other	270
Other aggregation set	271
Other aggregation deprecated	271
Access to other prov deprecated	272
Other report table set	272
Variants Other	272
Other Site - sample depreciated	273
Other Site - Site depreciated	273
Other Unit deprecated	274
Outbreak	274
Outbreak End	275
Outbreak - on-going	275
Outbreak set	275
Outbreak start	276
p100l delta-variant gene target	276
p2046l delta-variant gene target	276
p2287s delta-variant gene target	277
p3395h omicron-variant gene target { #p3395h }	277
p681r delta-variant gene target	277
Paired Layout	278
Pangolin nomenclature	278
Parent dataset ID	279
Parent sample ID	279
Parent Site ID	279
Part description	280
Part ID	280
Part instruction	281
Label	281
Parts Look-up table	282
Parts table column order	282
Parts header req	282
Part types	283
Primary Clarifier Effluent depreciated	283
Postal or Zip Code	284

Table of contents

PCR method	284
PCR quality set	284
PCR sequencingsSelection method	285
PCR Method set	285
Primary Clarifier Sludge depreciated	285
Primary clarifier effluent	286
Polyethyleneglycol (PEG) precipitation	286
Percent	286
Percent recovery	287
pH	287
Access to phac	287
PHAGE	288
PHAGE - notes	288
Phenol chloroform	289
Pseudomonas virus phi6	289
phi6 spike target	289
Contact phone	290
Total phosphorous	290
Total Phosphates	290
PHRED quality score	291
Physical property	291
Physical	291
Primary key	292
PMMoV- and flow-normalized mean	292
PMMoV-normalized mean	292
Population	293
Post-Grit depreciated	293
Polygon ID	294
Polygon table	294
Polygons table column order	294
Polygon header req	295
Polygon Population	295
Pooled	295
Type of date for case reporting depreciated	296

Table of contents

Case reporting date set	296
Population equivalents	297
Population group	297
Population Served	297
Population unit set	298
Poulatiion or site set	298
Population specimen set	298
Percent positive	299
parts per million	299
PMMoV-CP	299
Percent positivity rate	300
Percent primary sludge	300
Precipitation	301
Precipitation unit set	301
Pretreatment	301
Genetic primer	302
Genetic primer set	302
Private sector	302
Processing group	303
Program decription tables	303
promega automated tna kit	304
promega ht tna kit	304
promega manual tna kit	304
promega wastewater large volume tna capture kit	305
Proportion of total	305
Proportion of variant in sample	305
Protocol ID	306
Protocol ID container	306
Protocol ID object	306
Protocol ID subject	307
Protocol relationships table	307
Protocol relationships table column order	308
Protocol Relationship header req	308
Protocol Relationships set	309

Table of contents

Protocols table	309
Protocols table column order	309
Protocols header req	310
Protocol steps table	310
Protocol steps table column order	311
Protocol steps header req	311
Protocol table set	311
Protocol version	312
Access to prov ha	312
Provisional report	312
Primary clarifier sludge	313
Pumping station	313
Post-grit	314
Pretreatment description - depreciated	314
Number of positive tests	314
Access to public	315
Puro virus	315
Puro virus spike target	315
Purpose	316
Purpose category	316
q19e omicron-variant gene target { #q19e }	317
q493r omicron-variant gene target { #q493r }	317
q498r omicron-variant gene target { #q498r }	317
q954h omicron-variant gene target { #q954h }	318
Quaility assurance method	318
Quality concerns	318
qiagen allprep dna/rna kit	319
qiagen digital PCR	319
qiagen ez1 virus mini kit v2.0	319
qiagen allprep powerfecal dna/rna kit	320
qiagen allprep powerviral dna/rna kit	320
qiagen powerwater kit	320
qiagen qiaamp buffers with epoch columns	321
qiagen rneasy kit	321

Table of contents

qiagen rneasy powermicrobiome kit	321
Quantitative PCR	322
Qualitative	322
Quality control	322
Quality flag	323
Quality report	323
Quality indicators	323
Quality indicator sets	324
Quality reports table	324
Quality rerpots table column order	325
Quality report header req	325
Quality set	325
Quantitative	326
R squared	326
r203m delta-variant gene target	326
Raindance digital PCR	327
Rainy	327
Random sequencing selection method	327
Ratio	328
Raw wastewater	328
Raw sewage at site	329
Raw sewage downstream from a site	329
Raw sewage upstream from a site	329
IP2 rdrp sars-cov-2 gene target	330
IP4 rdrp sars-cov-2 gene target	330
Date sample recieved	330
Recovered patients	331
Reference link	331
Regular	332
Regular reports table set	332
Relative humidity	332
Relative humidity unit set	333
Representative grab sample	333
Replicate Type set	334

Table of contents

replicate number	334
Primary reporting authority	334
Secondary reporting authority	335
Reportable	336
Report date	336
Reporter ID depreciated	336
Reporter table depreciated	337
Replicate Type	337
Resuspend COSCa filter collection	338
Results tables	338
Retention pond	338
River, natural water body	339
Role of contact	339
Rate per 100,000	339
Raw post-pasteurized wastewater	340
Respiratory syncytial virus	340
s2083i omicron-variant gene target { #s2083i }	341
s371l omicron-variant gene target { #s371l }	341
s373p omicron-variant gene target { #s373p }	341
s375f omicron-variant gene target { #s375f }	342
s477n omicron-variant gene target { #s477n }	342
Sample	342
Sample material	343
Sample ID	343
Sample ID object	343
Sample ID subject	344
Sample material set	344
Sample quality set	345
Sample relationships table	345
Sample relationships table column order	345
Sample relationships header req	346
Sample relationship	346
Sample relationships set	347
Sample report table	347

Table of contents

Sample shed	347
Samples table column order	348
Sample header req	348
Collection other depreciated	349
SARS-CoV-2	349
Sample specimen set	349
Sparse coverage	350
Secondary clarifier effluent depreciated	350
School	351
Sequencing script version set	351
Genetic sequencing script version	351
Secondary clarifier sludge depreciated	352
Standard Deviation	352
Normalized Standard Deviation	352
Sea, natural water body	353
See unit (aggregation scales)	353
See unit (data type)	354
Secondary clarifier effluent	354
Access to self	354
Date sample was sent	355
Septic tank wastewater	355
Septic tank	355
Sequencing Layout	356
Sequencing Layout set	356
Sequencing quality set	357
Sequencing selection method	357
Sequencing selection method set	357
Sequencing Strategy	358
Sequencing Strategy set	358
Genetic sequences	358
Sewer Sediments Depreciated	359
Set	359
Sets look-up table	360
Sets table column order	360

Table of contents

Sets header req	360
Settled solids	361
Set type	361
Set Value	361
Severity indicator	362
Severity set	362
Sewer network file link depreciated	363
S sars-cov-2 gene target	363
Sampleshed set	363
Ship	364
Shipped on ice	364
Short names	364
Site	365
Single	365
Single Layout	365
Site or sample specimen set	366
Site specimen set	366
Site ID default	367
Site features group	367
Site ID	367
Site measure ID	368
Sites table	368
Sites table column order	369
Sites header req	369
Site Type	369
Site set	370
Skimmed milk flocculation	370
slope	370
Sewer network file link depreciated 2	371
Snowy	371
Solid fraction	371
Solid seperation	372
Solid separation set	372
Source Protocol ID	373

Table of contents

Protocol step source ID	373
Specific humidity	373
Specific humidity unit set	374
Specifies	374
Specimen	375
Specimens	375
Specimen sets	375
Spike material	376
Spike material set	376
Recovery efficiency spike target	376
Recovery efficiency spike target set	377
Sample spilled	377
Secondary clarifier sludge	378
Social services shelter	378
Stabilization pond	378
Standard concentrations	379
Standard curve	379
Departments, states, or provinces	379
State, Province, or Region	380
Status	380
Status set	381
Standard concentration unit set	381
Storage temp	381
Protocol Step ID	382
Step ID Object	382
Step ID Subject	382
Method step parent ID	383
Protocol step version	383
Sewer network file blob	384
Storage time	384
Sub-variant or lineage	385
Summary	385
Summary worksheet	385
Sunny	386

Table of contents

Surface compartment	386
Other surface	386
Surface compartment set	387
Surface and water compartment set	387
Surface swab	388
Sewage truck	388
Sewer sediment	388
Sewer catchment area	389
Synthetic sample	389
t19r delta-variant gene target	389
t3646a delta-variant gene target	390
t547k omicron-variant gene target { #t547k }	390
t9i omicron-variant gene target { #t9i }	390
Tables	391
Table column set	391
TaqPath N sars-cov-2 gene target	391
TaqPath S sars-cov-2 gene target	392
Sample temperature	392
Temperature	393
Temperature unit set	393
Nucleic acid template volume	393
Test date	394
Number of tests performed	394
Testing	394
thermo magmax microbiome ultra nucleic acid isolation kit	395
Time proportional sample	395
Total Kjeldahl Nitrogen	395
Total Nitrogen	396
Time Proportional 24hr sample	396
Translation look-up table	397
Translation table column order	397
Tranlation header req	397
trizol, zymo mag beads w/ zymo clean and concentrator	398
TRUE	398

Table of contents

Total solids concentration	398
Total suspended solids	399
Concentration of total suspended solids	399
Turbidity	400
Turbidity unit set	400
Tube weight empty	400
Tube weight full	401
Type other depreciate	401
Type other depreciate2	401
University campus	402
Undisclosed	402
Unique sample	403
Unit	403
Unitless measure	403
Unitless unit set	404
Units	404
Unit sets	404
Ultracentrifugation	405
Upper limit of a 95% confidence interval	405
Upstream sites	406
Micro-Siemens per centimetre	406
U site measure ID	406
v2930l delta-variant gene target	407
Validation study	407
Value	408
Variable character data type	408
Frequency of variants detected	408
Variants	409
Population with 1 dose of vaccine	409
Population with 2 doses of vaccine	409
Population with 3 doses of vaccine	410
ODM V1 category	410
ODM V1 location	410
ODM V1 table	411

Table of contents

ODM version 1 to 2 changes	411
ODM V1 variable	411
Miscellaneous viruses group	412
Size in volume	412
Volume proportional sample	413
Volume unit set	413
Volatile suspended solids	413
Volatile suspended solids - ignition	414
Water compartment	414
Water depreciated	414
Water compartment set	415
Wastewater Temperature	415
Weather	416
Weather set	416
Whole genome Sequencing	416
WHO nomenclature	417
Wide 95% interval	417
Wide name table	418
Wind Speed	418
Wind speed unit set	418
Measure identification (wide table)	419
Sample stored incorrectly	419
Sample stored at wrong temperature	419
World Sphere	420
W-SPHERE notes	420
Wastewater treatment plant	420
Back-up wastewater treatment site	421
Wastewater treatment plant designed capacity	421
Industrial wastewater treatment plant	422
Municipal wastewater treatment plant for combined sewage	422
Municipal wastewater treatment plant for sanitary dewage	
only	422
y505h omicron-variant gene target { #y505h }	423

Table of contents

zymo environ water rna kit/ zymo environ water rna kit (cat. r2042)	423
zymo quick-rna fungal/bacterial miniprep #r2014	424
Sets	425
Aggregations	425
Aggregation scales	427
Aggregation sets	428
Boolean aggregation set	428
Linear scale aggregation set	428
Logarithmic scale aggregation set	430
Aggregation set not applicable	431
Other aggregation set	431
Attributes	432
Boolean value set	448
Classes	449
Sample collection set	450
Compartments	451
Compartment sets	452
Air compartment set	452
Air and surface compartment set	453
Air and water compartment set	453
Any compartment set	453
Human compartment set	454
Compartment set not applicable	454
Surface compartment set	454
Surface and water compartment set	455
Water compartment set	455
Data types	456
Domains	457
Nucleic Acid Extraction set	457
Sample fraction set	460
Geographic set	461
Groups	461

Table of contents

Instrument set	462
Measures	463
Concentration method set	478
Methods	479
Missingness	481
Missingness sets	481
General missingness set	481
Not reported, Not applicable missing set	482
Nomenclatures	482
Organization type set	483
Sample Origin set	484
Outbreak set	484
PCR Method set	485
Genetic primer set	485
Protocol Relationships set	486
Purpose category	486
Quality indicators	487
Quality indicator sets	490
Generic quality flag set	490
Measures quality set	491
PCR quality set	493
Sample quality set	494
Sequencing quality set	495
Replicate Type set	496
Sample material set	497
Sample relationships set	499
Sequencing script version set	500
Sequencing Layout set	500
Sequencing selection method set	500
Sequencing Strategy set	501
Severity set	501
Sampleshed set	502
Site set	504
Solid separation set	506

Table of contents

Specimens	507
Specimen sets	508
Any specimen set	508
Specimen set not applicable	508
Poulatiion or site set	509
Population specimen set	509
Sample specimen set	509
Site or sample specimen set	510
Site specimen set	510
Spike material set	510
Recovery efficiency spike target set	511
Status set	512
Tables	513
Table column set	515
Units	516
Unit sets	522
Absolute humidity unit set	522
Allele unit set	522
Bacteria unit set	523
Capacity unit set	523
Conductivity unit set	524
Dissolved gas concentration unit set	524
Volume flow rate unit set	525
Genetics unit set	525
Unit set not applicable	526
Population unit set	526
Precipitation unit set	529
Relative humidity unit set	529
Specific humidity unit set	529
Standard concentration unit set	529
Temperature unit set	530
Turbidity unit set	530
Unitless unit set	530
Volume unit set	531

Table of contents

Wind speed unit set	531
Weather set	532
Tables	533
Table types	533
Header roles	534
Mandatory, optional, mandatoryIf	535
Table list	535
Contact table	536
Dataset table	537
Instrument table	539
Language Look-up table	540
Measure report table	542
Measure set report table	545
Organization table	546
Parts Look-up table	547
Polygon table	556
Protocol relationships table	557
Protocols table	559
Protocol steps table	560
Quality reports table	561
Sample relationships table	562
Sample report table	563
Sets look-up table	566
Sites table	567
Translation look-up table	569

Introduction

Make collecting and sharing high-quality, interoperable environmental surveillance data as easy as possible.

Warning

This documentation for the PHES-ODM Version 2.0.0 release candidate. Please be aware that changes will occur regularly in the next weeks. If there is something you'd like to see in the documentation that you cannot currently find, please please [create an issue on our GitHub repository](#) or post about it on [our Discourse page](#).

The Public Health and Environmental Surveillance Open Data Model (PHES-ODM or ODM) is an open-source data model and set of tools aimed at promoting global collaboration and use of environment and public health data. The ODM model facilitates interoperable environmental surveillance and research by improving quality reporting data, enhancing data management, and supporting data provenance and sharing.

This documentation is intended for those who want to understand how the ODM model works, its organizational structure, and how to use it effectively for their program. It is also a useful resource for those who currently use a different data dictionary and wish to share their data with a program that uses the ODM. If you need to share your data with another program that uses a different data dictionary, this documentation can help you compare and create a common cross-walk between dictionaries. The ODM starts as a conceptual model that represents real-world processes

Introduction

through data elements and data relationships, which are common to all environmental and public health surveillance systems.

How to use the documentation

There are three document sections, each with a different purpose.¹

Getting to know the ODM is oriented to understanding the ODM. This section explains the ODM, its structure, and how the model was designed and developed.

How to use the ODM are step-by-step guides to help you use different aspects of the ODM such as report templates.

Reference is a description of each part of the ODM. It serves as a reference for the other sections, and the definitive record of each part, table and set – along with their attributes.

Future documentation

Tutorials and videos are being develop. Tutorials will help newcomers get started. They will be lessons on how to:

- Describe a protocol using protocol tables.
- Make wide names and report templates using wide names.

¹[David Laing and Divio's Grand Unified Theory of Documentation.](#)

Quick start

You can get to know the ODM by imaging how to record a polymerase chain reaction (PCR) test result for SARS-CoV-2 from a wastewater sample.

1. Get oriented to the ODM tables where you record results

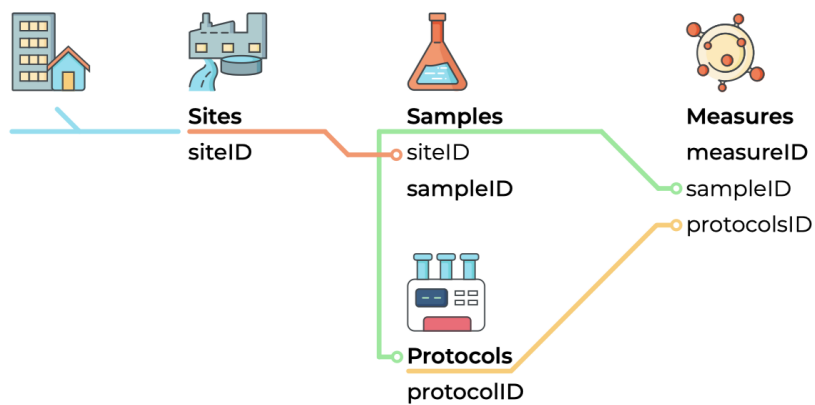


Figure 1.: **Figure 1 - Key ODM tables for recording data**

The **Sample** report table ([samples](#)) and the **Measure** report table ([measures](#)) are where you record information about samples and measures.

2. View mandatory and optional data fields in the tables reference guide

You can view mandatory and optional data field here: ([samples](#)) and ([measures](#)). ‘Mandatory’ means they are fields are needed to ensure interoperability.



What are ODM parts and how to find them?

In the ODM, tables are an example of an part. A part name can be identified by its text style like **this**. The part identifier often follows the part in parenthesis, with a link to the part in the reference section. So, **Sample report table** is the part name, and ([samples](#)) is the part ID and link to the documentation for that part.

You can find out more about parts on the [reference guide section](#).

3. Try entering example data in ODM Excel templates

Keep in mind, ODM templates are only examples of how to use the ODM dictionary to record data. Templates are available at [OSF.io](#). The ‘[How to use the ODM](#)’ section describes how to enter information into tables.

Storing sample data in an ODM template

A **sample** ([sampleID](#)) is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. The sample template stores one sample per row. The **Sample report table** has seven mandatory fields that are listed in the [reference[(tables.qmd#samples)] for the table.

<https://vimeo.com/802349692>

Storing measure data in an ODM template

A **measure** (**measure**) is an observation or measurement of any substance, including biological, physical, or chemical substances, obtained from a specimen such including a site, sample, person, or population. Each measure is stored in the **Measures** table template as a separate row, following the “long” table format, which is the primary method for storing measures. For example, if you recorded PCR results for two SARS-CoV-2 regions, covN1 and covN2, you would record each measurement on a separate line. Alternatively, you can use a “wide” table format to store multiple measures in a single row or create tables that combine both long and wide headers. The [European Union airport template](#) is a hybrid template that exemplifies this (found in the “Airports” folder).

(Video)

4. Next steps

In addition to recording measures from samples, you can also record measures for a specific site, such as wastewater flow rate (**flowRate**), or the population served by the site, such as the Covid-19 hospitalization rate for a wastewater treatment plant’s watershed or municipality. **Figure 2** illustrates the these three types of measurements: **sample**, **site**, and **population**.

The “who, where, and why” of these measurements is stored in contact tables, while the “how,” which refers to the measurement protocol or method, is stored in the protocols tables.

Getting help

We hope that the documentation serves its purpose and that you find it useful in figuring out how to best use the ODM and its suite of tools. If

Introduction

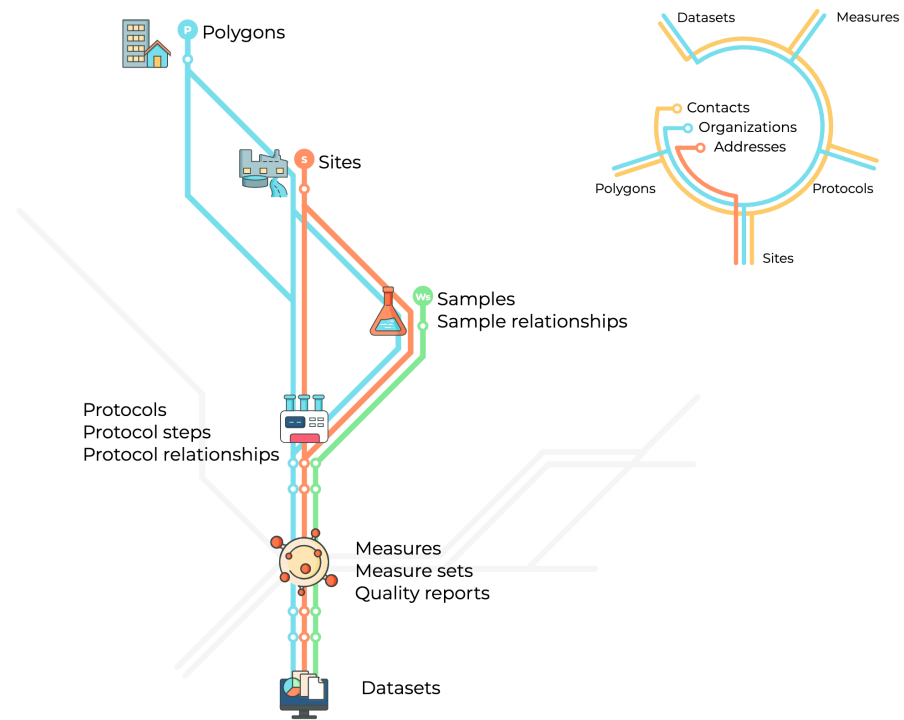


Figure 2.: **Figure 2 - ODM tables and 'subway' lines for storing measures**

Quick start

you do encounter any issues or questions that do not seem to have an easy answer, please [create an issue on our GitHub repository](#) or consult [our Discourse page](#).

Background

The ODM is designed to support a broad range of users in the field of environmental surveillance, including wastewater laboratories, municipalities, public health agencies, researchers, and citizen scientists. Its goal is to provide a comprehensive and extensible approach to data collection and management across the lifecycle of environmental surveillance, from sample collection and analysis to public reporting.

The ODM dictionary includes provisions for a wide range of biological, chemical, and physical properties, such as antimicrobial resistance, drugs, and toxins, as well as metadata for laboratory protocols that may affect measurement consistency across laboratories. The ODM supports laboratory information and management systems (LIMS) with measure to describe ensure data quality control and provenance.

To support the open science principle of FAIR data (Findable, Accessible, Interoperable, and Reusable), the ODM strives to make its data easily discoverable, accessible, and usable by all stakeholders.² It also adopts an open-source approach, allowing anyone to inspect, modify, and enhance the codebase in a transparent and collaborative manner on platforms like GitHub.

While the ODM is a complex dictionary, it includes tools to support its use, such as data validation, input templates, and SQL database definitions for data storage. Additionally, it supports collaborations among different organizations seeking to share and combine samples and measurements. As an open model, third parties can extend the model, such as

²[FAIR Principles](#).

Background

the European Union Digital European Exchange Platform (EU4S-DEEP) data templates for wastewater surveillance at airports, European Union [Digital European Exchange Platform \(EU4S-DEEP\)](#) data templates for wastewater surveillance at airports, the [CETO Epidemiologic platform](#), and [Ottawa Automatic Data Pipeline](#) to analyse qPCR data and transform it to ODM format.

Overall, the ODM aims to support the evolving field of environmental surveillance and science by providing a standardized, extensible approach to data collection and management.

Getting to know the ODM

What are “Parts” and “Part Types”?

Within the ODM, every component or element of the model and the dictionary is called a “part”. Parts are grouped into “part types”.

Because the term “parts” refers to every component of the ODM, “part types” are a way to differentiate between parts which have different functions and structures. The three most important part types are:

- **Measures:** A measurement or observation of any substance including a biological, physical or chemical substance.
- **Methods:** A procedure for collecting a sample or performing a measure.
- **Attributes:** A description of the who, where, when, and why of environmental surveillance.

See the [parts reference document](#) for more details.

There are additional part types to support the dictionary including:

- Aggregations
- Aggregation scales
- Categories
- Classes
- Compartments
- Dictionary support
- Domains

Getting to know the ODM

- Groups
- Missingness
- Nomenclature
- Quality indicators
- Specimens
- Tables
- Units

For information and details about all the part types, please see the [parts reference document](#).

“Sets” and the rationale behind sets

“Sets” are a way of grouping together different possible categorical inputs within the ODM. The purpose behind sets is to group options together so that users only need to check a few options from a drop-down menu, rather than scroll through a long list. Sets are also designed so that a single part can be used in multiple sets, which avoids the need to create multiple versions of this part for these different use contexts. For example, if you’re taking a measure of concentration, that measure will likely populate the unit choice drop-down list with units from the “standard concentration unit set”. Units in this set include Milligrams per litre, parts per million, and Percent Primary Sludge. Similarly, if one were measuring the amount of oxygen in a wastewater sample they would be offered units from the “dissolved gas unit set”. The units in this set include parts per million, and Milligrams per litre. You see here how Milligrams per litre is used twice across the sets, but exists still as a single part in the parts list. Unit sets are only one type of set, with the others being:

- Aggregation sets
- Compartment sets
- Quality sets
- Specimen sets

The “what” and “why” of Specimen IDs

Details about these sets are their similar - but unique - applications can be found in the [parts reference document](#). Importantly, “category sets” are not considered a “set” like these others. See the section below for an explanation of this difference.

Why category sets are treated differently

One of the reasons why “sets” are set up and organized the way that they are, is to allow for the re-use of parts in multiple sets. Category sets are, however, an exception to this rule because the categories that make up a category set are used in only one set. Furthermore, the general sets are needed for almost every measure or method (which require unit, aggregation, quality, and compartment information), while category sets are used only in specific circumstances for specific fields. These are things such as collection metadata (example: Sample collection category set), dictionary metadata and model architecture (example: Data Type Category Set), or more detailed types of certain methods (example: Nucleic Acid Extraction Category Set). Category sets are also a distinct part type, and more information about them can be found in the [parts reference document](#).

The “what” and “why” of Specimen IDs

Specimen IDs indicate the level at which a measure is being done. This can be: a site measure, ie. a measure at a site such as temperature or weather; a sample measure, ie. a laboratory measurement on a sample, such as gene copies of SARS-CoV-2 per mL; a person measure, ie. a measure of something at the level of an individual person, such as a blood pressure reading; or a population measure, ie. an aggregate measure of a population, such as the number of confirmed cases of a given illness. While there is currently nothing in version 2.0 of the ODM with uses the person specimen ID explicitly, it has been included to allow for the possibility of storing individual-level data in a future version. The reason

Getting to know the ODM

specimen IDs were created was so that there could be a single manner in which measures and methods are recorded, regardless of the level at which they were performed, while still maintaining the possibility to collect site, sample, and population measures as a distinct types of metadata.

The “what” and “why” of Groups and Classes

Groups and classes, similar to sets, are ways of grouping together different measures within the ODM. Given that the ODM aims to be as robust as possible and has a very long list of possible measures, groups and classes were designed to give shorter lists of measures in the drop down lists by specifying details about the kind of measure that a user is recording. Groups and classes can work together to further specify what kinds of measures are being reported. For example, the group `sarsCov2` contains many measures, but by specifying that the class is an `allele` or `variant`, the list of possible measures is pared down. Alternatively, class can be said to be non-applicable, paring down the measures in the `sarsCov2` group to only be unspecified measures of the quantity of the virus.

Tables and table types within the ODM

Within the ODM there are three types of tables: program description tables, results tables, and look-up tables. These table types exist to differentiate the function of these tables and to highlight these differences to users. The different tables take different types of inputs and maintenance from users, so understanding the differences can be important.

Measures, Methods, and Attributes: Key parts

Program description tables

Program description tables (represented in yellow in the ERD) are tables used to record metadata on the organizations, locations, methods, and appurtenance. These tables help to describe surveillance and testing programs, and are intended to be updated infrequently.

Results tables

Results tables (represented in blue in the ERD) are the tables used to record details on samples and measures. These tables record the main outcomes data and are updated daily, if not more frequently.

Look-up tables

Look-up tables (represented in green in the ERD) are the tables that are pre-programmed and pre-populated in the ODM. These hold information on sets, all parts, languages, and translation abilities. These are only updated by the ODM team in version updates.

Measures, Methods, and Attributes: Key parts

While there are many part types, there are three main part types users should be most familiar with: measures, methods, and attributes. These three have parallels with the three table types, and the differences and details of these parts are useful to understand.

Getting to know the ODM

Measures

Measures are actually types of measures that can be performed. These can range from temperature to the number of gene copies in a sample. The measure, or `measureID`, can be selected from a drop down in the templates and it specifies the kind of measure you intend to record. The actual value of the measure is then recorded in the `value` field of the `measures` table, with units and aggregation specified in the `unitID` and `aggregationID` fields.

Methods

Similar to measures, methods are types of methods that can be performed to accomplish a measurement. These can be diverse, ranging from incubation, qPCR, or nucleic acid extraction. The method itself (`methodID`) is selected from a drop down in the templates, specifying in general terms the kind of method the user wishes to record. From there, the `value` field of the `methodSteps` table can be populated by one of the inputs from the category set associated with that `methodID`. This provides a higher level of detail for a given method. For example, the `methodID` might be `solidSep` for solid separation of a sample. The value field might then be populated with `cent` for centrifugation, indicating more details about how the sample settling was done.

Attributes

Attributes are the largest category of part type, as these refer to any most of the fields in the ODM. Attributes are fields for metadata within the ODM and range greatly in the use. They comprise everything from collection dates, to sampling period, to names.

Data quality and reportability

Within the ODM there is a `qualityFlag` field present in both the `samples` and `measures` tables. This field serves to highlight whether or not there is a quality issues with the sample or the measure. The quality flag also allows a user to specify to type of quality issue. This is managed through the use of quality sets, which are the sets that contain the possible quality flags for a given measure or for samples. The idea is that this provides data on any issues with a sample and an indication of the nature of that issue. Having data about the nature and presence of a quality issue is often not sufficient for decision makers who are trying to use and interpret the data. As such, there is also a `reportable` field with is a Boolean indicator of whether or not data can or should be reported or included in final reports and decisions.

Time periods for samples and measures

For measures, there is a field `aDateStart` and `aDateEnd` which specify the date and time that an analysis was begun and finished. This allows for the recording of greater detail around timelines for especially long, multi-day analyses. For shorter analyses, which will likely make up the bulk of reported measures, the same date can input in both fields. The idea is to have the reporting tables for measures be as robust as possible to allow for various kinds of timelines.

For samples, there is `colldt`, `colldtStart`, and `colldtEnd`. The first field is for the collection date and time of a single grab sample, so start and end are not necessary and the other two fields can be left blank. For composite or pooled sample, the start and end date and time for collection is crucial information to know. When these two fields are populated, the `colldt` field can be left blank.

Dates within the ODM

There are a number of date fields within the ODM which all serve different purposes. Collection datetime (`collDT`) is the date a sample was collected, used only for grab samples. This field is left blank if the collection datetime start and end fields (`collDTstart` and `collDTend`) have been populated instead. Inversely, collection datetime start and end fields should be left blank if the collection datetime field has been populated. Similarly, analysis datetime start and analysis datetime end (`aDateStart` and `aDateEnd`) are used to report the date and or timeline of an analysis for a measure.

Date fields that are more related to data processing and laboratory infrastructure are: the last edited date (`lastEdited`) which indicates the last time a table, or a measure or sample details, were last edited or updated; the sent date (`sentDate`), or date that a sample was sent to the lab from the field; the received date (`recDate`), or the date the sample was received in the lab; and the report date (`repDate`), or the date that the analysis results or measures were reported. These give a greater indication of how up to date data is, but also on the speed and efficiency of the pipeline between sampling and reported results. Dates should be reported in day/month/year format to accommodate the most popular global convention.

Translation and language capabilities

The default language of the ODM is English, but French translations of all descriptive elements of the data model dictionary are also available. As additional nation states and partners adopt the ODM, we anticipate that these fields will be translated into other languages as well. The translation capacities are managed through the language look-up table (`languageLUs`), the translation look-up table (`translationLUs`), and the

Questions and online community

parts table (**partsLUs**). When a translation for a given part is not available, the dictionary will default to the English term. The language look-up table stores linguistic and classification codes for spoken human languages, with the most recent ISO639 code being the language ID (**languageID**). In the translation table, the language ID is paired with every part ID (**partID**), along with the label (**partLabel**), description (**partDesc**), and instructions (**partInstr**) translated to that language. The part IDs are linked to the full parts list which otherwise contains metadata that is coded in variables and requires no further translation.

Questions and online community

If any users have additional questions or issues with the ODM, we invite them to check out our [Discourse Page](#) for discussion boards and community support. For larger issues, and to ask the ODM team to add additional fields or variables into the model, we encourage users to visit the project's [GitHub repository](#) and to create an issue there so that a team member can respond.

How to use the ODM

How-To guides

1) How to use templates

ODM Excel report templates provide you with an easy way to enter your wastewater and other environmental data into the PHES-ODM format. The templates also serve as an interactive way to get to know the ODM. Below you will find a series of guides that will walk you through how to enter data into the report templates.

In these guides, “field” refers to a column in the template, and “entry” refers to a row of information that you entered. Guides also provide definitions and explanations for some terms; however, more complete information is found in the most can be found in the reference section of the ODM documentation.

The report templates can be found at the OSF.io [PHES-ODM reference file folder](#). More information about the tables can be found in the [getting-to-know](#) section.

Briefly,

- **samples** and **measures** tables are used day-to-day reporting. Information about a sample will be recorded in the **samples** table. Measurements are recorded in the **measures** table. A measure is

How to use the ODM

‘measurement or observation of any substance including a biological, physical or chemical substance’.

- There are three optional tables to support reporting samples and measures.
 - **sampleRelationships** is optional and used only when you have samples that are pooled or split.
 - **measureSets** is also an optional table that can be used to link measures together.
 - **qualityReports** is used for detailed quality assurance and control measures (you can—and should— record a **quality flag** attribute for each measure or sample in the **measures** and **samples** tables).
- Typically, there is one-time initial setup where information is entered the tables with the yellow tabs in the Excel template. These tables include information about **sites**, **addresses**, **organizations**, and **contacts**.
 - **sites** is a mandatory table records information about where testing is performed.
 - **addresses**, **organizations**, and **contacts** provide the contact information for sites and also other entities such as laboratories, or database custodians.
- There are optional supporting tables that include **protocols** and **instruments** where you record information about how samples are taken and measures are reported. These tables are usually completed once and then updated when there are changes to the protocol. **protocolSteps** and **protocolRelationships** tables are part used alongside the **protocols** table. Other supporting tables include **polygons** where information is store about the sample shed

2) How to report Measures and Measure Sets:

and **datasets** which includes information about the data custodian and owner of the data.

The templates contain four main types of fields, each with their own for instructional purposes:

- **Identifier, unique identifier and composite unique identifier (red):** In these fields you will create a value that will be used to identify something. For instance, **Sample** ID is used to identify each unique sample. An identifier can be repeated between entries, while an unique identifier cannot. There is usually only one unique identifier per template, and they correspond to the name of the table. e.g. **Measure** ID for the **measures** table, and **Contact** ID for the **contacts** table.
- **Drop-down and conditional drop-down menus (green):** These fields have a drop-down menu which contain the values that you can either select **or** type in. In conditional drop-down menus, the values depend on what you entered in a previous field.
- **Date fields (purple):** Dates are entered in these fields and must be in the form of the [ISO 8601](#) format: `yyyy-mm-dd`. All date fields also support entries for hours, minutes, and seconds: `2022-01-01T06:11:54` and `2022-01-01T06:11:54+13:30`
- **Free form fields (blue):** In these fields, you can enter your data. Any restrictions will be specified in the relevant sections. Keep in mind, different measures may have different data types that can not be validated in Excel. Data types each measure can be checked in the reference list for the measure. For example, the data type for **temperature** is a real number.

2) How to report Measures and Measure Sets:

In this guide you will learn how to record data in the **Measures** and **Measure Sets** table templates. In the **Measures** table, each row repre-

sents a measure. Along the same lines, in the **Measure Sets** template, each entry will be a collection or set of measures.

Quick Start

Below you will find the mandatory fields for the **measures** Table. You will also find definitions of terms and examples of entries of data, with all of the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed description.

Measures Template

1. Mandatory fields

- **Report ID**: Unique identifier for a Measure.
- **Sample ID**: Unique identifier for a Sample.
- **Site ID**: Unique identifier for the location where a Sample was taken.
- **Analysis Date End**: Date the Measure was completed.
- **Measure**: A measurement or observation of any substance including a biological, physical or chemical substances.
- **Value**: Value of the Measure.
- **Unit**: Units of a value.
- **Aggregation**: Statistical measures that the Measure represents.
- **Specimen**: Substance or thing upon which the observation was made. Specimens include water, air, and surface. The Specimen feild is mandatory if there is more than one specimen type that is collected in database.

Detailed description

2. Examples

Image will be added

Measure Sets Template

1. Mandatory fields

- **Report Set ID:** Unique identifier that links together a group of related Measures.

3. Examples

Image will be added.

Detailed description

Measures template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Report ID	Protocol ID	Sample ID	Purpose	Polygon ID	Site ID	Dataset ID	Report Set ID
2								
3								

i) Columns A-C, and E-H

- These are the identifier fields.

How to use the ODM

- **Report ID** (column A) is the unique identifier for this template, and it cannot be repeated. Each Report ID represents a unique measure.

The **Report ID** can be any text field up to 30 characters. If you don't have a method to generate a **Report ID** you can use the ODM recommended formula: {Mathew, can you add the formula here}.

- Other identifiers can be repeated between entries (if needed). For instance, if you enter two different measures from the same sample, then the **Sample ID** (column C) will be the same.
- You may have already created these identifiers in another template. For instance, you may have created **Sample ID** in the **samples** Template.

ii) Column D

- **Purpose** (optional) has a drop-down menu list. In the list, **Regular** is the default purpose.

2. Columns I to M

	I	J	K	L	M
1	Analysis date start	Analysis date end	Report Date	Specimen	Fraction analyzed
2					
3					

i) Columns I to K

- These are the date fields where you can enter the dates associated with your measure in the format yyyy-mm-dd. Including hours and minutes is optional, but all date-entries use the ISO 8601 format. 2022-01-01T06:11:54 or 2022-01-01T06:11:54+13:30

ii) Columns L and M

- **Specimen** and **Fraction analyzed** have drop-down menus. Both fields are optional, but both are recommended.
- **Specimen** is must be included your table has more than one type of specimen. When there is only one specimen, the database table record the specimen type. For example, if your **measures** table include only wastewater sample measures, you can indicate this information in the **databases** table.
- **Fraction analyzed** applies only for water and wastewater samples. You should record the fraction for all samples if they are component type water or wastewater.

3. Columns N to S

	N	O	P	Q	R	S
1	Group	Class	Measure	Value	Unit	Aggregation
2						
3						

i) Columns N and O

- **Group** and **Class** help organize **Measures** to make a specific measure easy to find. The fields have drop-down menus. Both fields are optional and used only if you need help finding a **Measure**.
- For example, to find the allele measures for SARS-CoV-2 you can start by selecting **SARS-CoV-2** in the **Group** column. Next, you select **Allele** in the **Class** column. Then, you'll see the drop-down in the **Measure** column that contains SARS-CoV-2 allele regions.
- Leave these fields **blank** if you do not use them.

ii) Column P

- **Measure** is where you select what is being measured. A measure in ODM is a 'measurement or observation of any

How to use the ODM

N	O	P	V
Group	Class	Measure	
SARS-CoV-2	Alleles class	<div> <div></div> <div> BRSV-N SARS-CoV-2-E SARS-CoV-2-N1 SARS-CoV-2-N2 SARS-CoV-2-N3 SARS-CoV-2-RdRp crAssphage-N PMMoV-CP S sars-cov-2 gene target Influenza virus A1 Influenza virus A2 Influenza virus B </div> </div>	

substance including a biological, physical or chemical substance.’ **Measure** has a conditional drop-down menu. See above section on **Group** and **Class**.

- If you chose to enter in a **Class** (column O), the drop down menu will contain the **Measures** found in that **Class**. If you did not enter a value in column O, the menu will contain all possible **Measures**.
- iii) Column Q - This where you will enter in the value of your **Measure**. - For instance, if your **Measure** is a temperature of 20 °C, you would enter “20”. - Note: Excel accepts any data type can be entered in the **Measure** column. However, each **Measure** has a specific data type. Data types can be found in the reference documentation for each measure. For example, the data type for temperature can be found [here](#).
- iv) Column R
- **Units** is where you enter the unit for your measure. There is a conditional drop-down list of the units that are valid for the measure in column P. For example, if the measure is [Storate temperature] (.../parts.qmd#sTemp), the only unit displayed is Celcsius.

v) Column S

- **Aggregations** is where you enter the aggregation of your value. There is a conditional drop-down menu field which contains the types of aggregations of your value. For example does your value represent a mean, median, etc.
- What appears in this menu list depends on what you entered in the Unit field (column R).

4. Columns T to AA

	T	U	V	W	X	Y	Z	AA
1	Nomenclature	Index	Measure License	Reportable data flag	Organization ID	Contact ID	Update Date - Last Edited	Notes
2								
3								

i) Columns U, V and AA

- **Index** (column U) is used if you have multiple entries with the exact same values in the other fields. If this is the case you could use “1”, “2”, etc. to differentiate them.
- **Measure license** (column V) refers to the access and use licensing of the measure that you are entering.

ii) Column T and W

- **Reportable** (column W) is where you can indicate if there are quality concerns with the measure and the measure should not be used for regular reporting. More details of the quality concerns can be recorded in the **quality sets** table. There is a drop down list this field.

iii) Column X and Y

- These are identifier fields that are used to indicate the **Organization** and **Contact** associated with the entry.

iv) Column Z

- **Last edited** field is where you indicate the date when the entry was updated. The field is used if you need to modify an entry (row) after your initial recording. You can leave

How to use the ODM

this field blank if the measure was entered with no updates.
Meaning, use this field only when you modify the entry (row).

Measure Sets Template

1. Columns A to G

i NOTE: this figure needs to be update with the correct labels.

	A	B	C	D	E	F
1	Report Set ID	Protocol ID	Name	Contact ID	Organization ID	Update - Last E
2						
3						

i) Columns A, B, D and E

- These are all Identifier fields.
- Report Set ID (column A) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Report Set ID represents a unique Measure Set.
- For the rest, you can repeat the Identifier between entries (if needed). For instance, if you are entering two different Measure Sets that came from the same Organization, then the Organization ID (column D) will be the same.
- You may have already created these Identifiers in another template. For instance you may have created Organization ID in the Organizations Template.

ii) Columns C and G

3) How to report Samples and Sample Relationships:

- These are Free Form fields in which you will enter the indicated information. Name (column C) refers to the name that you have given to the Measure Set.

iii) Column F

- This is a Date field in which you can enter in the date when the entry was updated.

You have now entered your data in the Measures and Measure Sets Template, congratulations!

3) How to report Samples and Sample Relationships:

In this guide you will learn how to enter Samples and Sample Relationships into their respective templates. In the Samples template, each entry represents a Sample. A Sample is the wastewater that you collected so that Measures can be made. Along the same lines, each entry in the Sample Relationships represents an interaction (or Relationship) between two Samples in the form, Subject - Relationship - Object. So, to specify that Sample A is a field sample replicate of Sample B, you would enter “Sample A” - “Field Sample Replicate” - “Sample B”

You can find additional information regarding terms and the colour-coding of columns at the beginning of this guide (1) [How to use templates](#)) and

Quick Start

Below you will find the fields that are mandatory for the Samples and Sample Relationships template. You will also find some examples of entries of data, with the all of the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed Description.

Samples Template

1. Mandatory fields

Sample ID	Site ID
Sample Material	Sample Collection Type
Collection Period	Collection Number
Collection Date Time	

2. Definitions

- Sample ID: Unique identifier for a Sample.
- Site ID: Unique identifier for the location where a wastewater sample was taken.
- Sample Material: Type of Sample.
- Sample Collection Type: Method used to collect the Sample.
- Collection Period: The time period over which the Sample was collected, in hours.
- Collection Number: The number of subsamples that were combined to create the Sample. Use NA for continuous, proportional or passive sampling.
- Collection Date Time: The date, time and time zone the sample was taken.

3. Examples

to be added

Sample Relationships Template

1. Mandatory fields

Sample ID Object	Relationship
Sample ID Subject	

2. Definitions

- Sample ID Object: The Object (or one of the samples) of a relationship between two Samples. This will always be a Sample ID that was previously created in the Samples Template.
- Relationship: Describes relationship between two Samples.
- Sample ID Subject: The Subject (or one of the samples) of a relationship between two Samples. This will always be a Sample ID that was previously created in the Samples Template.

3. Examples

Examples will be added

Detailed Description

Samples Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Sample ID	Protocol ID	Organization ID	Contact ID	Site ID	Purpose	Sample material
2							
3							

How to use the ODM

i) Columns A to E:

- These are Identifier fields.
- Sample ID (column A) is the Unique Identifier for this template and cannot be repeated between entries. You can think of each Sample ID value as representing a unique Sample.
- For the rest of the fields, you can repeat values between entries. For instance, if you are entering two different Samples from the same Site, then the Site ID (column E) will be the same.
- You may have already created these Identifiers in another template. For instance you may have created Site ID in the Sites Template.

ii) Columns F to G:

- These are Drop-down Menu fields where you can enter information regarding the Purpose and Sample Material of your Sample.

2. Columns H to N

	H	I	J	K	L	M	N
1	Dataset ID	Sample origin	Replicate Type	Sample collection type	Collection Period	Collection Number	Pooled
2							
3							

i) Columns H:

- These is an Identifier field for the Dataset that your Sample is associated with.

ii) Columns I to K, N:

- These are Drop-down Menu fields related to various characteristics of your Sample.

Detailed Description

- Sample Collection Type (column K) refers to the collection technique you used to obtain the Sample.
 - Pooled (column N) refers to whether or not the Sample that you are entering is made up of multiple child Samples.
- iii) Columns L and M:
- These are Free Form fields.
 - Collection Period (column L) refers to the number of hours that it took to collect your Sample.
 - Collection Number (column M) refers to the number of subsamples that were used to create Sample that you are entering.

3. Columns O to V

	O	P	Q	R	S	T	U	V
1	Collection date time	Collection date time start	Collection date time end	Date sample was sent	Date sample recieved	Reportable data flag	Update Date - Last Edited	Notes
2								
3								

- i) Columns O to S, U:
- These are Date fields that are related to various dates associated with your Sample. These dates will be in the format dd/mm/yyyy.
- ii) Column T:
- This is a Drop-down Menu field (Reportable Data Flag) where you can indicate whether or not your entry is of a good enough quality that it can be used by others.
- iii) Column V:
- This is a Free Form field where you can indicate anything of interest..

Sample Relationships Template

1. Columns A to E

How to use the ODM

	A	B	C	D	E
1	Sample ID Subject	Relationship	Sample ID Object	Update Date - Last Edited	Notes
2					
3					

i) Columns A and C:

- These are Identifier fields.
- Both of these Identifiers are Sample IDs and represent Samples. You would have created them previously in the Samples Template. Note that neither of them are Unique Identifiers and, thus, can be repeated between entries.

ii) Column B:

- This is a Drop-down Menu field, where you can select the type of relationship between two Samples. For instance, if Sample B was a child of Sample A, you would put the Sample ID of Sample A in the Sample ID Subject field (column A), and the Sample ID of Sample B in Sample ID Object field (column C). Then you would chose “Child relationship” from the menu in this column.

iii) Columns D and E:

- These are a Date field (column D) and a Free Form field (column E) where you can enter in the indicated information.

You have now entered your data in the Samples Template, congratulations!

4) How to report Protocols, Protocol Steps and Protocol Relationships:

4) How to report Protocols, Protocol Steps and Protocol Relationships:

In this guide you will learn how to enter Protocols, Protocol Steps and Protocol Relationships into their respective templates. As described in , a Protocol is “A procedure for collecting a sample or performing a measure”. Each entry in the Protocols Template represents a unique Protocol. A Protocol is made up of Protocol Steps. In the Protocol Steps Template, each entry is one of these steps. Finally, Protocols and Protocols can be linked to each other. Each entry in the Protocol Relationships Template represents one of these relationships in the form, Subject - Relationship - Object. So, for example, to specify that Protocol Step A needs to be done before Protocol Step B, you would enter “Protocol Step A” - “is Before” - “Protocol Step B”. .

You can find additional information regarding terms and the colour-coding of columns at the beginning of this guide ([1\) How to use templates](#)) and

Quick Start

Below you will find the fields that are mandatory for the Protocols, Protocol Steps and Protocol Relationships templates. You will also find some examples of entries of data, with all of the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed Description.

Protocol Steps Template

1. Mandatory fields

Protocol Step ID	Method (if entry is a Method)
Measure (if entry is a Measure)	Value (if entry is a Measure)
Aggregation (if entry is a Measure)	Unit (if entry is a Measure)

2. Definitions

- Protocol Step ID: Unique identifier for a Protocol Step.
- Measure: A measurement or observation of any substance including a biological, physical or chemical substance.
- Method: A procedure for collecting a sample or performing a Measure.
- Value: Value of the entry.
- Aggregation: Statistical measures used to report a Measure.
- Unit: The units of a Measure.

3. Examples

Image will be added

Protocols Template

1. Mandatory fields

Protocol ID

2. Definitions

- Protocol ID: A unique identifier for a given Protocol.

3. Examples

Image will be added

Protocol Relationships Template

1. Mandatory fields

Protocol ID Container	Relationship
Protocol ID Subject (if subject is a Protocol)	Protocol Step ID Subject (if subject is a Protocol Step)
Protocol ID Object (if object is a Protocol)	Protocol Step ID Object (if object is a Protocol Step)

2. Definitions

- Protocol ID Container: Unique identifier for a Protocol and the Steps and other Protocols that make it up.
- Relationship: Describes relationship between the Subject and Object.
- Protocol ID Subject: The Subject of a Relationship between a Protocol and a Protocol Step or Protocol.
- Protocol Step ID Subject: The Subject of a Relationship between a Protocol Step and a Protocol Step or Protocol.
- Protocol ID Object: The Object of a Relationship between a Protocol and a Protocol Step or Protocol.
- Protocol ID Step Object: The Object of a Relationship between a Protocol Step and a Protocol Step or Protocol.

3. Examples

Image will be added

Detailed Description

Protocol Steps Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Protocol Step ID	Method	Measure	Summary	Protocol step source	Protocol Step Version	Reference link
2							
3							

i) Column A:

- Protocol Step ID (column A) is the Unique Identifier field for this template and cannot be repeated between entries. You can think of each Protocol Step ID value as representing a unique Protocol Step.

ii) Columns B to C:

- These are Drop-down Menu fields.
- You only need to enter a value into Method ID (column B) or Measure ID (column C). The former is for when you are entering a Method as a Protocol Step. The latter is when you are entering a Measure as a Protocol Step. Leave the other field **blank**.

iii) Columns D to G:

- These are Free Form fields and you can enter the relevant information in any format that you wish.
- Summary (column D) is a short description of the Protocol Step that you are entering.
- Protocol Step Source (column E) specifies the Protocol Step which serves as a basis for the given Protocol Step.
- Protocol Step Version (column F) specifies the version of a given Protocol Step.

Detailed Description

	H	I	J	K	L	M	N	O
1	Organization ID	Contact ID	Instrument ID	Value	Unit	Aggregation	Update Date - Last Edited	Notes
2								
3								

2. Columns H to O

i) Columns H to J:

- These are Identifier fields.
- You can repeat values between entries. For instance, if you are entering two different Protocol Steps from the same Organization, then the Organization field (column H) will be the same.

ii) Column K:

- A Free Form field in which you will enter the Value of a Measure. If a Method is being entered, then you can put “NA” here.

iii) Columns L and M:

- These are Conditional Drop-down Menu fields.
- Unit (column L) depends on what you entered in Measure ID (column C). If you entered a Method, then you can put “NA” in this field.
- Aggregation (column M) depends on what you entered in Unit (column L). If you entered a Method, then you can put “NA” in this field.

iv) Columns N and O:

- These are Free Form fields where you can enter the indicated information.

Protocols Template

1. Columns A to F

How to use the ODM

	A	B	C	D	E	F
1	Source Protocol	Protocol ID	Dataset ID	Name	Summary	Reference link
2						
3						

i) Columns B and C:

- These are Identifier fields.
- Protocol ID (column B) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Protocol ID represents a unique Protocol.
- For Dataset ID, you can repeat between entries (if needed). For instance, if you are entering two different Protocols from the same Dataset, then the Dataset ID (column C) will be the same.
- You may have already created the Dataset IDs in the Datasets Template.

ii) Columns A, D to F:

- These are Free Form fields in which you can enter the indicated information.
- Source Protocol (column A) is where you can indicate the Protocol that served as the basis for the Protocol that is currently being entered.
- Name (column D) refers to the Name that you have decided to give your Protocol.

2. Columns G to K

i) Columns G, H:

- These are Identifier fields. You can repeat values between entries.

ii) Columns I, K:

Detailed Description

G	H	I	J	K
Organization ID	Contact ID	Protocol Version	Update Date - Last Edited	Notes

- These are Free Form fields and you can enter the relevant information in any format that you wish.
- You can use the Protocol Version field (column I) to indicate the version of the Protocol that you are entering.

iii) Columns J:

- This is a Date field where you can indicate when the entry was updated.

Protocol Relationships Template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Protocol ID Container	Protocol ID Object	Step ID Object	Relationship	Protocol ID Subject	Step ID Subject	Update Date - Last Edited	Notes
2								
3								

i) Columns A to C, E and F:

- These are Identifier fields.
- All of these Identifiers are either Protocol IDs or Protocol Step IDs and represent Protocols and Protocol Steps, respectively. This includes Protocol ID Container (column A), which represents the Protocol that contains the relationships that you are entering in this template. Note that none of these fields are Unique Identifiers and, thus, can be repeated between entries.

How to use the ODM

- You should only enter in a value for one of the two “Object” columns (column B and C); this depends on whether the Object of the relationship that you are entering is a Protocol or a Protocol Step. The same is true for the two “Subject” columns (column E and F).
- An example will be given in point (ii) below.

ii) Columns D:

- This is a Drop-down Menu field, where you can select the type of Relationship that is being entered.
- For example, say you wanted to enter in a Relationship stating that Protocol Step A takes place before Protocol Step B in a particular Protocol Container. You would enter the relevant Identifiers in the Step ID Object field (column B) and the Step ID Subject field (column F). You would then select “is Before” from the Relationship field. Finally, to identify this relationship, you would enter in the Identifier in the Protocol ID Container field (column A)

iii) Columns G and H:

- These are a Date field (column D) and a Free Form field (column E) where you can enter in the indicated information.

You have now entered your data in the Protocol Steps, Protocols and Protocol Relationships Template, congratulations!

5) How to report Metadata:

In this guide you will learn how to enter Metadata through a series of templates. Each entry in a template will usually represent what is found in the title of the template. For instance, in the Organizations template

each row will represent an Organization. In the context of the PHES-ODM, metadata is anything that gives information that is not found in the templates related to Measures, Protocols and Samples.

You can find additional information regarding terms and the colour-coding of columns at the beginning of this guide (1) [How to use templates](#)) and

Quick Start

Below you will find the fields that are mandatory for the metadata-related templates. You will also find some examples of entries of data, with the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed Description.

Addresses Template

1. Mandatory fields

Address ID	Dataset ID
Address Line 1	City
State, Province or Region	Country

- Address ID: A unique identifier for an Address.
- Dataset ID: The name of the Dataset that stores information for Measures, Samples and other reporting tables.
- Address Line 1: Line 1 (the street name, number and direction) for a given Address.
- City: The city where a Site or Organization is located; part of the Address.

How to use the ODM

- State, Province or Region: The state, province, or region where a Site or Organization is located; part of the Address.
- Country: The country where a Site or Organization is located; part of the Address.

2. Examples

Images will be added

Contacts Template

1. Mandatory fields

Contact ID	Organization ID
Contact email	

- Contact ID: A unique identifier for a given contact person.
- Organization ID: A unique identifier for the Organization to which the reporter is affiliated.
- Email: Contact e-mail address, for the labs.

2. Examples

Images will be added

Datasets Template

1. Mandatory fields

Dataset ID	License
Data Custodian	

- Dataset ID: A unique identifier for the Dataset that stores information.
- License: The license of the Dataset.
- Data Custodian: The data custodian(s) of the database.

2. Examples

Images will be added

Instruments Template

1. Mandatory fields

Instrument ID	Dataset ID
Model	Instrument Type

- Instrument ID: A unique identifier for an Instrument.
- Dataset ID: The name of the Dataset that stores information for Measures, Samples and other reporting tables.
- Model: Model number or version of the Instrument.
- Type of Instrument used to perform the measurement.

2. Examples

Images will be added

Organizations Template

1. Mandatory fields

Organization ID	Address ID
Organization Type	

- Organization ID: A unique identifier for the Organization to which the reporter is affiliated.
- Address ID: A unique identifier for an Address.
- Organization Type: Specifies the type or purpose of a given Organization.

2. Examples

Images will be added

Polygons Template

1. Mandatory fields

Polygon ID	Type of Geography
Well-known text	European Petroleum Survey Group Coordinates

- Polygon ID: Unique identifier for the Polygon.
- Type of Geography: Type of geography that is represented by the Polygon.
- Well-known text: Well-known text of the Polygon.
- European Petroleum Survey Group Coordinates: The unique EPSG code specifying a given geospatial area.

2. Examples

Images will be added

Quality Reports Template

1. Mandatory fields

Quality Report ID	Report ID (if entry is related to a Measure)
Sample ID (if entry is related to a Sample) Quality Flag	Report Set ID (if entry is related to a Measure in a Measure Set)

- Quality Report ID: A unique identifier for a given Quality Report.
- Report ID: Unique identifier for a Measure.
- Sample ID: Unique identifier for a Sample.
- Report Set ID: Unique identifier that links together a group of related Measures.
- Quality Flag: A field for reporting any quality concerns - of lack thereof - for a Sample or Measure.

2. Examples

Images will be added

Sites Template

1. Mandatory fields

Site ID	Site Type
Sample shed	Contact ID
Latitude	Longitude
European Petroleum Survey Group Coordinates	

- Site ID: Unique identifier for the location where wastewater sample was taken.

How to use the ODM

- Site Type: Type of Site where sample was taken.
- Sample shed: A geographic area, physical space, or structure. A Sample is taken from a sample shed for a representative measurement of a substance(s).
- Contact ID: A unique identifier for a given contact person.
- Latitude: Latitude in decimal coordinates of the Site.
- Longitude: Longitude in decimal coordinates of the Site.
- European Petroleum Survey Group Coordinates: The unique EPSG code specifying a given geospatial area.

2. Examples

Images will be added

Detailed Description

Addresses Template

1. Columns A to E

	A	B	C	D	E
1	Address ID	Dataset ID	Address Line 1	Address Line 2	City
2					
3					

i) Columns A and B:

- These are Identifier fields.
- Address ID (column A) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Address ID represents a unique Address.

Detailed Description

- For Dataset ID (column B), you can repeat the value between entries (if needed). For instance, if you are entering two different Addresses from the same Dataset, then this column will be the same.
- You may have already created the Dataset ID in the Datasets Template.

ii) Columns C to E:

- These are Free Form fields in which you can enter various information regarding the Address.

2. Columns F to J

F	G	H	I	J
Province, or Region	Postal or Zip Code	Country	Update Date - Last Edited	Notes

i) Columns F to J:

- Most of these fields are Free Form in which you can enter various information regarding the Contact.
- Update Date - Last Edited (column I) is a date field.

Contacts Template

1. Columns A to E

i) Columns A to C:

- These are Identifier fields.

How to use the ODM

	A	B	C	D	E
1	Contact ID	Dataset ID	Organization ID	First name of contact	Last Name of Contact
2					
3					

- Contact ID (column A) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Contact ID represents a unique Contact.
- For the rest, you can repeat the Identifier between entries (if needed). For instance, if you are entering two different Contacts that came from the same Organization, then the Organization ID (column C) will be the same.
- You may have already created these Identifiers in another template. For instance you may have created Organization ID in the Organizations Template.

ii) Columns D and E:

- These are Free Form fields in which you can enter the name of the Contact.

2. Columns F to J

	F	G	H	I	J
1	Contact email	Contact phone	Role of Contact	Update Date - Last Edited	Notes
2					
3					

i) Columns F to J:

- Most of these fields are Free Form in which you can enter various information regarding the Contact.

- Update Date - Last Edited (column J) is a date field.

Datasets Template

1. Columns A to G

	A	B	C	D	E	F	G
	Parent		Dataset				Reference
1	Dataset ID	Dataset ID	Creation Date	Name	License	Description	link
2							
3							

i) Columns A and B:

- These are Identifier fields.
- Dataset ID (column B) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Dataset ID represents a unique Dataset.
- For Parent Dataset ID (column A), you can repeat the value between entries (if needed). For instance, if you are entering two different Datasets that came from the same Parent Dataset, then the Parent Dataset ID will be the same.

ii) Column C

- This is a Date field where you can enter the date that the Dataset was created.

iii) Columns D to G:

- These are Free Form fields in which you can enter various information related to the Dataset.

2. Columns H to N

How to use the ODM

	H	I	J	K	L	M	N
1	Language ID	Funder Contact ID	Custody Contact ID	Funding agency	Data custodian	Update Date - Last Edited	Notes
2							
3							

i) Columns H to L:

- These are Identifier and Free Form fields.
- For information regarding Funders, you can enter their ID and Name in Funder Contact ID and Funding Agency, respectively (columns H and J).
- For information regarding Data Custodians, you can enter their ID and name in Custody Contact ID and Data Custodian, respectively (columns I and K).

ii) Columns M and N:

- These are a Date field (column L) and a Free Form field (column M) where you can enter in the indicated information.

Instruments Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Instrument ID	Dataset ID	Name	Model	Manufacturer	Contact ID	Organization ID
2							
3							

i) Columns A, B, F and G:

- These are Identifier fields.

Detailed Description

- Instrument ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Instrument ID represents a unique Instrument.
- For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Instruments from the same Dataset, then this column will be the same.
- You may have already created some of the Identifiers in other Templates. For instance, you may have already created Dataset ID in the Datasets template.

ii) Columns C to E

- These are Free Form fields where you can enter the Name (column C), Model (column D) and Manufacturer (column E) of the Instrument.

2. Columns H to N

	H	I	J	K	L	M	N
1	Description	Reference link	Instrument Type	Describe other instrument type, if applicable	Index	Update Date - Last Edited	Notes
2							
3							

i) Columns H, I, L to N:

- Most of these fields are Free Form in which you can enter relevant information regarding your Instrument
- You can use Index (column L) to differentiate two entries of data that are the same in the other fields.
- Update Date - Last Edited (column M) is a date field.

ii) Columns J and K:

- For Instrument Type (column J), you can select from the drop-down menu the type of Instrument. If you do not see your Instrument, you can enter “Other instrument type”

How to use the ODM

- If you entered “Other instrument type” in the Instrument Type field, you can use the Describe Other Instrument type, If Applicable field (column K) to describe your Instrument.

Organizations Template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Organization ID	Dataset ID	Address ID	Organization Type	Name	Description	Update Date - Last Edited	Notes
2								
3								

i) Columns A and C:

- These are Identifier fields.
- Organization ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Organization ID represents a unique Organization.
- For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Organizations from the same Address, then this column will be the same.
- You may have already created some of the Identifiers in other Templates. For instance, you may have already created the Address ID in the Addresses Template.

ii) Columns D to H:

- For Organization Type (column D), you can select from the menu the type of Organization that you are entering.
- The rest of the columns contain Free Form or Date fields in which you can enter additional information regarding the Organization.

Polygons Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Polygon ID	Dataset ID	Name	Description	Polygon Population	Type of geography	European Petroleum Survey Group Coordinates
2							
3							

i) Columns A and B:

- These are Identifier fields.
- Polygon ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. A Polygon is an area on a map. Each Polygon ID represents a unique Polygon.
- For Dataset ID (column B), you can repeat the value between entries (if needed). For instance, if you are entering two different Polygons from the same Dataset, then this column will have the same value.
- You may have already created the Dataset ID in the Datasets Template.

ii) Columns C to E, G:

- These are Free Form fields in which you can enter additional information regarding the Polygon.
- Name (column C) is the name of the Polygon.
- European Petroleum Survey Group Coordinates (column G) is a coordinate system that can be used to identify where your Polygon is.

iii) Columns F:

- For Type of geography (column F) of the Polygon, you can select one of the values from the drop-down menu.

2. Columns H to N

How to use the ODM

	H	I	J	K	L	M	N
1	Well-known text	File Location of Polygon	Reference link	Organization ID	Contact ID	Update Date - Last Edited	Notes
2							
3							

i) Columns H to J and N:

- These are Free Form or Date fields in which you can enter relevant information regarding the Polygon.
- Well-known text (column H) refers to the text markup language that can be used to represent the Polygon (

ii) Columns K and L:

- These are Identifier fields in which the Organization ID and the Contact ID associated with the Polygon are entered.

Quality Reports Template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Quality Report ID	Report ID	Sample ID	Report Set ID	Quality Flag	Severity Indicator	Update Date - Last Edited	Notes
2								
3								

i) Columns A to D:

- These are Identifier fields.
- Quality Report ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Quality Report ID represents a unique quality issue that you want to report.
- For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Quality Reports from the same Sample, then the Sample IDs will be the same.

Detailed Description

- You may have already created some of the Identifiers in other Templates. For instance, you may have already created Sample ID in the Samples template.
- ii) Columns E and F:
- Quality Flag (column E) is a Drop-down Menu field where you can select the type of quality issue that you would like to enter.
 - You can indicate the severity of the Quality Flag in the Severity Indicator field (column F).
- iii) Columns G and H:
- These are a Free Form field and Date field in which you can enter the indicated information.

Sites Template

1. Columns A to F

	A	B	C	D	E	F
1	Parent Site ID	Site ID	Dataset ID	Polygon ID	Site Type	Sample shed
2						
3						

- i) Columns A to C, and F:
- These are Identifier fields.
 - Site ID (column B) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Sites ID represents a unique Site.
 - For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Sites that are found in the same Polygon, then the Polygon IDs will be the same.

How to use the ODM

- Parent Site ID (column A) refers to the Site that is parent to the Site that is being entered. For example if the Site that is being entered is a room in a LTC facility, that the Parent Site ID would refer to the LTC facility.

ii) Columns D to E:

- These are Drop-down Menu fields in which you can specify the Site Type (column D) and the type of Sample shed (column E).

2. Columns G to M

	G	H	I	J	K	L	M
1	Address ID	Organization ID	Contact ID	Name	Description	Department or Region	Health Planning Region
2							
3							

i) Columns G to I:

- These are Identifier fields for additional information that are associated with the Site.

ii) Columns J to M:

- These are all Free Form fields in which you can enter in relevant information about your Site.
- Public Health Department or Region (column L) and Health Region (column M), refers to the Department or Region where your Site is located.

3. Columns N to S

	N	O	P	Q	R	S
1	Population Served	Latitude	Longitude	European Petroleum Survey Group Coordinates	Update Date - Last Edited	Notes
2						
3						

i) Columns N to S:

- These are Free Form fields and a Date field in which you can enter in additional information of your Site.

Example data for each of the 11 report tables.

- European Petroleum Survey Group Coordinates (column Q) is a coordinate system that can be used to identify where your Site is.

You have now entered in your Metadata, congratulations!

Example data for each of the 11 report tables.

Part I.

Reference guide

The reference guide is a technical description of the ODM **parts**. There are three sections to the reference guide.

- **parts** – Parts are the smallest way of describing anything in the ODM. An analogy to parts of a car can be used to illustrate this point. If you own a car, your garage can access a parts list that contains every part of your car, right down to every nut and bolt. There are over 900 parts, all listed in alphabetical order. Each part has up to 60 different attributes such as the part name and description.
- **tables** - The ODM database model has 21 tables in its full relational database. The full model is commonly referred to as “long” tables as it stores data with one measurement per row. The most important tables are the ‘report’ tables that are used for daily reporting of measures and samples.
- **sets** - Sets are lists or collections of parts. You can think of sets as the elements of drop down lists or categories.

💡 Did you know that each part ID has its own URL?

Every part, table, and set has its own URL, that is created using the partID. partIDs are usually parentheses that follow a part, table, or set name.

For example, ‘airport’ is a sample shed category which would appear as: **Airport** (**airport**). The permanent link to the airport partID is, <https://docs.phes-odm.org/parts.html#airport>. Airport also is listed in the **Sampleshed** category set (**shedCatSet**).

Parts

ODM version 2.0.0-rc.3.16.

Parts are the smallest way of describing anything in the ODM. An analogy to parts of a car can be used to illustrate this point. If you own a car, your garage can access a parts list that contains every part of your car, right down to every nut and bolt. Similarly, the ODM has a **part ID** ([partID](#)) and **part description** ([partDescription](#)) for every measure, method, and attribute. There are also **part IDs** for **units**, **aggregations**, and other **parts**. The part list is presented below.

💡 Did you know you can identify part IDs by their specific naming format?

Each part ID starts with a lowercase character and contains no special characters. Parts have less than 12 characters, while part types such as units and aggregations have a maximum of seven characters. The unique construction of part IDs serves several purposes in addition to generating valid URLs. Users can create wide table variable names (headers) using part IDs that are valid in all open software languages. ODM supports a standard method of generating wide variable names by concatenating partIDs with an underscore, `_`. For this reason, some part types have a limit of seven characters to ensure the total length of wide names remains less than 35 characters, which is the character limit length of some statistical language software. The wide variable naming method will be added to the documentation.

Parts

a1306s delta-variant gene target

(a1306s). a1306s delta-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScaleSet *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a1918v delta-variant gene target

(a1918v). a1918v delta-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a2710t omicron-variant gene target

(a2710t). a2710t omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a63t omicron-variant gene target

(a63t). a63t omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a67v omicron-variant gene target

(a67v). a67v omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Atline Analyzer

(aas). Atline analyzer with sampler. An atline analyzer with sampler. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Absolute humidity

(absHum). A measure of the total mass of water vapour present in the air per volume of air. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* humClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* absHumidUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Absolute humidity unit set

(absHumidUnitSet). Unit set for absolute humidity measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* absHumidUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Academic institution

(academ). The category of organization type used for academic institutions or research groups. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum*

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Number of active cases

(acti). Number of active cases. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Active

(active). Indicator that a part is in current use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Analysis date

(aDate). Date the measurement was performed in the lab. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date end

(aDateEnd). Date the measurement or analysis was completed. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date start

(aDateStart). Date the measurement or analysis was started. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Address Line 1

(addL1). Line 1 (the street name, number and direction) for a given address. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:*

naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Address Line 2

(addL2). Line 2 (the unit number) for a given address. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Address table

(addresses). The table that contains information about addresses. Addresses can be recorded for sites, organizations or contacts (individuals). The Sites, Organizations, and Contacts tables include a link to the addresses table through Address ID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: programDescr *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Addresses table column order

(addressesOrder). Specifies the order of the columns in the Addresses table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Address header req

(addressesRequired). Specifies the columns required in a Addresses table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Address ID

(addressID). A unique identifier for an address. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Administrative regions

(admRegLevel). Administrative regions *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: ISCU *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Air filter

(afu). Air filter as part of filtration or circulating unit. Typically the unit would have a fan or blower. *Part Type*: categories *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: airCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scale

(aggragationScale). A scale used for an aggregation. Only applicable for measures and units. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum*

Parts

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation

(aggregation). Statistical measures used to report a measure. Each aggregation has a corresponding value. The default aggregation is ‘single’, which corresponds to a single measurement. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Aggregations

(aggregations). Statistical measures used to report a measure. Each measure reported as a number should be reported with an aggregation. The default aggregation is ‘single’, which corresponds to a single measurement. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scales

(aggregationScales). The scale of an aggregation set. Aggregation scales include quantitative and qualitative. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation set

(aggregationSet). The aggregation set for a unit. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation sets

(aggregationSets). Sets of aggregations. Examples of aggregation sets include logarithm, linear and boolean. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0

Parts

Maximum Length: 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Air compartment

(air). A measure or observation made from a substance in the air. *Part Type:* compartmentSets *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* airCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Air compartment set

(airCompartmentSet). A compartment set for measures and methods in the air compartment. *Part Type:* compartmentSets *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* airCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Airplane

(airpln). Airplane *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* col-Grp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:*

NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Airport

(airport). Airport Sample Shed category type *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Airport lists worksheet

(airportLists). Parts lists used to generate airport and airplane data entry templates. *Part Type*: dictionarySupport *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Air and surface compartment set

(airSurfaceCompartmentSet). A compartment set for measures and methods in the air or surface compartments. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: airSurfaceCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*:

Parts

naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Environmental temperature

(airTemp). Environmental temperature. *Part Type*: measurements *Domain*: Phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: tempClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: temperatureUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Air and water compartment set

(airWaterCompartmentSet). A compartment set for measures and methods in the air or water compartments. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: airWaterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Alias deprecated

(aliasDep). ID of an assay that is the same or similar. A comma separated list. Deprecated as of version 2, please avoid use. Use notes instead.

Part Type: attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Alias ID deprecated

(aliasIDDep). Alias id ID of an assay that is the same or similar. a comma separated list. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

All domains

(allDo). Domain that specifies that it could apply to al domains; biological, chemical, and physical. *Part Type:* domains *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Alleles

(allele). Measures and methods related to alleles. *Part Type:* classes *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Allele unit set

(alleleUnitSet). Unit set for alleles. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* siSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* alleleUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Access to all org

(allOrgs). If this is 'no', this data will not be available to any partner organization. If missing, data will be available to the all organizations. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA

*Minimum Length: 0 Maximum Length: 1 Part Status: **deprecated**.* First released: 1.0.0. Last updated: 2.0.0.

Aloh3 precipitation

(aloh3). Aloh3 precipitation *Part Type: categories Domain: Che Speciment ID: naSpecimenSet Compartment Set: waterCompartmentSet Group: procGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Assay method id default deprecated

(amDefDep). Used as default when a new measurement is created for this lab. See ID in AssayMethod table. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: **deprecated**.* First released: 1.0.0. Last updated: 2.0.0.

Parts

Amicon ultrafiltration

(amiconUf). Amicon ultrafiltration *Part Type:* categories *Domain:* Chemistry *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Amicon filter, extract with MP96

(amiMP96). Nucleic acid extraction, usually used for the liquid fraction of a wastewater sample, using amicon filtration and using an MP96 for final extraction. *Part Type:* categories *Domain:* Bio *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Amplicon sequencing

(amp). Specifies the amplicon strategy for genetic sequencing *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:*

seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

AND Boolean aggregation

(andBoo). “AND” aggregation. If all values in the aggregation is “TRUE” then the AND aggregation is also “TRUE”. If all values in the aggregation is “FALSE” then the AND aggregation is also “TRUE”. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Any compartment set

(anyCompartmentSet). A compartment set for measures and methods in any compartment. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Any specimen set

(anySpecimenSet). A specimen set that inculdes any specimen. *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: anySpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass

Parts

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Area proportional sample

(areaPr). An area proportional sample. Used for surface testing. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Arrival temperature

(arTemp). The temperature of a sample upon arrival at the laboratory for analysis. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* tempClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Use of the articV3 primer

(articV3). Initial implementation of an ARTIC bioinformatics platform for nanopore sequencing of nCoV2019 novel coronavirus; Artic Foundation v3 primer. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Use of the articV4 primer

(articV4). Artic V4 sequencing primer for VOCs and VOIs. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Assay ID

(asID). Links with the AssayMethod used to perform the analysis. Use instrument.ID for measures that are not viral measures. Deprecated as of version 2, please avoid use. Refer to methodID instead. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Attributes

(attributes). An attribute of a domain, specimen, compartment, or table. Attributes are one of the three main components or ways to report environment surveillance data. The other two components are measures and methods. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Bacteria unit set

(bacteriaUnitSet). Unit set for bacteria-related measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* bacteriaUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous bacteria group

(bactMisc). A group of measures/methods related to miscellaneous bacteria. The miscallenous bacteria are often have only one measure or method.

When a bacteria has many measures, they will be given their own groupID in updated dictionary versions. *Part Type:* groups *Domain:* Bio *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* bact-Misc *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Bead beating, extract with MP96

(bbMP96). Nucleaic acid extraction, usually used for solid fraction of a wastewater sample, using bead beating and using an MP96 for final extraction. *Part Type:* categories *Domain:* Bio *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

bcov

(bcov). Measure of the amount of bovine coronavirus. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

bcov culture spike target

(bcovCul). Cultured bovine coronavirus is used as the recovery efficiency control target. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

bcov vaccine spike target

(bcovVac). The bovine coronavirus vaccine is used as the recovery efficiency control target. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Beef extract flocculation

(beeExtractFloc). Beef extract flocculation. *Part Type*: categories *Domain*: Che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Is Before

(before). Specifies that the object step or protocol occurs before the subject step or protocol in the overall protocol container. Use only for protocolRelationships table relationshipID. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Beta

(beta). B.1.351 *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* who *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Biologic

(bio). A living organism or biological substance. *Part Type:* domains *Domain:* Bio *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0

Parts

Maximum Length: 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Digital droplet emulsification PCR

(bioRadDdpcr). Describes a PCR analysis done using BioRad's digital droplet emulsification PCR technology. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Binary Large Object (BLOB) data type

(blob). The data type for Binary Large Object (BLOB) data. Only used for the dictionary entries of parts. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Breadth of coverage ($\geq 5\times$ depth)

(boc). Positions with read depth greater or equal to 5. Report as percentage of positions. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:*

naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: 100 *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

5-Day carbonaceous biochemical oxygen demand

(bod5c). The quantity of oxygen utilized for the biochemical degradation of organic matter under standard laboratory procedures in five (5) days in the presence of a nitrification inhibitor, expressed in milligrams per litre (mg/l). *Part Type*: measurements *Domain*: Bio *Speciment ID*: siSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: <https://doi.org/10.1002/wer.1541> *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

5-day total biochemical oxygen demand

(bod5t). 5 day total biochemical oxygen demand. *Part Type*: measurements *Domain*: Bio *Speciment ID*: siSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: <https://doi.org/10.1002/wer.1541> *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Boolean data type

(boolean). The data type for boolean/binary data. Encoded as 'TRUE' or 'FALSE' *Part Type*: dataTypes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Boolean aggregation set

(booleanAggrSet). Aggreagation set for boolean. *Part Type*: aggregationSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Boolean value set

(booleanSet). Set that contains the valid possible values for a boolean measure (TRUE or FALSE). *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: booleanSet *Unit Set*: NA *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Bovine respiratory syncytial virus group

(brsv). Bovine respiratory syncytial virus. *Part Type*: measurements *Domain*: Bio *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

brsv culture spike target

(brsvCul). Cultured bovine respiratory syncytial virus (BRSV) is used as the recovery efficiency control target. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

BRSV-N

(brsvN). bovine respiratory syncytial virus capsid protein gene region *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA

Parts

Reference Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

brsv vaccine spike target

(brsvVac). The bovine respiratory syncytial virus (BRSV) vaccine is used as the recovery efficiency control target. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Branch sewer pipeline

(bSwrPpl). Specifies a site type that is collection pipe that run lateral to other municipal sewer lines, allowing drainage into the main sewer. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Building cleanout

(buildCO). Specifies a site type that is a capped pipe that connects to a building's main sewer line. *Part Type:* categories *Domain:* naDomain

Speciment ID: naSpecimenSet *Compartment Set:* waterCompartmentSet
Group: colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology*
Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar
Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum*
Length: 30 *Part Status:* active. First released: 1.0.0. Last updated:
 2.0.0.

Capacity unit set

(capacityUnitSet). Unit set for capacity measurements. *Part Type:* unit-
 Set *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:*
 waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*
 naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*
 capacityUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Ref-*
erence Link: NA *Data Type:* categorical *Minimum Value:* NA *Maximum*
Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
 First released: 2.0.0. Last updated: 2.0.0.

Case report date

(caRepDate). Date that the numbers were reported publicly. Typically,
 reported data and this measure is most commonly reported and used. *Part*
Type: categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Com-*
partment Set: humanCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set*
ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Ref-*
erence Link: NA *Data Type:* datetime *Minimum Value:* NA *Maximum*
Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **de-**
preciated. First released: 1.0.0. Last updated: 2.0.0.

Parts

Categorical data type

(categorical). The data type for categorical data. Only used for the dictionary entries of parts. *Part Type*: dataTypes *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Categories

(categories). A discrete list of values that can be reported for a measure, method or attribute. A use-case example of categories would be how they are used to record the site where a measure is taken (partID = geoType). For the geoType attribute there are different type of sites where samples can be taken: a wastewater treatment plant, airplane holding tank, wastewater pumping station, etc. Each of these site types is recorded as a category that can be identified in categorySetID = sTypeCatSets *Part Type*: partType *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Long-term acute care hospital

(ccc). Acute care hospitals, or complex continuing care, that provide care for patients with average length of stay longer than 25 days. *Part Type*: cat-

egories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Child day care

(cdc). Child day care facility. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Degrees Celcius

(cel). Degrees Celsius. *Part Type:* units *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* -60 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Centrifugation

(cent). Describes solid separation from a wastewater sample via centrifugation. Likely connected to other method steps prior to analysis. *Part Type:* categories *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Centricon ultrafiltration

(centriconUf). Centricon ultrafiltration. *Part Type:* categories *Domain:* Che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ceres nanotrap

(ceres). Ceres nanotrap. *Part Type:* categories *Domain:* Che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

CFU per 100 ml

(cfu). Colony forming units per 100 ml of filtered sample. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Changes column

(changes). A column for recording the changes from a previous version. Use this column as a short-hand change log for dictionary development. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Chemical

(che). A chemical compound. *Part Type:* domains *Domain:* Che *Speciment ID:* anySpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Chemagic viral dna/rna 300 kit

(chemVir). Nucleic acid extraction performed using the chemagic viral dna/rna 300 kit. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Child relationship

(child). Indicated that this is a sample generated from (an)other sample(s) either because of pooling or sub-sampling. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

City

(city). The city where a site or organization is located; part of the address. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum*

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 1.0.0. Last updated: 2.0.0.

Clarified sample

(clari). Clarified sample. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Class

(class). A unique identifier for a class, which is akin to a sub-group; it's a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class. Currently, class is only used for biologics. For example, SARS-CoV-2 is a group of measures with the following classes of allele, variant, mutation, sequence, and protein. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Classes

(classes). A class is a collection of one or more related measures or methods within a group. A group can have one more more classes to describe different parts or the class. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Cloudy

(cloudy). Qualitative category for the weather measure, specifying an overcast day with no rain. *Part Type*: categories *Domain*: Phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Centimetres

(cm). Unit part for the SI unit of centimetres. *Part Type*: units *Domain*: Phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA

Maximum Length: NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Carbon Dioxide

(co2). A measure of an amount or concentraion of carbon dioxide. *Part Type:* measurements *Domain:* Che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* gasClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* dissGasUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Chemical Oxygen Demand

(cod). Chemical oxygen demand. *Part Type:* measurements *Domain:* Che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Collection group

(colGrp). A group of measurement-like attributes related to sample collection. *Part Type:* groups *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set*

Parts

ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Collection date time

(collIDT). For grab samples this is the date, time and timezone the sample was taken. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Collection date time end

(collDTEnd). For integrated time average samples this is the date, time and timezone the sample was finished being taken. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Collection date time start

(collDTStart). For integrated time averaged samples this is the date, time and timezone the sample was started being taken. *Part Type:* at-

tributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: datetime *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sample collection set

(collectSet). Methods for collection samples. Sample collection methods include water, air, and surface. See the attribute ‘Compartment set’ of the sample collection method. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: collectCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Collection number

(collNum). The number of subsamples that were combined to create the sample. Use NA for continuous, proportional or passive sampling. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Collection number and period

(collNumPer). Composite collection number.period. This a composite measure that combines ‘collection period’ (collPer) and ‘collection number’ (collNum). Collection number and period can be used to reduce the number of table headings during data collection. In data storage, ‘collNum’ and ‘collPeriod’ should be transformed to ‘collPer’ and ‘collNum’. For Grab sample, Surface swab, or Area proportional sample just enter 1

In case of Composite, Flow-proportional collection, the number of subsamples followed by a dot, followed by the period of the sampling, in hours. For example, for 4 composite subsamples covering a 8 hours period, the entry would be 4.2 For time-proportional 24 subsamples collected every hour, the entry would be 24.1 For volume-proportional, the entry for 24 subsamples collected over 24 hours should be 24.24, in this case the period part of the entry should be understood as the total collection time instead of the periodicity of the collection.

For a COSCA ball or Moore swab collecting for 24 hours, the entry would be 1.24 *Part Type:* attributes *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* 1:1 *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 6 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Collection period

(collPer). Collection period. The time period over which the sample was collected, in hours. Alternatively, use collectionStart and collectionEnd. Examples: 1 hour, 6 hours, 24 hours *Part Type:* attributes *Domain:*

naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample collection type

(collType). The type of collection. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: collectCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 12 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Co-located sample

(colocated). Second or multiples samples collected at same location but different time (water, air) or at a nearby location (soil, sediment). Sent blind to laboratory. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Column name depreciated

(columnNameDep). Name for the column Depreciated in ODM version 2. Look-up tables are now incorporated into parts and sets. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Composite sample - archival

(comp). A composite sample, usually generated by an autosampler. This is intended for use in updating and mapping archival data into the newest version of the ODM. Generally, autosamplers do time or flow-proportional sampling, and so those collection methods should be used instead. This is only for composite/autosampler samples where these additional details have not yet been specified. Use collectionPeriod to describe how many hours the sample was taken, and collectionNum to record the number of samples. *Part Type*: categories *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Composite grab sample of 3

(comp3). A grab-composite sample composed of 3 separate grab samples. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Composite grab sample of 3 depreciate

(comp3dep). A grab-composite sample composed of 3 separate grab samples. Deprecated in version 2, please do not use. *Part Type:* categories *Do-*
main: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* water-
CompartmentSet Group: colGrp *Class:* naClass *Nomenclature:* naNomen-
clature Ontology Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet
Aggregation Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:*
NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Max-*
imum Length: 30 *Part Status:* **deprecated**. First released: 1.0.0. Last
updated: 2.0.0.

Composite 3hr grab sample

(comp3h). A 3-hour composite with 3 grab samples each taken once per hour, generally manually performed. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCom-
partmentSet Group: colGrp *Class:* naClass *Nomenclature:* naNomencla-
ture Ontology Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet
Aggregation Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:*
varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0

Parts

Maximum Length: 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Composite 3hr grab sample depreciate

(comp3hdep). A 3-hour composite with 3 grab samples each taken once per hour, generally manually performed. Deprecated in version 2, please do not use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Composite 8hr grab sample depreciated

(comp8h). An 8-hour composite with 8 grab samples each taken once per hour, generally manually performed. Deprecated in version 2, please do not use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Composite 8hr grab sample deprecated

(comp8hDep). A 8-hour composite with 8 grab samples each taken once per hour, generally manually performed. *Part Type:* categories *Domain:*

naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Compartment

(compartment). The attribute identifying the substance from which where a sample was taken. For more information, see partID = compartment. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Compartments

(compartments). The substance from which a sample was taken. For example, wastewater has component of “water”. Compartments are attributes of measures, methods, units, and aggregations. *Part Type*: part-Type *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment sets

(compartmentSets). Sets of compartments. Compartment sets are used to identify when a measure can be recorded for more than one compartment. For example, SARS-CoV-2 can be measured in people (humans), water, surface, or air. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample concentrate

(conc). Sample concentrate. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Case report date

(conCase). Date that the numbers were reported publicly for a confirmed case. Typically, reported data and this measure is most commonly reported and used. *Part Type*: units *Domain*: Bio *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale

Quality Set: NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Is concurrent to

(concurrent). Specifies that the object step or protocol occurs at the same time as the subject step or protocol in the overall protocol container. Use only for protocolRelationships table relationshipID. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Water conductivity

(cond). Measurement of conductivity of sample or site. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* siteFeat *Class:* conductivityClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* conductivityUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Conductivity

(conductivity). Measures and methods related to conductivity. *Part Type:* classes *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* conductivityClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Conductivity unit set

(conductivityUnitSet). Unit set related to conductivity measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* conductivityUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Number of confirmed cases

(conf). Number of confirmed cases. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases episode date

(confEp). Episode date is the earliest of onset, test or reported date for a confirmed case. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Onset date

(confOn). Earliest that symptoms were reported for a confirmed case. This data is often not known and reported. In lieu, episode is used. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Contact ID

(contactID). A unique identifier for a given contact person. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Parts

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0.
Last updated: 2.0.0.

Contact name

(contactName). Contact person or group, for the lab. Depreciated in ODM version 2. see 'name' *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 100 Part Status: deprecated.* First released: 1.0.0. Last updated: 2.0.0.

Contact table

(contacts). The table that contains information about a contact; a person who is the contact of a site or laboratory. adapt from version 1 documentation *Part Type: tables Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: programDescr Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Contacts table column order

(contactsOrder). Specifies the order of the columns in a Contacts table. *Part Type: tableSupport Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class:*

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Contact header req

(contactsRequired). Specifies the columns required in a Contacts table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Contact table set

(contactTableSet). Tables are where measures, methods and attributes are recorded. Tables represent the main entities of the environmental surveillance. *Part Type:* dictSet *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Confirmed cases test date

(conTest). Date that the covid-19 test was performed for a confirmed case. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Cphd ID

(corFcil). Correctional facility *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

COSCa ball

(cosca). COSCa passive sampling device. Use collectionPeriod to describe how many hours the ball was used to collect the sample *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* <https://pubs.rsc.org/en/content/articlelanding/2021/ew/d1ew00207d> *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Country

(country). The country where a site or organization is located; part of the address. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Countries or sovereign states

(countryLevel). Countries or sovereign states *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: ISCU Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Districts, counties, regions

(countyLevel). Districts, counties, regions *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: ISCU Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation*

Parts

Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Covid-19

(cov). Covid-19 infection. Any form of Covid-19 human infection - including testing for Covid-19, symptomatic Covid-19, asymptomatic Covid-19, hospitalized Covid-19, long-Covid-19, etc. If needed, define the specific form Covid-19 with units. See populationUnits. *Part Type: measurements Domain: Bio Speciment ID: poSpecimenSet Compartment Set: humanCompartmentSet Group: sarsCov2 Class: disease Nomenclature: ICD Ontology Reference: NA Category Set ID: NA Unit Set: populationUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

SARS-CoV-2 measure

(cov2Me). Measure the amount of SARS-CoV-2 virus. This measure should only be used if there is no other SARS-CoV-2 measure. See groupID = sarsCov2. Use a note to describe the specific measure used. *Part Type: measurements Domain: Bio Speciment ID: anySpecimenSet Compartment Set: anyCompartmentSet Group: sarsCov2 Class: allele Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: geneticUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: seeUnitData Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

SARS-CoV-2-B.1.1.7

(covB117). Variant B.1.1.7 of the SARS-CoV-2 virus. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: pangolin *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-B.1.351

(covB135). Variant B.1.351 of the SARS-CoV-2 virus. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: pangolin *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-E

(covE). SARS-CoV-2 E gene. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N1

(covN1). SARS-CoV-2 nucleocapsid gene, allele 1. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: any-CompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N2

(covN2). SARS-CoV-2 nucleocapsid gene, allele 2. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: any-CompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N3

(covN3). SARS-CoV-2 nucleocapsid gene, allele 3. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: any-CompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-P.1

(covP1). Variant P.1 of the SARS-CoV-2 virus. *Part Type*: measurements *Domain*: Bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: pangolin *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-RdRp

(covRdrp). SARS-CoV-2 RdRp gene. *Part Type*: measurements *Domain*: Bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Covid-19 population measurement date

(cphDate). date of reporting for covid-19 measure. Depreciated in ODM version 2. The measure identification is now measureReportID with report dates as for any measure report. *Part Type*: attributes *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0

Parts

Maximum Length: 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Cphd ID

(cphid). Unique identifier for the table. Depreciated in ODM version 2. The measure identification is now measureReportID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

crAssphage-N

(cra). crAssphage virus capsid protein gene region *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Crosswal table set

(crosswalkTableSet). Tables used to translate to and from other environmental dictionaries and models. *Part Type:* dictSet *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet

Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Cycle threshold or quantification cycle (Ct or Cq)

(ct). Cycle thresholds in a PCR assay. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Cumulative count

(cuCo). Units for describing a population measure of a cumulative count of cases of a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Custody Contact ID

(custodyCont). A unique identifier for data custodians. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Data custodian ID

(custodyID). The data custodian of a database. Use Organization ID to populate this field, and the organizations table to describe contact information and other details for the data custodian. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d377y delta-variant gene target

(d377y). d377y delta-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference*

Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d3g omicron-variant gene target { #d3g }

(d3g). d3g omicron-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d63g delta-variant gene target

(d63g). d63g delta-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d796y omicron-variant gene target { #d796y }

(d796y). d796y omicron-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit*

Parts

Set: geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d950n delta-variant gene target

(d950n). d950n delta-variant gene target. *Part Type:* measurements *Domain:* Bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Daily count

(daiCo). Units for describing a population measure of a daily count of new cases of a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:* units *Domain:* Bio *Specimen ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dataset creation date

(datasetDate). Specifies the date a given dataset was created. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dataset ID

(datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables. Where possible, datasetID name should correspond the original data custodian responsible for generating environmental data. (suggestion for a default name convention.) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dataset table

(datasets). A report table for capturing details about data's parental data set and data custodians. Supplying attribution for data collectors. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

Parts

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Datasets table column order

(datasetsOrder). Specifies the order of the columns in a Datasets table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset header req

(datasetsRequired). Specifies the columns required in a Datasets table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Data types

(dataTypes). The data type for a part. Data types used in the ODM include: varchar, boolean, float, category, date, time, datetime, url, email. dataType corresponds the the entry or cell within a data table, most

commonly within a report table. If the data entry has a unit, then the `dataType` corresponds to the unit and the `dataType` refers to ‘unit’. If the data entry is a category, then the `dataType` refers to ‘category’. All categories are `varchar`. Otherwise the `dataType` is identified within the part entry. *TBA: dataType for dictionary. Part Type: partType Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 12 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Date

(date). Date date on which the `assayMethod` was created or updated (for version update). *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: datetime Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.*

Datetime data type

(datetime). The data type for date and time data. Only used for the dictionary entries of parts. *Part Type: dataTypes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum*

Parts

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ddcov_e sars-cov-2 gene target

(ddcovE). ddcov_e sars-cov-2 gene target. *Part Type:* measurements *Domain:* Bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ddcov_n sars-cov-2 gene target { #ddcovN }

(ddcovN). ddcov_n sars-cov-2 gene target. *Part Type:* measurements *Domain:* Bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Deaths

(death). Units for describing a population measure of patients who have died from a given cause. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:* units *Domain:* Bio *Specimen ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:*

naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

del143/145

(del143). 143 or 145 deletion omicron-variant gene target. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

del 157/158

(del157). 157 or 158 deletion delta-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

del2084/2084

(del2084). 2084 or 2084 deletion omicron-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

del212/212

(del212). 212 or 212 deletion omicron-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

del3674/3676

(del3674). 3674 or 3676 deletion omicron-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

del69/70

(del6970). 69 or 70 deletion omicron-variant gene target. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Delta

(delta). B.1.617.2 *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* who *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Deprecated

(deprecated). Indicator to say that a part is no longer in current use in the model. See partID = status *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

Parts

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Derived sample

(derived). Specifies a sample that is derived or made from another sample or material. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Description

(descr). A detailed description as an attribute for describing results or program elements. A description of the part that serves a clear presentation of the part to a wide audience including technical and not technical staff. A description should allow any person who generates environment surveillance data to know how to identify how to record their data elements in the ODM. partReference can be used to further describe a part. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 200 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Description of change

(descrChange). A description of change in a part from the previous version. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Desk or counter

(desk). Desk, table, countertop or other flat working surface. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dewatered solids

(deso). Dewatered solids. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Detected

(det). Substance detected or not detected. TRUE = detected, FALSE = not detected *Part Type*: units *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: booleanCatSet *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: boolean *Minimum Value*: 0 *Maximum Value*: 1 *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Access to details deprecated

(detailsDep). More details on the existing confidentiality requirements of this measurement. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: booleanSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: boolean *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 1 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Development

(development). Indicator that a part is under development use. See partID = status *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dictionary reference table set

(dictionaryRefTableSet). Reference or look-up tables. For example, the parts table describe all elements of the ODM, including tables, table headers, measures, methods, categories, and units. sets are collections of parts. For example, units can be grouped together in a unitSet. languages and translations support translations. *Part Type:* dictSet *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Additional dictionary sheets

(dictionarySupport). Additional sheets for the Excel version of the ODM dictionary. Sheets or tabs for the ODM Excel dictionary include: all ERD tables (look-up tables and data entry tables) and additional supportig tables. *Part Type:* partType *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Dictionary set type

(dictSet). Sets used to describe and group dictionary tables. *Part Type:* dictSet *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dilution factor

(dilFact). Specifies the extent to which a sample or aliquot was diluted prior to analysis. Dilution factor is reported as a unitless measure, where a value of 10 indicates a 10:1 dilution. *Part Type:* measurements *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* dillutionClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* unitlessUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Point dillutions

(dillute). Exact concentration or dillutions for generating a standard curve. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:*

NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Diseases (human)

(disease). Measure and methods related to disease or infection in humans. *Part Type*: classes *Domain*: Bio *Speciment ID*: poSpecimenSet *Compartment Set*: humanCompartmentSet *Group*: naGroup *Class*: disease *Nomenclature*: ICD *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Dissolved gas concentration unit set

(dissGasUnitSet). Unit set for carbon dioxide concentrations measurements in water or air. *Part Type*: unitSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: airWaterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: dissGasUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Depth of coverage

(doc). The sequencing read depth. Used to interpret the confidence in a presence or amount of a variant. *Part Type*: units *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*:

Parts

naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: quantAggScale Quality Set: NA Reference Link: NA Data Type: integer Minimum Value: 0 Maximum Value: 100 Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Domain

(domain). Domain is the highest level of describing of a measure. The domain ID that corresponds to a given part. Mostly applicable for measures, methods, units, and aggregations. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Domains

(domains). There are three domain types: biologic (i.e. Covid-19, chemical (i.e. nitrogen), physical measure (temperature). *Part Type: part-Type Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Higher education dormitory or residential building

(dorm). Higher education dormitory or residential building *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

e156g delta-variant gene target

(e156g). e156g delta-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

e484a omicron-variant gene target { #e484a }

(e484a). e484a omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Escherichia coli

(ecoli). Concentration of bacteria that are passed through the faecal excrement of humans, livestock and wildlife *Part Type*: measurements *Domain*: Bio *Speciment ID*: siSaSpecimenSet *Compartment Set*: surface-WaterCompartmentSet *Group*: bactMisc *Class*: bacteriaClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: bacteriaUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Education

(education). A measure or sample taken for education or training. Use this purpose, for example, when teaching how to use the PHES-ODM. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Efficiency

(efficient). The efficiency reported for a standard curve. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: standardCurve *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: NA *Maximum*

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 2.0.0. Last updated: 2.0.0.

Contact email

(email). Contact e-mail address, for the lab. *Part Type:* attributes
Domain: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:*
naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*
naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*
naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA
Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*
Length: 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0.
Last updated: 2.0.0.

European Nucleotide Association

(ena). ENA header *Part Type:* crosswalkTableSet *Domain:* naDomain
Speciment ID: naSpecimenSet *Compartment Set:* naCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology*
Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar
Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maxi-*
imum Length: 30 *Part Status:* active. First released: 2.0.0. Last updated:
2.0.0.

ENA - notes

(enaNotes). ENA notes *Part Type:* crosswalkTableSet *Domain:* naDomain
Speciment ID: naSpecimenSet *Compartment Set:* naCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology*
Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar

Parts

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Rainfall

(envRnF). Rainfall, i.e. amount of precipitation in the form of rain. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* precipClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* precipitationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Ground snow depth

(envSnwD). Total depth of snow on the ground. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* precipClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* precipitationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Snowfall

(envSnwF). Snowfall, i.e. amount of precipitation in the form of snow. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSpecimenSet

Compartment Set: anyCompartmentSet *Group:* siteFeat *Class:* precip-
Class Nomenclature: naNomenclature *Ontology Reference:* NA *Category*
Set ID: NA *Unit Set:* precipitationUnitSet *Aggregation Scale:* seeUnitAg-
gScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum*
Value: seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA
Maximum Length: NA *Part Status:* active. First released: 1.0.0. Last
 updated: 2.0.0.

Episode date

(epiDate). Episode date is the earliest of onset, test or reported date.
Part Type: categories *Domain:* naDomain *Speciment ID:* naSpecimenSet
Compartment Set: humanCompartmentSet *Group:* naGroup *Class:* na-
Class Nomenclature: naNomenclature *Ontology Reference:* NA *Category*
Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* datetime *Minimum Value:* NA *Maxi-*
imum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:*
deprecated. First released: 1.0.0. Last updated: 2.0.0.

ERD table set

(erdTableSet). All tables listed in the Entity Relationship Diagram. The
 full ODM model is commonly referred to as “long” tables as it stores
 data with one measurement per row. *Part Type:* dictSet *Domain:* naDo-
 main *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontol-*
ogy Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Min-*
imum Value: NA *Maximum Value:* NA *Minimum Length:* NA *Maximum*
Length: NA *Part Status:* active. First released: 2.0.0. Last updated:
 2.0.0.

Parts

Sarbecovirus-specific E sars-cov-2 gene target

(eSarbec). Sarbecovirus-specific E sars-cov-2 gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Estimated frequency of reads

(estFreqReads). Estimated frequency of reads in a sequencing assay *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGroup *Class*: sequence *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Estuary, natural water body

(estuary). Estuary, natural water body *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0

Maximum Length: 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Experiment Failed

(expFail). PCR experiment failed. No value reported. Value should be blank for a measure with this qualityFlag value. *Part Type:* qualityIndicators *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

4s method

(extract4s). Nucleic acid extraction performed using the 4s method. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* <https://www.protocols.io/view/v-4-direct-wastewater-rna-capture-and-purification-bpdfmi3n> *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nucleic acid extraction method

(extraction). Description of the xtraction method. Description of the method used to extract the sample *Part Type:* methods *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:*

Parts

procGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: extractCatSet Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 100 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.

Nucleic Acid Extraction set

(extractSet). set used for storing all the valid category values for the nucleic acid extraction method. *Part Type: mmaSets Domain: Bio Specimen ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: extractCatSet Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Extraction volume of sample

(exvol). Extraction volume of sample. Size of the sample that is analyzed. *Part Type: measurements Domain: Phy Specimen ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: procGrp Class: physical Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: volumeUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: float Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Fecal matter

(faeces). Fecal matter. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

FALSE

(false). Boolean data type = FALSE Use these values and their set for any boolean measure or attribute. Use only “FALSE” (case sensitive) *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 5 *Maximum Length:* 5 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Field sample

(field). Specifies a sample taken from the field; directly collected from an area for testing. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA

Parts

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Field sample replicate

(fieldReplicate). A sample divided into two or more homogeneous parts. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

File location of polygon

(fileLocation). The location of the file containing the geometry of the polygon. File path specified, or a URL *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* blob *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 65535 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Filtration

(filt). Describes solid separation from a wastewater sample via filtration. Proceeds further concentration or analysis of the liquid filtrate. *Part Type:* categories *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment*

Set: waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

First Nation

(fiNa). Used to categorize a sample that is a First Nation, or on reserve lands. Likely for internal use only, Indigenous data can and should not be shared without explicit consent of the nation in question. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

First name of contact

(firstName). Specifies the first name of a given contact. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

First released version

(firstReleased). The version in which a part was first released *Part Type*: partSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Foreign key

(fK). Foreign key for a table. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed

(flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution. The reported concentration estimate is the updated result after addressing inhibition. *Part Type*: qualityIndicators *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination

(flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. *Part Type: qualityIndicators Domain: Bio Specimen ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

FI - Inhibition present and unaddressed

(flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. *Part Type: qualityIndicators Domain: Bio Specimen ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

J - Weak signal extrapolation

(flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of the standard curve in the

Parts

non-linear region. *Part Type:* qualityIndicators *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ND - Non-detect

(flagND). No amplification occurred in the reaction; non-detect. For the value of a non-detected measure, report the actual value even if it is below the limit of detection. The flag here ensures that it's recorded as a non-detect regardless. In instances where the original data doesn't record a value, but only has the flag, please populate the value field with either a 0 if dealing with variant percentages, and a 1 for all other purposes to further indicate a null result. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

UJ - Trace signal extrapolation

(flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., "trace" signal observed); the reported value is based on the extrapolation of the standard curve in the non-linear region. *Part Type:* qualityIndicators *Domain:* Bio *Speciment ID:* saSpecimenSet

Compartment Set: anyCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
 First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable

(flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. *Part Type:* qualityIndicators *Domain:* Bio *Speciment ID:* saSpecimenSet
Compartment Set: anyCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA
Minimum Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Float data type

(float). The data type for float data. Only used for the dictionary entries of parts. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet
Compartment Set: anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA
Category Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA
Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Flow-normalized mean

(floMean). Mean measure normalized to wastewater flow. Mostly used for reporting specific alleles (N1, N2, E, etc.) *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Floor

(floor). Floor of a building or room. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: surfaceCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Flow rate

(floRate). Wastewater volumetric flow rate at the sample collection location over the 24-hr period during which the sample was collected. If only an instantaneous flow measurement is available, it may be reported in units of million gallons per day. *Part Type*: measurements *Domain*: Phy *Speciment ID*: siSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: siteFeat *Class*: flowClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: flowUnitSet *Aggregation*

Scale: seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Flow

(flow). Measures and methods related to flow. *Part Type:* classes *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* flowClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Flow proportional 24hr sample depreciated

(flow24hDep). A flow proportional 24-hour composite sample generally collected by an autosampler. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Flow proportional sample

(flowPr). A flow proportionalcomposite sample generally collected by an autosampler. Use collectionPeriod to describe how many hours the sam-

Parts

ple was taken. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Volume flow rate unit set

(flowUnitSet). Unit set for volume flow measurements. *Part Type:* unit-Set *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* flowUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Flow volume

(flowVol). Volume of influent. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* siteFeat *Class:* flowClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* flowUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus measure

(flu). General influenza virus measure This measure should only be used if there is no other influenza virus measure. See groupID = virusMisc. Use a note to describe the specific measure used. *Part Type*: measurements *Domain*: Bio *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: http://purl.obolibrary.org/obo/MONDO_0005812 *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus A1

(fluA1). Influenza virus A1 type *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

influenza virus A2

(fluA2). influenza virus A2 type *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*:

Parts

seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal
Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First
released: 2.0.0. Last updated: 2.0.0.

Influenza virus B

(fluB). Influenza virus B type *Part Type:* measurements *Domain:* Bio *Spec-*
iment ID: saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:*
virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Refer-*
ence: NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation*
Scale: seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:*
seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal
Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First
released: 2.0.0. Last updated: 2.0.0.

Fluidigm digital PCR

(fluidpcr). Describes a PCR analysis done using FluidIGM's digital PCR
technology. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpeci-
menSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:*
naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Cate-*
gory Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:*
NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Max-*
imum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:*
active. First released: 2.0.0. Last updated: 2.0.0.

Foggy

(foggy). Qualitative category for the weather measure, specifying a foggy
or hazy day. *Part Type:* categories *Domain:* Phy *Speciment ID:* siSpec-
imenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:*

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Fraction analyzed

(fraction). Fraction of the sample that is analyzed. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* fractionCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sample fraction set

(fractionSet). set for the fraction of the sample (solid, liquid, etc.). *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* fractionCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Freyja Script

(freyja). Freyja Script for sequencing. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet

Parts

Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* <https://github.com/andersen-lab/Freyja> *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

F-Specific RNA bacteriophages

(frna). A measure for amount of F-Specific RNA bacteriophages. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

F-Specific RNA bacteriophages, G2

(frnaG2). A measure for amount of G2 F-Specific RNA bacteriophages. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample frozen

(frozen). Sample was frozen before analysis or processing *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample temperature

(fst). Temperature that the sample is stored at while it is being sampled. This field is mainly relevant for composite samples which are either kept at ambient temperature or refrigerated while being sampled. *Part Type*: measurements *Domain*: Phy *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: tempClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: temperatureUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Full dictionary sheets set

(fullDictionarySheetSet). Worksheets in the full Excel dictionary. The full dictionary is 'ODM_full-dictionary.xlsx' *Part Type*: dictSet *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar

Parts

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Funder Contact ID

(funderCont). A unique identifier for funders. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Funding agency ID

(funderID). The funding agency of the dataset. Use Organization ID to populate this field, and the organizations table to describe contact information and other details for the funding agency. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

g215c delta-variant gene target

(g215c). g215c delta-variant gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartment

mentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

g339d omicron-variant gene target { #g339d }

(g339d). g339d omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

g496s omicron-variant gene target { #g496s }

(g496s). g496s omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

g662s delta-variant gene target

(g662s). g662s delta-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Gamma

(gam). P.1 *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per copy of crAssphage

(gcCrA). Gene or variant copies per copy of crAssphage. *Part Type*: units *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

gene copies per day per 100,000

(gcD100). The unit for measures reflecting the gene copies per day per 100,000 people in the population. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per gram solids

(gCGS). Gene or variant copies per gram solids. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per L

(gcL). Gene or variant copies per litre. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer

Parts

Minimum Value: 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per mL

(gcMl). Gene or variant copies per millilitre of solution. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per PMMoV copy

(gcPpmov). Gene or variant copies per copy of PMMoV. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Genetics unit set

(geneticUnitSet). Unit set for genetic-related measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

geneticUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

General missingness set

(genMissingnessSet). The general set for missingness values. *Part Type:* missingnessSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Generic quality flag set

(genQualitySet). A quality set to specify any generic quality concerns about a measure or sample. *Part Type:* qualityIndSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

European Petroleum Survey Group Coordinates

(geoEPSG). The unique EPSG code specifying a given geospatial area. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet

Parts

Compartment Set: naCompartmentSet *Group:* miscAttr *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 2.0.0. Last updated: 2.0.0.

Latitude

(geoLat). Geographical location, latitude in decimal coordinates, ie.: (45.424721) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* miscAttr *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* -90 *Maximum Value:* 90 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Longitude

(geoLong). Geographical location, longitude in decimal coordinates, ie.: (-75.695000) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* miscAttr *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* -180 *Maximum Value:* 180 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Type of geography

(geoType). Type of geography that is represented by the polygon. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* geoTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Geographic set

(geoTypeSet). set for different type of geographic components. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* geoTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Well-known text

(geoWKT). Well-known text of the polygon *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 63 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Gram per cubic metre

(gm3). Density unit. Used for absolute humidity and other measures. *Part Type*: units *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: airCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Geometric mean

(gmn). Geometric mean. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Government agency

(govt). The category of organization type used for government agencies, programs, or crown-owned bodies. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0

Maximum Length: 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Grams

(grams). A unit of mass or weight. *Part Type:* units *Domain:* Phy *Specimen ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Gravity settling

(graSet). Describes solid separation from a wastewater sample where the sample material is allowed to settle by gravity, and then separated. *Part Type:* categories *Domain:* allDo *Specimen ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Grab sample

(grb). A single large representative grab sample. If the sample was collected over a series of hours or is a composite sample, please use partID = comp *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:*

Parts

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gromstole 1.0 Script

(gromstole). Gromstole 1.0 Script for sequencing. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* <https://github.com/PoonLab/gromstole> *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Group

(group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations. A collection of related measures. For example, SARS-CoV-2 is a group. Within the SARS-CoV-2 group, there are measure classes tht include RNA alleles (N1, N2, E, etc.), mutations, (E484K0), an entire sequence, viral proteins, etc. Currently groups are used for measures only. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Groups

(groups). A collection of related measures. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Header

(header). Header for a table. Also known as a table variable or entity relationship ‘attribute’. Header is the top row or the variable name in a table. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Heat inactivated sars-cov-2 virus spike target

(heatInacSARS). Heat inactivated SARS-CoV-2 virus is used as the recovery efficiency control target. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*:

Parts

procGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

hep g armored rna

(hepGRna). Measure of the amount Hepatitis G Armored RNA. *Part Type: measurements Domain: Bio Speciment ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: virusMisc Class: allele Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: geneticUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: seeUnitData Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

hep g armored rna spike target

(hepGRnaMat). Hepatitis G armored RNA is used as the recovery efficiency control target. *Part Type: categories Domain: Bio Speciment ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: procGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

High severity

(high). Indicates a very sever quality issue, likely meaning the data should not be reported. *Part Type: categories Domain: naDomain Speciment ID:*

naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sewer Network Health Region

(hlthReg). Health region served by the sewer network *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Hand Measurement

(hma). Handheld measurement analyzer. A handheld measurement analyzer. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Hollow fiber dead end ultrafiltration

(hollowFiberUF). Hollow fiber dead end ultrafiltration *Part Type*: categories *Domain*: Che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

General Hospital Admissions

(hosa). Hospital admissions or patients newly admitted to hospital. *Part Type*: units *Domain*: Bio *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Hospital Census

(hosc). Hospital census or the number of people admitted with an ailment. *Part Type*: units *Domain*: Bio *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Hospital

(hosptl). Hospital *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Holding tank wastewater depreciated

(hoTaWa). Wastewater from a holding tank, such as from an airplane or ship Depreciated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Holding tank wastewater

(htSam). Wastewater sampled from a holding tank, such as from an airplane or ship *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA

Parts

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Holding tank

(htSite). Holding tank *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Human compartment

(hum). A measure or observation made about a human. *Part Type:* compartments *Domain:* Bio *Speciment ID:* anySpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Humidity

(human). Measures and methods related to humidity. *Part Type:* classes *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* airCompartmentSet *Group:* naGroup *Class:* humClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Human compartment set

(humanCompartmentSet). A compartment set for measures and methods in the human compartment. *Part Type:* compartmentSets *Domain:* Bio *Speciment ID:* anySpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

i1566v omicron-variant gene target { #i1566v }

(i1566v). i1566v omicron-variant gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

i3758v omicron-variant gene target { #i3758v }

(i3758v). i3758v omicron-variant gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:*

Parts

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

i82t delta-variant gene target

(i82t). i82t delta-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

International Classification of Diseases

(icd). Classification system for diseases in humans. *Part Type*: nomenclatures *Domain*: Bio *Speciment ID*: poSpecimenSet *Compartment Set*: humanCompartmentSet *Group*: naGroup *Class*: disease *Nomenclature*: ICD *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Intensive care unit patients

(icu). Units for describing a population measure of patients who are in intensive care due to a given cause. Unit for a population or disease

measure to explain population-level clinical effects. *Part Type*: units *Domain*: Bio *Speciment ID*: poSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Index

(index). Index number in case the measurement was taken multiple times. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 50 *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Inhibition measure

(inhibMe). Parameter to report whether or not inhibition was detected in the sample. Binary measure - detected or no. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: miscMeas *Class*: sequence *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: booleanUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Parts

Inhibition method

(inhibMeth). Description of the method used to evaluate molecular inhibition. Description of the inhibition parameters. *Part Type*: methods *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: measGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 100 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Innovaprep ultrafiltration

(innovaprepUF). Innovaprep ultrafiltration *Part Type*: categories *Domain*: Che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Input

(input). Input for a table. Indicates if a part can be used in a table. “Input” is used to indicate if a part can be used as an entry in a table (something with partType = table). For example, covN1 is a specific measure that can be entered in the measureID field of the measures table. Therefore, covN1 is has a value “input” in the measures column in the parts list. ODM users can generate their own custom template by modifying using the value “input” for only the parts that are applicable for their program.

Part Type: categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Instrument ID

(instrumentID). A unique identifier for an instrument. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Instrument table

(instruments). The table that contains information about instruments. adapt from version 1 documentation *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Instruments table column order

(instrumentsOrder). Specifies the order of the columns in the Instruments table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Instrutment header req

(instrumentsRequired). Specifies the columns required in the Instruments table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Other instrument

(instrumentTypeOther). Type of instrument other than those included in the PHES-ODM. An other type of measurement instrument. Add description to notes. See documentation for how to request new instruments added to the dictionary. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum*

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Instrument Type

(insType). Type of instrument used to perform the measurement. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* insTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Describe other instrument type, if applicable

(insTypeOth). Description of the instrument in case it is not listed in instrumentType. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Instrument set

(insTypeSet). List of instruments that are used for measures and methods *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* insTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality*

Parts

Set: NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Integer data type

(integer). The data type for integers. Only used for the dictionary entries of parts. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Intercept

(inter). Intercept value of the calibration curve. Used for storing calibration curve information, potentially also applicable for other curve types. *Part Type:* measurements *Domain:* allDo *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standard-Curve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Combined ip2 and ip4 sars-cov-2 gene target

(ip2ip4). Combined ip2 and ip4 sars-cov-2 gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:*

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-1

(iso6391). The first part of the ISO 639 series of international standards for language codes. Part 1 covers the registration of two-letter codes. There are 183 two-letter codes registered as of June 2021. The registered codes cover the world's major languages. Part of the metadata for a languageID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2B

(iso6392B). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'B' specifies the bibliographic code (B code). Part of the metadata for a languageID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA

Parts

Data Type: categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2T

(iso6392T). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The ‘T’ specifies the terminological code (T code). Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-3

(iso6393). A set of international standards that lists short codes for language names. ISO 639-3 extends the ISO 639-2 alpha-3 codes with an aim to cover all known natural languages. Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-6

(iso6396). A set of international standards that lists short codes for language names. ISO 639-6 builds off ISO639-3 with the use of four-letter codes, and allowing users to differentiate between variants of languages and language families, such as historical vs. revived versions of languages. Part of the metadata for a languageID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

k856r omicron-variant gene target { #k856r }

(k856r). k856r omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Kilogram per second

(kgS). Kilograms per second. *Part Type*: units *Domain*: Phy *Speciment ID*: siSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA

Parts

Category Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Kilolitres

(kl). Kilolitres of volume *Part Type:* units *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Litres

(l). Litres of volume *Part Type:* units *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

I452r delta-variant gene target

(I452r). I452r delta-variant gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet

Group: sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

1981f omicron-variant gene target { #1981f }

(1981f). 1981f omicron-variant gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Lab ID default deprecated

(labDefDep). Used as default when a new sample is created by this reporter. See ID in Lab table. Deprecated as of version 2, no defaults please specify. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Laboratory duplicate

(labDuplicate). Second (time or more) processing and analysis of sample. Usually for general chemistry or metals analyses. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Lab ID

(labID). Unique identifier for a laboratory. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Lagoon system

(lagoon). Logoon system for extensive wastewater treatment *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:*

NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Lake, natural water body

(lake). Lake, natural water body *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Lambda

(lambda). C.37 *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Language ID

(lang). Language code for translation purposes. Specifies the language for each translation, other than the default English. Follow the ISO639-3 codes for now. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA

Parts

Category Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Language family

(langFam). Specifies the language family of a given language for translation and language tracking purposes. Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Language name

(langName). Specifies the name of the language in roman alphabet characters for translation and language tracking purposes. Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Language Look-up table

(languages). Look up table for all languages, used to give structure to the translation table. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* lookup *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Language table column order

(languagesOrder). Specifies the order of the columns in the Languages table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Language header req

(languagesRequired). Required headers in the Languages table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Last edited

(lastEdited). The date the entry was last updated. Use lastEdited if an entry is updated. Leave lastEdited blank or 'NA' for the first entry. Updates can include additions or revisions of the entry for any reason. For example, a wastewater measure was repeated later with an improved method. You can revise the entry by updating the value and then adding the date the new measure was performed to 'lastEdited'. To ensure data provenance, the best practice in this example is to generate a new entry with the same measureID as the original measureID. The original and new measures are kept in the database with a date in the 'lastEdited' field for the new entry but not the original one. Some databases may choose a delete-and-replace approach where the original entry is deleted and replaced by a new entry. In this approach, the 'lastEdited' field indicates the delete and replace occurred. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Last name of contact

(lastName). Specifies the last name of a given contact. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Last updated version

(lastUpdated). The version in which the part was last updated. Any change to the part of list will result in a change in the lastUpdated field to the dictionary version where the update occurred. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

LaTeX expression

(latExp). LaTeX expression used to generate formulas, symbols, etc. Mainly relevant for units. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 20 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage

(IBC). The percentage of the genome covered by reads (the breadth) is low. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:*

Parts

NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory control sample duplicate

(lcsd). Known amounts of an analyte or representative compounds are added to a second “clean” matrix (lab water or clean sand) in laboratory. Duplicate of LCS *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage

(IDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Leaked sample

(leaked). Sample leaked, some volume and material was lost. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet

Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
 First released: 2.0.0. Last updated: 2.0.0.

Link

(li). Link to an external reference that describes the geometry of the polygon. Deprecated as of version 2, please avoid use. Use referenceLink instead. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

License

(license). The license of a dataset. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Linear scale aggregation set

(linearAggrSet). The aggregation set that contains all aggregations that exist on the linear scale. *Part Type:* aggregationSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Liquid fraction

(liq). Liquid fraction of a sample. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

List worksheet

(lists). Parts lists for template dropdowns and other documentation. Sheets or tabs for the ODM Excel dictionary include: all ERD tables (look-up tables and data entry tables) and additional supportig tables. *Part Type:* dictionarySupport *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality*

Set: NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

List set

(listSet). PartTypes that contain categories, sets, or lists. Used to generate documentation. *Part Type:* miscSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* airCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* windSpeedUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Access to local ha

(localHA). If this is 'no', the, data will not be available to local health authorities. If missing, data will be available to local health authorities. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Limit of detection (LOD)

(lod). Limit of detection. Limit of detection (LOD) for this method, if one exists. *Part Type:* measurements *Domain:* naDomain *Speciment ID:*

Parts

naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: pcrClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Limit of detection (LOD) - sewage

(lodSewa). Limit of detection for sewage samples. Limit of detection (LOD) for sewage *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: pcrClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Logarithmic scale aggregation set

(logAggrSet). The aggregation set that contains all aggregations that exist on the logarithmic scale (rather than linear or qualitative). *Part Type*: aggregationSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Dictionary tables

(lookup). Tables to describe and support the ODM dictionary. These tables hold information on sets, all parts, and support multiple languages and translations. The tables are updated by the ODM development team when a new ODM version is created. *Part Type*: classes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: lookup *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Limit of quantification (LOQ)

(loq). Loq Limit of quantification (LOQ) for this method if one exists. *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: pcrClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Low severity

(low). A marker for low severity *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar

Parts

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval

(lowerCI95). Specifies the the lower limit of a 95% confidence interval. Should typically be linked to an upper limit, and some indication of data spread (ie. standard deviation), and a mean/median measure. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Low-volume sample

(lowVol). Sample is low-volume, but was run regardless. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Long-term care facility

(ltcf). A residential healthcare facility that provides 24-medical care. These are also called skilled nursing facilities. *Part Type:* categories

Domain: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Long-term care - assisted living or retirement home

(ltcfAl). A residential facility that provides assistance with daily care but generally does not provide skilled nursing care. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other long-term care

(ltcfO). Other residential facilities that provide daily and/or medical care, but are not defined as nursing home/skilled nursing facilities or assisted living facilities. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Life technologies digital PCR

(ltDpcr). Describes a PCR analysis done using Life Technologies' digital PCR technology. *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Lysis buffer

(lysi). Lysis buffer. *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Cubic metres per day

(m3D). Cubic metres per day. *Part Type:* units *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per hour

(m3H). Cubic metres per hour. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per second

(m3S). Cubic metre per second. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Manufacturer

(manufacturer). Manufacturer of an instrument. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Mate pair layout

(mateLay). Specifies the mate pair layout for a sequencing method. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Maximum length

(maxLength). The maximum length of the value of a part or measure. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Maximum value part support

(maxValue). The maximum value of a part. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Maximum Value

(maxValueAgg). Highest value in a range of values for a measure. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Control-plate Flag - NTCs { #Maybe not something we're even reporting, so.. no include? }

(Maybe not something we're even reporting, so.. no include?). NA *Part Type:* NA *Domain:* NA *Speciment ID:* NA *Compartment Set:* NA *Group:* NA *Class:* NA *Nomenclature:* NA *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* NA *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Arithmetic mean

(me). Arithmetic mean. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Measurement group

(measGrp). A group of measures that cannot be otherwise categorized. *Part Type:* groups *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other Measure

(measOth). Other measure, not otherwise specified in measures. Add description to categoryOther. *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Measures quality set

(measQualitySet). A quality set for any measure, includes all the quality flag options for any measure. Used as an “any quality set” for measures, but excludes the sample quality flags to avoid confusion and data entry errors. *Part Type:* qualityIndSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA

Category Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measure

(measure). A measurement or observation of any substance including a biological, physical or chemical substance. An measurement is recorded from a specimen (i.e. a site, sample, person or population). Measures are organized into groups and classes (sub-groups). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Measures

(measurements). The attribute to describe a part type of measures. All measures have 'partType' = "measure". *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Report ID

(measureRepID). Unique identifier for a measurement. Report IDs cannot be repeated. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Measure report table

(measures). The table that contains information and details about a given measure adapt from version 1 documentation *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* results *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Measure set report table

(measureSets). The table that identifies sets of measures. Examples of measure sets include a set of replicates, dilutions (used to generate a Ct curve) or variants that are identified in a single sample. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* results *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measure sets table column order

(measureSetsOrder). Specifies the order of the columns in a Measure Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measure Set header req

(measureSetsRequired). Specifies the columns required in a Measure Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measures table column order

(measuresOrder). Specifies the order of the columns in a Measures table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:*

Parts

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measures header req

(measuresRequired). Specifies the columns required in a Measures table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Report set ID

(meaureSetRepID). Unique identifier that links together a group of related measures. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Median

(med). Median. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:*

naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized arithmetic mean

(menr). Arithmetic mean, normalized. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Other Measure Description

(meOthDe). Description of other measure (measOtherDesc). *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 60 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Measure license

(measureLic). Specifies the access and use licensing for a given single measurement. Populated by partID = licenseID *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Method

(method). A procedure for collecting a sample or performing a measure. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration method

(methodConc). Description of the method to concentrate a wastewater sample. Description of the method used to concentrate the sample *Part Type*: methods *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: methodConcCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*:

NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 100 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration method set

(methodConcSet). CatSeteogry set for concentration methods. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Methods

(methods). A procedure or step that is performed to make a measure. For example, methodExtract is a method tht describes the how a sample was extracted. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2

(mfAcidmgcl2). Membrane filtration with acidification and mgcl2 *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set*

Parts

ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2, membrane recombined with separated solids

(mfAcidmgcl2SS). Membrane filtration with acidification and mgcl2, membrane recombined with separated solids *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Membrane filtration with addition of mgcl2

(mfMgcl2). Membrane filtration with addition of mgcl2 *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Membrane filtration with mgcl2 and separated solids

(mfMgcl2SS). Membrane filtration with addition of mgcl2, membrane recombined with separated solids *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment

(mfNoAmend). Membrane filtration with no amendment *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment, membrane recombined with separated solids

(mfNoAmendSS). Membrane filtration with no amendment, membrane recombined with separated solids *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar

Parts

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with sample acidification

(mfSampleAcid). Membrane filtration with sample acidification *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with sample acidification, membrane recombined with separated solids

(mfSampleAcidSS). Membrane filtration with sample acidification, membrane recombined with separated solids *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Millions of gallons per day (MG/D)

(mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to mea-

sure flow. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre

(mgL). Milligrams per litre. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Murine Hepatitis Virus

(mhv). A measure for amount of murine hepatitis viurs. *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Multiple issues

(mI). Multiple issues have arisen in the sequencing process. When using the ‘multiple issues’ quality flag, please specify the issues and details in the notes section. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Mid-level severity

(mid). A marker for med-level severity, where there are some concerns. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Minimum length

(minLength). The maximum value of measure. *Part Type*: partSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: NA *Maximum Value*: NA *Minimum*

Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Minimum value part support

(minValue). The minimum value of part. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Minimum Value

(minValueAgg). Lowest value in a range of values for a measure. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Miscellaneous attribute group

(miscAttr). A group of miscellaneous measurement-like attributes. Examples of these would be longitude and latitude, or EPSG coordinates. *Part Type:* groups *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit*

Parts

Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous measure group

(miscMeas). A group of measures that cannot be otherwise categorized. *Part Type:* groups *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Missingness

(missingness). The part type for missingness values. Only used for the dictionary entries of parts. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Missingness sets

(missingnessSets). Missingness sets for measures, methods or attributes. Inputted data into an ODM field may have missing data. The missingness set lists valid missing categories for the measure, method, or attribute.

Part Type: partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Mixed/homogenized sample

(mix). Mixed or homogenized sample or fraction analyzed. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Millilitres

(ml). Millilitres of volume *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Megalitres per day (ML/d)

(mld). Megalitres per day. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres

(mm). Unit part for the SI unit of millimetres. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Measure, method, or attribute

(mmaSet). The set for a measure, method, or attribute. Only applicable for categorical fields, and the categories that populate them. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 12 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Measure, method, or attribute sets

(mmaSets). A set of categories. For example, site type has a list of different sites, such as a wastewater treatment plant, a septic tank, or a pumping station. Only used for the dictionary entries of parts. *Part Type:* part-Type *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Model

(model). Model number or version of the instrument. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Moore swab passive sample

(moorSw). Moore swab passive sample. Use collectionPeriod to describe how many hours the Moore swab was used to collect the sample *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference*

Parts

Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Metres per second

(ms). metres per second. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage

(ms2Col). Measure of the amount of ms2 coliphage. *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage spike target

(ms2ColMat). ms2 coliphage is used as the recovery efficiency control target. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:*

NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Matrix spike duplicate

(msd). A known amounts of an analyte or representative compounds are added in the laboratory to a second aliquot of the sample used for matrix spike. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Moore Swab Passive sample depreciated

(mspsDep). Moore swab passive sample. Deprecated in version 2, please do not use. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Major sewer pipeline

(mSwrPpl). Major sewer pipeline *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature

Parts

Ontology Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Mu

(mu). B.1.621 *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

murine coronavirus

(muCo). Measure of the amount of murine coronavirus. *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

murine coronavirus spike target

(muCoMat). Murine coronavirus is used as the recovery efficiency control target. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpec-

imenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple fraction

(mul). Multiple fractions were analyzed separately and aggregated together post-analysis. This fractionID is intended for use for archival data and should not be used for newly added new moving forward (2022). Please also specify in the notes section which fractions were used for the multiple fractions (ex. solid and liquid) to avoid loss of information. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple Purpose

(multiple). A measure or sample taken for multiple purposes, not easily captured by the other purpose categories. Useful for parent samples where various subsamples will have different purposes/ *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum*

Parts

Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Municipalities or communes

(municipalLevel). Municipalities or communes *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* ISCU *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Municipality

(municip). A complete municipality, this specifies an entire metropolitan area, either a city, town, etc. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Mutations

(mutation). Measures and methods related to mutations. *Part Type:* classes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

N sars-cov-2 gene target

(n). N sars-cov-2 gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Combined N1 and N2 sars-cov-2 gene target

(n1n2). Combined N1 and N2 sars-cov-2 gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n211i omicron-variant gene target { #n211i }

(n211i). n211i omicron-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*:

Parts

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n679k omicron-variant gene target { #n679k }

(n679k). n679k omicron-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n856k omicron-variant gene target { #n856k }

(n856k). n856k omicron-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n969k omicron-variant gene target { #n969k }

(n969k). n969k omicron-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*:

anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Not applicable

(NA). The field for which the expected value is not a property of the described object. Missing value indicator for missing data with no explanation. *Part Type*: missingness *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation not applicable

(naAggr). Not applicable for aggregations. Used for parts for which an aggregation value is not applicable or appropriate. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Aggregation scale not applicable

(naAggrScale). Not application for aggregation sets. Used for parts for which an aggregation scale value is not applicable or appropriate. *Part Type:* aggregationScales *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation set not applicable

(naAggrSet). Not application for aggregation sets. Used for parts for which an aggregation set value is not applicable or appropriate. *Part Type:* aggregationSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Class not applicable

(naClass). Class is not applicable. Use this for all parts that don't have classes (most non-measure or method parts) *Part Type:* classes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:*

varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment not applicable

(naCompartment). Compartment not applicable. *Part Type: compartments Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 13 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Compartment set not applicable

(naCompartmentSet). Compartment not applicable. *Part Type: compartmentSets Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 13 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Domain not applicable

(naDomain). Not applicable. Use ‘not applicable’ when there is no domain for the part *Part Type: domains Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA*

Parts

Category Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Group not applicable

(naGroup). Used when there is no group; ie. group is not applicable. Use 'not applicable' when there is no group for the part *Part Type:* groups *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Name

(name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity. Name of any entity. Use a prefix to distinguish names from different report tables. See fully specified name (FSN) list. i.e. Si-name = Site name, Po-name = Polygon name, Or-name = Organization name, Me-name = Method-name. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Not a number

(nan). The outcome of a measurement is not a valid number (plus or minus infinity, error, ...) Missing values indicator for when a numeric value is expected but not given. *Part Type*: missingness *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Nomenclature not applicable

(naNomenclature). Not applicable. Use 'not applicable' when there is no nomenclature for the part *Part Type*: nomenclatures *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen not applicable

(naSpecimen). Non applicable specimen. Use 'not applicable' when there is no specimen/set for the part *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*:

Parts

varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specimen set not applicable

(naSpecimenSet). A specimen set for when specimen/specimen set is not applicable. Use 'not applicable' when there is no specimen/set for the part *Part Type:* compartmentSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Native Name

(natName). The native name of the language, i.e. what the language is called by its speakers. Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Unit not applicable

(naUnit). Not appliable for units. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet

Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Unit set not applicable

(naUnitSet). Not applicable for unit sets. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

National Center for Biotechnology Information

(ncbi). NCBI header *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

NCBI - notes

(ncbiNotes). NCBI notes *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Neighbourhood

(neigh). A municipal neighbourhood, this specifies a sub-section of a larger municipality. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Local administrative units or neighborhoods

(neighborLevel). Local administrative units or neighborhoods *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* ISCU *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population of newly vaccinated persons

(newVax). Population of newly vaccinated persons *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

NextClade nomenclature

(nextclade). Specifies variant or genetic nomenclature as set out by NextClade. *Part Type*: nomenclatures *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Normalized geometric mean

(ngmn). Geometric mean, normalized. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Ammonium Nitrogen

(nH4N). Ammonium nitrogen concentration, as N. *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

niid_2019-ncov_n sars-cov-2 gene target

(niid19). niid_2019-ncov_n sars-cov-2 gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns

(noConcern). A flag to indicate there is are no quality concern about the measure or sample. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: 1 *Maximum Value*: 1 *Minimum Length*: NA *Maximum*

Length: NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sample not labelled

(noLabel). Sample had no label *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

No liquid concentration, liquid recombined with separated solids

(noliquid). No liquid concentration, liquid recombined with separated solids *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclature

(nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases. Currently only used for class = variant, mutation, or disease. Valid options are currently restricted to 'NextClade', 'Pangolin', 'WHO', or 'ICD'. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment*

Parts

Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclatures

(nomenclatures). A classification system to report the measure class. See partID = nomenclatureID. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

NORMAN

(norman). NORMAN header *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

NORMAN - notes

(normanNote). Norman notes *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Norovirus G1

(norog1). Norovirus genogroup 1 *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Norovirus G2

(norog2). Norovirus genogroup 2 *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Notes

(notes). A note used to describe details that are not captured in other attributes. There is no recommended structure of a note. Use notes as needed. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 200 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample missing time stamp

(noTime). Sample is missing the autosampler time; incomplete metadata on collection. *Part Type:* qualityIndicators *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Not reported

(nr). A value could have been recorded, however, it was not. Missing value indicator for data that is absent because it was not reported. *Part Type:* missingness *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:*

NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported, Not applicable missing set

(nrNAMissingnessSet). Not reported, Not applicable missing set. *Part Type*: missingnessSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sarbecovirus-specific N sars-cov-2 gene target

(nSarbec). Sarbecovirus-specific N sars-cov-2 gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Nephelometric turbidity unit

(ntu). Nephelometric Turbidity Units, a unit used in the measurement of turbidity (see “turbidity” under measureIDs). Nephelometric turbidity units (NTU) are based on white light (400–680nm) and 90° incident angle. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet

Parts

Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:*
NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference*
Link: <https://www.sciencedirect.com/topics/engineering/nephelometric-turbidity-unit> *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA
Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nuclisens automated magnetic bead extraction kit

(nucAuto). Nucleic acid extraction performed using the nuclisens automated magnetic bead extraction kit. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nuclisens manual magnetic bead extraction kit

(nucManu). Nucleic acid extraction performed using the nuclisens manual magnetic bead extraction kit. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Null

(null). A logical representation of a statement that is neither TRUE nor FALSE. Missing value indicator when the value is neither TRUE nor FALSE. *Part Type:* missingness *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

National Wastewater Surveillance System

(nwss). NWSS header *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

NWSS - notes

(nwssNotes). NWSS note *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

Parts

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

NWSS - process

(nwssProcess). NWSS process *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Water

(nww). Non-wastewater, coming from any kind of water body *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

coronavirus OC43

(oc43). Measure of the amount of coronavirus OC43. *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference*

Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

oc43 spike target

(oc43Mat). Human coronavirus OC43 is used as the recovery efficiency control target. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ocean, natural water body

(ocean). Ocean, natural water body *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Other Collection deprecated

(oColDep). Other type of collection method. Add description to collectionOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or

Parts

at phesd_odm@ohri.ca *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Lab Analysis

(ola). Offline laboratory analysis. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Omicron BA.1

(omicr1). Omicron B.1.1.529 *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* who *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2

(omicr2). Omicron BA.2 *Part Type*: measurements *Domain*: Bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2.75

(omicr275). Omicron BA.2.75 *Part Type*: measurements *Domain*: Bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.4

(omicr4). Omicron BA.4 *Part Type*: measurements *Domain*: Bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Omicron BA.5

(omicr5). Omicron BA.5 *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Onset date

(onsDate). Earliest that symptoms were reported for this case. This data is often not known and reported. In lieu, episode is used. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: datetime *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Online Sensor

(onse). Online sensor An online sensor *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0

Maximum Length: 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Ontology reference

(ontologyRef). Ontology reference for a part. This field contains a link to an existing ontology reference for the part. ODM content spans several existing ontologies. ODM does not strive to generate a new ontology, rather ODM seeks to harmonize and extend dictionaries for application in environmental and public health surveillance. *Part Type:* part-Support *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 200 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Orthophosphates

(ophos). Ortho-phosphate concentration. *Part Type:* measurements *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* water-CompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Other residential building

(orb). Individual residential buildings or institutions not captured in other categories *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

OR Boolean aggregation

(orBoo). “OR” aggregation. If any value in the aggregation is “TRUE” then the OR aggregation is also “TRUE”. If all values in the aggregation is “FALSE” then the OR aggregation is also “FALSE”. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ORF1a sars-cov-2 gene target

(orf1a). ORF1a sars-cov-2 gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal

Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ORF1ab sars-cov-2 gene target

(orf1ab). ORF1ab sars-cov-2 gene target *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ORF1b sars-cov-2 gene target

(orf1b). ORF1b sars-cov-2 gene target *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Organization ID

(organizationID). A unique identifier for the organization to which the reporter is affiliated. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:*

Parts

NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 100 *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Organization table

(organizations). The table that contains information about a laboratory. adapt from version 1 documentation *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: programDescr *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Organizations table column order

(organizationsOrder). Specifies the order of the columns in a Organizations table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization header req

(organizationsRequired). Specifies the columns required in a Organizations table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup

Class: naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Organization level

(orgLevel). The geographic level of an organization. There are six levels based on International Standard Classification of Administrative Units (ISCU). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* orgLevelCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Organization level set

(orgLevelSet). Categories of organization levels. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Organization sector

(orgSector). The sector of an organization *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: orgSectorCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization sector set

(orgSectorSet). Cateogries of organization sectors. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization Type

(orgType). Specifies the type or purpose of a given organization. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: orgTypeCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization type set

(orgTypeSet). The set for storing organization types. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample origin

(origin). An attribute of a sample specifying the origin. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample Origin set

(originSet). The set for storing the valid categorical values of sample origin. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Other Collection deprecated

(otco). Other type of collection method. Add description to collectionOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Other aggregation scale

(othAgg). Specifies an aggregation scale that can't be described as either "qualitative" or "quantitative". *Part Type:* aggregationScales *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other

(other). Other Other category. A categorical response for several variables in version 1. Also accompanied by a free text field. Deprecated in version 2, with common text fields added as new categories in version

2.0.0 *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation set

(otherAggrSet). Aggregation set used for unitless measures. *Part Type*: aggregationSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation deprecated

(otherDep). Other aggregation method. Add description to aggregationOther Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA

Parts

Maximum Length: NA *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Access to other prov deprecated

(otherProvDep). If this is ‘no’, this data will not be available to other data providers not listed before. If missing, data will be available to other data providers not listed before *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Other report table set

(otherReportTableSet). Additional tables that form the full data model. *Part Type:* dictSet *Domain:* naDomain *Specimen ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Variants Other

(otherV). Used for reporting variants detected for which a specific measureID doesn’t already exist. The ODM team will periodically migrate collected responses in this category into set measureIDs. Please write

in an official name for the new variant that was detected, and specify which nomenclature you're using with the nomenclatureID. *Part Type:* measurements *Domain:* Bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other Site - sample depreciated

(otsisaDep). Other type of site. Add description to typeOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at ph-esd_odm@ohri.ca *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Other Site - Site depreciated

(otsisiDep). Other site type. Add description to typeOther *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Other Unit depreciated

(otunDep). Other measurement of viral copies or wastewater treatment plant parameter. Add description to UnitOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Outbreak

(outb). Measure to indicate outbreak status. Given that when outbreaks occur, full data on case counts may not be available or completely accurate, using this measure indicates that case counts may be approximate. The datetime for the measure helps indicate that start and end periods for an outbreak. *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* humanCompartmentSet *Group:* miscMeas *Class:* outbreakClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* outbreakCatSet *Unit Set:* unit-lessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Data Type: categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak End

(outbEnd). Indicates an outbreak has ended. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak - on-going

(outbOngoing). Indicates an outbreak is on-going. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak set

(outbreakSet). set for the valid values of the outbreak measure. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* out-

Parts

breakCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak start

(outbStart). Indicates an outbreak began. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

p100l delta-variant gene target

(p100l). p100l delta-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

p2046l delta-variant gene target

(p2046l). p2046l delta-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*:

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

p2287s delta-variant gene target

(p2287s). p2287s delta-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

p3395h omicron-variant gene target { #p3395h }

(p3395h). p3395h omicron-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

p681r delta-variant gene target

(p681r). p681r delta-variant gene target *Part Type*: measurements *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*:

Parts

anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Paired Layout

(pairLay). Specifies the paired layout for a sequencing method *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Pangolin nomenclature

(pangolin). Specifies variant or genetic nomenclature as set out by the Phylogenetic Assignment of Named Global Outbreak Lineages (Pangolin). *Part Type*: nomenclatures *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parent dataset ID

(parDatasetID). The datasetID that is a parent to another datasetID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parent sample ID

(parent). If this sample has been pooled into one big sample for analysis this indicates the sample of the larger pooled sample. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.1.0. Last updated: 2.0.0.

Parent Site ID

(parSiteID). The siteID that is a parent to another siteID; usually the sub-site will be contained in the parentSite. An example would be surface testing in a LTC facility. The LTC facility would be the parentSiteID, and the various rooms/surfaces would be the sub-sites, or child-sites. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA

Parts

Unit Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Part description

(partDesc). The description of the part. Provides description of the part, used for translation. *Part Type:* partSupport *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Part ID

(partID). The unique identify of any entity within the dictionary. Every entity in the ODM has a partID. The partID must be unique with no duplicates. partIDs are short. partIDs that are constructed by different part attributes can be up to 28 characters (check max), but each part attribute should be 1 to 5 characters. partIDs are machine readable, but a knowledgeable user can understand what they represent. partID are used, for example, as a header in an Excel input table. Numbers are allowed but there are no other special characters other than ‘:’. ‘:’ is a reserved character that is used to generate a partID using a component or table attributes. A partID for a subComponent is automatically generated using attributes of a component. For example, the partID for the component ‘SARS-CoV-2’ is ‘covid’. ‘Covid’ can have many subComponents including different groups (Covid-19-al, Covid-19-va, covid-pr (allele, variant or protein) or units and aggregations (i.e. covid-al-N1-mean-cpl representing covid allele N1, mean

observation of viral copies per litre of effluent). A similar approach to naming subComponent is use for naming parts in a reprotTable. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 28 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Part instruction

(partInstr). Additional notes and instructions on how a part is used and/or defined. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 200 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Label

(partLabel). A human readable label of a part. Typically, a part label has no acrynomns (every word is spelled out). Equivalent to a LOINC common name. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Parts Look-up table

(parts). Look up table containing all parts in of the data model. Contains all parts, including self-referential parts. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: lookup *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts table column order

(partsOrder). Order of headers in the Parts table. *Part Type*: table-Support *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts header req

(partsRequired). Required headers in the Parts table. *Part Type*: table-Support *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Part types

(partType). Part types describe the purpose or use of the part. Environment data has three main part types: measure (i.e. covN1 viral measures, wastewater flow rate, temperature), method (e.g. how the measure was taken), and attribute (such as a site name). See description of each of the part types within categorySetID = partTypeCatSets. *Part Type: partType Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 15 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Primary Clarifier Effluent deprecated

(pceDep). Effluent obtained after primary clarifiers. Deprecated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: waterCompartmentSet Group: colGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: NA Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.*

Parts

Postal or Zip Code

(pCode). The zip code or postal code for a given address, specifying a specific geographic area. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

PCR method

(pcr). Description of the PCR method. Description of the PCR method used *Part Type:* methods *Domain:* Bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* pcrCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

PCR quality set

(pcrQualitySet). Quality set for PCR measures. *Part Type:* qualityIndSets *Domain:* naDomain *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

PCR sequencingsSelection method

(pcrSeq). Specifies the PCR selection method for sequencing, ie. that a pre-determined primer was used. Primer should be specified if this is selected. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

PCR Method set

(pcrSet). The set capturing containing all of the valid methods for the PCR methodID. *Part Type:* mmaSets *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* pcrCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Primary Clarifier Sludge deprecated

(pcs). Sludge produced by primary clarifiers. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* water-CompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:*

Parts

NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Primary clarifier effluent

(pEfflu). Effluent obtained after primary clarifiers *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Polyethyleneglycol (PEG) precipitation

(peg). Peg precipitation *Part Type*: categories *Domain*: Che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Percent

(perc). Percentage. *Part Type*: units *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum*

Value: 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Percent recovery

(percRec). Percent of the surrogate recovery for a recovery efficiency control assay. Use this for reporting the results of any recovery control assay. *Part Type:* units *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

pH

(ph). pH measurement *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* pHClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* unitlessUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 1 *Maximum Value:* 14 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Access to phac

(phac). If this is 'no', the data will not be available to employees of the Public Health Agency of Canada - PHAC. If missing, data will be available to employees of the Public Health Agency of Canada - PHAC. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment*

Parts

Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

PHAGE

(phage). PHAGE header *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

PHAGE - notes

(phageNotes). PHAGE note *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Phenol chloroform

(phenCl). Nucleic acid extraction performed using phenol chloroform. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Pseudomonas virus phi6

(phi6). Measure of the amount of pseudomonas virus phi6. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

phi6 spike target

(phi6Mat). Pseudomonas virus phi6is used as the recovery efficiency control target. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Contact phone

(phone). Contact phone number, for the lab. *Part Type:* attributes
Domain: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:*
naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*
naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*
naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA
Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*
Length: 10 *Maximum Length:* 12 *Part Status:* active. First released: 1.0.0.
Last updated: 2.0.0.

Total phosphorous

(phos). Total phosphorous *Part Type:* measurements *Domain:* Phy *Spec-*
iment ID: siSaSpecimenSet *Compartment Set:* waterCompartmentSet
Group: miscMeas *Class:* standardConc *Nomenclature:* naNomenclature
Ontology Reference: NA *Category Set ID:* NA *Unit Set:* stdConcentratio-
nUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference*
Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maxi-*
imum Value: seeUnitVal *Minimum Length:* NA *Maximum Length:* NA
Part Status: active. First released: 1.1.0. Last updated: 2.0.0.

Total Phosphates

(phostot). Total phosphates *Part Type:* measurements *Domain:* Phy
Speciment ID: siSaSpecimenSet *Compartment Set:* waterCompartmentSet
Group: miscMeas *Class:* standardConc *Nomenclature:* naNomenclature
Ontology Reference: NA *Category Set ID:* NA *Unit Set:* stdConcentratio-
nUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference*
Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maxi-*
imum Value: seeUnitVal *Minimum Length:* NA *Maximum Length:* NA
Part Status: active. First released: 1.1.0. Last updated: 2.0.0.

PHRED quality score

(phred). PHRED Quality Score For Sequencing *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGroup *Class*: sequence *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Physical property

(phy). A physical property or object not characterized by life or chemistry. *Part Type*: domains *Domain*: Phy *Speciment ID*: anySpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Physical

(physical). Measures and methods related to generic physical properties. *Part Type*: classes *Domain*: Phy *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: physical *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Primary key

(pK). Primary key for a table. All report tables have a primary key. Each row in a report table has unique identifier recorded in the primary key. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean

(pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Mostly used for reporting specific alleles (N1, N2, E, etc.) *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-normalized mean

(pmmovNorm). Mean measure normalized to amount of PMMoV. Mostly used for reporting specific alleles (N1, N2, E, etc.) *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*:

naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Population

(po). An measure or observation for a geographic area or population. Examples include human compartment measures such as case numbers, hospitalization rates for diseases or surface or air compartment measures such as snow coverage or ambient temperature. *Part Type*: specimens *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Post-Grit deprecated

(pogr). Raw wastewater after a treatment plant's headworks. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Polygon ID

(polygonID). Unique identifier for the polygon. A polygon is a geographic area representing the capment area of an environmental sample. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Polygon table

(polygons). The table that contains information about the geometry of a geographic area. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Polygons table column order

(polygonsOrder). Specifies the order of the columns in a Polygons table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum*

Value: NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Polygon header req

(polygonsRequired). Specifies the columns required in a Polygons table.
Part Type: tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Polygon Population

(polyPop). An attribute of a polygon, which specifies the population of that polygon. A rough estimate of the number of human residents.
Part Type: attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Pooled

(pooled). Is this a pooled sample, and therefore composed of multiple child samples obtained at different sites *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature

Parts

Ontology Reference: NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Type of date for case reporting depreciated

(popDateTypeDep). Type of date used for confirmed cases. Typically, report or episode are reported. Onset and test date is not usually reported within aggregate data. *Part Type:* attributes *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* popDateTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Case reporting date set

(popDateTypeSet). set for the date of case reporting. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* popDateTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Population equivalents

(popEq). A unit for measuring of design capacity for wastewater treatment plants, represented as the approximate amount of water a single person would use. *Part Type:* units *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Population group

(popGrp). A group of measures/methods related to population-level factors. Examples might be case numbers, hospitalization rates, etc. *Part Type:* groups *Domain:* allDo *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* popGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population Served

(popServ). An attribute of a site, which specifies the population of/population served by a given site. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:*

Parts

categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Population unit set

(populationUnitSet). Unit set for hospital-related measurements. *Part Type*: unitSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: humanCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: populationUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Poulatiion or site set

(poSiSpecimenSet). A specimen set that inculdes population or site specimen. *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: poSiSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Population specimen set

(poSpecimenSet). A specimen set that inculdes only a population specimen. *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*:

naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Percent positive

(pp). Percent positive of sample measures. Can use for Moore swabs, etc. *Part Type: units Domain: Bio Speciment ID: siSaSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: quantAggScale Quality Set: NA Reference Link: NA Data Type: float Minimum Value: 0 Maximum Value: 100 Minimum Length: NA Maximum Length: NA Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

parts per million

(ppm). Parts per million. *Part Type: units Domain: Phy Speciment ID: siSaSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: quantAggScale Quality Set: NA Reference Link: NA Data Type: float Minimum Value: 0 Maximum Value: NA Minimum Length: NA Maximum Length: NA Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

PMMoV-CP

(ppmv). Pepper mild mottle virus capsid protein gene region *Part Type: measurements Domain: Bio Speciment ID: saSpecimenSet Compartment*

Parts

Set: anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Percent positivity rate

(pprt). Percent positivity rate of tests conducted within a day in a given region. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Percent primary sludge

(pps). Percentage of total solids, for primary sludge. *Part Type:* units *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Precipitation

(precipitation). Measures and methods related to precipitation. *Part Type:* classes *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* precipClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Precipitation unit set

(precipitationUnitSet). Unit set for percipitation measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* precipitationUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Pretreatment

(pretreat). Was the sample chemically treated in anyway with the addition of stabilizers or other Pretrement method can be describe using a method. See methodID. *Part Type:* methods *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum*

Parts

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Genetic primer

(primer). Method ID used for indicating the primer used for a PCR or sequencing assay. *Part Type:* methods *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* primerCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Genetic primer set

(primerSet). The set for genetic primers used in sequencing or PCR assays. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* primerCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Private sector

(priv). The category of organization type used for private sector or non-profit groups that are not otherwise captured by the academic category. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set*

ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 10 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Processing group

(procGrp). A group of measures/methods related to processing samples for analysis. *Part Type: groups Domain: allDo Speciment ID: siSaSpecimenSet Compartment Set: anyCompartmentSet Group: procGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Program decription tables

(programDescr). Tables used to describe surveillance and testing programs. Program description tables record metadata on the organizations, locations, and protocols. Information include the name, location, and other contact information. These tables are usually completed once and they are updated only when contact information changes. *Part Type: classes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: programDescr Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

promega automated tna kit

(promAuto). Nucleic acid extraction performed using the promega automated tna kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

promega ht tna kit

(promHt). Nucleic acid extraction performed using the promega ht tna kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

promega manual tna kit

(promManu). Nucleic acid extraction performed using the promega manual tna kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

promega wastewater large volume tna capture kit

(promWW). Nucleic acid extraction performed using the promega wastewater large volume tna capture kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Proportion of total

(prop). Proportion as a percent of total. *Part Type:* units *Domain:* naDomain *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of variant in sample

(propV). Proportion of a variant as percent of total variants. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Parts

Protocol ID

(protocolID). A unique identifier for a given protocol. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID container

(protocolIDContainer). Unique identifier for the protocol and the steps and other protocols that make it up. Protocol ID Container is populated by a protocolID (partID = protocolID) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID object

(protocolIDObj). The object of the relationship between one protocol and another protocol or step, within a given protocol ID container. Populate this field using Protocol ID (partID = protocolID). Together with 'protocolIDSub' or 'stepIDSub' and 'relationshipID', the protocol relationship is defined as: [protocolIDSubject or stepIDSubject] [relationshipID]

[protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID subject

(protocolIDSub). The subject of the relationship between one protocol and another protocol or step, within a given protocol ID container. Populate this field using Protocol ID (partID = protocolID). Together with ‘protocolIDObj’ or ‘stepIDObj’ and ‘relationshipID’ the row of ‘protocolRelationships’ should say: [protocolIDSubject or stepIDSubject] [relationshipID] [protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol relationships table

(protocolRelationships). The table that contains the organizational information and details (such as order) about a given protocol, recorded as a series of protocol steps (protocolSteps). Protocols are a group of related protocol steps. For example a protocol step can describe a sample concentration, extraction, inhibition or a measurement. Protocols can include

Parts

both methodID (a procedure) and measureIDs (a measure). An example of a method is centrifugation; where the centrifugation rotation speed is described as a measure. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol relationships table column order

(protocolRelationshipsOrder). Specifies the order of the columns in the Protocol Organization table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol Relationship header req

(protocolRelationshipsRequired). Specifies the columns required in the Protocol Organization table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol Relationships set

(protocolRelSet). set for valid values of relationshipID in the protocolRelationships table. Use only for protocolRelationships table relationshipID. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: protocolRelSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocols table

(protocols). The table for protocols. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: programDescr *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocols table column order

(protocolsOrder). Specifies the order of the columns in the Protocols table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*:

Parts

naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocols header req

(protocolsRequired). Specifies the columns required in the Protocols table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol steps table

(protocolSteps). The table for collecting metadata on individual steps in a protocol, methodological process, or assay. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: programDescr *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Protocol steps table column order

(protocolStepsOrder). Specifies the order of the columns in a Protocol Steps table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol steps header req

(protocolStepsRequired). Specifies the columns required in a Protocol Steps table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol table set

(protocolTableSet). tables holding information on the methods used for sample collection or measurement. *Part Type*: dictSet *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum*

Parts

Length: NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol version

(protocolVersion). Specifies the version of a method set. Version of the method set. Semantic versioning is recommended. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Access to prov ha

(provHA). If this is ‘no’, this data will not be available to provincial health authorities. If missing, data will be available to provincial health authorities. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Provisional report

(provisional). A provisional or interm report. Use for measurement that are not yet finalized. Typically this purpose is not publicly reported. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet

Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
 First released: 2.0.0. Last updated: 2.0.0.

Primary clarifier sludge

(pSludge). Sludge produced by primary clarifiers *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Pumping station

(pStat). Pumping station *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Post-grit

(pstGrit). Raw wastewater after a treatment plant's headworks (post grit or removal of large solids at a treatment plant but prior to a primary clarifier) *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Pretreatment description - depreciated

(ptDescDep). If preTreatment then describe the treatment that was performed. Deprecated in version 2, please do not use. Use pretreatment methodID and specify further in summary and/or notes, if necessary. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 255 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Number of positive tests

(ptot). Number of positive tests conducted within a day in a given region. *Part Type*: units *Domain*: Bio *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference*

Link: NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Access to public

(public). If this is ‘no’, this data will not be available to the public. If missing, data will be available to the public. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Puro virus

(puro). Measure of the amount of puro virus. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Puro virus spike target

(puroMat). Puro Virus is used as the recovery efficiency control target. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass

Parts

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Purpose

(purpose). The reason the measure or sample was taken. See allowed categories: report, quality control, validation study, test *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* purposeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Purpose category

(purposeSet). Purpose set for a sample or measure. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* purposeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

q19e omicron-variant gene target { #q19e }

(q19e). q19e omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

q493r omicron-variant gene target { #q493r }

(q493r). q493r omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

q498r omicron-variant gene target { #q498r }

(q498r). q498r omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

q954h omicron-variant gene target { #q954h }

(q954h). q954h omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Quality assurance method

(qaqc). Quality ass qc Description of the quality control steps taken *Part Type*: methods *Domain*: naDomain *Speciment ID*: anySpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 100 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Quality concerns

(qf1). A flag to indicate there is a quality concern, not otherwise sepcified. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: 0 *Maximum Value*: 0 *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

qiagen allprep dna/rna kit

(qgDNARNA). Nucleic acid extraction performed using the qiagen allprep dna/rna kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen digital PCR

(qgDpcr). Describes a PCR analysis done using Qiagen's digital PCR technology. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen ez1 virus mini kit v2.0

(qgEz1). Nucleic acid extraction performed using the qiagen ez1 virus mini kit v2.0. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

qiagen allprep powerfecal dna/rna kit

(qgPwrFecal). Nucleic acid extraction performed using the qiagen allprep powerfecal dna/rna kit. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep powerviral dna/rna kit

(qgPwrViral). Nucleic acid extraction performed using the qiagen allprep powerviral dna/rna kit. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen powerwater kit

(qgPwrWtr). Nucleic acid extraction performed using the qiagen powerwater kit. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen qiaamp buffers with epoch columns

(qgQiAmp). Nucleic acid extraction performed using the qiagen qiaamp buffers with epoch columns. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy kit

(qgRneasy). Nucleic acid extraction performed using the qiagen rneasy kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy powermicrobiome kit

(qgRneasyPwr). Nucleic acid extraction performed using the qiagen rneasy powermicrobiome kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quantitative PCR

(qpcr). Real-time PCR, also called ‘quantAggScale’ PCR *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Qualitative

(qualAggScale). The qualitative aggregation scale. *Part Type:* aggregationScales *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* qualAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality control

(qualityControl). A measure or sample taken for the purpose of quality control. Measures with this category are take to assess quality control. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum*

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 2.0.0. Last updated: 2.0.0.

Quality flag

(qualityFlag). A field for reporting any quality concerns - of lack thereof - for a sample or measure. Types of quality flag vary on context; please see the qualityFlagSet to be confirm. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality report

(qualityID). A unique identifier for a given quality report. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicators

(qualityIndicators). A measure of the quality of a reported value or sample. Only used for the dictionary entries of parts. *Part Type:* part-*Type Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment*

Parts

Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicator sets

(qualityIndSets). Sets of quality indicators or measures. For example, PCR have a quality measure set. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality reports table

(qualityReports). The table for recording the various quality metrics and indicators for samples and measures. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* results *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality repots table column order

(qualityReportsOrder). Specifies the order of the columns in the Quality Reports table. *Part Type:* tableSupport *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality report header req

(qualityReportsRequired). Specifies the columns required in the Quality Reports table. *Part Type:* tableSupport *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality set

(qualitySetID). The quality set that corresponds to a given part. Only applicable for samples and measures. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0

Parts

Maximum Length: 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quantitative

(quantAggScale). The “quantAggScale” aggregation scale. *Part Type:* aggregationScales *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

R squared

(r2). R-squared value of the calibration curve. Used for storing calibration curve information, potentially also applicable for other curve types. *Part Type:* measurements *Domain:* allDo *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standard-Curve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

r203m delta-variant gene target

(r203m). r203m delta-variant gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* genetic-

UnitSet Aggregation Scale: seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Raindance digital PCR

(rainDpcr). Describes a PCR analysis done using Raindance's digital PCR technology. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Rainy

(rainy). Qualitative category for the weather measure, specifying a rainy day. *Part Type:* categories *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Random sequencing selection method

(ranSeq). Specifies the random selection method for sequencing No primer set should be used if this is selected. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet

Parts

Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ratio

(ratio). Ratio (unitless) Report as a real number or fraction. i.e 10:2 ratio is reported as 5. The specifics of the ratio (the numerator and the denominator) are described in the measure. *Part Type:* units *Domain:* naDomain *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Raw wastewater

(raws). Raw wastes water sample *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage at site

(rawWW). Wastewater without any form of treatment *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage downstream from a site

(rawWWdown). Downstream from a site. See partType = rawWWup. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Raw sewage upstream from a site

(rawWWup). Upstream from a site. Used when there is not direct access to a site's wastewater inconjuction with downstream sample. For example, if the site is a long-term care home, but ther is not access to the building or property cleanout. Use two sameple, one sample upstream and one sample downstream. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality*

Parts

Set: NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

IP2 rdrp sars-cov-2 gene target

(rdrpIP2). IP2 rdrp sars-cov-2 gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

IP4 rdrp sars-cov-2 gene target

(rdrpIP4). IP4 rdrp sars-cov-2 gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Date sample recieved

(recDate). The date the sample was received at the laboratory for analysis. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set*

ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: datetime Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Recovered patients

(recov). Units for describing a population measure of patients who have recovered from a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type: units Domain: Bio Speciment ID: poSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: quantAggScale Quality Set: NA Reference Link: NA Data Type: integer Minimum Value: 0 Maximum Value: NA Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Reference link

(refLink). Link to the reference material for a part. May link to literature on a method, measure, etc. Link to standard operating procedure for a measure or part. Part reference is a reference to a technical document that describes the part in detail. The reference could be a URL, DOI, or reference to a technical report. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 255 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Parts

Regular

(regular). A measure or sample taken for surveillance or epidemiology. A ‘regular’ measure is the default. Missing data will be recoded to this purpose. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Regular reports table set

(regularReportTableSet). Tables used for daily reporting of new measurements and information on sample collection. *Part Type*: dictSet *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Relative humidity

(relHum). The unit of relative humidity, or the air-water mixture. The ratio of the partial pressure of water vapor in the mixture to the equilibrium vapor pressure of water over a flat surface. of pure water at a given temperature. *Part Type*: measurements *Domain*: Phy *Speciment*

ID: siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* site-Feat *Class:* humClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* relHumidUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Relative humidity unit set

(relHumidUnitSet). Unit set for relative humidity measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* relHumidUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Representative grab sample

(repGrab). A single large representative grab sample. Deprecated in version 2, please do not use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Replicate Type set

(replicateSet). The set for storing valid categorical values of replicate type. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* replicateCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

replicate number

(repNum). The replicate number for a Ct or Cq value - used for specifying out a protocol standard curve. Only applies and used for aggregated standard curves, or curves that are used across multiple plates. Otherwise, these details can be captured in a measure report and specified by measure index. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* see-UnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Primary reporting authority

(repOrg1). The primary or most responsible authority for routine surveillance reports or findings, or where the site is located. Use the ‘organizations’ table to describe attributes such as the name and address of the

organization. Use 'orgType', 'orgLevel' and 'orgSector' to describe additional attributes. For example, whether the organization is a government agency ('orgType'), local ('orgLevel'), or public health ('orgSector'). If the attributes are not specified, the default is version 1.1.0 'publicHealthDepartment' defined as 'The public health department or region where the site is located. See also 'healthReg' if there is a separate regional health care delivery authority.' *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Secondary reporting authority

(repOrg2). The secondary, additional or alternative authority for routine surveillance reports or findings, or where the site is located. Use the 'organizations' table to describe attributes such as the name and address of the organization. Use 'orgType', 'orgLevel' and 'orgSector' to describe additional attributes. For example, whether the organization is a government agency ('orgType'), local ('orgLevel'), or public health ('orgSector'). If the attributes are not specified, the default is version 1.1.0 'healthRegion', The health planning authority where site is located. defined as 'The public health department or region where the site is located.' *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Parts

Reportable

(reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: booleanSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: boolean *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Report date

(reportDate). The date a measure was reported. This date is the first date the laboratory or entity reported the measure findings to the public health or other entity who ordered the measure or is responsible for action. Use the analyses end date if the measure is part of a research study or other purpose that does not have a reporting entity. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: datetime *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Reporter ID depreciated

(reporterIDDep). Unique identifier for the person or organization that is reporting the data. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*:

naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.1.0. Last updated: 2.0.0.

Reporter table deprecated

(reportersDep). The table that contains information about a reporter of a sample, method, measure, or attribute. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Replicate Type

(repType). Attribute of a sample, specifying whether the sample is unique, or a replicate. And if it is a replicate, what type of replicate is it? *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: replicateCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Resuspend COSCa filter collection

(resCosca). Nucleic acid extraction via resuspension of a sample collected using the COSCa ball method. See partID = cosca for additional details. *Part Type*: categories *Domain*: Bio *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Results tables

(results). Tables used to record samples, measures and quality reports. These tables are updated whenever a new sample is taken or a measure is reported. Quality reports are used to report quality assurance and quality control for samples and measures. *Part Type*: classes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: results *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Retention pond

(retPond). Retention pond *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation*

Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar
Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

River, natural water body

(river). River, natural water body *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Role of contact

(role). Specifies the organizational role of a given contact. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Rate per 100,000

(rP100). Units for describing a population measure of a rate of case incidence per 100,000 people of a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:*

Parts

units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Raw post-pasteurized wastewater

(rppw). Raw wastewater sample post-pasteurization *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Respiratory syncytial virus

(rsv). Respiratory syncytial virus. *Part Type:* measurements *Domain:* Bio *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* http://purl.obolibrary.org/obo/NCIT_C14267 *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

s2083i omicron-variant gene target { #s2083i }

(s2083i). s2083i omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

s371l omicron-variant gene target { #s371l }

(s371l). s371l omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

s373p omicron-variant gene target { #s373p }

(s373p). s373p omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

s375f omicron-variant gene target { #s375f }

(s375f). s375f omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

s477n omicron-variant gene target { #s477n }

(s477n). s477n omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample

(sa). A measure made on a compartment or property from a sample of a substance. *Part Type*: specimens *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample material

(saMaterial). Type of sample. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sampleMatCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sample ID

(sampleID). Unique identifier for a sample. Suggestion:siteID-date-index. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Sample ID object

(sampleIDObject). Populated by 'sampleID' - this specifies the object of a sample relationship. Use 'sampleID' for this field. Together with 'sampleIDSub' and 'relationshipID' the row of the 'sampleRelationships' should say: [sampleIDSubject] is a[relationshipID] of [sampleIDObject] *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA

Parts

Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample ID subject

(sampleIDSubject). Populated by 'sampleID' - this specifies the subject of a sample relationship. Use 'sampleID' for this field. Together with 'sampleIDSub' and 'relationshipID' the row of the 'sampleRelationships' should say: [sampleIDSubject] is a[relationshipID] of [sampleIDObject]
Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample material set

(sampleMatSet). set for the types of material that can be found in a sample.
Part Type: mmaSets Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: sampleMatCatSet Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample quality set

(sampleQualitySet). Quality set for a sample. *Part Type:* qualityIndSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationships table

(sampleRelationships). Table for recording the relationships between samples. Samples can be pooled or split. The sample relationships table holds information on parent-child relationships between samples, and allow for tracking sample lineage for single and pooled samples. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* results *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationships table column order

(sampleRelationshipsOrder). Specifies the order of the columns in a Sample Relationships table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*

Parts

Scale: NA Quality Set: NA Reference Link: NA Data Type: integer Minimum Value: 1 Maximum Value: NA Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationships header req

(sampleRelationshipsRequired). Specifies the columns required in a Sample Relationships table. *Part Type: tableSupport Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Sample relationship

(sampleRelID). Attribute for specifying the relationship between a sample subject and object. Together with ‘sampleIDSubject’ and ‘sampleRelID’, a row of the ‘sampleRelationships’ should say: [sampleIDSubject] is a [sampleID] of [sampleIDObject] *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: sampleRelSet, protocolRelSet Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Sample relationships set

(sampleRelSet). set for valid values of relationshipID in the sampleRelationships table. Use only for sampleRelationships table relationshipID. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sampleRelSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample report table

(samples). The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. Samples can be combined, split, stored and reused. The ‘Sample relationships’ (sampleRelationships) is used to describe the relationships between samples. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* results *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sample shed

(sampleShed). A geographic area, physical space, or structure. A sample is taken from a sampleShed for a representative measurement of a

Parts

substance(s). Examples: Airplane, Long-term care home, Neighbourhood, University campus, Municipality *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: shedCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Samples table column order

(samplesOrder). Specifies the order of the columns in a Samples table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample header req

(samplesRequired). Specifies the columns required in a Samples table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Collection other depreciated

(sampleTypeOther). Description for other type of method not listed in collection. *depreciated Part Type: methods Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: colGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 100 Part Status: **depreciated***. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2

(sarsCov2). A group of measures/methods related to the SARS-CoV-2 virus. SARS-CoV-2 measures have the following classes: proteins, alleles, variants, mutations, diseases. *Part Type: groups Domain: Bio Speciment ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: sarsCov2 Class: naClass Nomenclature: naNomenclature Ontology Reference: http://purl.obolibrary.org/obo/NCIT_C169076 Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active*. First released: 2.0.0. Last updated: 2.0.0.

Sample specimen set

(saSpecimenSet). A specimen set that inculdes only a sample specimen. *Part Type: specimenSets Domain: naDomain Speciment ID: saSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value:*

Parts

NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sparse coverage

(sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. If there are multiple issues to flag, use the multiple issues flag and explain the details in the notes section. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier effluent depreciated

(sce). Effluent obtained after secondary clarifiers. Depreciated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

School

(school). A school serving students in the kindergarten to 12th grade range
Part Type: categories *Domain:* naDomain *Speciment ID:* naSpecimenSet
Compartment Set: anyCompartmentSet *Group:* colGrp *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA
Unit Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA
Minimum Length: 0 *Maximum Length:* 30 *Part Status:* active.
 First released: 1.0.0. Last updated: 2.0.0.

Sequencing script version set

(scriptSet). The set for script versions *Part Type:* mmaSets *Domain:* Bio
Speciment ID: saSpecimenSet *Compartment Set:* anyCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA
Category Set ID: scriptCatSet *Unit Set:* naUnitSet
Aggregation Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData
Minimum Value: seeUnitVal *Maximum Value:* seeUnitVal
Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Genetic sequencing script version

(scriptVersion). Used to specify the script or script version used to extract summary information (i.e. coverage and mutation frequency statistics, etc.) from the sequencing output. *Part Type:* methods *Domain:* Bio *Speciment ID:* saSpecimenSet
Compartment Set: anyCompartmentSet *Group:* procGroup *Class:* sequence
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* scriptCatSet
Unit Set: geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA
Reference Link: NA *Data*

Parts

Type: seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier sludge deprecated

(scsDep). Sludge produced by secondary clarifiers. Deprecated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Standard Deviation

(sd). Standard deviation. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Normalized Standard Deviation

(sdn). Standard deviation, normalized. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature

Ontology Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sea, natural water body

(sea). Sea, natural water body *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

See unit (aggregation scales)

(seeUnitAggScale). A value for aggregation scale to be used when the aggregation scale depends on that specified in the parts list entry for the unit. When you encounter this value, see the unit being used for that entry and follow that aggregation scale entry. *Part Type:* aggregation-Scales *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

See unit (data type)

(seeUnitData). The data type entry when the data type for a given entry is specified by the unit used. Only used for measures, which is where the data type entry is used to validate the entry in the value field. For these cases, however, the data type is more dependant on the unit than the measure , and so this entry directs users to consult the data type of the unit. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier effluent

(sEfflu). Effluent obtained after secondary clarifiers *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Access to self

(self). If this is 'no', this data will not be shown on the portal when this reporter logs in. If missing, data will be available to this reporter. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Date sample was sent

(sentDate). The date the sample was sent for analyses at a laboratory. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Septic tank wastewater

(septage). Wastewater from within a septic tank *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Septic tank

(septTnk). Septic tank *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet

Parts

Group: colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sequencing Layout

(seqLay). The layout of genetic material used in sequencing *Part Type:* methods *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGroup *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* seqLayCatSet *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Layout set

(seqLaySet). The set for the layout of genetic material used in sequencing *Part Type:* mmaSets *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* seqLayCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing quality set

(seqQualitySet). Quality set for sequencing measures *Part Type*: qualityIndSets *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method

(seqSel). The primer sequence selection method used in sequencing *Part Type*: methods *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGroup *Class*: sequence *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: seqSelCatSet *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method set

(seqSelSet). The set for the selection method used in sequencing *Part Type*: mmaSets *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: seqSelCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum*

Parts

Value: seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy

(seqStrat). The sequencing strategy used for an analysis *Part Type:* methods *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGroup *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* seqStratCatSet *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy set

(seqStratSet). The set for the sequencing strategy used for an analysis *Part Type:* mmaSets *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* seqStratCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Genetic sequences

(sequence). Measures and methods related to genetic sequencing. *Part Type:* classes *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

naUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sewer Sediments Depreciated

(seseDep). Sediments obtained in sewer. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Set

(set). The unique identifier of a value set. A set is a group of units, attributes, categories, specimens, or compartments that can be used for a measure, method or attribute. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 28 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Sets look-up table

(sets). Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together. This table manages many-to-many relationships and category recycling within the ODM. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* lookup *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sets table column order

(setsOrder). Specifies the order of the columns in the Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sets header req

(setsRequired). Required headers in the Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Settled solids

(settsol). Amount of settled solids from a wastewater sample. May be a volume or a weight, see the unitID for a given measure row for details on a given sample. *Part Type: measurements Domain: Phy Speciment ID: saSpecimenSet Compartment Set: waterCompartmentSet Group: proc-Grp Class: physical Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: volumeUnitSet Aggregation Scale: quantAggScale Quality Set: NA Reference Link: NA Data Type: integer Minimum Value: 0 Maximum Value: NA Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Set type

(setType). The type of set. i.e. quality set, aggregation set, unit set *Part Type: partSupport Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: waterCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: NA Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Set Value

(setValue). The partID for the set value or category. Populated by any part where partType = aggregationSet, compartmentSet, qualitySet, specimenSet, or unitSet *Part Type: attributes Domain: naDomain Speciment*

Parts

ID: naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Severity indicator

(severity). An indicator of the severity or seriousness of a quality flag. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sevCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Severity set

(sevSet). A set for severity indicators. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sevCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sewer network file link deprecated

(sewerNetworkFileBLOB). Link to a file that has any detailed information about the sewer network associated with the site (any format). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

S sars-cov-2 gene target

(sGene). S sars-cov-2 gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sampleshed set

(shedSet). The set for all valid values of sampleshed. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* siSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA

Parts

Minimum Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ship

(ship). A cruise ship or other ship. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Shipped on ice

(shipOnIce). Was the sample kept cool while being shipped to the lab *Part Type:* methods *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Short names

(shortName). Shortened names for tables and other important parts for use in wide names. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:*

NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Site

(si). A measure made on a compartment or property at a site. An example is the wastewater flow rate, taken at a wastewater treatment plant. *Part Type*: specimens *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Single

(sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Single Layout

(sinLay). Specifies the single layout for a sequencing method *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*:

Parts

anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Site or sample specimen set

(siSaSpecimenSet). A specimen set that includes site or sample specimen. *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: siSaSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Site specimen set

(siSpecimenSet). A specimen set that includes only a site specimen. *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: siSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Site ID default

(siteDef). Used as default when a new sample is created by this reporter. See ID in Site table. Deprecated in version 2, please do not use. Just populate standard SiteID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Site features group

(siteFeat). A group of environmental measures/methods and those pertaining to the features of a site. *Part Type:* groups *Domain:* allDo *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Site ID

(siteID). Unique identifier for the location where a sample was taken. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum*

Parts

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 1.0.0. Last updated: 2.0.0.

Site measure ID

(siteMeasureID). Unique identifier for wide table only. Use when all measures are performed on a single sample. Depreciate in version 2. *siteMeasures Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Sites table

(sites). The table that contains information about a site; the location where an environmental sample was taken. The site of an enviromental sample. Information in the site table does not regularly change. Consider using the MeasureReport table if the infomation changes often. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sites table column order

(sitesOrder). Specifies the order of the columns in a Sites table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sites header req

(sitesRequired). Specifies the columns required in a Sites table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Site Type

(siteType). Type of site or institution where sample was taken. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* siteTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Site set

(siteTypeSet). set for the type of sampling site. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* siteTypeCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Skimmed milk flocculation

(skimMilkFloc). Skimmed milk flocculation *Part Type:* categories *Domain:* Che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

slope

(slope). Slope value of the calibration curve. Used for storing calibration curve information, potentially also applicable for other curve types. *Part Type:* measurements *Domain:* allDo *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sewer network file link depreciated 2

(snfl). Link to a file that has any detailed information about the sewer network associated with the site (any format). Depreciate in version 2. This attribute has been changed to referenceLink. Use po_referenceLink to specify polygon link to sewer network file. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Snowy

(snowy). Qualitative category for the weather measure, specifying a snowy day. *Part Type:* categories *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Solid fraction

(sol). Solid fraction of a sample. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet

Parts

Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Solid seperation

(solidSep). Process used to separate solid and liquid phases of the sample. Solid seperation is used either prior to or in the absence of the concentration method specified in 'methodConc'. *Part Type:* methods *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* nput *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Solid separation set

(solidSeparationSet). set for the separation of solids. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* solidSeparationCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Source Protocol ID

(sourceProtocol). A protocol that served as a basis for another protocol. This may be due to version changes, or a referencing a protocol that was later adapted. Populate with protocolID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step source ID

(sourceStep). Specifies the protocol step which serves as a basis for a given protocol step. Populate with a protocol step ID (stepID). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity

(specHum). Measure for specific humidity, the unit is the ratio of the mass of water vapour and the mass of total air. *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* humClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit*

Parts

Set: specHumidUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity unit set

(specHumidUnitSet). Unit set for specific humidity. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* specHumidUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specifies

(specifies). Specifies that the object step or protocol occurs specifies details for the subject step or protocol in the overall protocol container. Use only for protocolRelationships table relationshipID. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specimen

(specimen). The substance or thing upon which the observation was made. Specimens include: site, sample, person and population. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specimens

(specimens). Measures or observations are taken from three types of substances: populations (humans or geographic areas), samples, sites. See *partID = specimentID* *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specimen sets

(specimenSets). Sets of specimens. Specimen sets are used when a measure or attribute can be recorded for more than one specimen. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Spike material

(spikeMat). Material into which the recovery efficiency control target is spiked. *Part Type:* methods *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeMatCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 4 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Spike material set

(spikeMatSet). set for spikeMat (aterial into which the recovery efficiency control target is spiked). *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeMatCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Recovery efficiency spike target

(spikeTarget). Method specifying the recovery efficiency control target. This matches on to the the spike material method ID. Should typically be indexed as a step proceeding the spikeMaterial methodID step. *Part*

Type: methods *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeTargetCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Recovery efficiency spike target set

(spikeTargetSet). The set capturing containing all of the valid categories for the recovery efficiency control (or spike) target. *Part Type:* mmaSets *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeTargetCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample spilled

(spill). Sample contents spilled from container. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Secondary clarifier sludge

(sSludge). Sludge produced by secondary clarifiers *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Social services shelter

(sss). Other type of social services shelter *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Stabilization pond

(stabPnd). Specifies a site type which is a pond designed and built for wastewater treatment to reduce the organic content and remove pathogens from wastewater. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Standard concentrations

(standardConc). Measures and methods relating to standard concentrations. *Part Type:* classes *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Standard curve

(standardCurve). Measures and methods relating to generating or recording a standard curve. *Part Type:* classes *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Departments, states, or provinces

(stateProvLevel). Departments, states, or provinces *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* ISCU

Parts

Ontology Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

State, Province, or Region

(stateProvReg). The state, province, or region where a site or organization is located; part of the address. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Status

(status). Whether the part is still active and can be used in the most current ODM version. Values are 'active' or 'inactive'. The status of the part or list. 'active' means the part is used in the dictionary version. 'inactive' means the part has been retired and not longer used in the version. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* statusCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Status set

(statusSet). A set for partID = Status to indicate whether a part is in current use or not. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Standard concentration unit set

(stdConcentrationUnitSet). Unit set for concentration measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Storage temp

(sTemp). Temperature that the sample is stored at in Celsius. *Part Type:* measurements *Domain:* Phy *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* tempClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Protocol Step ID

(stepID). The unique identifier for a specific protocol step. Protocol Steps are the component parts of a larger protocols, the latter being linked to specific measure reports. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Step ID Object

(stepIDObj). The object of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Populate this field using Protocol Step ID (partID = stepID). Together with 'protocolIDSub' or 'stepIDSub' and 'relationshipID' the row of 'protocolRelationships' should say: [protocolIDSubject or stepIDSubject] [relationshipID] [protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Step ID Subject

(stepIDSub). The subject of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Populate

this field using Protocol ID (partID = stepID). Together with ‘protocolIDSub’ or ‘stepIDSub’ and ‘relationshipID’ the row of ‘protocolRelationships’ should say: [protocolIDSubject or stepIDSubject] [relationshipID] [protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Method step parent ID

(stepProvenanceID). A method step that served as a basis for another method step. This may be due to version changes, or a referencing a method set that was later adapted. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step version

(stepVer). Specifies the version of a given protocol step. Version of the method step. Semantic versioning is recommended. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

Parts

naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA
Data Type: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 50 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer network file blob

(storTempDef). A file blob that has any detailed information about the sewer network associated with the site (any format). Deprecated in version 2, please do not use. Refer to WKT and EPSG instead. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: blob *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 65535 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Storage time

(stoTim). Length of time that a sample was in storage. *Part Type*: measurements *Domain*: allDo *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: timeClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: timeUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sub-variant or lineage

(subVar). A unit used to report a specific genetic lineage, or to specify that a reported variant is a sub-variant of another. Combine in a measure report if necessary. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* variant *Nomenclature:* pangolin *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Summary

(summ). Short description of the assay and how it is different from the other assay methods. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 255 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Summary worksheet

(summary). Summary sheet of the dictionary.xlsx file. The summary sheet contains versioning and changelog. *Part Type:* dictionarySupport *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar

Parts

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sunny

(sunny). Qualitative category for the weather measure, specifying a clear and sunny day. *Part Type:* categories *Domain:* Phy *Specimen ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment

(surf). A measure or observation made from a substance on a surface. *Part Type:* compartments *Domain:* allDo *Specimen ID:* anySpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other surface

(surface). Surface other than floor or desk. See also categories for floor and desk. *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:*

naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment set

(surfaceCompartmentSet). A compartment set for measures and methods in the surface compartment. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: surfaceCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Surface and water compartment set

(surfaceWaterCompartmentSet). A compartment set for measures and methods in the surface or water compartments. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: surfaceWaterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Surface swab

(surfSw). Surface swab. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sewage truck

(swgTrck). Sewage truck *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sewer sediment

(swrSed). Sediments obtained in sewer *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sewer catchment area

(swrSet). Sewer catchment area. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Synthetic sample

(synthetic). Specifies a synthetic, or artificially derived sample. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

t19r delta-variant gene target

(t19r). t19r delta-variant gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

t3646a delta-variant gene target

(t3646a). t3646a delta-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

t547k omicron-variant gene target { #t547k }

(t547k). t547k omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

t9i omicron-variant gene target { #t9i }

(t9i). t9i omicron-variant gene target *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Tables

(tables). Tables are where measures, methods and attributes are recorded. Tables represent the main entities of the environmental surveillance. There An example of a table is ‘site’ a collection of attributes such as site name and address to describe where environment samples are taken. PHES-ODM has report tables, but users can create their own tables as well. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Table column set

(tableSet). The set for valid inputs in table columns. These categories are used in the parts list for columns named after tables to indicate the position and function of a given part in the entity relationship diagram of the data model. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* tableCatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

TaqPath N sars-cov-2 gene target

(taqpatN). TaqPath N sars-cov-2 gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:*

Parts

anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

TaqPath S sars-cov-2 gene target

(taqpatS). TaqPath N sars-cov-2 gene target *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample temperature

(temp). Temperature of the sample measured in degrees Celcius *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* tempClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Temperature

(temperature). Measures and methods related to temperatures. *Part Type:* classes *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* tempClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Temperature unit set

(temperatureUnitSet). Unit set for temperature measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nucleic acid template volume

(tempVol). The volume of DNA or RNA template used for PCR. *Part Type:* measurements *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Test date

(tesDate). Date that the covid-19 test was performed. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Number of tests performed

(test). Number of tests performed. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Testing

(testing). A measure or sample taken to test a new method. These measures are typically used for internal lab uses and not reported to external partners. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

thermo magmax microbiome ultra nucleic acid isolation kit

(thermMag). Nucleic acid extraction performed using the thermo magmax microbiome ultra nucleic acid isolation kit. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Time proportional sample

(timePr). A time proportional composite sample generally collected by an autosampler. Use collectionPeriod to describe how many hours the sample was taken. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Total Kjeldahl Nitrogen

(tkn). A measure of Total Kjeldahl Nitrogen; ie. the sum of nitrogen bound in organic substances, nitrogen in ammonia (NH₃-N) and in ammonium (NH₄⁺-N) in the chemical analysis of soil, water, or waste water. *Part Type:* measurements *Domain:* Che *Speciment ID:* siSaSpecimenSet

Parts

Compartment Set: waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Total Nitrogen

(tn). Total nitrogen concentration, as N. *Part Type:* measurements *Domain:* Che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Time Proportional 24hr sample

(tp24s). A time proportional 24-hour composite sample generally collected by an autosampler. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* col-Grp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Translation look-up table

(translations). Look up table for translations of the description, label, and instruction for all parts. The default language if a translation is not specified is English. *Part Type*: table *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: lookup *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Translation table column order

(translationsOrder). Specifies the order of the columns in the Translation table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Tranalation header req

(translationsRequired). Required headers in the Translations table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum*

Parts

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 2.0.0. Last updated: 2.0.0.

trizol, zymo mag beads w/ zymo clean and concentrator

(trizol). Nucleic acid extraction performed using the trizol, zymo mag beads w/ zymo clean and concentrator method. *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

TRUE

(true). Boolean data type = TRUE Use only “TRUE” (case sensitive) *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 4 *Maximum Length:* 4 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Total solids concentration

(ts). Total solids concentration *Part Type:* measurements *Domain:* Che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature

Ontology Reference: NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Total suspended solids

(tss). Total suspended solids *Part Type:* measurements *Domain:* Che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Concentration of total suspended solids

(tssConc). Total suspended solids concentration of the wastewater. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Turbidity

(turb). A measure for the turbidity of water, or a liquid sample, quantifying the relative clarity or opcaity of a liquid. Highly indicative of water quality. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* turbidityClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* turbidityUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Turbidity unit set

(turbidityUnitSet). Unit set for turbidity measurements. *Part Type:* unitSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* turbidityUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Tube weight empty

(tWeighE). The weight of the tube used for analysis while full of analyte. *Part Type:* measurements *Domain:* Phy *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0

Maximum Value: NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Tube weight full

(tWeighF). The weight of the tube used for analysis while empty. *Part Type:* measurements *Domain:* Phy *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Type other depreciate

(tyOtDep). Description for other type of sample not listed in *Part Type:* methods *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Type other depreciate2

(tyOtDep2). Description of the site when the site is not listed. See site-Type. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:*

Parts

NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 100 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

University campus

(uCampus). University or college campus - comprising an entire campus or part of a campus. See also 'Higher Education Dormitory'. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Undisclosed

(undisc). A value has been recorded, however it was not included in the dataset. Missing value indicator for data that has not been shared or disclosed with outside parties. *Part Type*: missingness *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Unique sample

(unique). A unique sample, not duplicated. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Unit

(unit). The units of a measurement. Different units that are used to describe measurement values. Units are combined into UnitSetIDs. Please requests new units by creating an issue in the ODM repository. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Unitless measure

(unitless). A unit for unitless measures, like pH. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

Parts

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Unitless unit set

(unitlessUnitSet). Unit set for measurements that does not have units. *Part Type*: unitSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: unitlessUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Units

(units). The unit of the measurement. Every measurement must have a unit. Meaning, ‘units’ are a mandatory field whenever a measurement is recorded. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Unit sets

(unitSets). Sets of units. Contains units that are associated with part types such as measures. Sets are lists of input values. Each input value within a set is a part (has a partID) that can be reused between sets. For example, a set A could contain the values of ‘yes’ and ‘no’; and set B

could contain the values of ‘yes’, ‘no’, and ‘maybe.’ partTypes with sets include aggregations, compartments, missingness, quality flags, specimens, and units. Measures with categories and dataTypes have their own specific sets (catSetID) because their input values are unique and are not reused. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ultracentrifugation

(untracent). An idenfication of a set of units that can be used for a measure or method. *Part Type:* categories *Domain:* Che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval

(upperCI95). Ultracentrifugation Should typically be linked to a lower limit, and some indication of data spread (ie. standard deviation), and a mean/median measure. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:*

Parts

varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal
Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Upstream sites

(upstream). Specifies the the upper limit of a 95% confidence interval. A catch-all category for upstream wastewater sampling sites, including major sewer pipelines and pumping stations. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Micro-Siemens per centimetre

(uSCm). A general site type for upstream wastewater sampling sites. *Part Type:* units *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

U site measure ID

(uSiteMeasureID). Micro-siemens per centimetre. Deprecated in version 2, please do not use. MeasureID will be associated with site, with site spec-

ified as the specimen as well. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 12 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

v2930l delta-variant gene target

(v2930l). Unique identifier for each measurement for a site. *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Validation study

(validationStudy). v2930l delta-variant gene target *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Value

(value). A measure or sample taken for a validation study. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 28 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Variable character data type

(varchar). Value of a measure, observation or attribute. Only used for the dictionary entries of parts. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Frequency of variants detected

(varFreq). The data type for variable character data. Measured by sequencing. *Part Type:* units *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum*

Value: 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Variants

(variant). A description of the frequency of a variants detected. *Part Type:* classes *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* variant *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population with 1 dose of vaccine

(vax1). Measures and methods relating to variants. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population with 2 doses of vaccine

(vax2). A measure of the population with a single dose of vaccine. *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population with 3 doses of vaccine

(vax3). Population with 2 doses of vaccine *Part Type:* units *Domain:* Bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ODM V1 category

(version1Category). Population with 3 doses of vaccine *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ODM V1 location

(version1Location). ODM V1 category *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ODM V1 table

(version1Table). ODM V1 location *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ODM version 1 to 2 changes

(version1to2Changes). ODM V1 table *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ODM V1 variable

(version1Variable). Changes for part between ODM v1 and V2 *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass

Parts

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous viruses group

(virusMisc). ODM V1 variable The miscallenous viruses are often measured for normalization or standardization and have only one measure or method. When a virus has many measures, they will be given their own groupID in updated dictionary versions. *Part Type:* groups *Domain:* Bio *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Size in volume

(vol). A group of measures/methods related to miscellaneous viruses. *Part Type:* measurements *Domain:* Phy *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Volume proportional sample

(volPr). Total volume of water or sludge sampled. Used for water or air testing *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: airWaterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Volume unit set

(volumeUnitSet). A volume proportional sample generally collected by an autosampler. *Part Type*: unitSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: volumeUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Volatile suspended solids

(vss). Unit set of volume measurements. *Part Type*: measurements *Domain*: Che *Speciment ID*: siSaSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum*

Parts

Length: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Volatile suspended solids - ignition

(vssIg). Volatile suspended solids *Part Type*: measurements *Domain*: Phy *Speciment ID*: saSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: measGrp *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment

(wat). A water quality measure, captured by via the loss on ignition of the mass of measured total suspended solids. This ignition generally takes place in an oven at a temperature of 550 °C to 600 °C. It represents the amount of volatile matter present in the solid fraction of the measured solution. *Part Type*: compartments *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 1 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Water depreciated

(water). A measure or observation made from a substance in the water, including wastewater. *Part Type*: categories *Domain*: naDomain *Speciment*

ID: naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* col-Grp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Water compartment set

(waterCompartmentSet). Non-wastewater, coming from any kind of water body. *Part Type:* compartmentSets *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* waterCompartmentSet *Group:* na-Group *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Wastewater Temperature

(watTemp). A compartment set for measures and methods in the water compartment. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* tempClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Weather

(weath). Temperature of the wastewater. *Part Type*: measurements *Domain*: Phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: weatherClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: weathCatSet *Unit Set*: unitlessUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Weather set

(weathSet). A measure for weather captured through qualitative categories; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: weathCatSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Whole genome Sequencing

(wgs). A set of the valid qualitative categories for the qualitative weather measure. *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*:

NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

WHO nomenclature

(who). Specifies the whole genome sequencing (WGS) strategy for genetic sequencing *Part Type*: nomenclatures *Domain*: Bio *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Wide 95% interval

(wI). Specifies variant or genetic nomenclature as set out by the World Health Organization. If there are multiple issues to flage, use the multiple issues flag and explain the details in the notes section. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Wide name table

(wideNames). The 95% interval is too wide, low confidence in results. *Part Type:* tables *Domain:* naDomain *Specimen ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Wind Speed

(wind). The table for wide names. *Part Type:* measurements *Domain:* Phy *Specimen ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* speedClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* windSpeedUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Wind speed unit set

(windSpeedUnitSet). A measure for wind speed; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources. *Part Type:* unitSet *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* airCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* windSpeedUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum*

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measure identification (wide table)

(wMeasureID). Unit set for wind speed measurements. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Sample stored incorrectly

(wrongStorage). Unique identifier for wide table only. Use when all measures are performed on a single sample at the same time and same laboratory. Suggestion: siteID_sampleID_LabID_reportDate_ID. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored at wrong temperature

(wrongTemp). Sample was stored inappropriately. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*

Parts

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

World Sphere

(wSphere). Sample was stored at an inappropriate temperature. *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

W-SPHERE notes

(wSphereNotes). World Sphere header *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Wastewater treatment plant

(wwtp). W-Shere notes Used as a catch-all category for industrial or municipal wastewater treatment plants that carry combined or sanitary sewage,

or lagoons. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Back-up wastewater treatment site

(wwtpBack). A general site type for wastewater treatment plants. Sites of this type should ideally should be paired via parent site ID with the general site information. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Wastewater treatment plant designed capacity

(wwtpCap). Indicates a wastewater treatment plant where multiple labs are or were sampling, and where secondary samples are being taken taken. Maps to units of either million gallons per day, or population equivalents. *Part Type:* measurements *Domain:* Phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* siteFeat *Class:* capacity-Class *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* capacityUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA

Parts

Maximum Length: NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Industrial wastewater treatment plant

(wwtpInd). A measure for the designed capacity of a wastewater treatment plant, in terms of how much water it was designed to be able to hold and process. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for combined sewage

(wwtpMuC). Industrial wastewater treatment plant *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for sanitary dewage only

(wwtpMuS). Municipal wastewater treatment plant for combined sewage *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

y505h omicron-variant gene target { #y505h }

(y505h). Municipal wastewater treatment plant for sanitary sewage only *Part Type*: measurements *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

zymo environ water rna kit/ zymo environ water rna kit (cat. r2042)

(zymoEnv). y505h omicron-variant gene target *Part Type*: categories *Domain*: Bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

zymo quick-rna fungal/bacterial miniprep #r2014

(zymoQuick). Nucleic acid extraction performed using the zymo environ water rna kit/ zymo environ water rna kit (cat. r2042). *Part Type:* categories *Domain:* Bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sets

ODM version 2.0.0-rc.3.16.

Sets are lists of categories and other parts. For example, `Collection category set` ([collectionCatSet](#)) contains a list of the different collection types including: `Moore swab`, `Grab sample`, `Flow proportional sample`, etc.

Sets can have parts that are reused. For example, there are different quality measures for sample collection (e.g. `Sample not collected`), PCR measures (`Inhibition preset`), and sequencing measures (`Low depth of coverage`), but all quality sets share the `No quality concerns` indicator.

Aggregations

(aggregations). Statistical measures used to report a measure. Each measure reported as a number should be reported with an aggregation. The default aggregation is `single`, which corresponds to a single measurement.

AND Boolean aggregation (`andBoo`). “AND” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (`floMean`). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Geometric mean (gmn). Geometric mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum Value (maxValueAgg). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arithmetic mean (me). Arithmetic mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Median (med). Median. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized arithmetic mean (menr). Arithmetic mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum Value (minValueAgg). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation not applicable (naAggr). Not applicable for aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Normalized geometric mean (ngmn). Geometric mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

OR Boolean aggregation (orBoo). “OR” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation deprecated (otherDep). Other aggregation method. Add description to aggregationOther Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean (pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scales

PMMoV-normalized mean (pmmovNorm). Mean measure normalized to amount of PMMoV. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard Deviation (sd). Standard deviation. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized Standard Deviation (sdn). Standard deviation, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Ultracentrifugation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scales

(aggregationScales). The scale of an aggregation set. Aggregation scales include quantitative and qualitative.

Aggregation scale not applicable (naAggrScale). Not application for aggregation sets. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation scale (othAgg). Specifies an aggregation scale that can't be described as either "qualitative" or "quantitative". Status: active. First released: 2.0.0. Last updated: 2.0.0.

Qualitative (qualAggScale). The qualitative aggregation scale. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quantitative (quantAggScale). The "quantAggScale" aggregation scale. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

See unit (aggregation scales) (seeUnitAggScale). A value for aggregation scale to be used when the aggregation scale depends on that specified in the parts list entry for the unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation sets

(aggregationSets). Sets of aggregations. Examples of aggregation sets include logarithm, linear and boolean.

Boolean aggregation set

(booleanAggrSet)Aggreagation set for boolean. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

AND Boolean aggregation (andBoo). “AND” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

OR Boolean aggregation (orBoo). “OR” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Linear scale aggregation set

(linearAggrSet)The aggregation set that contains all aggregations that exist on the linear scale. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (floMean). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum Value (maxValueAgg). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arithmetic mean (me). Arithmetic mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Median (med). Median. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized arithmetic mean (menr). Arithmetic mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum Value (minValueAgg). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean (pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-normalized mean (pmmovNorm). Mean measure normalized to amount of PMMoV. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard Deviation (sd). Standard deviation. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized Standard Deviation (sdn). Standard deviation, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Ultracentrifugation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Logarithmic scale aggregation set

(logAggrSet) The aggregation set that contains all aggregations that exist on the logarithmic scale (rather than linear or qualitative). Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (floMean). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Geometric mean (gm). Geometric mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum Value (maxValueAgg). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum Value (minValueAgg). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized geometric mean (ngm). Geometric mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean (pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-normalized mean (pmmovNorm). Mean measure normalized to amount of PMMoV. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Ultracentrifugation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation set not applicable

(naAggrSet)Not application for aggregation sets. Used for parts for which an aggregation set value is not applicable or appropriate. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Aggregation not applicable (naAggr). Not applicable for aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation set

(otherAggrSet)Aggregation set used for unitless measures. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (floMean). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Geometric mean (gmn). Geometric mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum Value (maxValueAgg). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arithmetic mean (me). Arithmetic mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Median (med). Median. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Normalized arithmetic mean (menr). Arithmetic mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum Value (minValueAgg). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation not applicable (naAggr). Not applicable for aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Normalized geometric mean (ngmn). Geometric mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Standard Deviation (sd). Standard deviation. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized Standard Deviation (sdn). Standard deviation, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Ultracentrifugation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Attributes

(attributes). An attribute of a domain, specimen, compartment, or table. Attributes are one of the three main components or ways to report environment surveillance data. The other two components are measures and methods.

Analysis date (aDate). Date the measurement was performed in the lab. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date end (aDateEnd). Date the measurement or analysis was completed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date start (aDateStart). Date the measurement or analysis was started. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Address Line 1 (addL1). Line 1 (the street name, number and direction) for a given address. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Address Line 2 (addL2). Line 2 (the unit number) for a given address. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Address table (addresses). The table that contains information about addresses. Addresses can be recorded for sites, organizations or contacts (individuals). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Address ID (addressID). A unique identifier for an address. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Administrative regions (admRegLevel). Administrative regions Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation (aggregation). Statistical measures used to report a measure. Each aggregation has a corresponding value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation scale (aggregationScale). A scale used for an aggregation. Only applicable for measures and units. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation set (aggregationSet). The aggregation set for a unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Alias depreciated (aliasDep). ID of an assay that is the same or similar. A comma separated list. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Alias ID depreciated (aliasIDDep). Alias id Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Sets

Access to all org (allOrgs). If this is ‘no’, this data will not be available to any partner organization. If missing, data will be available to the all organizations. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Assay method id default depreciated (amDefDep). Used as default when a new measurement is created for this lab. See ID in AssayMethod table. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Assay ID (asID). Links with the AssayMethod used to perform the analysis. Use instrument.ID for measures that are not viral measures. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Measure, method, or attribute (mmaSet). The set for a measure, method, or attribute. Only applicable for categorical fields, and the categories that populate them. Status: active. First released: 2.0.0. Last updated: 2.0.0.

City (city). The city where a site or organization is located; part of the address. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Class (class). A unique identifier for a class, which is akin to a sub-group; it’s a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection date time (collDT). For grab samples this is the date, time and timezone the sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Collection date time end (collDTEnd). For integrated time average samples this is the date, time and timezone the sample was finished being taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Collection date time start (collDTStart). For integrated time averaged samples this is the date, time and timezone the sample was started being taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Collection number (collNum). The number of subsamples that were combined to create the sample. Use NA for continuous, proportional or passive sampling. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection number and period (collNumPer). Composite collection number.period. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection period (collPer). Collection period. The time period over which the sample was collected, in hours. Alternatively, use collectionStart and collectionEnd. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample collection type (collType). The type of collection. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Column name depreciated (columnNameDep). Name for the column Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Compartment (compartment). The attribute identifying the substance from which where a sample was taken. For more information, see partID = compartment. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Contact ID (contactID). A unique identifier for a given contact person. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Contact name (contactName). Contact person or group, for the lab. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Country (country). The country where a site or organization is located; part of the address. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Countries or sovereign states (countryLevel). Countries or sovereign states Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Districts, counties, regions (countyLevel). Districts, counties, regions
Status: active. First released: 2.0.0. Last updated: 2.0.0.

Covid-19 population measurement date (cphDate). date of reporting
for covid-19 measure. Status: depreciated. First released: 1.0.0. Last
updated: 2.0.0.

Cphd ID (cphid). Unique identifier for the table. Status: depreciated.
First released: 1.0.0. Last updated: 2.0.0.

Custody Contact ID (custodyCont). A unique identifier for data cus-
todians. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Data custodian ID (custodyID). The data custodian of a database. Sta-
tus: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset creation date (datasetDate). Specifies the date a given dataset
was created. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset ID (datasetID). The name of the dataset that stores information
for MeasureReport, SampleReport and other reporting tables. Status:
active. First released: 2.0.0. Last updated: 2.0.0.

Date (date). Date Status: depreciated. First released: 1.0.0. Last up-
dated: 2.0.0.

Description (descr). A detailed description as an attribute for describing
results or program elements. Status: active. First released: 2.0.0. Last
updated: 2.0.0.

Access to details deprecated (detailsDep). More details on the existing
confidentiality requirements of this measurement. Status: depreciated.
First released: 1.0.0. Last updated: 2.0.0.

Domain (domain). Domain is the highest level of describing of a measure.
Status: active. First released: 2.0.0. Last updated: 2.0.0.

Contact email (email). Contact e-mail address, for the lab. Status:
active. First released: 1.0.0. Last updated: 2.0.0.

File location of polygon (fileLocation). The location of the file containing the geometry of the polygon. Status: active. First released: 1.0.0. Last updated: 2.0.0.

First name of contact (firstName). Specifies the first name of a given contact. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Fraction analyzed (fraction). Fraction of the sample that is analyzed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Funder Contact ID (funderCont). A unique identifier for funders. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Funding agency ID (funderID). The funding agency of the dataset. Status: active. First released: 2.0.0. Last updated: 2.0.0.

European Petroleum Survey Group Coordinates (geoEPSG). The unique EPSG code specifying a given geospatial area. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Latitude (geoLat). Geographical location, latitude in decimal coordinates, ie.: (45.424721) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Longitude (geoLong). Geographical location, longitude in decimal coordinates, ie.: (-75.695000) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Type of geography (geoType). Type of geography that is represented by the polygon. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Well-known text (geoWKT). Well-known text of the polygon Status: active. First released: 1.0.0. Last updated: 2.0.0.

Group (group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Index (index). Index number in case the measurement was taken multiple times. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Instrument ID (instrumentID). A unique identifier for an instrument. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Instrument Type (insType). Type of instrument used to perform the measurement. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Describe other instrument type, if applicable (insTypeOth). Description of the instrument in case it is not listed in instrumentType. Status: active. First released: 1.0.0. Last updated: 2.0.0.

ISO639-1 (iso6391). The first part of the ISO 639 series of international standards for language codes. Part 1 covers the registration of two-letter codes. There are 183 two-letter codes registered as of June 2021. The registered codes cover the world's major languages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2B (iso6392B). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'B' specifies the bibliographic code (B code). Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2T (iso6392T). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'T' specifies the terminological code (T code). Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-3 (iso6393). A set of international standards that lists short codes for language names. ISO 639-3 extends the ISO 639-2 alpha-3 codes with an aim to cover all known natural languages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-6 (iso6396). A set of international standards that lists short codes for language names. ISO 639-6 builds off ISO639-3 with the use of four-letter codes, and allowing users to differentiate between variants of languages and language families, such as historical vs. revived versions of languages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Lab ID default deprecated (labDefDep). Used as default when a new sample is created by this reporter. See ID in Lab table. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Lab ID (labID). Unique identifier for a laboratory. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Language family (langFam). Specifies the language family of a given language for translation and language tracking purposes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Language name (langName). Specifies the name of the language in roman alphabet characters for translation and language tracking purposes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Last edited (lastEdited). The date the entry was last updated. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Last name of contact (lastName). Specifies the last name of a given contact. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Link (li). Link to an external reference that describes the geometry of the polygon. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

License (license). The license of a dataset. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Access to local ha (localHA). If this is ‘no’, the, data will not be available to local health authorities. If missing, data will be available to local health authorities. Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

Manufacturer (manufacturer). Manufacturer of an instrument. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure (measure). A measurement or observation of any substance including a biological, physical or chemical substance. An measurement is recorded from a specimen (i.e. a site, sample, person or population). Status: active. First released: 1.0.0. Last updated: 2.0.0.

Report ID (measureRepID). Unique identifier for a measurement. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Report set ID (measureSetRepID). Unique identifier that links together a group of related measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure license (measureLic). Specifies the access and use licensing for a given single measurement. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Method (method). A procedure for collecting a sample or performing a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Model (model). Model number or version of the instrument. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Municipalities or communes (municipalLevel). Municipalities or communes Status: active. First released: 2.0.0. Last updated: 2.0.0.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Native Name (natName). The native name of the language, i.e. what the language is called by its speakers. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Local administrative units or neighborhoods (neighborLevel). Local administrative units or neighborhoods Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclature (nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Notes (notes). A note used to describe details that are not captured in other attributes Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Organization level (orgLevel). The geographic level of an organization. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization sector (orgSector). The sector of an organization Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization Type (orgType). Specifies the type or purpose of a given organization. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample origin (origin). An attribute of a sample specifying the origin. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to other prov deprecated (otherProvDep). If this is 'no', this data will not be available to other data providers not listed before. If missing, data will be available to other data providers not listed before Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

Sets

Parent dataset ID (parDatasetID). The datasetID that is a parent to another datasetID. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Parent sample ID (parent). If this sample has been pooled into one big sample for analysis this indicates the sample of the larger pooled sample. Status: depreciated. First released: 1.1.0. Last updated: 2.0.0.

Parent Site ID (parSiteID). The siteID that is a parent to another siteID; usually the sub-site will be contained in the parentSite. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Postal or Zip Code (pCode). The zip code or postal code for a given address, specifying a specific geographic area. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to phac (phac). If this is 'no', the data will not be available to employees of the Public Health Agency of Canada - PHAC. If missing, data will be available to employees of the Public Health Agency of Canada - PHAC. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Contact phone (phone). Contact phone number, for the lab. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Polygon ID (polygonID). Unique identifier for the polygon. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Polygon Population (polyPop). An attribute of a polygon, which specifies the population of that polygon. A rough estimate of the number of human residents. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Pooled (pooled). Is this a pooled sample, and therefore composed of multiple child samples obtained at different sites Status: active. First released: 1.0.0. Last updated: 2.0.0.

Type of date for case reporting depreciated (popDateTypeDep). Type of date used for confirmed cases. Typically, report or episode are

reported. Onset and test date is not usually reported within aggregate data. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Population Served (popServ). An attribute of a site, which specifies the population of/population served by a given site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID (protocolID). A unique identifier for a given protocol. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID container (protocolIDContainer). Unique identifier for the protocol and the steps and other protocols that make it up. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID object (protocolIDObj). The object of the relationship between one protocol and another protocol or step, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID subject (protocolIDSub). The subject of the relationship between one protocol and another protocol or step, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol version (protocolVersion). Specifies the version of a method set. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to prov ha (provHA). If this is ‘no’, this data will not be available to provincial health authorities. If missing, data will be available to provincial health authorities. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Pretreatment description - depreciated (ptDescDep). If preTreatment then describe the treatment that was performed. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Access to public (public). If this is ‘no’, this data will not be available to the public. If missing, data will be available to the public. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Sets

Purpose (purpose). The reason the measure or sample was taken. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality flag (qualityFlag). A field for reporting any quality concerns - of lack thereof - for a sample or measure. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality report (qualityID). A unique identifier for a given quality report. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality set (qualitySetID). The quality set that corresponds to a given part. Only applicable for samples and measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Date sample recieved (recDate). The date the sample was received at the laboratory for analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Primary reporting authority (repOrg1). The primary or most responsible authority for routine surveillance reports or findings, or where the site is located. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Secondary reporting authority (repOrg2). The secondary, additional or alternative authority for routine surveillance reports or findings, or where the site is located. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Reportable (reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Report date (reportDate). The date a measure was reported. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Reporter ID deprecated (reporterIDDep). Unique identifier for the person or organization that is reporting the data. Status: deprecated. First released: 1.1.0. Last updated: 2.0.0.

Replicate Type (repType). Attribute of a sample, specifying whether the sample is unique, or a replicate. And if it is a replicate, what type of replicate is it? Status: active. First released: 2.0.0. Last updated: 2.0.0.

Role of contact (role). Specifies the organizational role of a given contact. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample material (saMaterial). Type of sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample ID (sampleID). Unique identifier for a sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sample ID object (sampleIDObject). Populated by **sampleID** - this specifies the object of a sample relationship. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample ID subject (sampleIDSubject). Populated by **sampleID** - this specifies the subject of a sample relationship. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationship (sampleRelID). Attribute for specifying the relationship between a sample subject and object. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample shed (sampleShed). A geographic area, physical space, or structure. A sample is taken from a sampleshed for a representative measurement of a substance(s). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to self (self). If this is 'no', this data will not be shown on the portal when this reporter logs in. If missing, data will be available to

Sets

this reporter. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Date sample was sent (sentDate). The date the sample was sent for analyses at a laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Set (set). The unique identifier of a value set. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Set Value (setValue). The partID for the set value or category. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Severity indicator (severity). An indicator of the severity or seriousness of a quality flag. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sewer network file link depreciated (sewerNetworkFileBLOB). Link to a file that has any detailed information about the sewer network associated with the site (any format). Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Site ID default (siteDef). Used as default when a new sample is created by this reporter. See ID in Site table. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Site ID (siteID). Unique identifier for the location where a sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Site measure ID (siteMeasureID). Unique identifier for wide table only. Use when all measures are performed on a single sample. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Site Type (siteType). Type of site or institution where sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer network file link depreciated 2 (snfl). Link to a file that has any detailed information about the sewer network associated with the site

(any format). Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

Source Protocol ID (sourceProtocol). A protocol that served as a basis for another protocol. This may be due to version changes, or a referencing a protocol that was later adapted. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step source ID (sourceStep). Specifies the protocol step which serves as a basis for a given protocol step. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen (specimen). The substance or thing upon which the observation was made. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Departments, states, or provinces (stateProvLevel). Departments, states, or provinces Status: active. First released: 2.0.0. Last updated: 2.0.0.

State, Province, or Region (stateProvReg). The state, province, or region where a site or organization is located; part of the address. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Protocol Step ID (stepID). The unique identifier for a specific protocol step. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Step ID Object (stepIDObj). The object of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Step ID Subject (stepIDSub). The subject of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Method step parent ID (stepProvenanceID). A method step that served as a basis for another method step. This may be due to version changes,

Sets

or a referencing a method set that was later adapted. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step version (stepVer). Specifies the version of a given protocol step. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer network file blob (storTempDef). A file blob that has any detailed information about the sewer network associated with the site (any format). Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Summary (summ). Short description of the assay and how it is different from the other assay methods. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Type other depreciate2 (tyOtDep2). Description of the site when the site is not listed. See siteType. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Unit (unit). The units of a measurement. Status: active. First released: 1.0.0. Last updated: 2.0.0.

U site measure ID (uSiteMeasureID). Micro-siemens per centimetre. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Value (value). A measure or sample taken for a validation study. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure identification (wide table) (wMeasureID). Unit set for wind speed measurements. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Boolean value set

(booleanSet). Set that contains the valid possible values for a boolean measure (TRUE or FALSE).

TRUE (true). Boolean data type = TRUE Status: active. First released: 1.1.0. Last updated: 2.0.0.

FALSE (false). Boolean data type = FALSE Status: active. First released: 1.1.0. Last updated: 2.0.0.

Classes

(classes). A class is a collection of one or more related measures or methods within a group. A group can have one more more classes to describe different parts or the class.

Alleles (allele). Measures and methods related to alleles. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Conductivity (conductivity). Measures and methods related to conductivity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Diseases (human) (disease). Measure and methods related to disease or infection in humans. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Flow (flow). Measures and methods related to flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Humidity (human). Measures and methods related to humidity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dictionary tables (lookup). Tables to describe and support the ODM dictionary. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Mutations (mutation). Measures and methods related to mutations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Class not applicable (naClass). Class is not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Physical (physical). Measures and methods related to generic physical properties. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Precipitation (precipitation). Measures and methods related to precipitation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Program decription tables (programDescr). Tables used to describe surveillance and testing programs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Results tables (results). Tables used to record samples, measures and quality reports. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Genetic sequences (sequence). Measures and methods related to genetic sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard concentrations (standardConc). Measures and methods relating to standard concentrations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard curve (standardCurve). Measures and methods relating to generating or recording a standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Temperature (temperature). Measures and methods related to temperatures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Variants (variant). A description of the frequency of a variants detected. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample collection set

(collectSet). Methods for collection samples. Sample collection methods include water, air, and surface. See the attribute **Compartment set** of the sample collection method.

Area proportional sample (areaPr). An area proportional sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

COSCa ball (cosca). COSCa passive sampling device. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Flow proportional sample (flowPr). A flow proportional composite sample generally collected by an autosampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Grab sample (grb). A single large representative grab sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Moore swab passive sample (moorSw). Moore swab passive sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Surface swab (surfSw). Surface swab. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Time proportional sample (timePr). A time proportional composite sample generally collected by an autosampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Volume proportional sample (volPr). Total volume of water or sludge sampled. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Composite sample - archival (comp). A composite sample, usually generated by an autosampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Compartments

(compartments). The substance from which a sample was taken. For example, wastewater has component of “water”. Compartments are attributes of measures, methods, units, and aggregations.

Sets

Human compartment (hum). A measure or observation made about a human. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment not applicable (naCompartment). Compartment not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment (wat). A water quality measure, captured by via the loss on ignition of the mass of measured total suspended solids. This ignition generally takes place in an oven at a temperature of 550 °C to 600 °C. It represents the amount of volatile matter present in the solid fraction of the measured solution. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment sets

(compartmentSets). Sets of compartments. Compartment sets are used to identify when a measure can be recorded for more than one compartment. For example, SARS-CoV-2 can be measured in people (humans), water, surface, or air.

Air compartment set

(airCompartmentSet)A compartment set for measures and methods in the air compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Air and surface compartment set

(airSurfaceCompartmentSet)A compartment set for measures and methods in the air or surface compartments. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Air and water compartment set

(airWaterCompartmentSet)A compartment set for measures and methods in the air or water compartments. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Any compartment set

(anyCompartmentSet)A compartment set for measures and methods in any compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment (wat). A water quality measure, captured by via the loss on ignition of the mass of measured total suspended solids. This ignition generally takes place in an oven at a temperature of 550 °C to 600 °C. It represents the amount of volatile matter present in the solid fraction of the measured solution. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Human compartment set

(humanCompartmentSet)A compartment set for measures and methods in the human compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Human compartment (hum). A measure or observation made about a human. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment set not applicable

(naCompartmentSet)Compartment not applicable. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Compartment not applicable (naCompartment). Compartment not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment set

(surfaceCompartmentSet)A compartment set for measures and methods in the surface compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface and water compartment set

(surfaceWaterCompartmentSet)A compartment set for measures and methods in the surface or water compartments. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment (wat). A water quality measure, captured by via the loss on ignition of the mass of measured total suspended solids. This ignition generally takes place in an oven at a temperature of 550 °C to 600 °C. It represents the amount of volatile matter present in the solid fraction of the measured solution. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment set

(waterCompartmentSet)Non-wastewater, coming from any kind of water body. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Water compartment (wat). A water quality measure, captured by via the loss on ignition of the mass of measured total suspended solids. This ignition generally takes place in an oven at a temperature of 550 °C to 600 °C. It represents the amount of volatile matter present in the solid fraction of the measured solution. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Data types

(dataTypes). The data type for a part. Data types used in the ODM include: varchar, boolean, float, category, date, time, datetime, url, email. dataType corresponds to the entry or cell within a data table, most commonly within a report table. If the data entry has a unit, then the dataType corresponds to the unit and the dataType refers to 'unit'. If the data entry is a category, then the dataType refers to 'category'. All categories are varchar. Otherwise the dataType is identified within the part entry. TBA: dataType for dictionary.

Binary Large Object (BLOB) data type (blob). The data type for Binary Large Object (BLOB) data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Boolean data type (boolean). The data type for boolean/binary data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Categorical data type (categorical). The data type for categorical data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Datetime data type (datetime). The data type for date and time data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Float data type (float). The data type for float data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Integer data type (integer). The data type for integers. Status: active. First released: 2.0.0. Last updated: 2.0.0.

See unit (data type) (seeUnitData). The data type entry when the data type for a given entry is specified by the unit used. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Variable character data type (varchar). Value of a measure, observation or attribute. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Domains

(domains). There are three domain types: biologic (i.e. Covid-19, chemical (i.e. nitrogen), physical measure (temperature).

All domains (allDo). Domain that specifies that it could apply to all domains; biological, chemical, and physical. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Biologic (bio). A living organism or biological substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Chemical (che). A chemical compound. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Domain not applicable (naDomain). Not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Physical property (phy). A physical property or object not characterized by life or chemistry. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nucleic Acid Extraction set

(extractSet). set used for storing all the valid category values for the nucleic acid extraction method.

4s method (extract4s). Nucleic acid extraction performed using the 4s method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Amicon filter, extract with MP96 (amiMP96). Nucleic acid extraction, usually used for the liquid fraction of a wastewater sample, using

Sets

amicon filtration and using an MP96 for final extraction. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Bead beating, extract with MP96 (bbMP96). Nucleic acid extraction, usually used for solid fraction of a wastewater sample, using bead beating and using an MP96 for final extraction. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Chemagic viral dna/rna 300 kit (chemVir). Nucleic acid extraction performed using the chemagic viral dna/rna 300 kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nuclisens automated magnetic bead extraction kit (nucAuto). Nucleic acid extraction performed using the nuclisens automated magnetic bead extraction kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nuclisens manual magnetic bead extraction kit (nucManu). Nucleic acid extraction performed using the nuclisens manual magnetic bead extraction kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Phenol chloroform (phenCl). Nucleic acid extraction performed using phenol chloroform. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega automated tna kit (promAuto). Nucleic acid extraction performed using the promega automated tna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega ht tna kit (promHt). Nucleic acid extraction performed using the promega ht tna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega manual tna kit (promManu). Nucleic acid extraction performed using the promega manual tna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega wastewater large volume tna capture kit (promWW). Nucleic acid extraction performed using the promega wastewater large volume tna capture kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep dna/rna kit (qgDNARNA). Nucleic acid extraction performed using the qiagen allprep dna/rna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen ez1 virus mini kit v2.0 (qgEz1). Nucleic acid extraction performed using the qiagen ez1 virus mini kit v2.0. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep powerfecal dna/rna kit (qgPwrFecal). Nucleic acid extraction performed using the qiagen allprep powerfecal dna/rna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep powerviral dna/rna kit (qgPwrViral). Nucleic acid extraction performed using the qiagen allprep powerviral dna/rna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen powerwater kit (qgPwrWtr). Nucleic acid extraction performed using the qiagen powerwater kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen qiaamp buffers with epoch columns (qgQiAmp). Nucleic acid extraction performed using the qiagen qiaamp buffers with epoch columns. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy kit (qgRneasy). Nucleic acid extraction performed using the qiagen rneasy kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy powermicrobiome kit (qgRneasyPwr). Nucleic acid extraction performed using the qiagen rneasy powermicrobiome kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Resuspend COSCa filter collection (resCosca). Nucleic acid extraction via resuspension of a sample collected using the COSCa ball method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

thermo magmax microbiome ultra nucleic acid isolation kit (thermMag). Nucleic acid extraction performed using the thermo magmax microbiome ultra nucleic acid isolation kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

trizol, zymo mag beads w/ zymo clean and concentrator (trizol). Nucleic acid extraction performed using the trizol, zymo mag beads w/ zymo clean and concentrator method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

zymo environ water rna kit/ zymo environ water rna kit (cat. r2042) (zymoEnv). y505h omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

zymo quick-rna fungal/bacterial miniprep #r2014 (zymoQuick). Nucleic acid extraction performed using the zymo environ water rna kit/ zymo environ water rna kit (cat. r2042). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample fraction set

(fractionSet). set for the fraction of the sample (solid, liquid, etc.).

Liquid fraction (liq). Liquid fraction of a sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Mixed/homogenized sample (mix). Mixed or homogenized sample or fraction analyzed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Solid fraction (sol). Solid fraction of a sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Geographic set

(geoTypeSet). set for different type of geographic components.

Sewer Network Health Region (hlthReg). Health region served by the sewer network Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer catchment area (swrSet). Sewer catchment area. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Groups

(groups). A collection of related measures.

Miscellaneous bacteria group (bactMisc). A group of measures/methods related to miscellaneous bacteria. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection group (colGrp). A group of measurement-like attributes related to sample collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measurement group (measGrp). A group of measures that cannot be otherwise categorized. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Miscellaneous attribute group (miscAttr). A group of miscellaneous measurement-like attributes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous measure group (miscMeas). A group of measures that cannot be otherwise categorized. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Group not applicable (naGroup). Used when there is no group; ie. group is not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population group (popGrp). A group of measures/methods related to population-level factors. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Processing group (procGrp). A group of measures/methods related to processing samples for analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

SARS-CoV-2 (sarsCov2). A group of measures/methods related to the SARS-CoV-2 virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site features group (siteFeat). A group of environmental measures/methods and those pertaining to the features of a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous viruses group (virusMisc). ODM V1 variable Status: active. First released: 2.0.0. Last updated: 2.0.0.

Instrument set

(insTypeSet). List of instruments that are used for measures and methods

Atline Analyzer (aas). Atline analyzer with sampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Hand Measurement (hma). Handheld measurement analyzer. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Other instrument (instrumentTypeOther). Type of instrument other than those included in the PHES-ODM. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lab Analysis (ola). Offline laboratory analysis. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Online Sensor (onse). Online sensor Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measures

(measurements). The attribute to describe a part type of measures. All measures have **partType** = “measure”.

a1306s delta-variant gene target (a1306s). a1306s delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a1918v delta-variant gene target (a1918v). a1918v delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a2710t omicron-variant gene target (a2710t). a2710t omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a63t omicron-variant gene target (a63t). a63t omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a67v omicron-variant gene target (a67v). a67v omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Absolute humidity (absHum). A measure of the total mass of water vapour present in the air per volume of air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Environmental temperature (airTemp). Environmental temperature. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arrival temperature (arTemp). The temperature of a sample upon arrival at the laboratory for analysis. Status: active. First released: 1.0.0. Last updated: 2.0.0.

bcov (bcov). Measure of the amount of bovine coronavirus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Beta (beta). B.1.351 Status: active. First released: 2.0.0. Last updated: 2.0.0.

5-Day carbonaceous biochemical oxygen demand (bod5c). The quantity of oxygen utilized for the biochemical degradation of organic matter under standard laboratory procedures in five (5) days in the presence of a nitrification inhibitor, expressed in milligrams per litre (mg/l). Status: active. First released: 1.1.0. Last updated: 2.0.0.

5-day total biochemical oxygen demand (bod5t). 5 day total biochemical oxygen demand. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Bovine respiratory syncytial virus group (brsv). Bovine respiratory syncytial virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

BRSV-N (brsvN). bovine respiratory syncytial virus capsid protein gene region Status: active. First released: 1.0.0. Last updated: 2.0.0.

Carbon Dioxide (co2). A measure of an amount or concentraion of carbon dioxide. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Chemical Oxygen Demand (cod). Chemical oxygen demand. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Water conductivity (cond). Measurement of conductivity of sample or site. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Covid-19 (cov). Covid-19 infection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

SARS-CoV-2 measure (cov2Me). Measure the amount of SARS-CoV-2 virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

SARS-CoV-2-B.1.1.7 (covB117). Variant B.1.1.7 of the SARS-CoV-2 virus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-B.1.351 (covB135). Variant B.1.351 of the SARS-CoV-2 virus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-E (covE). SARS-CoV-2 E gene. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N1 (covN1). SARS-CoV-2 nucleocapsid gene, allele 1. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N2 (covN2). SARS-CoV-2 nucleocapsid gene, allele 2. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N3 (covN3). SARS-CoV-2 nucleocapsid gene, allele 3. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-P.1 (covP1). Variant P.1 of the SARS-CoV-2 virus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-RdRp (covRdrp). SARS-CoV-2 RdRp gene. Status: active. First released: 1.0.0. Last updated: 2.0.0.

crAssphage-N (cra). crAssphage virus capsid protein gene region Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cycle threshold or quantification cycle (Ct or Cq) (ct). Cycle thresholds in a PCR assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

d377y delta-variant gene target (d377y). d377y delta-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

d3g omicron-variant gene target (d3g). d3g omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

d63g delta-variant gene target (d63g). d63g delta-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

d796y omicron-variant gene target (d796y). d796y omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

d950n delta-variant gene target (d950n). d950n delta-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ddcov__e sars-cov-2 gene target (ddcovE). ddcov__e sars-cov-2 gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ddcov__n sars-cov-2 gene target (ddcovN). ddcov__n sars-cov-2 gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del143/145 (del143). 143 or 145 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del 157/158 (del157). 157 or 158 deletion delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

del2084/2084 (del2084). 2084 or 2084 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del212/212 (del212). 212 or 212 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del3674/3676 (del3674). 3674 or 3676 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del69/70 (del6970). 69 or 70 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Delta (delta). B.1.617.2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dilution factor (dilFact). Specifies the extent to which a sample or aliquot was diluted prior to analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Point dillutions (dillute). Exact concentration or dillutions for generating a standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

e156g delta-variant gene target (e156g). e156g delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

e484a omicron-variant gene target (e484a). e484a omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Escherichia coli (ecoli). Concentration of bacteria that are passed through the faecal excrement of humans, livestock and wildlife Status: active. First released: 2.0.0. Last updated: 2.0.0.

Efficiency (efficient). The efficiency reported for a standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rainfall (envRnF). Rainfall, i.e. amount of precipitation in the form of rain. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ground snow depth (envSnwD). Total depth of snow on the ground. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Snowfall (envSnwF). Snowfall, i.e. amount of precipitation in the form of snow. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sarbecovirus-specific E sars-cov-2 gene target (eSarbec). Sarbecovirus-specific E sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Estimated frequency of reads (estFreqReads). Estimated frequency of reads in a sequencing assay Status: active. First released: 2.0.0. Last updated: 2.0.0.

Extraction volume of sample (exvol). Extraction volume of sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Flow rate (floRate). Wastewater volumetric flow rate at the sample collection location over the 24-hr period during which the sample was collected. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Flow volume (flowVol). Volume of influent. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus measure (flu). General influenza virus measure Status: active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus A1 (fluA1). Influenza virus A1 type Status: active. First released: 2.0.0. Last updated: 2.0.0.

influenza virus A2 (fluA2). influenza virus A2 type Status: active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus B (fluB). Influenza virus B type Status: active. First released: 2.0.0. Last updated: 2.0.0.

F-Specific RNA bacteriophages (frna). A measure for amount of F-Specific RNA bacteriophages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

F-Specific RNA bacteriophages, G2 (frnaG2). A measure for amount of G2 F-Specific RNA bacteriophages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample temperature (fst). Temperature that the sample is stored at while it is being sampled. This field is mainly relevant for composite samples which are either kept at ambient temperature or refrigerated while being sampled. Status: active. First released: 1.0.0. Last updated: 2.0.0.

g215c delta-variant gene target (g215c). g215c delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

g339d omicron-variant gene target (g339d). g339d omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

g496s omicron-variant gene target (g496s). g496s omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

g662s delta-variant gene target (g662s). g662s delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gamma (gam). P.1 Status: active. First released: 2.0.0. Last updated: 2.0.0.

hep g armored rna (hepGRna). Measure of the amount Hepatitis G Armored RNA. Status: active. First released: 2.0.0. Last updated: 2.0.0.

i1566v omicron-variant gene target (i1566v). i1566v omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

i3758v omicron-variant gene target (i3758v). i3758v omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

i82t delta-variant gene target (i82t). i82t delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Inhibition measure (inhibMe). Parameter to report whether or not inhibition was detected in the sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Intercept (inter). Intercept value of the calibration curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Combined ip2 and ip4 sars-cov-2 gene target (ip2ip4). Combined ip2 and ip4 sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

k856r omicron-variant gene target (k856r). k856r omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

l452r delta-variant gene target (l452r). l452r delta-variant gene target
Status: active. First released: 2.0.0. Last updated: 2.0.0.

l981f omicron-variant gene target (l981f). l981f omicron-variant gene
target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Lambda (lamba). C.37 Status: active. First released: 2.0.0. Last up-
dated: 2.0.0.

Limit of detection (LOD) (lod). Limit of detection. Status: active.
First released: 1.1.0. Last updated: 2.0.0.

Limit of detection (LOD) - sewage (lodSewa). Limit of detection
for sewage samples. Status: active. First released: 2.0.0. Last updated:
2.0.0.

Limit of quantification (LOQ) (loq). Loq Status: active. First re-
leased: 1.1.0. Last updated: 2.0.0.

Other Measure (measOth). Other measure, not otherwise specified
in measures. Status: depreciated. First released: 1.0.0. Last updated:
2.0.0.

Other Measure Description (meOthDe). Description of other mea-
sure (measOtherDesc). Status: depreciated. First released: 1.0.0. Last
updated: 2.0.0.

Murine Hepatitis Virus (mhv). A measure for amount of murine hep-
atitis viurs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage (ms2Col). Measure of the amount of ms2 coliphage. Sta-
tus: active. First released: 2.0.0. Last updated: 2.0.0.

Mu (mu). B.1.621 Status: active. First released: 2.0.0. Last updated:
2.0.0.

murine coronavirus (muCo). Measure of the amount of murine coron-
avirus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

N sars-cov-2 gene target (n). N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Combined N1 and N2 sars-cov-2 gene target (n1n2). Combined N1 and N2 sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n211i omicron-variant gene target (n211i). n211i omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n679k omicron-variant gene target (n679k). n679k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n856k omicron-variant gene target (n856k). n856k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n969k omicron-variant gene target (n969k). n969k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Ammonium Nitrogen (nH4N). Ammonium nitrogen concentration, as N. Status: active. First released: 1.0.0. Last updated: 2.0.0.

niid_2019-ncov_n sars-cov-2 gene target (niid19). niid_2019-ncov_n sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Norovirus G1 (noroG1). Norovirus genogroup 1 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Norovirus G2 (noroG2). Norovirus genogroup 2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sarbecovirus-specific N sars-cov-2 gene target (nSarbec). Sarbecovirus-specific N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

coronavirus OC43 (oc43). Measure of the amount of coronavirus OC43. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Omicron BA.1 (omicr1). Omicron B.1.1.529 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2 (omicr2). Omicron BA.2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2.75 (omicr275). Omicron BA.2.75 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.4 (omicr4). Omicron BA.4 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.5 (omicr5). Omicron BA.5 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Orthophosphates (ophos). Ortho-phosphate concentration. Status: active. First released: 1.0.0. Last updated: 2.0.0.

ORF1a sars-cov-2 gene target (orf1a). ORF1a sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

ORF1ab sars-cov-2 gene target (orf1ab). ORF1ab sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

ORF1b sars-cov-2 gene target (orf1b). ORF1b sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Variants Other (otherV). Used for reporting variants detected for which a specific measureID doesn't already exist. The ODM team will periodically migrate collected responses in this category into set measureIDs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak (outb). Measure to indicate outbreak status. Given that when outbreaks occur, full data on case counts may not be available or completely accurate, using this measure indicates that case counts may be approximate. Status: active. First released: 2.0.0. Last updated: 2.0.0.

p100l delta-variant gene target (p100l). p100l delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p2046l delta-variant gene target (p2046l). p2046l delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p2287s delta-variant gene target (p2287s). p2287s delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p3395h omicron-variant gene target (p3395h). p3395h omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p681r delta-variant gene target (p681r). p681r delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

pH (ph). pH measurement Status: active. First released: 1.1.0. Last updated: 2.0.0.

Pseudomonas virus phi6 (phi6). Measure of the amount of pseudomonas virus phi6. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Total phosphorous (phos). Total phosphorous Status: active. First released: 1.1.0. Last updated: 2.0.0.

Total Phosphates (phostot). Total phosphates Status: active. First released: 1.1.0. Last updated: 2.0.0.

PHRED quality score (phred). PHRED Quality Score For Sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-CP (ppmv). Pepper mild mottle virus capsid protein gene region Status: active. First released: 1.0.0. Last updated: 2.0.0.

Puro virus (puro). Measure of the amount of puro virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

q19e omicron-variant gene target (q19e). q19e omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

q493r omicron-variant gene target (q493r). q493r omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

q498r omicron-variant gene target (q498r). q498r omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

q954h omicron-variant gene target (q954h). q954h omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

R squared (r2). R-squared value of the calibration curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

r203m delta-variant gene target (r203m). r203m delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

IP2 rdrp sars-cov-2 gene target (rdrpIP2). IP2 rdrp sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

IP4 rdrp sars-cov-2 gene target (rdrpIP4). IP4 rdrp sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Relative humidity (relHum). The unit of relative humidity, or the air-water mixture. The ratio of the partial pressure of water vapor in the mixture to the equilibrium vapor pressure of water over a flat surface. of pure water at a given temperature. Status: active. First released: 2.0.0. Last updated: 2.0.0.

replicate number (repNum). The replicate number for a Ct or Cq value - used for specifying out a protocol standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Respiratory syncytial virus (rsv). Respiratory syncytial virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

s2083i omicron-variant gene target (s2083i). s2083i omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s371l omicron-variant gene target (s371l). s371l omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s373p omicron-variant gene target (s373p). s373p omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s375f omicron-variant gene target (s375f). s375f omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s477n omicron-variant gene target (s477n). s477n omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Settled solids (settsol). Amount of settled solids from a wastewater sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

S sars-cov-2 gene target (sGene). S sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

slope (slope). Slope value of the calibration curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity (specHum). Measure for specific humidity, the unit is the ratio of the mass of water vapour and the mass of total air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Storage temp (sTemp). Temperature that the sample is stored at in Celsius. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Storage time (stoTim). Length of time that a sample was in storage. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sub-variant or lineage (subVar). A unit used to report a specific genetic lineage, or to specify that a reported variant is a sub-variant of another. Status: active. First released: 2.0.0. Last updated: 2.0.0.

t19r delta-variant gene target (t19r). t19r delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

t3646a delta-variant gene target (t3646a). t3646a delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

t547k omicron-variant gene target (t547k). t547k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

t9i omicron-variant gene target (t9i). t9i omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

TaqPath N sars-cov-2 gene target (taqpatN). TaqPath N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

TaqPath S sars-cov-2 gene target (taqpatS). TaqPath N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample temperature (temp). Temperature of the sample measured in degrees Celcius Status: active. First released: 1.0.0. Last updated: 2.0.0.

Nucleic acid template volume (tempVol). The volume of DNA or RNA template used for PCR. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Total Kjeldahl Nitrogen (tkn). A measure of Total Kjeldahl Nitrogen; ie. the sum of nitrogen bound in organic substances, nitrogen in ammonia (NH₃-N) and in ammonium (NH₄⁺-N) in the chemical analysis of soil, water, or waste water. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Total Nitrogen (tn). Total nitrogen concentration, as N. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Total solids concentration (ts). Total solids concentration Status: active. First released: 1.0.0. Last updated: 2.0.0.

Total suspended solids (tss). Total suspended solids Status: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration of total suspended solids (tssConc). Total suspended solids concentration of the wastewater. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Turbidity (turb). A measure for the turbidity of water, or a liquid sample, quantifying the relative clarity or opcaity of a liquid. Highly indicative of water quality. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Tube weight empty (tWeighE). The weight of the tube used for analysis while full of analyte. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Tube weight full (tWeighF). The weight of the tube used for analysis while empty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

v2930l delta-variant gene target (v2930l). Unique identifier for each measurement for a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Size in volume (vol). A group of measures/methods related to miscellaneous viruses. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Volatile suspended solids (vss). Unit set of volume measurements. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Volatile suspended solids - ignition (vssIg). Volatile suspended solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wastewater Temperature (watTemp). A compartment set for measures and methods in the water compartment. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Weather (weath). Temperature of the wastewater. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wind Speed (wind). The table for wide names. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wastewater treatment plant designed capacity (wwtpCap). Indicates a wastewater treatment plant where multiple labs are or were sampling, and where secondary samples are being taken taken. Status: active. First released: 2.0.0. Last updated: 2.0.0.

y505h omicron-variant gene target (y505h). Municipal wastewater treatment plant for sanitary sewage only Status: active. First released: 2.0.0. Last updated: 2.0.0.

Concentration method set

(methodConcSet). CatSeteogry set for concentration methods.

Aloh3 precipitation (aloh3). Aloh3 precipitation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Amicon ultrafiltration (amiconUf). Amicon ultrafiltration Status: active. First released: 2.0.0. Last updated: 2.0.0.

Beef extract flocculation (beeExtractFloc). Beef extract flocculation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Centricon ultrafiltration (centriconUf). Centricon ultrafiltration. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Ceres nanotrap (ceres). Ceres nanotrap. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Hollow fiber dead end ultrafiltration (hollowFiberUF). Hollow fiber dead end ultrafiltration Status: active. First released: 2.0.0. Last updated: 2.0.0.

Innovaprep ultrafiltration (innovaprepUF). Innovaprep ultrafiltration Status: active. First released: 2.0.0. Last updated: 2.0.0.

Polyethyleneglycol (PEG) precipitation (peg). Peg precipitation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Skimmed milk flocculation (skimMilkFloc). Skimmed milk flocculation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Ultracentrifugation (untracent). An idenfication of a set of units that can be used for a measure or method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Methods

(methods). A procedure or step that is performed to make a measure. For example, methodExtract is a method tht describes the how a sample was extracted.

Nucleic acid extraction method (extraction). Description of the xtraction method. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Inhibition method (inhibMeth). Description of the method used to evaluate molecular inhibition. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration method (methodConc). Description of the method to concentrate a wastewater sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

PCR method (pcr). Description of the PCR method. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Pretreatment (pretreat). Was the sample chemically treated in anyway with the addition of stabilizers or other Status: active. First released: 1.0.0. Last updated: 2.0.0.

Genetic primer (primer). Method ID used for indicating the primer used for a PCR or sequencing assay. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quaility assurance method (qaqc). Quality ass qc Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Collection other depreciated (sampleTypeOther). Description for other type of method not listed in collection. depreciated Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Sets

Genetic sequencing script version (scriptVersion). Used to specify the script or script version used to extract summary information (i.e. coverage and mutation frequency statistics, etc.) from the sequencing output. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Layout (seqLay). The layout of genetic material used in sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method (seqSel). The primer sequence selection method used in sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy (seqStrat). The sequencing strategy used for an analysis Status: active. First released: 2.0.0. Last updated: 2.0.0.

Shipped on ice (shipOnIce). Was the sample kept cool while being shipped to the lab Status: active. First released: 1.0.0. Last updated: 2.0.0.

Solid seperation (solidSep). Process used to separate solid and liquid phases of the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Spike material (spikeMat). Material into which the recovery efficiency control target is spiked. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Recovery efficiency spike target (spikeTarget). Method specifying the recovery efficiency control target. This matches on to the the spike material method ID. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Type other depreciate (tyOtDep). Description for other type of sample not listed in Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Missingness

(missingness). The part type for missingness values. Only used for the dictionary entries of parts.

Not applicable (NA). The field for which the expected value is not a property of the described object. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not a number (nan). The outcome of a measurement is not a valid number (plus or minus infinity, error, ...) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported (nr). A value could have been recorded, however, it was not. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Null (null). A logical representation of a statement that is neither TRUE nor FALSE. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Undisclosed (undisc). A value has been recorded, however it was not included in the dataset. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Missingness sets

(missingnessSets). Missingness sets for measures, methods or attributes. Inputted data into an ODM field may have missing data. The missingness set lists valid missing categories for the measure, method, or attribute.

General missingness set

(genMissingnessSet)The general set for missingness values. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sets

Not applicable (NA). The field for which the expected value is not a property of the described object. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not a number (nan). The outcome of a measurement is not a valid number (plus or minus infinity, error, ...) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported (nr). A value could have been recorded, however, it was not. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Null (null). A logical representation of a statement that is neither TRUE nor FALSE. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Undisclosed (undisc). A value has been recorded, however it was not included in the dataset. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported, Not applicable missing set

(nrNAMissingnessSet)Not reported, Not applicable missing set. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Not applicable (NA). The field for which the expected value is not a property of the described object. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Nomenclatures

(nomenclatures). A classification system to report the measure class. See partID = nomenclatureID.

International Classification of Diseases (icd). Classification system for diseases in humans. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization type set

Nomenclature not applicable (naNomenclature). Not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

NextClade nomenclature (nextclade). Specifies variant or genetic nomenclature as set out by NextClade. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Pangolin nomenclature (pangolin). Specifies variant or genetic nomenclature as set out by the Phylogenetic Assignment of Named Global Outbreak Lineages (Pangolin). Status: active. First released: 2.0.0. Last updated: 2.0.0.

WHO nomenclature (who). Specifies the whole genome sequencing (WGS) strategy for genetic sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization type set

(orgTypeSet). The set for storing organization types.

Academic institution (academ). The category of organization type used for academic institutions or research groups. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Government agency (govt). The category of organization type used for government agencies, programs, or crown-owned bodies. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Private sector (priv). The category of organization type used for private sector or non-profit groups that are not otherwise captured by the academic category. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Sample Origin set

(originSet). The set for storing the valid categorical values of sample origin.

Derived sample (derived). Specifies a sample that is derived or made from another sample or material. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample (field). Specifies a sample taken from the field; directly collected from an area for testing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Synthetic sample (synthetic). Specifies a synthetic, or artificially derived sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak set

(outbreakSet). set for the valid values of the outbreak measure.

Outbreak End (outbEnd). Indicates an outbreak has ended. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak - on-going (outbOngoing). Indicates an outbreak is on-going. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak start (outbStart). Indicates an outbreak began. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PCR Method set

(pcrSet). The set capturing containing all of the valid methods for the PCR methodID.

Digital droplet emulsification PCR (bioRadDdpcr). Describes a PCR analysis done using BioRad's digital droplet emulsification PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Fluidigm digital PCR (fluidpcr). Describes a PCR analysis done using FluidIGM's digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Life technologies digital PCR (ltDpcr). Describes a PCR analysis done using Life Technologies' digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen digital PCR (qgDpcr). Describes a PCR analysis done using Qiagen's digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quantitative PCR (qpcr). Real-time PCR, also called 'quantAggScale' PCR Status: active. First released: 2.0.0. Last updated: 2.0.0.

Raindance digital PCR (rainDpcr). Describes a PCR analysis done using Raindance's digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Genetic primer set

(primerSet). The set for genetic primers used in sequencing or PCR assays.

Sets

Use of the articV3 primer (articV3). Initial implementation of an ARTIC bioinformatics platform for nanopore sequencing of nCoV2019 novel coronavirus; Artic Foundation v3 primer. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Use of the articV4 primer (articV4). Artic V4 sequencing primer for VOCs and VOIs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol Relationships set

(protocolRelSet). set for valid values of relationshipID in the protocolRelationships table. Use only for protocolRelationships table relationshipID.

Is Before (before). Specifies that the object step or protocol occurs before the subject step or protocol in the overall protocol container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specifies (specifies). Specifies that the object step or protocol occurs specifies details for the subject step or protocol in the overall protocol container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Is concurrent to (concurrent). Specifies that the object step or protocol occurs at the same time as the subject step or protocol in the overall protocol container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Purpose category

(purposeSet). Purpose set for a sample or measure.

Quality indicators

Education (education). A measure or sample taken for education or training. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple Purpose (multiple). A measure or sample taken for multiple purposes, not easily captured by the other purpose categories. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Provisional report (provisional). A provisional or interm report. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality control (qualityControl). A measure or sample taken for the purpose of quality control. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Regular (regular). A measure or sample taken for surveillance or epidemiology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Testing (testing). A measure or sample taken to test a new method. These measures are typically used for internal lab uses and not reported to external partners. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Validation study (validationStudy). v2930l delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicators

(qualityIndicators). A measure of the quality of a reported value or sample. Only used for the dictionary entries of parts.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed (flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution.

Sets

The reported concentration estimate is the updated result after addressing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination (flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

FI - Inhibition present and unaddressed (flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

J - Weak signal extrapolation (flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ND - Non-detect (flagND). No amplification occurred in the reaction; non-detect. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UJ - Trace signal extrapolation (flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., “trace” signal observed); the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable (flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample frozen (frozen). Sample was frozen before analysis or processing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage (IBC). The percentage of the genome covered by reads (the breadth) is low. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage (IDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Leaked sample (leaked). Sample leaked, some volume and material was lost. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low-volume sample (lowVol). Sample is low-volume, but was run regardless. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple issues (mI). Multiple issues have arisen in the sequencing process. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there is are no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample not labelled (noLabel). Sample had no label Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality concerns (qf1). A flag to indicate there is a quality concern, not otherwise sepcified. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sparse coverage (sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample spilled (spill). Sample contents spilled from container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Wide 95% interval (wI). Specifies variant or genetic nomenclature as set out by the World Health Organization. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored incorrectly (wrongStorage). Unique identifier for wide table only. Use when all measures are performed on a single sample at the same time and same laboratory. Suggestion: siteID_sampleID_LabID_reportDate_ID. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored at wrong temperature (wrongTemp). Sample was stored inappropriately. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicator sets

(qualityIndSets). Sets of quality indicators or measures. For example, PCR have a quality measure set.

Generic quality flag set

(genQualitySet) A quality set to specify any generic quality concerns about a measure or sample. Status: active. First Released: 1.0.0. Last updated: 2.0.0.

Multiple issues (mI). Multiple issues have arisen in the sequencing process. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there are no quality concerns about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Quality concerns (qf1). A flag to indicate there is a quality concern, not otherwise specified. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measures quality set

(measQualitySet)A quality set for any measure, includes all the quality flag options for any measure. Used as an “any quality set” for measures, but excludes the sample quality flags to avoid confusion and data entry errors. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed (flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution. The reported concentration estimate is the updated result after addressing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination (flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

FI - Inhibition present and unaddressed (flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

J - Weak signal extrapolation (flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

ND - Non-detect (flagND). No amplification occurred in the reaction; non-detect. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UJ - Trace signal extrapolation (flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., “trace” signal observed); the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable (flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage (IBC). The percentage of the genome covered by reads (the breadth) is low. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage (IDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple issues (mI). Multiple issues have arisen in the sequencing process. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there are no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality concerns (qf1). A flag to indicate there is a quality concern, not otherwise specified. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sparse coverage (sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PCR quality set

(pcrQualitySet)Quality set for PCR measures. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed (flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution. The reported concentration estimate is the updated result after addressing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination (flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

FI - Inhibition present and unaddressed (flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

J - Weak signal extrapolation (flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ND - Non-detect (flagND). No amplification occurred in the reaction; non-detect. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

UJ - Trace signal extrapolation (flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., “trace” signal observed); the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable (flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there is are no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample quality set

(sampleQualitySet)Quality set for a sample. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sample frozen (frozen). Sample was frozen before analysis or processing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Leaked sample (leaked). Sample leaked, some volume and material was lost. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low-volume sample (lowVol). Sample is low-volume, but was run regardless. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there is no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample not labelled (noLabel). Sample had no label Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample spilled (spill). Sample contents spilled from container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored incorrectly (wrongStorage). Unique identifier for wide table only. Use when all measures are performed on a single sample at the same time and same laboratory. Suggestion: siteID_sampleID_LabID_reportDate_ID. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored at wrong temperature (wrongTemp). Sample was stored inappropriately. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing quality set

(seqQualitySet)Quality set for sequencing measures Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage (IBC). The percentage of the genome covered by reads (the breadth) is low. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage (IDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there is no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sparse coverage (sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Replicate Type set

(replicateSet). The set for storing valid categorical values of replicate type.

Co-located sample (colocated). Second or multiples samples collected at same location but different time (water, air) or at a nearby location (soil, sediment). Sent blind to laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample replicate (fieldReplicate). A sample divided into two or more homogeneous parts. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory duplicate (labDuplicate). Second (time or more) processing and analysis of sample. Usually for general chemistry or metals analyses. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory control sample duplicate (lcsd). Known amounts of an analyte or representative compounds are added to a second “clean” matrix (lab water or clean sand) in laboratory. Duplicate of LCS Status: active. First released: 2.0.0. Last updated: 2.0.0.

Matrix spike duplicate (msd). A known amounts of an analyte or representative compounds are added in the laboratory to a second aliquot

Sample material set

of the sample used for matrix spike. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unique sample (unique). A unique sample, not duplicated. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample material set

(sampleMatSet). set for the types of material that can be found in a sample.

Air filter (afu). Air filter as part of filtration or circulating unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Desk or counter (desk). Desk, table, countertop or other flat working surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Fecal matter (faeces). Fecal matter. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Floor (floor). Floor of a building or room. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Holding tank wastewater (htSam). Wastewater sampled from a holding tank, such as from an airplane or ship Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water (nww). Non-wastewater, coming from any kind of water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Primary clarifier effluent (pEfflu). Effluent obtained after primary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Primary clarifier sludge (pSludge). Sludge produced by primary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Post-grit (pstGrit). Raw wastewater after a treatment plant's headworks (post grit or removal of large solids at a treatment plant but prior to a primary clarifier) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage at site (rawWW). Wastewater without any form of treatment Status: active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage downstream from a site (rawWWdown). Downstream from a site. See partType = rawWWup. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Raw sewage upstream from a site (rawWWup). Upstream from a site. Used when there is not direct access to a site's wastewater inconjuction with downstream sample. For example, if the site is a long-term care home, but there is not access to the building or property cleanout. Use two sample, one sample upstream and one sample downstream. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier effluent (sEfflu). Effluent obtained after secondary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Septic tank wastewater (septage). Wastewater from within a septic tank Status: active. First released: 1.0.0. Last updated: 2.0.0.

Secondary clarifier sludge (sSludge). Sludge produced by secondary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Other surface (surface). Surface other than floor or desk. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sewer sediment (swrSed). Sediments obtained in sewer Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample relationships set

(sampleRelSet). set for valid values of relationshipID in the sampleRelationships table. Use only for sampleRelationships table relationshipID.

Child relationship (child). Indicated that this is a sample generated from (an)other sample(s) either because of pooling or sub-sampling. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Co-located sample (colocated). Second or multiples samples collected at same location but different time (water, air) or at a nearby location (soil, sediment). Sent blind to laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample replicate (fieldReplicate). A sample divided into two or more homogeneous parts. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Matrix spike duplicate (msd). A known amounts of an analyte or representative compounds are added in the laboratory to a second aliquot of the sample used for matrix spike. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory duplicate (labDuplicate). Second (time or more) processing and analysis of sample. Usually for general chemistry or metals analyses. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory control sample duplicate (lcsd). Known amounts of an analyte or representative compounds are added to a second “clean” matrix (lab water or clean sand) in laboratory. Duplicate of LCS Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing script version set

(scriptSet). The set for script versions

Gromstole 1.0 Script (gromstole). Gromstole 1.0 Script for sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Freyja Script (freyja). Freyja Script for sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Layout set

(seqLaySet). The set for the layout of genetic material used in sequencing

Single Layout (sinLay). Specifies the single layout for a sequencing method Status: active. First released: 2.0.0. Last updated: 2.0.0.

Paired Layout (pairLay). Specifies the paired layout for a sequencing method Status: active. First released: 2.0.0. Last updated: 2.0.0.

Mate pair layout (mateLay). Specifies the mate pair layout for a sequencing method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method set

(seqSelSet). The set for the selection method used in sequencing

Random sequencing selection method (ranSeq). Specifies the random selection method for sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

PCR sequencingsSelection method (pcrSeq). Specifies the PCR selection method for sequencing, ie. that a pre-determined primer was used. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy set

(seqStratSet). The set for the sequencing strategy used for an analysis

Amplicon sequencing (amp). Specifies the amplicon strategy for genetic sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Whole genome Sequencing (wgs). A set of the valid qualitative categories for the qualitative weather measure. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Severity set

(sevSet). A set for severity indicators.

Low severity (low). A marker for low severity Status: active. First released: 2.0.0. Last updated: 2.0.0.

Mid-level severity (mid). A marker for med-level severity, where there are some concerns. Status: active. First released: 2.0.0. Last updated: 2.0.0.

High severity (high). Indicates a very sever quality issue, likely meaning the data should not be reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sampled set

(shedSet). The set for all valid values of sampled.

Airplane (airpln). Airplane Status: active. First released: 1.0.0. Last updated: 2.0.0.

Long-term acute care hospital (ccc). Acute care hospitals, or complex continuing care, that provide care for patients with average length of stay longer than 25 days. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Child day care (cdc). Child day care facility. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Cphd ID (corFcil). Correctional facility Status: active. First released: 1.0.0. Last updated: 2.0.0.

Higher education dormitory or residential building (dorm). Higher education dormitory or residential building Status: active. First released: 2.0.0. Last updated: 2.0.0.

First Nation (fiNa). Used to categorize a sampled that is a First Nation, or on reserve lands. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Hospital (hosptl). Hospital Status: active. First released: 1.0.0. Last updated: 2.0.0.

Long-term care facility (ltcf). A residential healthcare facility that provides 24-medical care. These are also called skilled nursing facilities. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Long-term care - assisted living or retirement home (ltcfAl). A residential facility that provides assistance with daily care but generally does not provide skilled nursing care. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other long-term care (ltcfO). Other residential facilities that provide daily and/or medical care, but are not defined as nursing home/skilled nursing facilities or assisted living facilities. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Municipality (municp). A complete municipality, this specifies an entire metropolitan area, either a city, town, etc. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Neighbourhood (neigh). A municipal neighbourhood, this specifies a sub-section of a larger municipality. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other residential building (orb). Individual residential buildings or institutions not captured in other categories Status: active. First released: 2.0.0. Last updated: 2.0.0.

School (school). A school serving students in the kindergarten to 12th grade range Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ship (ship). A cruise ship or other ship. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Social services shelter (sss). Other type of social services shelter Status: active. First released: 2.0.0. Last updated: 2.0.0.

University campus (uCampus). Universityor college campus - comprising an entire campus or part of a campus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Airport (airport). Airport Sample Shed category type Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site set

(siteTypeSet). set for the type of sampling site.

Branch sewer pipeline (bSwrPpl). Specifies a site type that is collection pipe that run lateral to other municipal sewer lines, allowing drainage into the main sewer. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Building cleanout (buildCO). Specifies a site type that is a capped pipe that connects to a building's main sewer line. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Estuary, natural water body (estuary). Estuary, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lake, natural water body (lake). Lake, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Holding tank (htSite). Holding tank Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lagoon system (lagoon). Logoon system for extensive wastewater treatment Status: active. First released: 1.0.0. Last updated: 2.0.0.

Major sewer pipeline (mSwrPpl). Major sewer pipeline Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ocean, natural water body (ocean). Ocean, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Pumping station (pStat). Pumping station Status: active. First released: 1.0.0. Last updated: 2.0.0.

Retention pond (retPond). Retention pond Status: active. First released: 1.0.0. Last updated: 2.0.0.

River, natural water body (river). River, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sea, natural water body (sea). Sea, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Septic tank (septTnk). Septic tank Status: active. First released: 1.0.0. Last updated: 2.0.0.

Stabilization pond (stabPnd). Specifies a site type which is a pond designed and built for wastewater treatment to reduce the organic content and remove pathogens from wastewater. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewage truck (swgTrck). Sewage truck Status: active. First released: 1.0.0. Last updated: 2.0.0.

Industrial wastewater treatment plant (wwtpInd). A measure for the designed capacity of a wastewater treatment plant, in terms of how much water it was designed to be able to hold and process. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for combined sewage (wwtpMuC). Industrial wastewater treatment plant Status: active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for sanitary dewage only (wwtpMuS). Municipal wastewater treatment plant for combined sewage Status: active. First released: 1.0.0. Last updated: 2.0.0.

Wastewater treatment plant (wwtp). W-Shere notes Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upstream sites (upstream). Specifies the the upper limit of a 95% confidence interval. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Back-up wastewater treatment site (wwtpBack). A general site type for wastewater treatment plants. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Solid separation set

(solidSeparationSet). set for the separation of solids.

Centrifugation (cent). Describes solid separation from a wastewater sample via centrifugation. Likely connected to other method steps prior to analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Filtration (filt). Describes solid separation from a wastewater sample via filtration. Proceeds further concentration or analysis of the liquid filtrate. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gravity settling (graSet). Describes solid separation from a wastewater sample where the sample material is allowed to settle by gravity, and then separated. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2 (mfAcidmgcl2). Membrane filtration with acidification and mgcl2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2, membrane recombined with separated solids (mfAcidmgcl2SS). Membrane filtration with acidification and mgcl2, membrane recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with addition of mgcl2 (mfMgcl2). Membrane filtration with addition of mgcl2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with mgcl2 and separated solids (mfMgcl2SS). Membrane filtration with addition of mgcl2, membrane recombined with

Specimens

separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment (mfNoAmend). Membrane filtration with no amendment Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment, membrane recombined with separated solids (mfNoAmendSS). Membrane filtration with no amendment, membrane recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with sample acidification (mfSampleAcid). Membrane filtration with sample acidification Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with sample acidification, membrane recombined with separated solids (mfSampleAcidSS). Membrane filtration with sample acidification, membrane recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

No liquid concentration, liquid recombined with separated solids (noliquid). No liquid concentration, liquid recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimens

(specimens). Measures or observations are taken from three types of substances: populations (humans or geographic areas), samples, sites. See partID = specimenID

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen sets

(specimenSets). Sets of specimens. Specimen sets are used when a measure or attribute can be recorded for more than one specimen.

Any specimen set

(anySpecimenSet) A specimen set that includes any specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen set not applicable

(naSpecimenSet) A specimen set for when specimen/specimen set is not applicable. Use 'not applicable' when there is no specimen/set for the part Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Specimen not applicable (naSpecimen). Non applicable specimen. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Poulatiion or site set

(poSiSpecimenSet)A specimen set that inculdes population or site specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population specimen set

(poSpecimenSet)A specimen set that inculdes only a population specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample specimen set

(saSpecimenSet)A specimen set that inculdes only a sample specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Site or sample specimen set

(siSaSpecimenSet) A specimen set that includes site or sample specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site specimen set

(siSpecimenSet) A specimen set that includes only a site specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Spike material set

(spikeMatSet). set for spikeMat (material into which the recovery efficiency control target is spiked).

Clarified sample (clari). Clarified sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample concentrate (conc). Sample concentrate. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Dewatered solids (deso). Dewatered solids. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Recovery efficiency spike target set

Lysis buffer (lysi). Lysis buffer. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Raw wastewater (raws). Raw wastes water sample Status: active. First released: 1.0.0. Last updated: 2.0.0.

Raw post-pasteurized wastewater (rppw). Raw wastewater sample post-pasteurization Status: active. First released: 1.0.0. Last updated: 2.0.0.

Recovery efficiency spike target set

(spikeTargetSet). The set capturing containing all of the valid categories for the recovery efficiency control (or spike) target.

bcov culture spike target (bcovCul). Cultured bovine coronavirus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

bcov vaccine spike target (bcovVac). The bovine coronavirus vaccine is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

brsv culture spike target (brsvCul). Cultured bovine respiratory syncytial virus (BRSV) is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

brsv vaccine spike target (brsvVac). The bovine respiratory syncytial virus (BRSV) vaccine is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Heat inactivated sars-cov-2 virus spike target (heatInacSARS). Heat inactivated SARS-CoV-2 virus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

hep g armored rna spike target (hepGRnaMat). Hepatitis G armored RNA is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage spike target (ms2ColMat). ms2 coliphage is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

murine coronavirus spike target (muCoMat). Murine coronavirus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

oc43 spike target (oc43Mat). Human coronavirus OC43 is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

phi6 spike target (phi6Mat). Pseudomonas virus phi6 is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Puro virus spike target (puroMat). Puro Virus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Status set

(statusSet). A set for partID = Status to indicate whether a part is in current use or not.

Active (active). Indicator that a part is in current use. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Development (development). Indicator that a part is under development use. See partID = status Status: active. First released: 2.0.0. Last updated: 2.0.0.

Deprecated (deprecated). Indicator to say that a part is no longer in current use in the model. See `partID = status` Status: active. First released: 2.0.0. Last updated: 2.0.0.

Tables

(tables). Tables are where measures, methods and attributes are recorded. Tables represent the main entities of the environmental surveillance. There An example of a table is **site** a collection of attributes such as site name and address to describe where environment samples are taken. PHES-ODM has report tables, but users can create their own tables as well.

Contact table (contacts). The table that contains information about a contact; a person who is the contact of a site or laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset table (datasets). A report table for capturing details about data's parental data set and data custodians. Supplying attribution for data collectors. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Instrument table (instruments). The table that contains information about instruments. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Language Look-up table (languages). Look up table for all languages, used to give structure to the translation table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure report table (measures). The table that contains information and details about a given measure Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measure set report table (measureSets). The table that identifies sets of measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Organization table (organizations). The table that contains information about a laboratory. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Parts Look-up table (parts). Look up table containing all parts in of the data model. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Polygon table (polygons). The table that contains information about the geometry of a geographic area. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Protocol relationships table (protocolRelationships). The table that contains the organizational information and details (such as order) about a given protocol, recorded as a series of protocol steps (protocolSteps). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocols table (protocols). The table for protocols. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol steps table (protocolSteps). The table for collecting metadata on individual steps in a protocol, methodological process, or assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Quality reports table (qualityReports). The table for recording the various quality metrics and indicators for samples and measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Reporter table deprecated (reportersDep). The table that contains information about a reporter of a sample, method, measure, or attribute. Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

Sample relationships table (sampleRelationships). Table for recording the relationships between samples. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample report table (samples). The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets look-up table (sets). Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sites table (sites). The table that contains information about a site; the location where an environmental sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Translation look-up table (translations). Look up table for translations of the description, label, and instruction for all parts. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wide name table (wideNames). The 95% interval is too wide, low confidence in results. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Table column set

(tableSet). The set for valid inputs in table columns. These categories are used in the parts list for columns named after tables to indicate the position and function of a given part in the entity relationship diagram of the data model.

Foreign key (fK). Foreign key for a table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Header (header). Header for a table. Also known as a table variable or entity relationship ‘attribute’. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Input (input). Input for a table. Indicates if a part can be used in a table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Primary key (pK). Primary key for a table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Units

(units). The unit of the measurement. Every measurement must have a unit. Meaning, **units** are a mandatory field whenever a measurement is recorded.

Number of active cases (acti). Number of active cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Breadth of coverage (≥ 5 x depth) (boc). Positions with read depth greater or equal to 5. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Degrees Celcius (cel). Degrees Celsius. Status: active. First released: 1.1.0. Last updated: 2.0.0.

CFU per 100 ml (cfu). Colony forming units per 100 ml of filtered sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Centimetres (cm). Unit part for the SI unit of centimetres. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Confirmed Cases Case report date (conCase). Date that the numbers were reported publicly for a confirmed case. Typically, reported data and this measure is most commonly reported and used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Number of confirmed cases (conf). Number of confirmed cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases episode date (confEp). Episode date is the earliest of onset, test or reported date for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Onset date (confOn). Earliest that symptoms were reported for a confirmed case. This data is often not known and reported. In lieu, episode is used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases test date (conTest). Date that the covid-19 test was performed for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cumulative count (cuCo). Units for describing a population measure of a cumulative count of cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Daily count (daiCo). Units for describing a population measure of a daily count of new cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Deaths (death). Units for describing a population measure of patients who have died from a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Detected (det). Substance detected or not detected. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Depth of coverage (doc). The sequencing read depth. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per copy of crAssphage (gcCrA). Gene or variant copies per copy of crAssphage. Status: active. First released: 1.0.0. Last updated: 2.0.0.

gene copies per day per 100,000 (gcD100). The unit for measures reflecting the gene copies per day per 100,000 people in the population. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Gene copies per gram solids (gCGS). Gene or variant copies per gram solids. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per L (gcL). Gene or variant copies per litre. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per mL (gcMl). Gene or variant copies per millilitre of solution. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per PMMoV copy (gcPpmov). Gene or variant copies per copy of PMMoV. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gram per cubic metre (gm3). Density unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Grams (grams). A unit of mass or weight. Status: active. First released: 2.0.0. Last updated: 2.0.0.

General Hospital Admissions (hosa). Hospital admissions or patients newly admitted to hospital. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Hospital Census (hosc). Hospital census or the number of people admitted with an ailment. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Intensive care unit patients (icu). Units for describing a population measure of patients who are in intensive care due to a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Kilogram per second (kgS). Kilograms per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Kilolitres (kl). Kilolitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Litres (l). Litres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Cubic metres per day (m3D). Cubic metres per day. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per hour (m3H). Cubic metres per hour. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per second (m3S). Cubic metre per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Millions of gallons per day (MG/D) (mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to measure flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre (mgL). Milligrams per litre. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Millilitres (ml). Millilitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Megalitres per day (ML/d) (mld). Megalitres per day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres (mm). Unit part for the SI unit of millimetres. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Metres per second (ms). metres per second. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unit not applicable (naUnit). Not applicable for units. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population of newly vaccinated persons (newVax). Population of newly vaccinated persons Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nephelometric turbidity unit (ntu). Nephelometric Turbidity Units, a unit used in the measurement of turbidity (see “turbidity” under measureIDs). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Other Unit deprecated (otunDep). Other measurement of viral copies or wastewater treatment plant parameter. Add description to UnitOther. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Percent (perc). Percentage. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Percent recovery (percRec). Percent of the surrogate recovery for a recovery efficiency control assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Population equivalents (popEq). A unit for measuring of design capacity for wastewater treatment plants, represented as the approximate amount of water a single person would use. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent positive (pp). Percent positive of sample measures. Status: active. First released: 1.0.0. Last updated: 2.0.0.

parts per million (ppm). Parts per million. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent positivity rate (pprt). Percent positivity rate of tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent primary sludge (pps). Percentage of total solids, for primary sludge. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of total (prop). Proportion as a percent of total. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of variant in sample (propV). Proportion of a variant as percent of total variants. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Number of positive tests (ptot). Number of positive tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ratio (ratio). Ratio (unitless) Status: active. First released: 2.0.0. Last updated: 2.0.0.

Recovered patients (recov). Units for describing a population measure of patients who have recovered from a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rate per 100,000 (rP100). Units for describing a population measure of a rate of case incidence per 100,000 people of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Number of tests performed (test). Number of tests performed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Unitless measure (unitless). A unit for unitless measures, like pH. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Micro-Siemens per centimetre (uSCm). A general site type for upstream wastewater sampling sites. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Frequency of variants detected (varFreq). The data type for variable character data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 1 dose of vaccine (vax1). Measures and methods relating to variants. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 2 doses of vaccine (vax2). A measure of the population with a single dose of vaccine. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 3 doses of vaccine (vax3). Population with 2 doses of vaccine Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unit sets

(unitSets). Sets of units. Contains units that are associated with part types such as measures. Sets are lists of input values. Each input value within a set is a part (has a partID) that can be reused between sets. For example, a set A could contain the values of 'yes' and 'no'; and set B could contain the values of 'yes', 'no', and 'maybe.' partTypes with sets include aggregations, compartments, missingness, quality flags, specimens, and units. Measures with categories and dataTypes have their own specific sets (catSetID) because their input values are unique and are not reused.

Absolute humidity unit set

(absHumidUnitSet)Unit set for absolute humidity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Gram per cubic metre (gm3). Density unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Allele unit set

(alleleUnitSet)Unit set for alleles. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Cycle threshold or quantification cycle (Ct or Cq) (ct). Cycle thresholds in a PCR assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Detected (det). Substance detected or not detected. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Depth of coverage (doc). The sequencing read depth. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per copy of crAssphage (gcCrA). Gene or variant copies per copy of crAssphage. Status: active. First released: 1.0.0. Last updated: 2.0.0.

gene copies per day per 100,000 (gcD100). The unit for measures reflecting the gene copies per day per 100,000 people in the population. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per gram solids (gCGS). Gene or variant copies per gram solids. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per mL (gcMl). Gene or variant copies per millilitre of solution. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per PMMoV copy (gcPpmov). Gene or variant copies per copy of PMMoV. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Bacteria unit set

(bacteriaUnitSet)Unit set for bacteria-related measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

CFU per 100 ml (cfu). Colony forming units per 100 ml of filtered sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Capacity unit set

(capacityUnitSet)Unit set for capacity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Millions of gallons per day (MG/D) (mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to measure flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Megalitres per day (ML/d) (mld). Megalitres per day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population equivalents (popEq). A unit for measuring of design capacity for wastewater treatment plants, represented as the approximate amount of water a single person would use. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Conductivity unit set

(conductivityUnitSet)Unit set related to conductivity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Micro-Siemens per centimetre (uSCm). A general site type for upstream wastewater sampling sites. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Dissolved gas concentration unit set

(dissGasUnitSet)Unit set for carbon dioxide concentrations measurements in water or air. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre (mgL). Milligrams per litre. Status: active. First released: 1.1.0. Last updated: 2.0.0.

parts per million (ppm). Parts per million. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Volume flow rate unit set

(flowUnitSet)Unit set for volume flow measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Kilogram per second (kgS). Kilograms per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cubic metres per day (m3D). Cubic metres per day. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per hour (m3H). Cubic metres per hour. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per second (m3S). Cubic metre per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Millions of gallons per day (MG/D) (mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to measure flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Megalitres per day (ML/d) (mld). Megalitres per day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Genetics unit set

(geneticUnitSet)Unit set for genetic-related measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Breadth of coverage ($\geq 5\times$ depth) (boc). Positions with read depth greater or equal to 5. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Detected (det). Substance detected or not detected. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sets

Depth of coverage (doc). The sequencing read depth. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per L (gcL). Gene or variant copies per litre. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent recovery (percRec). Percent of the surrogate recovery for a recovery efficiency control assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of total (prop). Proportion as a percent of total. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of variant in sample (propV). Proportion of a variant as percent of total variants. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Unit set not applicable

(naUnitSet)Not applicable for unit sets. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Unit not applicable (naUnit). Not applicable for units. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population unit set

(populationUnitSet)Unit set for hospital-related measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Number of active cases (acti). Number of active cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Case report date (conCase). Date that the numbers were reported publicly for a confirmed case. Typically, reported data and

this measure is most commonly reported and used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Number of confirmed cases (conf). Number of confirmed cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases episode date (confEp). Episode date is the earliest of onset, test or reported date for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Onset date (confOn). Earliest that symptoms were reported for a confirmed case. This data is often not known and reported. In lieu, episode is used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases test date (conTest). Date that the covid-19 test was performed for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cumulative count (cuCo). Units for describing a population measure of a cumulative count of cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Daily count (daiCo). Units for describing a population measure of a daily count of new cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Deaths (death). Units for describing a population measure of patients who have died from a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

General Hospital Admissions (hosa). Hospital admissions or patients newly admitted to hospital. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Hospital Census (hosc). Hospital census or the number of people admitted with an ailment. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Intensive care unit patients (icu). Units for describing a population measure of patients who are in intensive care due to a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population of newly vaccinated persons (newVax). Population of newly vaccinated persons Status: active. First released: 2.0.0. Last updated: 2.0.0.

Percent positivity rate (pprt). Percent positivity rate of tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Number of positive tests (ptot). Number of positive tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Recovered patients (recov). Units for describing a population measure of patients who have recovered from a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rate per 100,000 (rP100). Units for describing a population measure of a rate of case incidence per 100,000 people of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Number of tests performed (test). Number of tests performed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Population with 1 dose of vaccine (vax1). Measures and methods relating to variants. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 2 doses of vaccine (vax2). A measure of the population with a single dose of vaccine. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 3 doses of vaccine (vax3). Population with 2 doses of vaccine Status: active. First released: 2.0.0. Last updated: 2.0.0.

Precipitation unit set

(precipitationUnitSet)Unit set for precipitation measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Centimetres (cm). Unit part for the SI unit of centimetres. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres (mm). Unit part for the SI unit of millimetres. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Relative humidity unit set

(relHumidUnitSet)Unit set for relative humidity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Percent (perc). Percentage. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity unit set

(specHumidUnitSet)Unit set for specific humidity. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Ratio (ratio). Ratio (unitless) Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard concentration unit set

(stdConcentrationUnitSet)Unit set for concentration measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre (mgL). Milligrams per litre. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sets

parts per million (ppm). Parts per million. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent primary sludge (pps). Percentage of total solids, for primary sludge. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Temperature unit set

(temperatureUnitSet)Unit set for temperature measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Degrees Celcius (cel). Degrees Celsius. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Turbidity unit set

(turbidityUnitSet)Unit set for turbidity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Nephelometric turbidity unit (ntu). Nephelometric Turbidity Units, a unit used in the measurement of turbidity (see “turbidity” under measureIDs). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unitless unit set

(unitlessUnitSet)Unit set for measurements that does not have units. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Unitless measure (unitless). A unit for unitless measures, like pH. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Volume unit set

(volumeUnitSet) A volume proportional sample generally collected by an autosampler. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Grams (grams). A unit of mass or weight. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Kilolitres (kl). Kilolitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Litres (l). Litres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millilitres (ml). Millilitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres (mm). Unit part for the SI unit of millimetres. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Percent positive (pp). Percent positive of sample measures. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sub-variant or lineage (subVar). A unit used to report a specific genetic lineage, or to specify that a reported variant is a sub-variant of another. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Frequency of variants detected (varFreq). The data type for variable character data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wind speed unit set

(windSpeedUnitSet) A measure for wind speed; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sets

Metres per second (ms). metres per second. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Weather set

(weathSet). A measure for weather captured through qualitative categories; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources.

Cloudy (cloudy). Qualitative category for the weather measure, specifying an overcast day with no rain. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Foggy (foggy). Qualitative category for the weather measure, specifying a foggy or hazy day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rainy (rainy). Qualitative category for the weather measure, specifying a rainy day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Snowy (snowy). Qualitative category for the weather measure, specifying a snowy day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sunny (sunny). Qualitative category for the weather measure, specifying a clear and sunny day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Tables

ODM version 2.0.0-rc.3.16.

The ODM database model has 21 tables in its full relational database. The full model is commonly referred to as “long” tables as it stores data with one measurement per row. Each table has mandatory and optional headers.

i Data storage

Store all ODM data in the 21 recommended tables to ensure optimal interoperability and support by the [ODM validation toolkit](#) and other ODM libraries.

Users can, however, create custom tables or ‘views’ by merging fields from any table. A common format for data entry format are “wide” table format, with *one day per row* and multiple measures or attributes per row. It is important to note that this guide does not provide instructions on generating wide tables.

Table types

There are five table types:

1) **Regular report tables.** The **measures** and **samples** tables are used for daily reporting of new measurements and information on sample collection. There are three tables that support regular reporting:- **measureSets**

Tables

(optional) is used to link measures, `sampleRelationships` (optional) is used when samples are pooled or split, and `qualityReports` is used for detailed quality assurance and control measures.

2) **Contact information tables.** The `sites`, `organizations`, `contacts`, and `addresses` tables store contact information for the testing site and the person who performed the testing.

3) **Protocol or methods tables.** The `protocols` and `instruments` tables hold information on the methods used for sample collection or measurement. `protocolSteps` and `protocolRelationships` tables are used in conjunction with the `protocols` table.

4) **Other report tables.** The `polygons` table stores information on the sample area and the `datasets` table includes information on the data custodian and owner.

5) **Dictionary reference tables.** The `parts` table describes all elements of the ODM, including tables, table headers, measures, methods, categories, and units. The `sets` table stores collections of parts, such as units grouped together in a `unitSet`. The `languages` and `translations` tables support translations.

Header roles

Each table in the ODM has column headers, also known as table variables, fields, or entity relationship attributes. The header serves as the top row and contains the variable name.

Primary key (PK) - With the exception of the dictionary tables, all tables have one primary key. The primary key serves as a unique identifier for each row in the table, meaning that a primary key value cannot be repeated.

Mandatory, optional, mandatoryIf

Foreign key (FK) - Tables may have one or more foreign keys. A foreign key connects an entry to another table with a primary key. This allows for relationships to be established between tables, making it possible to link data between different tables.

Header (header) - This label is used in a table description when a header is not a primary or foreign key.

Mandatory, optional, mandatoryIf

Each header may be classified as either **mandatory** or **optional**. In the **measures** table, examples of mandatory headers are the **measureID** and **measurement date**. The ODM validation toolkit can be used to validate data, and will return an error if a table lacks a mandatory header or if there are missing row entries for mandatory fields.

Additionally, there are special headers called **mandatoryIf** headers. These headers are mandatory if certain conditions apply. For example, the **protocol steps** table has row entries that are either **measures** or **methods**. If the row entry is a **measure**, then a **unit** becomes **mandatory** since all measures must have a unit. Therefore, units are **mandatoryIf** in the protocol steps table.

Table list

Contact table

Dataset table

Instrument table

Language Look-up table

Measure report table

Tables

Measure set report table

Organization table

Parts Look-up table

Polygon table

Protocol relationships table

Protocols table

Protocol steps table

Quality reports table

Sample relationships table

Sample report table

Sets look-up table

Sites table

Translation look-up table

Contact table

partID: contacts. The table that contains information about a contact; a person who is the contact of a site or laboratory. adapt from version 1 documentation

Contact ID (contactID). A unique identifier for a given contact person.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

First name of contact (firstName). Specifies the first name of a given contact.Role: Header.Requirement: optional.Data type: varchar.

Last name of contact (lastName). Specifies the last name of a given contact.Role: Header.Requirement: optional.Data type: varchar.

Contact email (email). Contact e-mail address, for the lab.Role: Header.Requirement: mandatory.Data type: varchar.

Contact phone (phone). Contact phone number, for the lab.Role: Header.Requirement: optional.Data type: varchar.

Role of contact (role). Specifies the organizational role of a given contact.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Dataset table

partID: datasets. A report table for capturing details about data's parental data set and data custodians. Supplying attribution for data collectors.

Parent dataset ID (parDatasetID). The datasetID that is a parent to another datasetID.Role: Foreign Key.Requirement: optional.Data type: varchar.

Tables

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset creation date (datasetDate). Specifies the date a given dataset was created.Role: Header.Requirement: recommended.Data type: date-time.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

License (license). The license of a dataset.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English.Role: Foreign Key.Requirement: optional.Data type: varchar.

Funder Contact ID (funderCont). A unique identifier for funders.Role: Foreign Key.Requirement: recommended.Data type: varchar.

Custody Contact ID (custodyCont). A unique identifier for data custodians.Role: Foreign Key.Requirement: recommended.Data type: varchar.

Funding agency ID (funderID). The funding agency of the dataset.Role: Foreign Key.Requirement: optional.Data type: varchar.

Data custodian ID (custodyID). The data custodian of a database.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Instrument table

partID: instruments. The table that contains information about instruments. adapt from version 1 documentation

Instrument ID (instrumentID). A unique identifier for an instrument.Role: Primary Key.Requirement: mandatory.Data type: categorical.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Model (model). Model number or version of the instrument.Role: Header.Requirement: mandatory.Data type: varchar.

Manufacturer (manufacturer). Manufacturer of an instrument.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Tables

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Instrument Type (insType). Type of instrument used to perform the measurement.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: insTypeCatSet

Describe other instrument type, if applicable (insTypeOth). Description of the instrument in case it is not listed in instrumentType.Role: Header.Requirement: optional.Data type: varchar.

Index (index). Index number in case the measurement was taken multiple times.Role: Header.Requirement: optional.Data type: integer.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Language Look-up table

partID: languages. Look up table for all languages, used to give structure to the translation table.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Language family (langFam). Specifies the language family of a given language for translation and language tracking purposes.Role: Header.Requirement: mandatory.Data type: varchar.

Language name (langName). Specifies the name of the language in roman alphabet characters for translation and language tracking purposes.Role: Header.Requirement: mandatory.Data type: varchar.

Native Name (natName). The native name of the language, i.e. what the language is called by its speakers.Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-1 (iso6391). The first part of the ISO 639 series of international standards for language codes. Part 1 covers the registration of two-letter codes. There are 183 two-letter codes registered as of June 2021. The registered codes cover the world's major languages.Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-2B (iso6392B). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'B' specifies the bibliographic code (B code).Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-2T (iso6392T). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'T' specifies the terminological code (T code).Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-3 (iso6393). A set of international standards that lists short codes for language names. ISO 639-3 extends the ISO 639-2 alpha-3 codes with an aim to cover all known natural languages.Role: Header.Requirement: mandatory.Data type: categorical.

Tables

ISO639-6 (iso6396). A set of international standards that lists short codes for language names. ISO 639-6 builds off ISO639-3 with the use of four-letter codes, and allowing users to differentiate between variants of languages and language families, such as historical vs. revived versions of languages.Role: Header.Requirement: mandatory.Data type: categorical.

First released version (firstReleased). The version in which a part was first released.Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: optional.Data type: varchar.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Measure report table

partID: measures. The table that contains information and details about a given measure adapt from version 1 documentation

Report ID (measureRepID). Unique identifier for a measurement.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Foreign Key.Requirement: optional.Data type: varchar.

Sample ID (sampleID). Unique identifier for a sample.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Purpose (purpose). The reason the measure or sample was taken.Role: Foreign Key.Requirement: optional.Data type: categorical.Measure, method or attribute set: purposeCatSet

Polygon ID (polygonID). Unique identifier for the polygon.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site ID (siteID). Unique identifier for the location where a sample was taken.Role: Foreign Key.Requirement: recommended.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Report set ID (measureSetRepID). Unique identifier that links together a group of related measures.Role: Foreign Key.Requirement: optional.Data type: varchar.

Analysis date start (aDateStart). Date the measurement or analysis was started.Role: Header.Requirement: optional.Data type: datetime.

Analysis date end (aDateEnd). Date the measurement or analysis was completed.Role: Header.Requirement: mandatory.Data type: datetime.

Report date (reportDate). The date a measure was reported.Role: Header.Requirement: optional.Data type: datetime.

Specimen (specimen). The substance or thing upon which the observation was made.Role: Foreign Key.Requirement: mandatory.Data type: categorical.

Fraction analyzed (fraction). Fraction of the sample that is analyzed.Role: Foreign Key.Requirement: recommended.Data type: categorical.Measure, method or attribute set: fractionCatSet

Group (group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations.Role: Foreign Key.Requirement: optional.Data type: categorical.

Tables

Class (class). A unique identifier for a class, which is akin to a sub-group; it's a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class.Role: Foreign Key.Requirement: optional.Data type: categorical.

Measure (measure). A measurement or observation of any substance including a biological, physical or chemical substance. An measurement is recorded from a specimen (i.e. a site, sample, person or population).Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Value (value). A measure or sample taken for a validation study.Role: Header.Requirement: mandatory.

Unit (unit). The units of a measurement.Role: Foreign Key.Requirement: mandatory.Data type: categorical.

Aggregation (aggregation). Statistical measures used to report a measure. Each aggregation has a corresponding value.Role: Foreign Key.Requirement: mandatory.Data type: categorical.

Nomenclature (nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases.Role: Foreign Key.Requirement: optional.Data type: categorical.

Index (index). Index number in case the measurement was taken multiple times.Role: Header.Requirement: optional.Data type: integer.

Measure license (measureLic). Specifies the access and use licensing for a given single measurement.Role: Header.Requirement: optional.Data type: varchar.

Reportable (reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied.Role: Header.Requirement: optional.Data type: boolean.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Measure set report table

partID: measureSets. The table that identifies sets of measures. Examples of measure sets include a set of replicates, dilutions (used to generate a Ct curve) or variants that are identified in a single sample.

Report set ID (measureSetRepID). Unique identifier that links together a group of related measures.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Foreign Key.Requirement: optional.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Tables

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Organization table

partID: organizations. The table that contains information about a laboratory. adapt from version 1 documentation

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Address ID (addressID). A unique identifier for an address.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization Type (orgType). Specifies the type or purpose of a given organization.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: orgTypeCatSet

Organization level (orgLevel). The geographic level of an organization.Role: Header.Requirement: recommended.Data type: categorical.Measure, method or attribute set: orgLevelCatSet

Organization sector (orgSector). The sector of an organizationRole: Header.Requirement: recommended.Data type: categorical.Measure, method or attribute set: orgSectorCatSet

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Parts Look-up table

partID: parts. Look up table containing all parts in of the data model. Contains all parts, including self-referential parts.

Part ID (partID). The unique identify of any entity within the dictionary.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Label (partLabel). A human readable label of a part.Role: Header.Requirement: mandatory.Data type: varchar.

Part types (partType). Part types describe the purpose or use of the part.Role: Header.Requirement: mandatory.Data type: varchar.

Short names (shortName). Shortened names for tables and other important parts for use in wide names.Role: Header.Requirement: mandatory.Data type: varchar.

Part description (partDesc). The description of the part.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Part instruction (partInstr). Additional notes and instructions on how a part is used and/or defined.Role: Header.Requirement: mandatory.Data type: varchar.

Domain (domain). Domain is the highest level of describing of a measure.Role: Header.Requirement: mandatory.Data type: categorical.

Specimen sets (specimenSets). Sets of specimens.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Group (group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations.Role: Header.Requirement: mandatory.Data type: categorical.

Class (class). A unique identifier for a class, which is akin to a subgroup; it's a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class.Role: Header.Requirement: mandatory.Data type: categorical.

Nomenclature (nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases.Role: Header.Requirement: mandatory.Data type: categorical.

Ontology reference (ontologyRef). Ontology reference for a part.Role: Header.Requirement: mandatory.Data type: varchar.

LaTeX expression (latExp). LaTeX expression used to generate formulas, symbols, etc. Mainly relevant for units.Role: Header.Requirement: mandatory.Data type: varchar.

Measure, method, or attribute (mmaSet). The set for a measure, method, or attribute. Only applicable for categorical fields, and the categories that populate them.Role: Header.Requirement: mandatory.Data type: categorical.

Aggregation scale (aggregationScale). A scale used for an aggregation. Only applicable for measures and units.Role: Header.Requirement: mandatory.Data type: categorical.

Aggregation set (aggregationSet). The aggregation set for a unit.Role: Header.Requirement: mandatory.Data type: categorical.

Quality set (qualitySetID). The quality set that corresponds to a given part. Only applicable for samples and measures.Role: Header.Requirement: mandatory.Data type: varchar.

Missingness sets (missingnessSets). Missingness sets for measures, methods or attributes.Role: Header.Requirement: mandatory.Data type: varchar.

Status (status). Whether the part is still active and can be used in the most current ODM version. Values are 'active' or 'inactive'.Role: Header.Requirement: mandatory.Data type: varchar.

First released version (firstReleased). The version in which a part was first released.Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: mandatory.Data type: varchar.

Protocol steps table (protocolSteps). The table for collecting metadata on individual steps in a protocol, methodological process, or assay.Role: Header.Requirement: mandatory.Data type: varchar.

Protocol Relationship header req (protocolRelationshipsRequired). Specifies the columns required in the Protocol Organization table.Role: Header.Requirement: mandatory.Data type: varchar.

Protocol steps header req (protocolStepsRequired). Specifies the columns required in a Protocol Steps table.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Protocol steps table column order (protocolStepsOrder). Specifies the order of the columns in a Protocol Steps table.Role: Header.Requirement: mandatory.Data type: integer.

Protocol relationships table (protocolRelationships). The table that contains the organizational information and details (such as order) about a given protocol, recorded as a series of protocol steps (protocolSteps).Role: Header.Requirement: mandatory.Data type: varchar.

Protocol relationships table column order (protocolRelationshipsOrder). Specifies the order of the columns in the Protocol Organization table.Role: Header.Requirement: mandatory.Data type: integer.

Measure report table (measures). The table that contains information and details about a given measureRole: Header.Requirement: mandatory.Data type: varchar.

Measures header req (measuresRequired). Specifies the columns required in a Measures table.Role: Header.Requirement: mandatory.Data type: varchar.

Measures table column order (measuresOrder). Specifies the order of the columns in a Measures table.Role: Header.Requirement: mandatory.Data type: integer.

Measure set report table (measureSets). The table that identifies sets of measures.Role: Header.Requirement: mandatory.Data type: varchar.

Measure sets table column order (measureSetsOrder). Specifies the order of the columns in a Measure Sets table.Role: Header.Requirement: mandatory.Data type: integer.

Measure Set header req (measureSetsRequired). Specifies the columns required in a Measure Sets table.Role: Header.Requirement: mandatory.Data type: varchar.

Dataset table (datasets). A report table for capturing details about data's parental data set and data custodians. Supplying attribution for

data collectors.Role: Header.Requirement: mandatory.Data type: varchar.

Dataset header req (datasetsRequired). Specifies the columns required in a Datasets table.Role: Header.Requirement: mandatory.Data type: varchar.

Datasets table column order (datasetsOrder). Specifies the order of the columns in a Datasets table.Role: Header.Requirement: mandatory.Data type: integer.

Sites table (sites). The table that contains information about a site; the location where an environmental sample was taken.Role: Header.Requirement: mandatory.Data type: varchar.

Sites header req (sitesRequired). Specifies the columns required in a Sites table.Role: Header.Requirement: mandatory.Data type: varchar.

Sites table column order (sitesOrder). Specifies the order of the columns in a Sites table.Role: Header.Requirement: mandatory.Data type: integer.

Sample report table (samples). The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site.Role: Header.Requirement: mandatory.Data type: varchar.

Sample header req (samplesRequired). Specifies the columns required in a Samples table.Role: Header.Requirement: mandatory.Data type: varchar.

Samples table column order (samplesOrder). Specifies the order of the columns in a Samples table.Role: Header.Requirement: mandatory.Data type: integer.

Address table (addresses). The table that contains information about addresses. Addresses can recorded for sites, organizations or contacts (individuals).Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Address header req (addressesRequired). Specifies the columns required in a Addresses table.Role: Header.Requirement: mandatory.Data type: varchar.

Addresses table column order (addressesOrder). Specifies the order of the columns in the Addresses table.Role: Header.Requirement: mandatory.Data type: integer.

Contact table (contacts). The table that contains information about a contact; a person who is the contact of a site or laboratory.Role: Header.Requirement: mandatory.Data type: varchar.

Contact header req (contactsRequired). Specifies the columns required in a Contacts table.Role: Header.Requirement: mandatory.Data type: varchar.

Contacts table column order (contactsOrder). Specifies the order of the columns in a Contacts table.Role: Header.Requirement: mandatory.Data type: integer.

Organization table (organizations). The table that contains information about a laboratory.Role: Header.Requirement: mandatory.Data type: varchar.

Organziation header req (organizationsRequired). Specifies the columns required in a Organizations table.Role: Header.Requirement: mandatory.Data type: varchar.

Organizations table column order (organizationsOrder). Specifies the order of the columns in a Organizations table.Role: Header.Requirement: mandatory.Data type: integer.

Instrument table (instruments). The table that contains information about instruments.Role: Header.Requirement: mandatory.Data type: varchar.

Instrutment header req (instrumentsRequired). Specifies the columns required in the Instruments table.Role: Header.Requirement: mandatory.Data type: varchar.

Instruments table column order (instrumentsOrder). Specifies the order of the columns in the Instruments table.Role: Header.Requirement: mandatory.Data type: integer.

Polygon table (polygons). The table that contains information about the geometry of a geographic area.Role: Header.Requirement: mandatory.Data type: varchar.

Polygon header req (polygonsRequired). Specifies the columns required in a Polygons table.Role: Header.Requirement: mandatory.Data type: varchar.

Polygons table column order (polygonsOrder). Specifies the order of the columns in a Polygons table.Role: Header.Requirement: mandatory.Data type: integer.

Language Look-up table (languages). Look up table for all languages, used to give structure to the translation table.Role: Header.Requirement: mandatory.Data type: varchar.

Language header req (languagesRequired). Required headers in the Languages table.Role: Header.Requirement: mandatory.Data type: varchar.

Language table column order (languagesOrder). Specifies the order of the columns in the Languages table.Role: Header.Requirement: mandatory.Data type: integer.

Translation look-up table (translations). Look up table for translations of the description, label, and instruction for all parts.Role: Header.Requirement: mandatory.Data type: varchar.

Parts Look-up table (parts). Look up table containing all parts in of the data model.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Tranalation header req (translationsRequired). Required headers in the Translations table.Role: Header.Requirement: mandatory.Data type: varchar.

Parts table column order (partsOrder). Order of headers in the Parts table.Role: Header.Requirement: mandatory.Data type: varchar.

Translation table column order (translationsOrder). Specifies the order of the columns in the Translation table.Role: Header.Requirement: mandatory.Data type: integer.

Parts header req (partsRequired). Required headers in the Parts table.Role: Header.Requirement: mandatory.Data type: varchar.

Sets look-up table (sets). Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together.Role: Header.Requirement: mandatory.Data type: varchar.

Sets table column order (setsOrder). Specifies the order of the columns in the Sets table.Role: Header.Requirement: mandatory.Data type: integer.

Sets header req (setsRequired). Required headers in the Sets table.Role: Header.Requirement: mandatory.Data type: varchar.

Quality reports table (qualityReports). The table for recording the various quality metrics and indicators for samples and measures.Role: Header.Requirement: mandatory.Data type: varchar.

Quality report header req (qualityReportsRequired). Specifies the columns required in the Quality Reports table.Role: Header.Requirement: mandatory.Data type: varchar.

Quality rerpots table column order (qualityReportsOrder). Specifies the order of the columns in the Quality Reports table.Role: Header.Requirement: mandatory.Data type: integer.

Sample relationships table (sampleRelationships). Table for recording the relationships between samples.Role: Header.Requirement: mandatory.Data type: varchar.

Sample relationships header req (sampleRelationshipsRequired). Specifies the columns required in a Sample Relationships table.Role: Header.Requirement: mandatory.Data type: varchar.

Sample relationships table column order (sampleRelationshipsOrder). Specifies the order of the columns in a Sample Relationships table.Role: Header.Requirement: mandatory.Data type: integer.

Protocols table (protocols). The table for protocols.Role: Header.Requirement: mandatory.Data type: varchar.

Protocols header req (protocolsRequired). Specifies the columns required in the Protocols table.Role: Header.Requirement: mandatory.Data type: varchar.

Protocols table column order (protocolsOrder). Specifies the order of the columns in the Protocols table.Role: Header.Requirement: mandatory.Data type: integer.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: mandatory.Data type: varchar.

Data types (dataTypes). The data type for a part.Role: Header.Requirement: mandatory.Data type: varchar.

Minimum value part support (minValue). The minimum value of part.Role: Header.Requirement: mandatory.Data type: varchar.

Maximum value part support (maxValue). The maximum value of a part.Role: Header.Requirement: mandatory.Data type: varchar.

Minimum length (minLength). The maximum value of measure.Role: Header.Requirement: mandatory.Data type: integer.

Tables

Maximum length (maxLength). The maximum length of the value of a part or measure.Role: Header.Requirement: mandatory.Data type: integer.

Polygon table

partID: polygons. The table that contains information about the geometry of a geographic area.

Polygon ID (polygonID). Unique identifier for the polygon.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Polygon Population (polyPop). An attribute of a polygon, which specifies the population of that polygon. A rough estimate of the number of human residents.Role: Header.Requirement: optional.Data type: integer.

Type of geography (geoType). Type of geography that is represented by the polygon.Role: Header.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: geoTypeCatSet

European Petroleum Survey Group Coordinates (geoEPSG). The unique EPSG code specifying a given geospatial area.Role: Header.Requirement: mandatory.Data type: float.

Well-known text (geoWKT). Well-known text of the polygonRole: Header.Requirement: mandatory.Data type: varchar.

File location of polygon (fileLocation). The location of the file containing the geometry of the polygon.Role: Header.Requirement: optional.Data type: blob.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Protocol relationships table

partID: protocolRelationships. The table that contains the organizational information and details (such as order) about a given protocol, recorded as a series of protocol steps (protocolSteps). Protocols are a group of related protocol steps. For example a protocol step can describe a sample concentration, extraction, inhibition or a measurement. Protocols can include both methodID (a procedure) and measureIDs (a measure). An example of a method is centrifugation; where the centrifugation rotation speed is described as a measure.

Tables

Protocol ID container (protocolIDContainer). Unique identifier for the protocol and the steps and other protocols that make it up.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Protocol ID object (protocolIDObj). The object of the relationship between one protocol and another protocol or step, within a given protocol ID container.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Step ID Object (stepIDObj). The object of the relationship between one protocol step and another step or protocol, within a given protocol ID container.Role: Foreign Key.Requirement: mandatoryIf.Data type: integer.

Sample relationship (sampleRelID). Attribute for specifying the relationship between a sample subject and object.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Protocol ID subject (protocolIDSub). The subject of the relationship between one protocol and another protocol or step, within a given protocol ID container.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Step ID Subject (stepIDSub). The subject of the relationship between one protocol step and another step or protocol, within a given protocol ID container.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Protocols table

partID: protocols. The table for protocols.

Source Protocol ID (sourceProtocol). A protocol that served as a basis for another protocol. This may be due to version changes, or a referencing a protocol that was later adapted.Role: Foreign Key.Requirement: optional.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: optional.Data type: varchar.

Summary (summ). Short description of the assay and how it is different from the other assay methods.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Protocol version (protocolVersion). Specifies the version of a method set.Role: Header.Requirement: optional.Data type: integer.

Tables

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Protocol steps table

partID: protocolSteps. The table for collecting metadata on individual steps in a protocol, methodological process, or assay.

Protocol Step ID (stepID). The unique identifier for a specific protocol step.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Method (method). A procedure for collecting a sample or performing a measure.Role: Foreign Key.Requirement: mandatoryIf.Data type: categorical.

Measure (measure). A measurement or observation of any substance including a biological, physical or chemical substance. An measurement is recorded from a specimen (i.e. a site, sample, person or population).Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Summary (summ). Short description of the assay and how it is different from the other assay methods.Role: Header.Requirement: optional.Data type: varchar.

Protocol step source ID (sourceStep). Specifies the protocol step which serves as a basis for a given protocol step.Role: Foreign Key.Requirement: optional.Data type: varchar.

Protocol step version (stepVer). Specifies the version of a given protocol step.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Instrument ID (instrumentID). A unique identifier for an instrument.Role: Foreign Key.Requirement: optional.Data type: categorical.

Value (value). A measure or sample taken for a validation study.Role: Header.Requirement: mandatoryIf.

Unit (unit). The units of a measurement.Role: Foreign Key.Requirement: mandatoryIf.Data type: categorical.

Aggregation (aggreation). Statistical measures used to report a measure. Each aggregation has a corresponding value.Role: Foreign Key.Requirement: mandatoryIf.Data type: categorical.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Quality reports table

partID: qualityReports. The table for recording the various quality metrics and indicators for samples and measures.

Tables

Quality report (qualityID). A unique identifier for a given quality report.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Report ID (measureRepID). Unique identifier for a measurement.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Sample ID (sampleID). Unique identifier for a sample.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Report set ID (meaureSetRepID). Unique identifier that links together a group of related measures.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Quality flag (qualityFlag). A field for reporting any quality concerns - of lack thereof - for a sample or measure.Role: Header.Requirement: mandatory.Data type: categorical.

Severity indicator (severity). An indicator of the severity or seriousness of a quality flag.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Sample relationships table

partID: sampleRelationships. Table for recording the relationships between samples. Samples can be pooled or split. The sample relationships table holds information on parent-child relationships between samples, and allow for tracking sample lineage for single and pooled samples.

Sample ID subject (sampleIDSubject). Populated by **sampleID** - this specifies the subject of a sample relationship.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Sample relationship (sampleRelID). Attribute for specifying the relationship between a sample subject and object.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Sample ID object (sampleIDObject). Populated by **sampleID** - this specifies the object of a sample relationship.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Sample report table

partID: samples. The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. Samples can be combined, split, stored and reused. The **Sample relationships** (sampleRelationships) is used to describe the relationships between samples.

Sample ID (sampleID). Unique identifier for a sample.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Tables

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site ID (siteID). Unique identifier for the location where a sample was taken.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Purpose (purpose). The reason the measure or sample was taken.Role: Foreign Key.Requirement: optional.Data type: categorical.Measure, method or attribute set: purposeCatSet

Sample material (saMaterial). Type of sample.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: sampleMatCatSet

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Sample origin (origin). An attribute of a sample specifying the origin.Role: Header.Requirement: optional.Data type: categorical.

Replicate Type (repType). Attribute of a sample, specifying whether the sample is unique, or a replicate. And if it is a replicate, what type of replicate is it?Role: Foreign Key.Requirement: optional.Data type: categorical.Measure, method or attribute set: replicateCatSet

Sample collection type (collType). The type of collection.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: collectCatSet

Collection period (collPer). Collection period. The time period over which the sample was collected, in hours. Alternatively, use collectionStart and collectionEnd.Role: Header.Requirement: mandatory.Data type: float.

Collection number (collNum). The number of subsamples that were combined to create the sample. Use NA for continuous, proportional or

passive sampling.Role: Header.Requirement: mandatory.Data type: integer.

Collection number and period (collNumPer). Composite collection number.period.Role: Header.Requirement: mandatory.Data type: varchar.

Pooled (pooled). Is this a pooled sample, and therefore composed of multiple child samples obtained at different sitesRole: Header.Requirement: optional.Data type: boolean.

Collection date time (collDT). For grab samples this is the date, time and timezone the sample was taken.Role: Header.Requirement: mandatory.Data type: datetime.

Collection date time start (collDTStart). For integrated time averaged samples this is the date, time and timezone the sample was started being taken.Role: Header.Requirement: optional.Data type: datetime.

Collection date time end (collDTEnd). For integrated time average samples this is the date, time and timezone the sample was finished being taken.Role: Header.Requirement: optional.Data type: datetime.

Date sample was sent (sentDate). The date the sample was sent for analyses at a laboratory.Role: Header.Requirement: optional.Data type: datetime.

Date sample recieved (recDate). The date the sample was received at the laboratory for analysis.Role: Header.Requirement: optional.Data type: datetime.

Reportable (reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied.Role: Header.Requirement: optional.Data type: boolean.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Tables

Notes (notes). A note used to describe details that are not captured in other attributes.
Role: Header.Requirement: optional.Data type: varchar.

Sets look-up table

partID: sets. Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together. This table manages many-to-many relationships and category recycling within the ODM.

Set (set). The unique identifier of a value set.
Role: Primary Key.Requirement: mandatory.Data type: varchar.

Set type (setType). The type of set. i.e. quality set, aggregation set, unit set.
Role: Header.Requirement: mandatory.

Set Value (setValue). The partID for the set value or category.
Role: Foreign Key.Requirement: mandatory.

Label (partLabel). A human readable label of a part.
Role: Header.Requirement: mandatory.Data type: varchar.

Part instruction (partInstr). Additional notes and instructions on how a part is used and/or defined.
Role: Header.Requirement: mandatory.Data type: varchar.

Status (status). Whether the part is still active and can be used in the most current ODM version. Values are 'active' or 'inactive'.
Role: Header.Requirement: mandatory.Data type: varchar.

First released version (firstReleased). The version in which a part was first released.
Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.
Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: optional.Data type: varchar.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Sites table

partID: sites. The table that contains information about a site; the location where an environmental sample was taken. The site of an environmental sample. Information in the site table does not regularly change. Consider using the MeasureReport table if the information changes often.

Parent Site ID (parSiteID). The siteID that is a parent to another siteID; usually the sub-site will be contained in the parentSite.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site ID (siteID). Unique identifier for the location where a sample was taken.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Polygon ID (polygonID). Unique identifier for the polygon.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site Type (siteType). Type of site or institution where sample was taken.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: siteTypeCatSet

Tables

Sample shed (sampleShed). A geographic area, physical space, or structure. A sample is taken from a sampleshed for a representative measurement of a substance(s).Role: Header.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: shedCatSet

Address ID (addressID). A unique identifier for an address.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Primary reporting authority (repOrg1). The primary or most responsible authority for routine surveillance reports or findings, or where the site is located.Role: Header.Requirement: optional.Data type: varchar.

Secondary reporting authority (repOrg2). The secondary, additional or alternative authority for routine surveillance reports or findings, or where the site is located.Role: Header.Requirement: optional.Data type: varchar.

Population Served (popServ). An attribute of a site, which specifies the population of/population served by a given site.Role: Header.Requirement: optional.Data type: categorical.

Latitude (geoLat). Geographical location, latitude in decimal coordinates, ie.: (45.424721)Role: Header.Requirement: mandatoryIf.Data type: float.

Longitude (geoLong). Geographical location, longitude in decimal coordinates, ie.: (-75.695000)Role: Header.Requirement: mandatoryIf.Data type: float.

European Petroleum Survey Group Coordinates (geoEPSG). The unique EPSG code specifying a given geospatial area.Role: Header.Requirement: mandatoryIf.Data type: float.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Translation look-up table

partID: translations. Look up table for translations of the description, label, and instruction for all parts. The default language if a translation is not specified is English.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Part ID (partID). The unique identify of any entity within the dictionary.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Label (partLabel). A human readable label of a part.Role: Header.Requirement: mandatory.Data type: varchar.

Part description (partDesc). The description of the part.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Part instruction (partInstr). Additional notes and instructions on how a part is used and/or defined.Role: Header.Requirement: optional.Data type: varchar.

First released version (firstReleased). The version in which a part was first released.Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: optional.Data type: varchar.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Page created on Thu Mar 9 00:00:00 2023.