

Public Health Environmental Surveillance Open Data Model (PHES-ODM) Documentation

Version 2.0.0-rc.3.18

2023-03-17

Table of contents

Introduction	1
How to use the documentation	2
Future documentation	2
Quick start	3
1. Get oriented to the ODM tables where you record results	3
2. View mandatory and optional data fields in the tables	
reference guide	4
3. Try entering example data in ODM Excel templates . . .	4
4. Next steps	5
Getting help	5
Background	9
Getting to know the ODM	11
What are “Parts” and “Part Types”?	11
“Sets” and the rationale behind sets	12
Why category sets are treated differently	13
The “what” and “why” of Specimen IDs	13
The “what” and “why” of Groups and Classes	14
Tables and table types within the ODM	14
Program description tables	15
Results tables	15
Look-up tables	15
Measures, Methods, and Attributes: Key parts	15
Measures	16
Methods	16

Table of contents

Attributes	16
Data quality and reportability	17
Time periods for samples and measures	17
Dates within the ODM	18
Translation and language capabilities	18
Questions and online community	19
How to use the ODM	21
How-To guides	21
1) How to use templates	21
2) How to report Measures and Measure Sets:	23
Quick Start	24
Measures Template	24
Measure Sets Template	25
Detailed description	25
Measures template	25
Measure Sets Template	30
3) How to report Samples and Sample Relationships:	31
Quick Start	31
Samples Template	32
Sample Relationships Template	33
Detailed Description	33
Samples Template	33
Sample Relationships Template	35
4) How to report Protocols, Protocol Steps and Protocol	
Relationships:	37
Quick Start	37
Protocol Steps Template	37
Protocols Template	38
Protocol Relationships Template	39
Detailed Description	40
Protocol Steps Template	40
Protocols Template	41
Protocol Relationships Template	43

Table of contents

5) How to report Metadata:	44
Quick Start	45
Addresses Template	45
Contacts Template	46
Datasets Template	46
Instruments Template	47
Organizations Template	48
Polygons Template	48
Quality Reports Template	49
Sites Template	49
Detailed Description	50
Addresses Template	50
Contacts Template	51
Datasets Template	53
Instruments Template	54
Organizations Template	56
Polygons Template	57
Quality Reports Template	58
Sites Template	59
Example data for each of the 11 report tables.	61
 I. Reference guide	 63
 Parts	 67
a1306s delta-variant gene target	68
a1918v delta-variant gene target	68
a2710t omicron-variant gene target	68
a63t omicron-variant gene target	69
a67v omicron-variant gene target	69
Atline Analyzer	69
Absolute humidity	70
Absolute humidity unit set	70
Academic institution	70

Table of contents

Number of active cases	71
Active	71
Analysis date	71
Analysis date end	72
Analysis date start	72
Address Line 1	72
Address Line 2	73
Address table	73
Addresses table column order	74
Address table required headers	74
Address ID	74
Administrative regions	75
Air filter	75
Aggregation scale	75
Aggregation	76
Aggregations	76
Aggregation scales	77
Aggregation set	77
Aggregation sets	77
Air compartment	78
Air compartment set	78
Airplane	78
Airport	79
Airport lists worksheet	79
Air and surface compartment set	79
Environmental temperature	80
Air and water compartment set	80
Alias depreciated	81
Alias ID depreciated	81
All domains	81
Alleles class	82
Allele unit set	82
Access to all org	82
Aloh3 precipitation	83

Table of contents

Assay method id default depreciated	83
Amicon ultrafiltration	84
Amicon filter, extract with MP96	84
Amplicon sequencing	84
AND Boolean aggregation	85
Any compartment set	85
Any specimen set	85
Area proportional sample	86
Arrival temperature	86
Use of the articV3 primer	87
Use of the articV4 primer	87
Assay ID	87
Attributes	88
Bacteria Class	88
Bacteria unit set	89
Miscellaneous bacteria group	89
Bead beating, extract with MP96	89
bcov	90
bcov culture spike target	90
bcov vaccine spike target	90
Beef extract flocculation	91
Is Before	91
Beta	92
Biologic	92
Digital droplet emulsification PCR	92
Binary Large Object (BLOB) data type	93
Breadth of coverage ($\geq 5\times$ depth)	93
5-Day carbonaceous biochemical oxygen demand	93
5-day total biochemical oxygen demand	94
Boolean data type	94
Boolean aggregation set	95
Boolean value set	95
Bovine respiratory syncytial virus group	95
brsv culture spike target	96

Table of contents

BRSV-N	96
brsv vaccine spike target	96
Branch sewer pipeline	97
Building cleanout	97
Capacity class	97
Capacity unit set	98
Case report date	98
Categorical data type	99
Categories	99
Long-term acute care hospital	99
Child day care	100
Degrees Celcius	100
Centrifugation	101
Centricon ultrafiltration	101
Ceres nanotrap	101
CFU per 100 ml	102
Changes column	102
Chemical	102
Chemagic viral dna/rna 300 kit	103
Child relationship	103
City	103
Clarified sample	104
Class	104
Classes	105
Cloudy	105
Centimetres	105
Carbon Dioxide	106
Chemical Oxygen Demand	106
Collection group	106
Collection date time	107
Collection date time end	107
Collection date time start	107
Sample collection set	108
Collection number	108

Table of contents

Collection number and period	109
Collection period	109
Sample collection type	110
Co-located sample	110
Column name depreciated	111
Composite sample - archival	111
Composite grab sample of 3	112
Composite grab sample of 3 depreciate	112
Composite 3hr grab sample	112
Composite 3hr grab sample depreciate	113
Composite 8hr grab sample depreciated	113
Composite 8hr grab sample deprecated	113
Compartment	114
Compartments	114
Compartment sets	115
Sample concentrate	115
Confirmed Cases Case report date	115
Is concurrent to	116
Water conductivity	116
Conductivity class	117
Conductivity unit set	117
Number of confirmed cases	117
Confirmed cases episode date	118
Confirmed Cases Onset date	118
Contact ID	118
Contact name	119
Contact table	119
Contacts table column order	119
Contact table required headers	120
Contact table set	120
Confirmed cases test date	121
Cphd ID	121
COSCa ball	121
Country	122

Table of contents

Countries or sovereign states	122
Districts, counties, regions	122
Covid-19	123
SARS-CoV-2 measure	123
SARS-CoV-2-B.1.1.7	124
SARS-CoV-2-B.1.351	124
SARS-CoV-2-E	124
SARS-CoV-2-N1	125
SARS-CoV-2-N2	125
SARS-CoV-2-N3	125
SARS-CoV-2-P.1	126
SARS-CoV-2-RdRp	126
Covid-19 population measurement date	126
Cphd ID	127
crAssphage-N	127
Crosswalk table set	127
Cycle threshold or quantification cycle (Ct or Cq)	128
Cumulative count	128
Custody Contact ID	129
Data custodian ID	129
d377y delta-variant gene target	129
d3g omicron-variant gene target { #d3g }	130
d63g delta-variant gene target	130
d796y omicron-variant gene target { #d796y }	130
d950n delta-variant gene target	131
Daily count	131
Dataset creation date	132
Dataset ID	132
Dataset table	132
Datasets table column order	133
Dataset table required headers	133
Data types	133
Date	134
Datetime data type	134

Table of contents

Days	135
ddcov_e sars-cov-2 gene target	135
ddcov_n sars-cov-2 gene target { #ddcovN }	135
Deaths	136
del143/145	136
del 157/158	136
del2084/2084	137
del212/212	137
del3674/3676	138
del69/70	138
Delta	138
Depreciated	139
Derived sample	139
Description	139
Description of change	140
Desk or counter	140
Dewatered solids	141
Detected	141
Access to details deprecated	141
Development	142
Dictionary reference table set	142
Additional dictionary sheets	143
Dictionary set type	143
Dilution factor	143
Point dillutions	144
Dilution Class	144
Diseases (human) class	144
Dissolved gas concentration unit set	145
Depth of coverage	145
Domain	146
Domains	146
Higher education domitory or residential building	146
e156g delta-variant gene target	147
e484a omicron-variant gene target { #e484a }	147

Table of contents

Escherichia coli	147
Education	148
Efficiency	148
Contact email	149
European Nucleotide Association	149
European Nucleotide Association (ENA) - notes	149
Rainfall	150
Ground snow depth	150
Snowfall	150
Episode date	151
ERD table set	151
Sarbecovirus-specific E sars-cov-2 gene target	151
Estimated frequency of reads	152
Estuary, natural water body	152
Experiment Failed	153
4s method	153
Nucleic acid extraction method	153
Nucleic Acid Extraction set	154
Extraction volume of sample	154
Fecal matter	154
FALSE	155
Field sample	155
Field sample replicate	156
File location of polygon	156
Filtration	156
First Nation	157
First name of contact	157
First released version	157
Foreign key	158
AI - Inhibition present but addressed	158
B - Trace levels of contamination	159
FI - Inhibition present and unaddressed	159
J - Weak signal extrapolation	159
ND - Non-detect	160

Table of contents

UJ - Trace signal extrapolation	160
UQ - Unquantifiable	161
Float data type	161
Flow-normalized mean	161
Floor	162
Flow rate	162
Flow class	163
Flow proportional 24hr sample depreciated	163
Flow proportional sample	163
Volume flow rate unit set	164
Flow volume	164
Influenza virus measure	164
Influenza virus A1	165
influenza virus A2	165
Influenza virus B	166
Fluidigm digital PCR	166
Foggy	166
Fraction analyzed	167
Sample fraction set	167
Freyja Script	167
F-Specific RNA bacteriophages	168
F-Specific RNA bacteriophages, G2	168
Sample frozen	168
Field sample temperature	169
Full dictionary sheets set	169
Funder Contact ID	169
Funding agency ID	170
g215c delta-variant gene target	170
g339d omicron-variant gene target { #g339d }	171
g496s omicron-variant gene target { #g496s }	171
g662s delta-variant gene target	171
Gamma	172
Gas class	172
Gene copies per copy of crAssphage	172

Table of contents

gene copies per day per 100,000	173
Gene copies per gram solids	173
Gene copies per L	173
Gene copies per mL	174
Gene copies per PMMoV copy	174
Genetics unit set	174
General missingness set	175
Generic quality flag set	175
European Petroleum Survey Group Coordinates	175
Latitude	176
Longitude	176
Type of geography	177
Geographic set	177
Well-known text	177
Gram per cubic metre	178
Geometric mean	178
Government agency	178
Grams	179
Gravity settling	179
Grab sample	179
Gromstole 1.0 Script	180
Group	180
Groups	181
Header	181
Health administration or planning agency	181
Heat inactivated sars-cov-2 virus spike target	182
hep g armored rna	182
hep g armored rna spike target	182
High severity	183
Sewer Network Health Region	183
Hand Measurement	183
Hollow fiber dead end ultrafiltration	184
General Hospital Admissions	184
Hospital Census	184

Table of contents

Hospital	185
Holding tank wastewater depreciated	185
Hours	186
Holding tank wastewater	186
Holding tank	186
Human compartment	187
Human compartment set	187
Humidity class	187
i1566v omicron-variant gene target { #i1566v }	188
i3758v omicron-variant gene target { #i3758v }	188
i82t delta-variant gene target	188
International Classification of Diseases	189
Intensive care unit patients	189
Index	189
Inhibition set	190
Inhibition measure	190
Inhibition method	191
Innovaprep ultrafiltration	191
Input	191
Instrument ID	192
Instrument table	192
Instruments table column order	193
Instrutment table required headers	193
Other instrument	193
Instrument Type	194
Describe other instrument type, if applicable	194
Instrument set	194
Integer data type	195
Intercept	195
Combined ip2 and ip4 sars-cov-2 gene target	195
ISO639-1	196
ISO639-2B	196
ISO639-2T	197
ISO639-3	197

Table of contents

ISO639-6	198
k856r omicron-variant gene target { #k856r }	198
Kilogram per second	198
Kilolitres	199
Litres	199
l452r delta-variant gene target	199
l981f omicron-variant gene target { #l981f }	200
Laboratory	200
Lab ID default deprecated	201
Laboratory duplicate	201
Lab ID	201
Lagoon system	202
Lake, natural water body	202
Lambda	202
Language ID	203
Language family	203
Language name	203
Language Look-up table	204
Language table column order	204
Language table required headers	205
Last edited	205
Last name of contact	206
Last updated version	206
LaTeX expression	206
Low breadth of coverage	207
Laboratory control sample duplicate	207
Low depth of coverage	207
Leaked sample	208
Link	208
License	209
Linear scale aggregation set	209
Liquid fraction	209
List worksheet	210
List set	210

Table of contents

Access to local ha	210
Limit of detection (LOD)	211
Limit of detection (LOD) - sewage	211
Logarithmic scale aggregation set	212
Dictionary tables	212
Limit of quantification (LOQ)	212
Low severity	213
Lower limit of a 95% confidence interval	213
Low-volume sample	214
Long-term care facility	214
Long-term care - assisted living or retirement home	214
Other long-term care	215
Life technologies digital PCR	215
Lysis buffer	215
Cubic metres per day	216
Cubic metres per hour	216
Cubic metres per second	216
Manufacturer	217
Mate pair layout	217
Maximum length	217
Maximum value	218
Maximum value part support	218
Arithmetic mean	218
Measurement group	219
Other Measure	219
Measures quality set	220
Measure	220
Measures	220
Report ID	221
Measure report table	221
Measure set report table	222
Measure sets table column order	222
Measure Set table required headers	222
Measures table column order	223

Table of contents

Measures table required headers	223
Report set ID	223
Median	224
Normalized arithmetic mean	224
Other Measure Description	224
Measure license	225
Method	225
Concentration method	225
Concentration method set	226
Methods	226
Membrane filtration with acidification and mgcl2	227
Membrane filtration with acidification and mgcl2, mem- brane recombined with separated solids	227
Membrane filtration with addition of mgcl2	227
Membrane filtration with mgcl2 and separated solids	228
Membrane filtration with no amendment	228
Membrane filtration with no amendment, membrane recom- bined with separated solids	228
Membrane filtration with sample acidification	229
Membrane filtration with sample acidification, membrane recombined with separated solids	229
Millions of gallons per day (MG/D)	230
Milligrams per litre	230
Murine Hepatitis Virus	230
Multiple issues	231
Mid-level severity	231
Minimum length	231
Minutes	232
Minimum value	232
Minimum value part support	232
Miscellaneous attribute group	233
Miscellaneous measure group	233
Missingness	233
Missingness sets	234

Table of contents

Mixed/homogenized sample	234
Millilitres	235
Megalitres per day (ML/d)	235
Millimetres	235
Measure, method, or attribute	236
Measure, method, or attribute sets	236
Model	236
Months	237
Moore swab passive sample	237
Metres per second	237
ms2 coliphage	238
ms2 coliphage spike target	238
Matrix spike duplicate	238
Moore swab passive sample depreciated	239
Major sewer pipeline	239
Mu	240
murine coronavirus	240
murine coronavirus spike target	240
Multiple fraction	241
Multiple Purpose	241
Municipalities or communes	241
Municipality	242
Mutations class	242
N sars-cov-2 gene target	243
Combined N1 and N2 sars-cov-2 gene target	243
n211i omicron-variant gene target { #n211i }	243
n679k omicron-variant gene target { #n679k }	244
n856k omicron-variant gene target { #n856k }	244
n969k omicron-variant gene target { #n969k }	244
Not applicable	245
Aggregation not applicable	245
Aggregation scale not applicable	245
Aggregation set not applicable	246
Class not applicable	246

Table of contents

Compartment not applicable	246
Compartment set not applicable	247
Domain not applicable	247
Group not applicable	248
Name	248
Not a number	248
Nomenclature not applicable	249
Specimen not applicable	249
Specimen set not applicable	250
Native Name	250
Unit not applicable	250
Unit set not applicable	251
National Center for Biotechnology Information	251
NCBI - notes	251
Neighbourhood	252
Local administrative units or neighborhoods	252
Population of newly vaccinated persons	252
NextClade nomenclature	253
Normalized geometric mean	253
Ammonium Nitrogen	253
niid_2019-ncov_n sars-cov-2 gene target	254
No quality concerns	254
Sample not labelled	255
No liquid concentration, liquid recombined with separated solids	255
Nomenclature	255
Nomenclatures	256
NORMAN	256
NORMAN - notes	256
Norovirus G1	257
Norovirus G2	257
Notes	257
Sample missing time stamp	258
Not reported	258

Table of contents

Not reported, Not applicable missing set	259
Sarbecovirus-specific N sars-cov-2 gene target	259
Nephelometric turbidity unit	259
Nuclisens automated magnetic bead extraction kit	260
Nuclisens manual magnetic bead extraction kit	260
Null	261
National Wastewater Surveillance System	261
NWSS - notes	261
NWSS - process	262
Water	262
coronavirus OC43	262
oc43 spike target	263
Ocean, natural water body	263
Other Collection depreciated	263
Lab Analysis	264
Omicron BA.1	264
Omicron BA.2	265
Omicron BA.2.75	265
Omicron BA.4	265
Omicron BA.5	266
Onset date	266
Online Sensor	266
Ontology reference	267
Orthophosphates	267
Other residential building	268
OR Boolean aggregation	268
ORF1a sars-cov-2 gene target	268
ORF1ab sars-cov-2 gene target	269
ORF1b sars-cov-2 gene target	269
Organization ID	269
Organization table	270
Organizations table column order	270
Organziation table required headers	270
Organization level	271

Table of contents

Organization level set	271
Organization sector	272
Organization sector set	272
Organization Type	272
Organization type set	273
Sample origin	273
Sample Origin set	273
Other Collection depreciated	274
Other aggregation scale	274
Other	274
Other aggregation set	275
Other aggregation deprecated	275
Access to other prov deprecated	276
Other report table set	276
Variants Other	276
Other Site - sample depreciated	277
Other Site - Site depreciated	277
Other Unit deprecated	278
Outbreak	278
Outbreak End	279
Outbreak - on-going	279
Outbreak class	279
Outbreak set	280
Outbreak start	280
p100l delta-variant gene target	280
p2046l delta-variant gene target	281
p2287s delta-variant gene target	281
p3395h omicron-variant gene target { #p3395h }	282
p681r delta-variant gene target	282
Paired Layout	282
Pangolin nomenclature	283
Parent dataset ID	283
Parent sample ID	283
Parent Site ID	284

Table of contents

Part description	284
Part identifier	284
Part instruction	285
Label	286
Parts Look-up table	286
Parts table column order	286
Parts table required headers	287
Part types	287
Primary Clarifier Effluent depreciated	288
Postal or Zip Code	288
qPCR Class	288
PCR method	289
PCR quality set	289
PCR sequencingsSelection method	289
PCR Method set	290
Primary Clarifier Sludge depreciated	290
Primary clarifier effluent	290
Polyethyleneglycol (PEG) precipitation	291
Percent	291
Percent recovery	291
pH	292
Access to phac	292
PHAGE	293
PHAGE - notes	293
pH class	293
Phenol chloroform	294
Pseudomonas virus phi6	294
phi6 spike target	294
Contact phone	295
Total phosphorous	295
Total Phosphates	295
PHRED quality score	296
Physical property	296
Physical class	296

Table of contents

Primary key	297
PMMoV- and flow-normalized mean	297
PMMoV-normalized mean	297
Population	298
Post-Grit depreciated	298
Polygon ID	299
Polygon table	299
Polygons table column order	299
Polygon table required headers	300
Polygon Population	300
Pooled	300
Type of date for case reporting depreciated	301
Case reporting date set	301
Population equivalents	302
Population group	302
Population Served	302
Population unit set	303
Population or sample set	303
Population or site set	303
Population specimen set	304
Percent positive	304
parts per million	304
PMMoV-CP	305
Percent positivity rate	305
Percent primary sludge	305
Precipitation class	306
Precipitation unit set	306
Pretreatment	307
Genetic primer	307
Genetic primer set	307
Private sector	308
Processing group	308
Program decription tables	308
promega automated tna kit	309

Table of contents

promega ht tna kit	309
promega manual tna kit	310
promega wastewater large volume tna capture kit	310
Proportion of total	310
Proportion of variant in sample	311
Protocol ID	311
Protocol ID container	311
Protocol ID object	312
Protocol ID subject	312
Protocol relationships table	313
Protocol relationships table column order	313
Protocol Relationship table required headers	314
Protocol Relationships set	314
Protocols table	314
Protocols table column order	315
Protocols table required headers	315
Protocol steps table	315
Protocol steps table column order	316
Protocol steps table required headers	316
Protocol table set	316
Protocol version	317
Access to prov ha	317
Provisional report	318
Primary clarifier sludge	318
Pumping station	318
Post-grit	319
Pretreatment description - depreciated	319
Number of positive tests	320
Public healith agency	320
Access to public	320
Puro virus	321
Puro virus spike target	321
Purpose	321
Purpose set	322

Table of contents

q19e omicron-variant gene target { #q19e }	322
q493r omicron-variant gene target { #q493r }	322
q498r omicron-variant gene target { #q498r }	323
q954h omicron-variant gene target { #q954h }	323
Quaility assurance method	323
Quality concerns	324
qiagen allprep dna/rna kit	324
qiagen digital PCR	325
qiagen ez1 virus mini kit v2.0	325
qiagen allprep powerfecal dna/rna kit	325
qiagen allprep powerviral dna/rna kit	326
qiagen powerwater kit	326
qiagen qiaamp buffers with epoch columns	326
qiagen rneasy kit	327
qiagen rneasy powermicrobiome kit	327
Quantitative PCR	327
Qualitative	328
Quality control	328
Quality flag	328
Quality report ID	329
Quality indicators	329
Quality indicator sets	329
Quality reports table	330
Quality rerpots table column order	330
Quality report table required headers	331
Quality set	331
Quantitative	331
R squared	332
r203m delta-variant gene target	332
Raindance digital PCR	332
Rainy	333
Random sequencing selection method	333
Ratio	334
Raw wastewater	334

Table of contents

Raw sewage at site	334
Raw sewage downstream from a site	335
Raw sewage upstream from a site	335
IP2 rdrp sars-cov-2 gene target	335
IP4 rdrp sars-cov-2 gene target	336
Date sample recieved	336
Recovered patients	337
Reference link	337
Regular	337
Regular reports table set	338
Relative humidity	338
Relative humidity unit set	339
Representative grab sample	339
Replicate Type set	339
replicate number	340
Primary reporting authority	340
Secondary reporting authority	341
Reportable	341
Report date	342
Reporter ID depreciated	342
Reporter table depreciated	343
Replicate Type	343
Resuspend COSCa filter collection	343
Results tables	344
Retention pond	344
River, natural water body	345
Role of contact	345
Rate per 100,000	345
Raw post-pasteurized wastewater	346
Respiratory syncytial virus	346
s2083i omicron-variant gene target { #s2083i }	346
s371l omicron-variant gene target { #s371l }	347
s373p omicron-variant gene target { #s373p }	347
s375f omicron-variant gene target { #s375f }	347

Table of contents

s477n omicron-variant gene target { #s477n }	348
Sample	348
Sample material	348
Sample ID	349
Sample ID object	349
Sample ID subject	350
Sample material set	350
Sample quality set	350
Sample relationships table	351
Sample relationships table column order	351
Sample relationships table required headers	352
Sample relationship	352
Sample relationships set	352
Sample report table	353
Sample shed	353
Samples table column order	354
Sample table required headers	354
Collection other depreciated	354
SARS-CoV-2	355
Sample specimen set	355
Sparse coverage	355
Secondary clarifier effluent depreciated	356
Standard curve frequency	356
School	357
Sequencing script version set	357
Genetic sequencing script version	357
Secondary clarifier sludge depreciated	358
Standard deviation	358
Normalized standard deviation	359
Sea, natural water body	359
See unit (aggregation scales)	359
See unit (data type)	360
See unit (value)	360
Secondary clarifier effluent	361

Table of contents

Access to self	361
Date sample was sent	361
Septic tank wastewater	362
Septic tank	362
Sequencing Layout	362
Sequencing Layout set	363
Sequencing quality set	363
Sequencing selection method	363
Sequencing selection method set	364
Sequencing Strategy	364
Sequencing Strategy set	364
Genetic sequences class	365
Serial dilution	365
Sewer Sediments Depreciated	366
Set ID	366
Sets look-up table	366
Sets table column order	367
Sets table required headers	367
Settled solids	367
Set type	368
Set Value	368
Severity indicator	369
Severity set	369
Sewer network file link depreciated	369
S sars-cov-2 gene target	370
Sampleshed set	370
Ship	370
Shipped on ice	371
Short names	371
Site	371
Single	372
Single Layout	372
Site or sample specimen set	373
Site specimen set	373

Table of contents

Site ID default	373
Site features group	374
Site ID	374
Site measure ID	374
Sites table	375
Sites table column order	375
Sites table required headers	376
Site Type	376
Site set	376
Skimmed milk flocculation	377
slope	377
Sewer network file link depreciated 2	377
Snowy	378
Solid fraction	378
Solid seperation	378
Solid separation set	379
Source Protocol ID	379
Protocol step source ID	380
Specific humidity	380
Specific humidity unit set	380
Specifies	381
Specimen	381
Specimens	381
Specimen sets	382
Speed class	382
Spike matrix and recovery	383
Spike material	383
Spike material set	383
Recovery efficiency spike target	384
Recovery efficiency spike target set	384
Sample spilled	385
Secondary clarifier sludge	385
Social services shelter	385
Stabilization pond	386

Table of contents

Standard concentrations class	386
Standard curve class	386
Departments, states, or provinces	387
State, Province, or Region	387
Status	387
Status set	388
Standard concentration unit set	388
Storage temp	389
Protocol Step ID	389
Step ID Object	389
Step ID Subject	390
Method step parent ID	390
Protocol step version	391
Sewer network file blob	391
Storage time	392
Sub-variant or lineage	392
Summary	392
Summary worksheet	393
Sunny	393
Surface compartment	393
Other surface	394
Surface compartment set	394
Surface and water compartment set	394
Surface swab	395
Sewage truck	395
Sewer sediment	396
Sewer catchment area	396
Synthetic sample	396
t19r delta-variant gene target	397
t3646a delta-variant gene target	397
t547k omicron-variant gene target { #t547k }	397
t9i omicron-variant gene target { #t9i }	398
Tables	398
Table column set	398

Table of contents

TaqPath N sars-cov-2 gene target	399
TaqPath S sars-cov-2 gene target	399
Sample temperature	400
Temperature class	400
Temperature unit set	400
Nucleic acid template volume	401
Airport terminal	401
Test date	401
Number of tests performed	402
Testing	402
thermo magmax microbiome ultra nucleic acid isolation kit	402
Time class	403
Time proportional sample	403
Time unit set	404
Total Kjeldahl Nitrogen	404
Total Nitrogen	404
Time Proportional 24hr sample	405
Translation look-up table	405
Translation table column order	405
Translation table required headers	406
trizol, zymo mag beads w/ zymo clean and concentrator . .	406
TRUE	407
Total solids concentration	407
Total suspended solids	407
Concentration of total suspended solids	408
Turbidity	408
Turbidity class	408
Turbidity unit set	409
Tube weight empty	409
Tube weight full	409
Type other depreciate	410
Type other depreciate2	410
University campus	410
Undisclosed	411

Table of contents

Unique sample	411
Unit	412
Unitless measure	412
Unitless unit set	412
Units	413
Unit set	413
Unit sets	413
Ultracentrifugation	414
Upper limit of a 95% confidence interval	414
Upstream sites	415
Micro-Siemens per centimetre	415
U site measure ID	416
v2930l delta-variant gene target	416
Validation study	416
Value	417
Variable character data type	417
Frequency of variants detected	417
Variants class	418
Population with 1 dose of vaccine	418
Population with 2 doses of vaccine	418
Population with 3 doses of vaccine	419
ODM V1 category	419
ODM V1 location	419
ODM V1 table	420
ODM version 1 to 2 changes	420
ODM V1 variable	420
Miscellaneous viruses group	421
Size in volume	421
Volume proportional sample	421
Volume unit set	422
Volatile suspended solids	422
Volatile suspended solids - ignition	423
Water compartment	423
Water depreciated	423

Table of contents

Water compartment set	424
Wastewater temperature	424
Weather	425
Weather class	425
Weather set	425
Whole genome sequencing	426
WHO nomenclature	426
Wide 95% interval	427
Wide name table	427
Wind Speed	427
Wind speed unit set	428
Measure identification (wide table)	428
Sample stored incorrectly	428
Sample stored at wrong temperature	429
World Sphere (W-SPHERE)	429
World Sphere (W-SPHERE) notes	429
Wastewater treatment plant	430
Back-up wastewater treatment site	430
Wastewater treatment plant designed capacity	431
Industrial wastewater treatment plant	431
Municipal wastewater treatment plant for combined sewage	431
Municipal wastewater treatment plant for sanitary dewage	
only	432
y505h omicron-variant gene target { #y505h }	432
Years	432
zyzo environ water rna kit/ zyzo environ water rna kit	
(cat. r2042)	433
zyzo quick-rna fungal/bacterial miniprep #r2014	433

Sets	435
Aggregations	435
Aggregation scales	437
Aggregation sets	438
Boolean aggregation set	438

Table of contents

Linear scale aggregation set	438
Logarithmic scale aggregation set	440
Aggregation set not applicable	441
Other aggregation set	441
Attributes	443
Classes	459
Compartments	462
Compartment sets	462
Air compartment set	462
Air and surface compartment set	463
Air and water compartment set	463
Any compartment set	463
Human compartment set	464
Compartment set not applicable	464
Surface compartment set	465
Surface and water compartment set	465
Water compartment set	465
Data types	466
Domains	467
Groups	467
Measures	469
Methods	483
Missingness	486
Missingness sets	486
General missingness set	486
Not reported, Not applicable missing set	487
Measure, method, or attribute sets	487
Boolean value set	488
Sample collection set	488
Nucleic Acid Extraction set	489
Sample fraction set	492
Geographic set	492
Inhibition set	493
Instrument set	493

Table of contents

Concentration method set	494
Organization level set	495
Organization sector set	495
Organization type set	497
Sample Origin set	497
Outbreak set	498
PCR Method set	498
Genetic primer set	499
Protocol Relationships set	499
Purpose set	500
Replicate Type set	501
Sample material set	501
Sample relationships set	503
Sequencing script version set	504
Sequencing Layout set	504
Sequencing selection method set	505
Sequencing Strategy set	505
Severity set	506
Sampleshed set	506
Site set	508
Solid separation set	510
Spike material set	512
Recovery efficiency spike target set	513
Status set	514
Table column set	515
Weather set	515
Nomenclatures	516
Quality indicators	517
Quality indicator sets	519
Generic quality flag set	519
Measures quality set	520
PCR quality set	522
Sample quality set	523
Sequencing quality set	524

Table of contents

Specimens	525
Specimen sets	525
Any specimen set	526
Specimen set not applicable	526
Population or sample set	526
Population or site set	527
Population specimen set	527
Sample specimen set	527
Site or sample specimen set	528
Site specimen set	528
Tables	528
Units	531
Unit sets	537
Absolute humidity unit set	537
Allele unit set	537
Bacteria unit set	538
Capacity unit set	539
Conductivity unit set	539
Dissolved gas concentration unit set	539
Volume flow rate unit set	540
Genetics unit set	540
Unit set not applicable	541
Population unit set	541
Precipitation unit set	544
Relative humidity unit set	544
Specific humidity unit set	544
Standard concentration unit set	545
Temperature unit set	545
Time unit set	545
Turbidity unit set	546
Unitless unit set	546
Volume unit set	546
Wind speed unit set	547

Table of contents

Tables	549
Table types	549
Header roles	550
Mandatory, optional, mandatoryIf	551
Table list	551
Address table	552
Contact table	554
Dataset table	555
Instrument table	556
Language Look-up table	558
Measure report table	559
Measure set report table	562
Organization table	563
Parts Look-up table	564
Polygon table	573
Protocol relationships table	575
Protocols table	576
Protocol steps table	577
Quality reports table	579
Sample relationships table	580
Sample report table	580
Sets look-up table	583
Sites table	584
Translation look-up table	586

Introduction

Make collecting and sharing high-quality, interoperable environmental surveillance data as easy as possible.

Warning

This documentation for the PHES-ODM Version 2.0.0 release candidate. Please be aware that changes will occur regularly in the next weeks. If there is something you'd like to see in the documentation that you cannot currently find, please please [create an issue on our GitHub repository](#) or post about it on [our Discourse page](#).

The Public Health and Environmental Surveillance Open Data Model (PHES-ODM or ODM) is an open-source data model and set of tools aimed at promoting global collaboration and use of environment and public health data. The ODM model facilitates interoperable environmental surveillance and research by improving quality reporting data, enhancing data management, and supporting data provenance and sharing.

This documentation is intended for those who want to understand how the ODM model works, its organizational structure, and how to use it effectively for their program. It is also a useful resource for those who currently use a different data dictionary and wish to share their data with a program that uses the ODM. If you need to share your data with another program that uses a different data dictionary, this documentation can help you compare and create a common cross-walk between dictionaries. The ODM starts as a conceptual model that represents real-world processes

Introduction

through data elements and data relationships, which are common to all environmental and public health surveillance systems.

How to use the documentation

There are three document sections, each with a different purpose.¹

Getting to know the ODM is oriented to understanding the ODM. This section explains the ODM, its structure, and how the model was designed and developed.

How to use the ODM are step-by-step guides to help you use different aspects of the ODM such as report templates.

Reference is a description of each part of the ODM. It serves as a reference for the other sections, and the definitive record of each part, table and set – along with their attributes.

Future documentation

Tutorials and videos are being develop. Tutorials will help newcomers get started. They will be lessons on how to:

- Describe a protocol using protocol tables.
- Make wide names and report templates using wide names.

¹[David Laing and Divio's Grand Unified Theory of Documentation.](#)

Quick start

You can get to know the ODM by imaging how to record a polymerase chain reaction (PCR) test result for SARS-CoV-2 from a wastewater sample.

1. Get oriented to the ODM tables where you record results

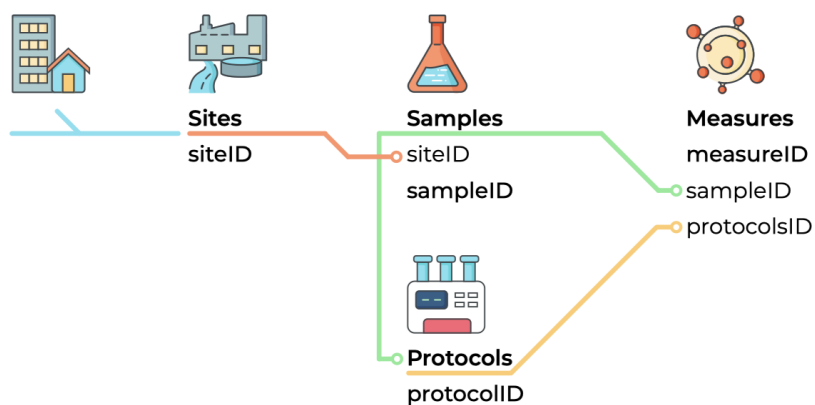


Figure 1.: **Figure 1 - Key ODM tables for recording data**

The **Sample** report table ([samples](#)) and the **Measure** report table ([measures](#)) are where you record information about samples and measures.

2. View mandatory and optional data fields in the tables reference guide

You can view mandatory and optional data field here: ([samples](#)) and ([measures](#)). ‘Mandatory’ means they are fields are needed to ensure interoperability.



What are ODM parts and how to find them?

In the ODM, tables are an example of an part. A part name can be identified by its text style like **this**. The part identifier often follows the part in parenthesis, with a link to the part in the reference section. So, **Sample report table** is the part name, and ([samples](#)) is the part ID and link to the documentation for that part.

You can find out more about parts on the [reference guide section](#).

3. Try entering example data in ODM Excel templates

Keep in mind, ODM templates are only examples of how to use the ODM dictionary to record data. Templates are available at [OSF.io](#). The ‘[How to use the ODM](#)’ section describes how to enter information into tables.

Storing sample data in an ODM template

A **sample** ([sampleID](#)) is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. The sample template stores one sample per row. The **Sample report table** has seven mandatory fields that are listed in the [reference[(tables.qmd#samples)] for the table.

<https://vimeo.com/802349692>

Storing measure data in an ODM template

A **measure** ([measure](#)) is an observation or measurement of any substance, including biological, physical, or chemical substances, obtained from a specimen such including a site, sample, person, or population. Each measure is stored in the **Measures** table template as a separate row, following the “long” table format, which is the primary method for storing measures. For example, if you recorded PCR results for two SARS-CoV-2 regions, covN1 and covN2, you would record each measurement on a separate line. Alternatively, you can use a “wide” table format to store multiple measures in a single row or create tables that combine both long and wide headers. The [European Union airport template](#) is a hybrid template that exemplifies this (found in the “Airports” folder).

(Video)

4. Next steps

In addition to recording measures from samples, you can also record measures for a specific site, such as wastewater flow rate ([flowRate](#)), or the population served by the site, such as the Covid-19 hospitalization rate for a wastewater treatment plant’s watershed or municipality. **Figure 2** illustrates the these three types of measurements: **sample**, **site**, and **population**.

The “who, where, and why” of these measurements is stored in contact tables, while the “how,” which refers to the measurement protocol or method, is stored in the protocols tables.

Getting help

We hope that the documentation serves its purpose and that you find it useful in figuring out how to best use the ODM and its suite of tools. If

Introduction

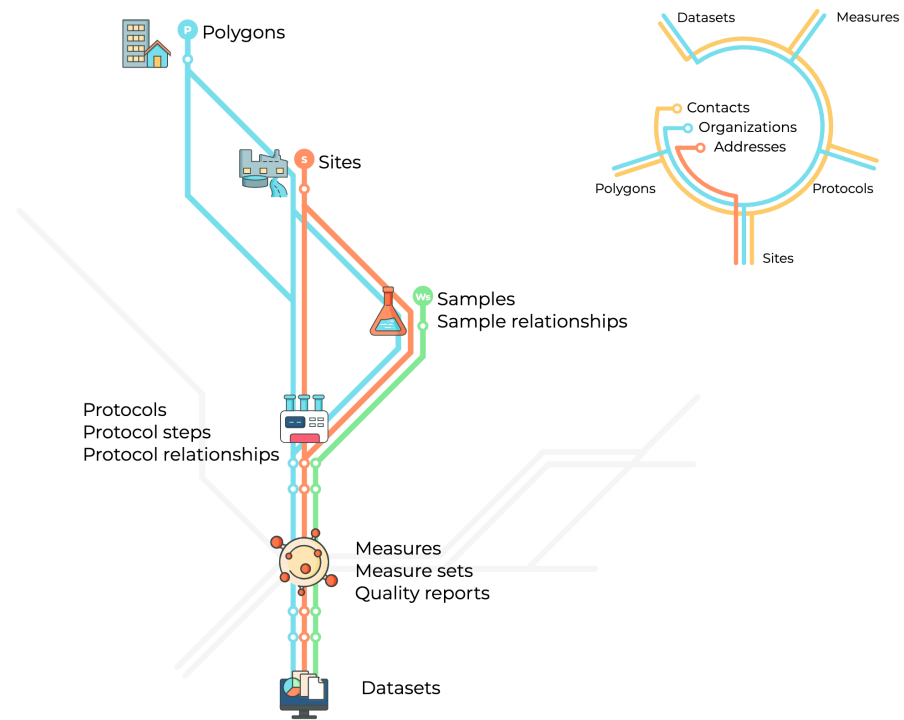


Figure 2.: **Figure 2 - ODM tables and 'subway' lines for storing measures**

Quick start

you do encounter any issues or questions that do not seem to have an easy answer, please [create an issue on our GitHub repository](#) or consult [our Discourse page](#).

Background

The ODM is designed to support a broad range of users in the field of environmental surveillance, including wastewater laboratories, municipalities, public health agencies, researchers, and citizen scientists. Its goal is to provide a comprehensive and extensible approach to data collection and management across the lifecycle of environmental surveillance, from sample collection and analysis to public reporting.

The ODM dictionary includes provisions for a wide range of biological, chemical, and physical properties, such as antimicrobial resistance, drugs, and toxins, as well as metadata for laboratory protocols that may affect measurement consistency across laboratories. The ODM supports laboratory information and management systems (LIMS) with measure to describe ensure data quality control and provenance.

To support the open science principle of FAIR data (Findable, Accessible, Interoperable, and Reusable), the ODM strives to make its data easily discoverable, accessible, and usable by all stakeholders.² It also adopts an open-source approach, allowing anyone to inspect, modify, and enhance the codebase in a transparent and collaborative manner on platforms like GitHub.

While the ODM is a complex dictionary, it includes tools to support its use, such as data validation, input templates, and SQL database definitions for data storage. Additionally, it supports collaborations among different organizations seeking to share and combine samples and measurements. As an open model, third parties can extend the model, such as

²[FAIR Principles](#).

Background

the European Union Digital European Exchange Platform (EU4S-DEEP) data templates for wastewater surveillance at airports, European Union [Digital European Exchange Platform \(EU4S-DEEP\)](#) data templates for wastewater surveillance at airports, the [CETO Epidemiologic platform](#), and [Ottawa Automatic Data Pipeline](#) to analyse qPCR data and transform it to ODM format.

Overall, the ODM aims to support the evolving field of environmental surveillance and science by providing a standardized, extensible approach to data collection and management.

Getting to know the ODM

What are “Parts” and “Part Types”?

Within the ODM, every component or element of the model and the dictionary is called a “part”. Parts are grouped into “part types”.

Because the term “parts” refers to every component of the ODM, “part types” are a way to differentiate between parts which have different functions and structures. The three most important part types are:

- **Measures:** A measurement or observation of any substance including a biological, physical or chemical substance.
- **Methods:** A procedure for collecting a sample or performing a measure.
- **Attributes:** A description of the who, where, when, and why of environmental surveillance.

See the [parts reference document](#) for more details.

There are additional part types to support the dictionary including:

- Aggregations
- Aggregation scales
- Categories
- Classes
- Compartments
- Dictionary support
- Domains

Getting to know the ODM

- Groups
- Missingness
- Nomenclature
- Quality indicators
- Specimens
- Tables
- Units

For information and details about all the part types, please see the [parts reference document](#).

“Sets” and the rationale behind sets

“Sets” are a way of grouping together different possible categorical inputs within the ODM. The purpose behind sets is to group options together so that users only need to check a few options from a drop-down menu, rather than scroll through a long list. Sets are also designed so that a single part can be used in multiple sets, which avoids the need to create multiple versions of this part for these different use contexts. For example, if you’re taking a measure of concentration, that measure will likely populate the unit choice drop-down list with units from the “standard concentration unit set”. Units in this set include Milligrams per litre, parts per million, and Percent Primary Sludge. Similarly, if one were measuring the amount of oxygen in a wastewater sample they would be offered units from the “dissolved gas unit set”. The units in this set include parts per million, and Milligrams per litre. You see here how Milligrams per litre is used twice across the sets, but exists still as a single part in the parts list. Unit sets are only one type of set, with the others being:

- Aggregation sets
- Compartment sets
- Quality sets
- Specimen sets

The “what” and “why” of Specimen IDs

Details about these sets are their similar - but unique - applications can be found in the [parts reference document](#). Importantly, “category sets” are not considered a “set” like these others. See the section below for an explanation of this difference.

Why category sets are treated differently

One of the reasons why “sets” are set up and organized the way that they are, is to allow for the re-use of parts in multiple sets. Category sets are, however, an exception to this rule because the categories that make up a category set are used in only one set. Furthermore, the general sets are needed for almost every measure or method (which require unit, aggregation, quality, and compartment information), while category sets are used only in specific circumstances for specific fields. These are things such as collection metadata (example: Sample collection category set), dictionary metadata and model architecture (example: Data Type Category Set), or more detailed types of certain methods (example: Nucleic Acid Extraction Category Set). Category sets are also a distinct part type, and more information about them can be found in the [parts reference document](#).

The “what” and “why” of Specimen IDs

Specimen IDs indicate the level at which a measure is being done. This can be: a site measure, ie. a measure at a site such as temperature or weather; a sample measure, ie. a laboratory measurement on a sample, such as gene copies of SARS-CoV-2 per mL; a person measure, ie. a measure of something at the level of an individual person, such as a blood pressure reading; or a population measure, ie. an aggregate measure of a population, such as the number of confirmed cases of a given illness. While there is currently nothing in version 2.0 of the ODM with uses the person specimen ID explicitly, it has been included to allow for the possibility of storing individual-level data in a future version. The reason

Getting to know the ODM

specimen IDs were created was so that there could be a single manner in which measures and methods are recorded, regardless of the level at which they were performed, while still maintaining the possibility to collect site, sample, and population measures as a distinct types of metadata.

The “what” and “why” of Groups and Classes

Groups and classes, similar to sets, are ways of grouping together different measures within the ODM. Given that the ODM aims to be as robust as possible and has a very long list of possible measures, groups and classes were designed to give shorter lists of measures in the drop down lists by specifying details about the kind of measure that a user is recording. Groups and classes can work together to further specify what kinds of measures are being reported. For example, the group `sarsCov2` contains many measures, but by specifying that the class is an `allele` or `variant`, the list of possible measures is pared down. Alternatively, class can be said to be non-applicable, paring down the measures in the `sarsCov2` group to only be unspecified measures of the quantity of the virus.

Tables and table types within the ODM

Within the ODM there are three types of tables: program description tables, results tables, and look-up tables. These table types exist to differentiate the function of these tables and to highlight these differences to users. The different tables take different types of inputs and maintenance from users, so understanding the differences can be important.

Program description tables

Program description tables (represented in yellow in the ERD) are tables used to record metadata on the organizations, locations, methods, and appurtenance. These tables help to describe surveillance and testing programs, and are intended to be updated infrequently.

Results tables

Results tables (represented in blue in the ERD) are the tables used to record details on samples and measures. These tables record the main outcomes data and are updated daily, if not more frequently.

Look-up tables

Look-up tables (represented in green in the ERD) are the tables that are pre-programmed and pre-populated in the ODM. These hold information on sets, all parts, languages, and translation abilities. These are only updated by the ODM team in version updates.

Measures, Methods, and Attributes: Key parts

While there are many part types, there are three main part types users should be most familiar with: measures, methods, and attributes. These three have parallels with the three table types, and the differences and details of these parts are useful to understand.

Getting to know the ODM

Measures

Measures are actually types of measures that can be performed. These can range from temperature to the number of gene copies in a sample. The measure, or `measureID`, can be selected from a drop down in the templates and it specifies the kind of measure you intend to record. The actual value of the measure is then recorded in the `value` field of the `measures` table, with units and aggregation specified in the `unitID` and `aggregationID` fields.

Methods

Similar to measures, methods are types of methods that can be performed to accomplish a measurement. These can be diverse, ranging from incubation, qPCR, or nucleic acid extraction. The method itself (`methodID`) is selected from a drop down in the templates, specifying in general terms the kind of method the user wishes to record. From there, the `value` field of the `methodSteps` table can be populated by one of the inputs from the category set associated with that `methodID`. This provides a higher level of detail for a given method. For example, the `methodID` might be `solidSep` for solid separation of a sample. The value field might then be populated with `cent` for centrifugation, indicating more details about how the sample settling was done.

Attributes

Attributes are the largest category of part type, as these refer to any most of the fields in the ODM. Attributes are fields for metadata within the ODM and range greatly in the use. They comprise everything from collection dates, to sampling period, to names.

Data quality and reportability

Within the ODM there is a `qualityFlag` field present in both the `samples` and `measures` tables. This field serves to highlight whether or not there is a quality issues with the sample or the measure. The quality flag also allows a user to specify to type of quality issue. This is managed through the use of quality sets, which are the sets that contain the possible quality flags for a given measure or for samples. The idea is that this provides data on any issues with a sample and an indication of the nature of that issue. Having data about the nature and presence of a quality issue is often not sufficient for decision makers who are trying to use and interpret the data. As such, there is also a `reportable` field with is a Boolean indicator of whether or not data can or should be reported or included in final reports and decisions.

Time periods for samples and measures

For measures, there is a field `aDateStart` and `aDateEnd` which specify the date and time that an analysis was begun and finished. This allows for the recording of greater detail around timelines for especially long, multi-day analyses. For shorter analyses, which will likely make up the bulk of reported measures, the same date can input in both fields. The idea is to have the reporting tables for measures be as robust as possible to allow for various kinds of timelines.

For samples, there is `col1DT`, `col1DTStart`, and `col1DTEnd`. The first field is for the collection date and time of a single grab sample, so start and end are not necessary and the other two fields can be left blank. For composite or pooled sample, the start and end date and time for collection is crucial information to know. When these two fields are populated, the `col1DT` field can be left blank.

Dates within the ODM

There are a number of date fields within the ODM which all serve different purposes. Collection datetime (`collDT`) is the date a sample was collected, used only for grab samples. This field is left blank if the collection datetime start and end fields (`collDTstart` and `collDTend`) have been populated instead. Inversely, collection datetime start and end fields should be left blank if the collection datetime field has been populated. Similarly, analysis datetime start and analysis datetime end (`aDateStart` and `aDateEnd`) are used to report the date and or timeline of an analysis for a measure.

Date fields that are more related to data processing and laboratory infrastructure are: the last edited date (`lastEdited`) which indicates the last time a table, or a measure or sample details, were last edited or updated; the sent date (`sentDate`), or date that a sample was sent to the lab from the field; the received date (`recDate`), or the date the sample was received in the lab; and the report date (`repDate`), or the date that the analysis results or measures were reported. These give a greater indication of how up to date data is, but also on the speed and efficiency of the pipeline between sampling and reported results. Dates should be reported in day/month/year format to accommodate the most popular global convention.

Translation and language capabilities

The default language of the ODM is English, but French translations of all descriptive elements of the data model dictionary are also available. As additional nation states and partners adopt the ODM, we anticipate that these fields will be translated into other languages as well. The translation capacities are managed through the language look-up table (`languageLUs`), the translation look-up table (`translationLUs`), and the

Questions and online community

parts table (**partsLUs**). When a translation for a given part is not available, the dictionary will default to the English term. The language look-up table stores linguistic and classification codes for spoken human languages, with the most recent ISO639 code being the language ID (**languageID**). In the translation table, the language ID is paired with every part ID (**partID**), along with the label (**partLabel**), description (**partDesc**), and instructions (**partInstr**) translated to that language. The part IDs are linked to the full parts list which otherwise contains metadata that is coded in variables and requires no further translation.

Questions and online community

If any users have additional questions or issues with the ODM, we invite them to check out our [Discourse Page](#) for discussion boards and community support. For larger issues, and to ask the ODM team to add additional fields or variables into the model, we encourage users to visit the project's [GitHub repository](#) and to create an issue there so that a team member can respond.

How to use the ODM

How-To guides

1) How to use templates

ODM Excel report templates provide you with an easy way to enter your wastewater and other environmental data into the PHES-ODM format. The templates also serve as an interactive way to get to know the ODM. Below you will find a series of guides that will walk you through how to enter data into the report templates.

In these guides, “field” refers to a column in the template, and “entry” refers to a row of information that you entered. Guides also provide definitions and explanations for some terms; however, more complete information is found in the most can be found in the reference section of the ODM documentation.

The report templates can be found at the OSF.io [PHES-ODM reference file folder](#). More information about the tables can be found in the [getting-to-know](#) section.

Briefly,

- **samples** and **measures** tables are used day-to-day reporting. Information about a sample will be recorded in the **samples** table. Measurements are recorded in the **measures** table. A measure is

How to use the ODM

‘measurement or observation of any substance including a biological, physical or chemical substance’.

- There are three optional tables to support reporting samples and measures.
 - **sampleRelationships** is optional and used only when you have samples that are pooled or split.
 - **measureSets** is also an optional table that can be used to link measures together.
 - **qualityReports** is used for detailed quality assurance and control measures (you can—and should— record a **quality flag** attribute for each measure or sample in the **measures** and **samples** tables).
- Typically, there is one-time initial setup where information is entered the tables with the yellow tabs in the Excel template. These tables include information about **sites**, **addresses**, **organizations**, and **contacts**.
 - **sites** is a mandatory table records information about where testing is performed.
 - **addresses**, **organizations**, and **contacts** provide the contact information for sites and also other entities such as laboratories, or database custodians.
- There are optional supporting tables that include **protocols** and **instruments** where you record information about how samples are taken and measures are reported. These tables are usually completed once and then updated when there are changes to the protocol. **protocolSteps** and **protocolRelationships** tables are part used alongside the **protocols** table. Other supporting tables include **polygons** where information is store about the sample shed

2) How to report Measures and Measure Sets:

and **datasets** which includes information about the data custodian and owner of the data.

The templates contain four main types of fields, each with their own for instructional purposes:

- **Identifier, unique identifier and composite unique identifier (red):** In these fields you will create a value that will be used to identify something. For instance, **Sample** ID is used to identify each unique sample. An identifier can be repeated between entries, while an unique identifier cannot. There is usually only one unique identifier per template, and they correspond to the name of the table. e.g. **Measure** ID for the **measures** table, and **Contact** ID for the **contacts** table.
- **Drop-down and conditional drop-down menus (green):** These fields have a drop-down menu which contain the values that you can either select **or** type in. In conditional drop-down menus, the values depend on what you entered in a previous field.
- **Date fields (purple):** Dates are entered in these fields and must be in the form of the [ISO 8601](#) format: `yyyy-mm-dd`. All date fields also support entries for hours, minutes, and seconds: `2022-01-01T06:11:54` and `2022-01-01T06:11:54+13:30`
- **Free form fields (blue):** In these fields, you can enter your data. Any restrictions will be specified in the relevant sections. Keep in mind, different measures may have different data types that can not be validated in Excel. Data types each measure can be checked in the reference list for the measure. For example, the data type for **temperature** is a real number.

2) How to report Measures and Measure Sets:

In this guide you will learn how to record data in the **Measures** and **Measure Sets** table templates. In the **Measures** table, each row repre-

sents a measure. Along the same lines, in the **Measure Sets** template, each entry will be a collection or set of measures.

Quick Start

Below you will find the mandatory fields for the **measures** Table. You will also find definitions of terms and examples of entries of data, with all of the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed description.

Measures Template

1. Mandatory fields

- **Report ID**: Unique identifier for a Measure.
- **Sample ID**: Unique identifier for a Sample.
- **Site ID**: Unique identifier for the location where a Sample was taken.
- **Analysis Date End**: Date the Measure was completed.
- **Measure**: A measurement or observation of any substance including a biological, physical or chemical substances.
- **Value**: Value of the Measure.
- **Unit**: Units of a value.
- **Aggregation**: Statistical measures that the Measure represents.
- **Specimen**: Substance or thing upon which the observation was made. Specimens include water, air, and surface. The Specimen feild is mandatory if there is more than one specimen type that is collected in database.

Detailed description

2. Examples

Image will be added

Measure Sets Template

1. Mandatory fields

- **Report Set ID:** Unique identifier that links together a group of related Measures.

3. Examples

Image will be added.

Detailed description

Measures template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Report ID	Protocol ID	Sample ID	Purpose	Polygon ID	Site ID	Dataset ID	Report Set ID
2								
3								

i) Columns A-C, and E-H

- These are the identifier fields.

How to use the ODM

- **Report ID** (column A) is the unique identifier for this template, and it cannot be repeated. Each Report ID represents a unique measure.

The **Report ID** can be any text field up to 30 characters. If you don't have a method to generate a **Report ID** you can use the ODM recommended formula: {Mathew, can you add the formula here}.

- Other identifiers can be repeated between entries (if needed). For instance, if you enter two different measures from the same sample, then the **Sample ID** (column C) will be the same.
- You may have already created these identifiers in another template. For instance, you may have created **Sample ID** in the **samples** Template.

ii) Column D

- **Purpose** (optional) has a drop-down menu list. In the list, **Regular** is the default purpose.

2. Columns I to M

	I	J	K	L	M
1	Analysis date start	Analysis date end	Report Date	Specimen	Fraction analyzed
2					
3					

i) Columns I to K

- These are the date fields where you can enter the dates associated with your measure in the format yyyy-mm-dd. Including hours and minutes is optional, but all date-entries use the ISO 8601 format. 2022-01-01T06:11:54 or 2022-01-01T06:11:54+13:30

ii) Columns L and M

- **Specimen** and **Fraction analyzed** have drop-down menus. Both fields are optional, but both are recommended.
- **Specimen** is must be included your table has more than one type of specimen. When there is only one specimen, the database table record the specimen type. For example, if your **measures** table include only wastewater sample measures, you can indicate this information in the **databases** table.
- **Fraction analyzed** applies only for water and wastewater samples. You should record the fraction for all samples if they are component type water or wastewater.

3. Columns N to S

	N	O	P	Q	R	S
1	Group	Class	Measure	Value	Unit	Aggregation
2						
3						

i) Columns N and O

- **Group** and **Class** help organize **Measures** to make a specific measure easy to find. The fields have drop-down menus. Both fields are optional and used only if you need help finding a **Measure**.
- For example, to find the allele measures for SARS-CoV-2 you can start by selecting **SARS-CoV-2** in the **Group** column. Next, you select **Allele** in the **Class** column. Then, you'll see the drop-down in the **Measure** column that contains SARS-CoV-2 allele regions.
- Leave these fields **blank** if you do not use them.

ii) Column P

- **Measure** is where you select what is being measured. A measure in ODM is a 'measurement or observation of any

How to use the ODM

N	O	P	V
Group	Class	Measure	
SARS-CoV-2	Alleles class	<div> <div></div> <div> BRSV-N SARS-CoV-2-E SARS-CoV-2-N1 SARS-CoV-2-N2 SARS-CoV-2-N3 SARS-CoV-2-RdRp crAssphage-N PMMoV-CP S sars-cov-2 gene target Influenza virus A1 Influenza virus A2 Influenza virus B </div> </div>	

substance including a biological, physical or chemical substance.’ **Measure** has a conditional drop-down menu. See above section on **Group** and **Class**.

- If you chose to enter in a **Class** (column O), the drop down menu will contain the **Measures** found in that **Class**. If you did not enter a value in column O, the menu will contain all possible **Measures**.
- iii) Column Q - This where you will enter in the value of your **Measure**. - For instance, if your **Measure** is a temperature of 20 °C, you would enter “20”. - Note: Excel accepts any data type can be entered in the **Measure** column. However, each **Measure** has a specific data type. Data types can be found in the reference documentation for each measure. For example, the data type for temperature can be found [here](#).
- iv) Column R
- **Units** is where you enter the unit for your measure. There is a conditional drop-down list of the units that are valid for the measure in column P. For example, if the measure is [Storate temperature] (.../parts.qmd#sTemp), the only unit displayed is Celcsius.

v) Column S

- **Aggregations** is where you enter the aggregation of your value. There is a conditional drop-down menu field which contains the types of aggregations of your value. For example does your value represent a mean, median, etc.
- What appears in this menu list depends on what you entered in the Unit field (column R).

4. Columns T to AA

	T	U	V	W	X	Y	Z	AA
1	Nomenclature	Index	Measure License	Reportable data flag	Organization ID	Contact ID	Update Date - Last Edited	Notes
2								
3								

i) Columns U, V and AA

- **Index** (column U) is used if you have multiple entries with the exact same values in the other fields. If this is the case you could use “1”, “2”, etc. to differentiate them.
- **Measure license** (column V) refers to the access and use licensing of the measure that you are entering.

ii) Column T and W

- **Reportable** (column W) is where you can indicate if there are quality concerns with the measure and the measure should not be used for regular reporting. More details of the quality concerns can be recorded in the **quality sets** table. There is a drop down list this field.

iii) Column X and Y

- These are identifier fields that are used to indicate the **Organization** and **Contact** associated with the entry.

iv) Column Z

- **Last edited** field is where you indicate the date when the entry was updated. The field is used if you need to modify an entry (row) after your initial recording. You can leave

How to use the ODM

this field blank if the measure was entered with no updates.
Meaning, use this field only when you modify the entry (row).

Measure Sets Template

1. Columns A to G

i NOTE: this figure needs to be update with the correct labels.

	A	B	C	D	E	F
1	Report Set ID	Protocol ID	Name	Contact ID	Organization ID	Update - Last E
2						
3						

i) Columns A, B, D and E

- These are all Identifier fields.
- Report Set ID (column A) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Report Set ID represents a unique Measure Set.
- For the rest, you can repeat the Identifier between entries (if needed). For instance, if you are entering two different Measure Sets that came from the same Organization, then the Organization ID (column D) will be the same.
- You may have already created these Identifiers in another template. For instance you may have created Organization ID in the Organizations Template.

ii) Columns C and G

3) How to report Samples and Sample Relationships:

- These are Free Form fields in which you will enter the indicated information. Name (column C) refers to the name that you have given to the Measure Set.

iii) Column F

- This is a Date field in which you can enter in the date when the entry was updated.

You have now entered your data in the Measures and Measure Sets Template, congratulations!

3) How to report Samples and Sample Relationships:

In this guide you will learn how to enter Samples and Sample Relationships into their respective templates. In the Samples template, each entry represents a Sample. A Sample is the wastewater that you collected so that Measures can be made. Along the same lines, each entry in the Sample Relationships represents an interaction (or Relationship) between two Samples in the form, Subject - Relationship - Object. So, to specify that Sample A is a field sample replicate of Sample B, you would enter “Sample A” - “Field Sample Replicate” - “Sample B”

You can find additional information regarding terms and the colour-coding of columns at the beginning of this guide (1) [How to use templates](#)) and

Quick Start

Below you will find the fields that are mandatory for the Samples and Sample Relationships template. You will also find some examples of entries of data, with the all of the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed Description.

Samples Template

1. Mandatory fields

Sample ID	Site ID
Sample Material	Sample Collection Type
Collection Period	Collection Number
Collection Date Time	

2. Definitions

- Sample ID: Unique identifier for a Sample.
- Site ID: Unique identifier for the location where a wastewater sample was taken.
- Sample Material: Type of Sample.
- Sample Collection Type: Method used to collect the Sample.
- Collection Period: The time period over which the Sample was collected, in hours.
- Collection Number: The number of subsamples that were combined to create the Sample. Use NA for continuous, proportional or passive sampling.
- Collection Date Time: The date, time and time zone the sample was taken.

3. Examples

to be added

Sample Relationships Template

1. Mandatory fields

Sample ID Object	Relationship
Sample ID Subject	

2. Definitions

- **Sample ID Object:** The Object (or one of the samples) of a relationship between two Samples. This will always be a Sample ID that was previously created in the Samples Template.
- **Relationship:** Describes relationship between two Samples.
- **Sample ID Subject:** The Subject (or one of the samples) of a relationship between two Samples. This will always be a Sample ID that was previously created in the Samples Template.

3. Examples

Examples will be added

Detailed Description

Samples Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Sample ID	Protocol ID	Organization ID	Contact ID	Site ID	Purpose	Sample material
2							
3							

How to use the ODM

i) Columns A to E:

- These are Identifier fields.
- Sample ID (column A) is the Unique Identifier for this template and cannot be repeated between entries. You can think of each Sample ID value as representing a unique Sample.
- For the rest of the fields, you can repeat values between entries. For instance, if you are entering two different Samples from the same Site, then the Site ID (column E) will be the same.
- You may have already created these Identifiers in another template. For instance you may have created Site ID in the Sites Template.

ii) Columns F to G:

- These are Drop-down Menu fields where you can enter information regarding the Purpose and Sample Material of your Sample.

2. Columns H to N

	H	I	J	K	L	M	N
1	Dataset ID	Sample origin	Replicate Type	Sample collection type	Collection Period	Collection Number	Pooled
2							
3							

i) Columns H:

- These is an Identifier field for the Dataset that your Sample is associated with.

ii) Columns I to K, N:

- These are Drop-down Menu fields related to various characteristics of your Sample.

Detailed Description

- Sample Collection Type (column K) refers to the collection technique you used to obtain the Sample.
 - Pooled (column N) refers to whether or not the Sample that you are entering is made up of multiple child Samples.
- iii) Columns L and M:
- These are Free Form fields.
 - Collection Period (column L) refers to the number of hours that it took to collect your Sample.
 - Collection Number (column M) refers to the number of subsamples that were used to create Sample that you are entering.

3. Columns O to V

	O	P	Q	R	S	T	U	V
1	Collection date time	Collection date time start	Collection date time end	Date sample was sent	Date sample recieved	Reportable data flag	Update Date - Last Edited	Notes
2								
3								

- i) Columns O to S, U:
- These are Date fields that are related to various dates associated with your Sample. These dates will be in the format dd/mm/yyyy.
- ii) Column T:
- This is a Drop-down Menu field (Reportable Data Flag) where you can indicate whether or not your entry is of a good enough quality that it can be used by others.
- iii) Column V:
- This is a Free Form field where you can indicate anything of interest..

Sample Relationships Template

1. Columns A to E

How to use the ODM

	A	B	C	D	E
1	Sample ID Subject	Relationship	Sample ID Object	Update Date - Last Edited	Notes
2					
3					

i) Columns A and C:

- These are Identifier fields.
- Both of these Identifiers are Sample IDs and represent Samples. You would have created them previously in the Samples Template. Note that neither of them are Unique Identifiers and, thus, can be repeated between entries.

ii) Column B:

- This is a Drop-down Menu field, where you can select the type of relationship between two Samples. For instance, if Sample B was a child of Sample A, you would put the Sample ID of Sample A in the Sample ID Subject field (column A), and the Sample ID of Sample B in Sample ID Object field (column C). Then you would chose “Child relationship” from the menu in this column.

iii) Columns D and E:

- These are a Date field (column D) and a Free Form field (column E) where you can enter in the indicated information.

You have now entered your data in the Samples Template, congratulations!

4) How to report Protocols, Protocol Steps and Protocol Relationships:

4) How to report Protocols, Protocol Steps and Protocol Relationships:

In this guide you will learn how to enter Protocols, Protocol Steps and Protocol Relationships into their respective templates. As described in , a Protocol is “A procedure for collecting a sample or performing a measure”. Each entry in the Protocols Template represents a unique Protocol. A Protocol is made up of Protocol Steps. In the Protocol Steps Template, each entry is one of these steps. Finally, Protocols and Protocols can be linked to each other. Each entry in the Protocol Relationships Template represents one of these relationships in the form, Subject - Relationship - Object. So, for example, to specify that Protocol Step A needs to be done before Protocol Step B, you would enter “Protocol Step A” - “is Before” - “Protocol Step B”. .

You can find additional information regarding terms and the colour-coding of columns at the beginning of this guide ([1\) How to use templates](#)) and

Quick Start

Below you will find the fields that are mandatory for the Protocols, Protocol Steps and Protocol Relationships templates. You will also find some examples of entries of data, with all of the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed Description.

Protocol Steps Template

1. Mandatory fields

Protocol Step ID	Method (if entry is a Method)
Measure (if entry is a Measure)	Value (if entry is a Measure)
Aggregation (if entry is a Measure)	Unit (if entry is a Measure)

2. Definitions

- Protocol Step ID: Unique identifier for a Protocol Step.
- Measure: A measurement or observation of any substance including a biological, physical or chemical substance.
- Method: A procedure for collecting a sample or performing a Measure.
- Value: Value of the entry.
- Aggregation: Statistical measures used to report a Measure.
- Unit: The units of a Measure.

3. Examples

Image will be added

Protocols Template

1. Mandatory fields

Protocol ID

2. Definitions

- Protocol ID: A unique identifier for a given Protocol.

3. Examples

Image will be added

Protocol Relationships Template

1. Mandatory fields

Protocol ID Container	Relationship
Protocol ID Subject (if subject is a Protocol)	Protocol Step ID Subject (if subject is a Protocol Step)
Protocol ID Object (if object is a Protocol)	Protocol Step ID Object (if object is a Protocol Step)

2. Definitions

- Protocol ID Container: Unique identifier for a Protocol and the Steps and other Protocols that make it up.
- Relationship: Describes relationship between the Subject and Object.
- Protocol ID Subject: The Subject of a Relationship between a Protocol and a Protocol Step or Protocol.
- Protocol Step ID Subject: The Subject of a Relationship between a Protocol Step and a Protocol Step or Protocol.
- Protocol ID Object: The Object of a Relationship between a Protocol and a Protocol Step or Protocol.
- Protocol ID Step Object: The Object of a Relationship between a Protocol Step and a Protocol Step or Protocol.

3. Examples

Image will be added

Detailed Description

Protocol Steps Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Protocol Step ID	Method	Measure	Summary	Protocol step source	Protocol Step Version	Reference link
2							
3							

i) Column A:

- Protocol Step ID (column A) is the Unique Identifier field for this template and cannot be repeated between entries. You can think of each Protocol Step ID value as representing a unique Protocol Step.

ii) Columns B to C:

- These are Drop-down Menu fields.
- You only need to enter a value into Method ID (column B) or Measure ID (column C). The former is for when you are entering a Method as a Protocol Step. The latter is when you are entering a Measure as a Protocol Step. Leave the other field **blank**.

iii) Columns D to G:

- These are Free Form fields and you can enter the relevant information in any format that you wish.
- Summary (column D) is a short description of the Protocol Step that you are entering.
- Protocol Step Source (column E) specifies the Protocol Step which serves as a basis for the given Protocol Step.
- Protocol Step Version (column F) specifies the version of a given Protocol Step.

Detailed Description

	H	I	J	K	L	M	N	O
1	Organization ID	Contact ID	Instrument ID	Value	Unit	Aggregation	Update Date - Last Edited	Notes
2								
3								

2. Columns H to O

i) Columns H to J:

- These are Identifier fields.
- You can repeat values between entries. For instance, if you are entering two different Protocol Steps from the same Organization, then the Organization field (column H) will be the same.

ii) Column K:

- A Free Form field in which you will enter the Value of a Measure. If a Method is being entered, then you can put “NA” here.

iii) Columns L and M:

- These are Conditional Drop-down Menu fields.
- Unit (column L) depends on what you entered in Measure ID (column C). If you entered a Method, then you can put “NA” in this field.
- Aggregation (column M) depends on what you entered in Unit (column L). If you entered a Method, then you can put “NA” in this field.

iv) Columns N and O:

- These are Free Form fields where you can enter the indicated information.

Protocols Template

1. Columns A to F

How to use the ODM

	A	B	C	D	E	F
1	Source Protocol	Protocol ID	Dataset ID	Name	Summary	Reference link
2						
3						

i) Columns B and C:

- These are Identifier fields.
- Protocol ID (column B) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Protocol ID represents a unique Protocol.
- For Dataset ID, you can repeat between entries (if needed). For instance, if you are entering two different Protocols from the same Dataset, then the Dataset ID (column C) will be the same.
- You may have already created the Dataset IDs in the Datasets Template.

ii) Columns A, D to F:

- These are Free Form fields in which you can enter the indicated information.
- Source Protocol (column A) is where you can indicate the Protocol that served as the basis for the Protocol that is currently being entered.
- Name (column D) refers to the Name that you have decided to give your Protocol.

2. Columns G to K

i) Columns G, H:

- These are Identifier fields. You can repeat values between entries.

ii) Columns I, K:

Detailed Description

G	H	I	J	K
Organization ID	Contact ID	Protocol Version	Update Date - Last Edited	Notes

- These are Free Form fields and you can enter the relevant information in any format that you wish.
- You can use the Protocol Version field (column I) to indicate the version of the Protocol that you are entering.

iii) Columns J:

- This is a Date field where you can indicate when the entry was updated.

Protocol Relationships Template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Protocol ID Container	Protocol ID Object	Step ID Object	Relationship	Protocol ID Subject	Step ID Subject	Update Date - Last Edited	Notes
2								
3								

i) Columns A to C, E and F:

- These are Identifier fields.
- All of these Identifiers are either Protocol IDs or Protocol Step IDs and represent Protocols and Protocol Steps, respectively. This includes Protocol ID Container (column A), which represents the Protocol that contains the relationships that you are entering in this template. Note that none of these fields are Unique Identifiers and, thus, can be repeated between entries.

How to use the ODM

- You should only enter in a value for one of the two “Object” columns (column B and C); this depends on whether the Object of the relationship that you are entering is a Protocol or a Protocol Step. The same is true for the two “Subject” columns (column E and F).
- An example will be given in point (ii) below.

ii) Columns D:

- This is a Drop-down Menu field, where you can select the type of Relationship that is being entered.
- For example, say you wanted to enter in a Relationship stating that Protocol Step A takes place before Protocol Step B in a particular Protocol Container. You would enter the relevant Identifiers in the Step ID Object field (column B) and the Step ID Subject field (column F). You would then select “is Before” from the Relationship field. Finally, to identify this relationship, you would enter in the Identifier in the Protocol ID Container field (column A)

iii) Columns G and H:

- These are a Date field (column D) and a Free Form field (column E) where you can enter in the indicated information.

You have now entered your data in the Protocol Steps, Protocols and Protocol Relationships Template, congratulations!

5) How to report Metadata:

In this guide you will learn how to enter Metadata through a series of templates. Each entry in a template will usually represent what is found in the title of the template. For instance, in the Organizations template

each row will represent an Organization. In the context of the PHES-ODM, metadata is anything that gives information that is not found in the templates related to Measures, Protocols and Samples.

You can find additional information regarding terms and the colour-coding of columns at the beginning of this guide (1) [How to use templates](#)) and

Quick Start

Below you will find the fields that are mandatory for the metadata-related templates. You will also find some examples of entries of data, with the fields filled and colour-coded. A more detailed description of some of the key concepts can be found in the following section, Detailed Description.

Addresses Template

1. Mandatory fields

Address ID	Dataset ID
Address Line 1	City
State, Province or Region	Country

- Address ID: A unique identifier for an Address.
- Dataset ID: The name of the Dataset that stores information for Measures, Samples and other reporting tables.
- Address Line 1: Line 1 (the street name, number and direction) for a given Address.
- City: The city where a Site or Organization is located; part of the Address.

How to use the ODM

- State, Province or Region: The state, province, or region where a Site or Organization is located; part of the Address.
- Country: The country where a Site or Organization is located; part of the Address.

2. Examples

Images will be added

Contacts Template

1. Mandatory fields

Contact ID	Organization ID
Contact email	

- Contact ID: A unique identifier for a given contact person.
- Organization ID: A unique identifier for the Organization to which the reporter is affiliated.
- Email: Contact e-mail address, for the labs.

2. Examples

Images will be added

Datasets Template

1. Mandatory fields

Dataset ID	License
Data Custodian	

- Dataset ID: A unique identifier for the Dataset that stores information.
- License: The license of the Dataset.
- Data Custodian: The data custodian(s) of the database.

2. Examples

Images will be added

Instruments Template

1. Mandatory fields

Instrument ID	Dataset ID
Model	Instrument Type

- Instrument ID: A unique identifier for an Instrument.
- Dataset ID: The name of the Dataset that stores information for Measures, Samples and other reporting tables.
- Model: Model number or version of the Instrument.
- Type of Instrument used to perform the measurement.

2. Examples

Images will be added

Organizations Template

1. Mandatory fields

Organization ID	Address ID
Organization Type	

- Organization ID: A unique identifier for the Organization to which the reporter is affiliated.
- Address ID: A unique identifier for an Address.
- Organization Type: Specifies the type or purpose of a given Organization.

2. Examples

Images will be added

Polygons Template

1. Mandatory fields

Polygon ID	Type of Geography
Well-known text	European Petroleum Survey Group Coordinates

- Polygon ID: Unique identifier for the Polygon.
- Type of Geography: Type of geography that is represented by the Polygon.
- Well-known text: Well-known text of the Polygon.
- European Petroleum Survey Group Coordinates: The unique EPSG code specifying a given geospatial area.

2. Examples

Images will be added

Quality Reports Template

1. Mandatory fields

Quality Report ID	Report ID (if entry is related to a Measure)
Sample ID (if entry is related to a Sample) Quality Flag	Report Set ID (if entry is related to a Measure in a Measure Set)

- Quality Report ID: A unique identifier for a given Quality Report.
- Report ID: Unique identifier for a Measure.
- Sample ID: Unique identifier for a Sample.
- Report Set ID: Unique identifier that links together a group of related Measures.
- Quality Flag: A field for reporting any quality concerns - of lack thereof - for a Sample or Measure.

2. Examples

Images will be added

Sites Template

1. Mandatory fields

Site ID	Site Type
Sample shed	Contact ID
Latitude	Longitude
European Petroleum Survey Group Coordinates	

- Site ID: Unique identifier for the location where wastewater sample was taken.

How to use the ODM

- Site Type: Type of Site where sample was taken.
- Sample shed: A geographic area, physical space, or structure. A Sample is taken from a sample shed for a representative measurement of a substance(s).
- Contact ID: A unique identifier for a given contact person.
- Latitude: Latitude in decimal coordinates of the Site.
- Longitude: Longitude in decimal coordinates of the Site.
- European Petroleum Survey Group Coordinates: The unique EPSG code specifying a given geospatial area.

2. Examples

Images will be added

Detailed Description

Addresses Template

1. Columns A to E

	A	B	C	D	E
1	Address ID	Dataset ID	Address Line 1	Address Line 2	City
2					
3					

i) Columns A and B:

- These are Identifier fields.
- Address ID (column A) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Address ID represents a unique Address.

Detailed Description

- For Dataset ID (column B), you can repeat the value between entries (if needed). For instance, if you are entering two different Addresses from the same Dataset, then this column will be the same.
- You may have already created the Dataset ID in the Datasets Template.

ii) Columns C to E:

- These are Free Form fields in which you can enter various information regarding the Address.

2. Columns F to J

F	G	H	I	J
Province, or Region	Postal or Zip Code	Country	Update Date - Last Edited	Notes

i) Columns F to J:

- Most of these fields are Free Form in which you can enter various information regarding the Contact.
- Update Date - Last Edited (column I) is a date field.

Contacts Template

1. Columns A to E

i) Columns A to C:

- These are Identifier fields.

How to use the ODM

	A	B	C	D	E
1	Contact ID	Dataset ID	Organization ID	First name of contact	Last Name of Contact
2					
3					

- Contact ID (column A) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Contact ID represents a unique Contact.
- For the rest, you can repeat the Identifier between entries (if needed). For instance, if you are entering two different Contacts that came from the same Organization, then the Organization ID (column C) will be the same.
- You may have already created these Identifiers in another template. For instance you may have created Organization ID in the Organizations Template.

ii) Columns D and E:

- These are Free Form fields in which you can enter the name of the Contact.

2. Columns F to J

	F	G	H	I	J
1	Contact email	Contact phone	Role of Contact	Update Date - Last Edited	Notes
2					
3					

i) Columns F to J:

- Most of these fields are Free Form in which you can enter various information regarding the Contact.

- Update Date - Last Edited (column J) is a date field.

Datasets Template

1. Columns A to G

	A	B	C	D	E	F	G
	Parent		Dataset				Reference
1	Dataset ID	Dataset ID	Creation Date	Name	License	Description	link
2							
3							

i) Columns A and B:

- These are Identifier fields.
- Dataset ID (column B) is the Unique Identifier for this template, and cannot be repeated between entries. In essence, each value of Dataset ID represents a unique Dataset.
- For Parent Dataset ID (column A), you can repeat the value between entries (if needed). For instance, if you are entering two different Datasets that came from the same Parent Dataset, then the Parent Dataset ID will be the same.

ii) Column C

- This is a Date field where you can enter the date that the Dataset was created.

iii) Columns D to G:

- These are Free Form fields in which you can enter various information related to the Dataset.

2. Columns H to N

How to use the ODM

	H	I	J	K	L	M	N
1	Language ID	Funder Contact ID	Custody Contact ID	Funding agency	Data custodian	Update Date - Last Edited	Notes
2							
3							

i) Columns H to L:

- These are Identifier and Free Form fields.
- For information regarding Funders, you can enter their ID and Name in Funder Contact ID and Funding Agency, respectively (columns H and J).
- For information regarding Data Custodians, you can enter their ID and name in Custody Contact ID and Data Custodian, respectively (columns I and K).

ii) Columns M and N:

- These are a Date field (column L) and a Free Form field (column M) where you can enter in the indicated information.

Instruments Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Instrument ID	Dataset ID	Name	Model	Manufacturer	Contact ID	Organization ID
2							
3							

i) Columns A, B, F and G:

- These are Identifier fields.

Detailed Description

- Instrument ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Instrument ID represents a unique Instrument.
- For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Instruments from the same Dataset, then this column will be the same.
- You may have already created some of the Identifiers in other Templates. For instance, you may have already created Dataset ID in the Datasets template.

ii) Columns C to E

- These are Free Form fields where you can enter the Name (column C), Model (column D) and Manufacturer (column E) of the Instrument.

2. Columns H to N

	H	I	J	K	L	M	N
1	Description	Reference link	Instrument Type	Describe other instrument type, if applicable	Index	Update Date - Last Edited	Notes
2							
3							

i) Columns H, I, L to N:

- Most of these fields are Free Form in which you can enter relevant information regarding your Instrument
- You can use Index (column L) to differentiate two entries of data that are the same in the other fields.
- Update Date - Last Edited (column M) is a date field.

ii) Columns J and K:

- For Instrument Type (column J), you can select from the drop-down menu the type of Instrument. If you do not see your Instrument, you can enter “Other instrument type”

How to use the ODM

- If you entered “Other instrument type” in the Instrument Type field, you can use the Describe Other Instrument type, If Applicable field (column K) to describe your Instrument.

Organizations Template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Organization ID	Dataset ID	Address ID	Organization Type	Name	Description	Update Date - Last Edited	Notes
2								
3								

i) Columns A and C:

- These are Identifier fields.
- Organization ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Organization ID represents a unique Organization.
- For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Organizations from the same Address, then this column will be the same.
- You may have already created some of the Identifiers in other Templates. For instance, you may have already created the Address ID in the Addresses Template.

ii) Columns D to H:

- For Organization Type (column D), you can select from the menu the type of Organization that you are entering.
- The rest of the columns contain Free Form or Date fields in which you can enter additional information regarding the Organization.

Polygons Template

1. Columns A to G

	A	B	C	D	E	F	G
1	Polygon ID	Dataset ID	Name	Description	Polygon Population	Type of geography	European Petroleum Survey Group Coordinates
2							
3							

i) Columns A and B:

- These are Identifier fields.
- Polygon ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. A Polygon is an area on a map. Each Polygon ID represents a unique Polygon.
- For Dataset ID (column B), you can repeat the value between entries (if needed). For instance, if you are entering two different Polygons from the same Dataset, then this column will have the same value.
- You may have already created the Dataset ID in the Datasets Template.

ii) Columns C to E, G:

- These are Free Form fields in which you can enter additional information regarding the Polygon.
- Name (column C) is the name of the Polygon.
- European Petroleum Survey Group Coordinates (column G) is a coordinate system that can be used to identify where your Polygon is.

iii) Columns F:

- For Type of geography (column F) of the Polygon, you can select one of the values from the drop-down menu.

2. Columns H to N

How to use the ODM

	H	I	J	K	L	M	N
1	Well-known text	File Location of Polygon	Reference link	Organization ID	Contact ID	Update Date - Last Edited	Notes
2							
3							

i) Columns H to J and N:

- These are Free Form or Date fields in which you can enter relevant information regarding the Polygon.
- Well-known text (column H) refers to the text markup language that can be used to represent the Polygon (

ii) Columns K and L:

- These are Identifier fields in which the Organization ID and the Contact ID associated with the Polygon are entered.

Quality Reports Template

1. Columns A to H

	A	B	C	D	E	F	G	H
1	Quality Report ID	Report ID	Sample ID	Report Set ID	Quality Flag	Severity Indicator	Update Date - Last Edited	Notes
2								
3								

i) Columns A to D:

- These are Identifier fields.
- Quality Report ID (column A) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Quality Report ID represents a unique quality issue that you want to report.
- For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Quality Reports from the same Sample, then the Sample IDs will be the same.

Detailed Description

- You may have already created some of the Identifiers in other Templates. For instance, you may have already created Sample ID in the Samples template.
- ii) Columns E and F:
- Quality Flag (column E) is a Drop-down Menu field where you can select the type of quality issue that you would like to enter.
 - You can indicate the severity of the Quality Flag in the Severity Indicator field (column F).
- iii) Columns G and H:
- These are a Free Form field and Date field in which you can enter the indicated information.

Sites Template

1. Columns A to F

	A	B	C	D	E	F
1	Parent Site ID	Site ID	Dataset ID	Polygon ID	Site Type	Sample shed
2						
3						

- i) Columns A to C, and F:
- These are Identifier fields.
 - Site ID (column B) is the Unique identifier for this template, and cannot be repeated between entries. In essence, each Sites ID represents a unique Site.
 - For the rest of the columns, you can repeat the value between entries (if needed). For instance, if you are entering two different Sites that are found in the same Polygon, then the Polygon IDs will be the same.

How to use the ODM

- Parent Site ID (column A) refers to the Site that is parent to the Site that is being entered. For example if the Site that is being entered is a room in a LTC facility, that the Parent Site ID would refer to the LTC facility.

ii) Columns D to E:

- These are Drop-down Menu fields in which you can specify the Site Type (column D) and the type of Sample shed (column E).

2. Columns G to M

	G	H	I	J	K	L	M
1	Address ID	Organization ID	Contact ID	Name	Description	Department or Region	Health Planning Region
2							
3							

i) Columns G to I:

- These are Identifier fields for additional information that are associated with the Site.

ii) Columns J to M:

- These are all Free Form fields in which you can enter in relevant information about your Site.
- Public Health Department or Region (column L) and Health Region (column M), refers to the Department or Region where your Site is located.

3. Columns N to S

	N	O	P	Q	R	S
1	Population Served	Latitude	Longitude	European Petroleum Survey Group Coordinates	Update Date - Last Edited	Notes
2						
3						

i) Columns N to S:

- These are Free Form fields and a Date field in which you can enter in additional information of your Site.

Example data for each of the 11 report tables.

- European Petroleum Survey Group Coordinates (column Q) is a coordinate system that can be used to identify where your Site is.

You have now entered in your Metadata, congratulations!

Example data for each of the 11 report tables.

Part I.

Reference guide

The reference guide is a technical description of the ODM **parts**. There are three sections to the reference guide.

- **parts** – Parts are the smallest way of describing anything in the ODM. An analogy to parts of a car can be used to illustrate this point. If you own a car, your garage can access a parts list that contains every part of your car, right down to every nut and bolt. There are over 900 parts, all listed in alphabetical order. Each part has up to 60 different attributes such as the part name and description.
- **tables** - The ODM database model has 21 tables in its full relational database. The full model is commonly referred to as “long” tables as it stores data with one measurement per row. The most important tables are the ‘report’ tables that are used for daily reporting of measures and samples.
- **sets** - Sets are lists or collections of parts. You can think of sets as the elements of drop down lists or categories.

💡 Did you know that each part ID has its own URL?

Every part, table, and set has its own URL, that is created using the partID. partIDs are usually parentheses that follow a part, table, or set name.

For example, ‘airport’ is a sample shed category which would appear as: **Airport** (**airport**). The permanent link to the airport partID is, <https://docs.phes-odm.org/parts.html#airport>. Airport also is listed in the **Sampleshed** category set (**shedCatSet**).

Parts

ODM version v2.0.0.

Parts are the smallest way of describing anything in the ODM. An analogy to parts of a car can be used to illustrate this point. If you own a car, your garage can access a parts list that contains every part of your car, right down to every nut and bolt. Similarly, the ODM has a **part ID** ([partID](#)) and **part description** ([partDescription](#)) for every measure, method, and attribute. There are also **part IDs** for **units**, **aggregations**, and other **parts**. The part list is presented below.

💡 Did you know you can identify part IDs by their specific naming format?

Each part ID starts with a lowercase character and contains no special characters. Parts have less than 12 characters, while part types such as units and aggregations have a maximum of seven characters. The unique construction of part IDs serves several purposes in addition to generating valid URLs. Users can create wide table variable names (headers) using part IDs that are valid in all open software languages. ODM supports a standard method of generating wide variable names by concatenating partIDs with an underscore, `_`. For this reason, some part types have a limit of seven characters to ensure the total length of wide names remains less than 35 characters, which is the character limit length of some statistical language software. The wide variable naming method will be added to the documentation.

Parts

a1306s delta-variant gene target

(a1306s). a1306s delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a1918v delta-variant gene target

(a1918v). a1918v delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a2710t omicron-variant gene target

(a2710t). a2710t omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a63t omicron-variant gene target

(a63t). a63t omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

a67v omicron-variant gene target

(a67v). a67v omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Atline Analyzer

(aas). Atline analyzer with sampler. An atline analyzer with sampler. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Absolute humidity

(absHum). A measure of the total mass of water vapour present in the air per volume of air. *Part Type*: measurements *Domain*: phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: humid *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: absHumidUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Absolute humidity unit set

(absHumidUnitSet). Unit set for absolute humidity measurements. *Part Type*: unitSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: absHumidUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Academic institution

(academ). The category of organization type used for academic institutions or research groups. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum*

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Number of active cases

(acti). Number of active cases. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Active

(active). Indicator that a part is in current use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Analysis date

(aDate). Date the measurement was performed in the lab. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date end

(aDateEnd). Date the measurement or analysis was completed. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date start

(aDateStart). Date the measurement or analysis was started. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Address Line 1

(addL1). Line 1 (the street name, number and direction) for a given address. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:*

naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Address Line 2

(addL2). Line 2 (the unit number) for a given address. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Address table

(addresses). The table that contains information about addresses. Addresses can be recorded for sites, organizations or contacts (individuals). The Sites, Organizations, and Contacts tables include a link to the addresses table through Address ID. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: programDescr *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Addresses table column order

(addressesOrder). Specifies the order of the columns in the Addresses table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Address table required headers

(addressesRequired). Specifies the columns required in a Addresses table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Address ID

(addressID). A unique identifier for an address. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Administrative regions

(admRegLevel). Administrative regions *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: ISCU *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Air filter

(afu). Air filter as part of filtration or circulating unit. Typically the unit would have a fan or blower. *Part Type*: categories *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: airCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scale

(aggragationScale). A scale used for an aggregation. Only applicable for measures and units. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum*

Parts

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation

(aggregation). Statistical measures used to report a measure. Each aggregation has a corresponding value. The default aggregation is ‘single’, which corresponds to a single measurement. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Aggregations

(aggregations). Statistical measures used to report a measure. Each measure reported as a number should be reported with an aggregation. The default aggregation is ‘single’, which corresponds to a single measurement. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scales

(aggregationScales). The scale of an aggregation set. Aggregation scales include quantitative and qualitative. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation set

(aggregationSet). The aggregation set for a unit. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation sets

(aggregationSets). Sets of aggregations. Examples of aggregation sets include logarithm, linear and boolean. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

Parts

0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Air compartment

(air). A measure or observation made from a substance in the air. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: airCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 1 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Air compartment set

(airCompartmentSet). A compartment set for measures and methods in the air compartment. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: airCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Airplane

(airpln). Airplane sample shed category type *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*:

varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Airport

(airport). Airport sample shed category type *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Airport lists worksheet

(airportLists). Parts lists used to generate airport and airplane data entry templates. *Part Type:* dictionarySupport *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Air and surface compartment set

(airSurfaceCompartmentSet). A compartment set for measures and methods in the air or surface compartments. *Part Type:* compartmentSets *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* airSurfaceCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*

Parts

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Environmental temperature

(airTemp). Environmental temperature. *Part Type*: measurements *Domain*: phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: temperature *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: temperatureUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Air and water compartment set

(airWaterCompartmentSet). A compartment set for measures and methods in the air or water compartments. *Part Type*: compartmentSets *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: airWaterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Alias deprecated

(aliasDep). ID of an assay that is the same or similar. A comma separated list. Deprecated as of version 2, please avoid use. Use notes instead. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Alias ID deprecated

(aliasIDDep). Alias id ID of an assay that is the same or similar. a comma separated list. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

All domains

(allDo). Domain that specifies that it could apply to al domains; biological, chemical, and physical. *Part Type:* domains *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA

Parts

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Alleles class

(allele). Measures and methods related to alleles. *Part Type:* classes *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Allele unit set

(alleleUnitSet). Unit set for alleles. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* siSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* alleleUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Access to all org

(allOrgs). If this is 'no', this data will not be available to any partner organization. If missing, data will be available to the all organizations. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: boolean *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 1 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Aloh3 precipitation

(aloh3). Aloh3 precipitation *Part Type*: categories *Domain*: che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Assay method id default depreciated

(amDefDep). Used as default when a new measurement is created for this lab. See ID in AssayMethod table. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Amicon ultrafiltration

(amiconUf). Amicon ultrafiltration *Part Type:* categories *Domain:* che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Amicon filter, extract with MP96

(amiMP96). Nucleic acid extraction, usually used for the liquid fraction of a wastewater sample, using amicon filtration and using an MP96 for final extraction. *Part Type:* categories *Domain:* bio *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Amplicon sequencing

(amp). Specifies the amplicon strategy for genetic sequencing *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 Maximum Length: 10 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

AND Boolean aggregation

(andBoo). “AND” aggregation. If all values in the aggregation is “TRUE” then the AND aggregation is also “TRUE”. If all values in the aggregation is “FALSE” then the AND aggregation is also “TRUE”. *Part Type: aggregations Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: NA Minimum Value: NA Maximum Value: NA Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Any compartment set

(anyCompartmentSet). A compartment set for measures and methods in any compartment. *Part Type: compartmentSets Domain: allDo Speciment ID: anySpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 15 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Any specimen set

(anySpecimenSet). A specimen set that inculdes any specimen. *Part Type: specimenSets Domain: naDomain Speciment ID: anySpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass*

Parts

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Area proportional sample

(areaPr). An area proportional sample. Used for surface testing. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Arrival temperature

(arTemp). The temperature of a sample upon arrival at the laboratory for analysis. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* temperature *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Use of the articV3 primer

(articV3). Initial implementation of an ARTIC bioinformatics platform for nanopore sequencing of nCoV2019 novel coronavirus; Artic Foundation v3 primer. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Use of the articV4 primer

(articV4). Artic V4 sequencing primer for VOCs and VOIs. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Assay ID

(asID). Links with the AssayMethod used to perform the analysis. Use instrument.ID for measures that are not viral measures. Deprecated as of version 2, please avoid use. Refer to methodID instead. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Attributes

(attributes). Attributes describe the who, where, when, and why of environmental surveillance. There are attributes for domain, specimen, compartment, or table. Attributes are one of the three main components or ways to report environment surveillance data. The other two components are measures and methods. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Bacteria Class

(bacteria). Measures and methods relating to bacteria. *Part Type:* classes *Domain:* bio *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* bacteria *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* bacteriaUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Bacteria unit set

(bacteriaUnitSet). Unit set for bacteria-related measurements. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* bacteriaUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous bacteria group

(bactMisc). A group of measures/methods related to miscellaneous bacteria. The miscallenous bacteria are often have only one measure or method. When a bacteria has many measures, they will be given their own groupID in updated dictionary versions. *Part Type:* groups *Domain:* bio *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* bact-Misc *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Bead beating, extract with MP96

(bbMP96). Nucleaic acid extraction, usually used for solid fraction of a wastewater sample, using bead beating and using an MP96 for final extraction. *Part Type:* categories *Domain:* bio *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

Parts

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

bcov

(bcov). Measure of the amount of bovine coronavirus. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

bcov culture spike target

(bcovCul). Cultured bovine coronavirus is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

bcov vaccine spike target

(bcovVac). The bovine coronavirus vaccine is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp

Class: naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Beef extract flocculation

(beeExtractFloc). Beef extract flocculation. *Part Type:* categories *Domain:* che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Is Before

(before). Specifies that the object step or protocol occurs before the subject step or protocol in the overall protocol container. Use only for protocolRelationships table relationshipID. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Beta

(beta). B.1.351 *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Biologic

(bio). A living organism or biological substance. *Part Type*: domains *Domain*: bio *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Digital droplet emulsification PCR

(bioRadDdpcr). Describes a PCR analysis done using BioRad's digital droplet emulsification PCR technology. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Binary Large Object (BLOB) data type

(blob). The data type for Binary Large Object (BLOB) data. Only used for the dictionary entries of parts. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Breadth of coverage ($\geq 5\times$ depth)

(boc). Positions with read depth greater or equal to 5. Report as percentage of positions. *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

5-Day carbonaceous kiochemical oxygen demand

(bod5c). The quantity of oxygen utilized for the biochemical degradation of organic matter under standard laboratory procedures in five (5) days in the presence of a nitrification inhibitor, expressed in milligrams

Parts

per litre (mg/l). *Part Type*: measurements *Domain*: bio *Speciment ID*: siSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: <https://doi.org/10.1002/wer.1541> *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

5-day total biochemical oxygen demand

(bod5t). 5 day total biochemical oxygen demand. *Part Type*: measurements *Domain*: bio *Speciment ID*: siSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: <https://doi.org/10.1002/wer.1541> *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Boolean data type

(boolean). The data type for boolean/binary data. Encoded as 'TRUE' or 'FALSE' *Part Type*: dataTypes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Boolean aggregation set

(booleanAggrSet). Aggreagation set for boolean. *Part Type*: aggregationSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Boolean value set

(booleanSet). Set that contains the valid possible values for a boolean measure (TRUE or FALSE). *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: booleanSet *Unit Set*: NA *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Bovine respiratory syncytial virus group

(brsv). Bovine respiratory syncytial virus. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

brsv culture spike target

(brsvCul). Cultured bovine respiratory syncytial virus (BRSV) is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

BRSV-N

(brsvN). bovine respiratory syncytial virus capsid protein gene region *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

brsv vaccine spike target

(brsvVac). The bovine respiratory syncytial virus (BRSV) vaccine is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Branch sewer pipeline

(bSwrPpl). Specifies a site type that is collection pipe that run lateral to other municipal sewer lines, allowing drainage into the main sewer. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Building cleanout

(buildCO). Specifies a site type that is a capped pipe that connects to a building's main sewer line. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Capacity class

(capacity). Measures and methods relating to capacity. *Part Type:* classes *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:*

Parts

anyCompartmentSet *Group:* naGroup *Class:* capacity *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* capacityUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Capacity unit set

(capacityUnitSet). Unit set for capacity measurements. *Part Type:* unit-Sets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* capacityUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Case report date

(caRepDate). Date that the numbers were reported publicly. Typically, reported data and this measure is most commonly reported and used. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Categorical data type

(categorical). The data type for categorical data. Only used for the dictionary entries of parts. *Part Type*: dataTypes *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Categories

(categories). A discrete list of values that can be reported for a measure, method or attribute. A use-case example of categories would be how they are used to record the site where a measure is taken (*partID* = *geoType*). For the *geoType* attribute there are different type of sites where samples can be taken: a wastewater treatment plant, airplane holding tank, wastewater pumping station, etc. Each of these site types is recorded as a category that can be identified in *categorySetID* = *sTypeCatSets* *Part Type*: part-Type *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Long-term acute care hospital

(ccc). Acute care hospitals, or complex continuing care, that provide care for patients with average length of stay longer than 25 days. *Part Type*: cat-

Parts

egories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Child day care

(cdc). Child day care facility. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Degrees Celcius

(cel). Degrees Celsius. *Part Type:* units *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* -60 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Centrifugation

(cent). Describes solid separation from a wastewater sample via centrifugation. Likely connected to other method steps prior to analysis. *Part Type:* categories *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Centricon ultrafiltration

(centriconUf). Centricon ultrafiltration. *Part Type:* categories *Domain:* che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ceres nanotrap

(ceres). Ceres nanotrap. *Part Type:* categories *Domain:* che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

CFU per 100 ml

(cfu). Colony forming units per 100 ml of filtered sample. *Part Type:* units *Domain:* bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Changes column

(changes). A column for recording the changes from a previous version. Use this column as a short-hand change log for dictionary development. *Part Type:* tableSupport *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Chemical

(che). A chemical compound. *Part Type:* domains *Domain:* che *Specimen ID:* anySpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Chemagic viral dna/rna 300 kit

(chemVir). Nucleic acid extraction performed using the chemagic viral dna/rna 300 kit. *Part Type:* categories *Domain:* bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Child relationship

(child). Indicated that this is a sample generated from (an)other sample(s) either because of pooling or sub-sampling. *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

City

(city). The city where a site or organization is located; part of the address. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum*

Parts

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 1.0.0. Last updated: 2.0.0.

Clarified sample

(clari). Clarified sample. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Class

(class). A unique identifier for a class, which is akin to a sub-group; it's a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class. Currently, class is only used for biologics. For example, SARS-CoV-2 is a group of measures with the following classes of allele, variant, mutation, sequence, and protein. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Classes

(classes). A class is a collection of one or more related measures or methods within a group. A group can have one more more classes to describe different parts or the class. *Part Type*: partType *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Cloudy

(cloudy). Qualitative category for the weather measure, specifying an overcast day with no rain. *Part Type*: categories *Domain*: phy *Specimen ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Centimetres

(cm). Unit part for the SI unit of centimetres. *Part Type*: units *Domain*: phy *Specimen ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA

Parts

Maximum Length: NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Carbon Dioxide

(co2). A measure of an amount or concentraion of carbon dioxide. *Part Type:* measurements *Domain:* che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* gas *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* dissGasUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Chemical Oxygen Demand

(cod). Chemical oxygen demand. *Part Type:* measurements *Domain:* che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Collection group

(colGrp). A group of measurement-like attributes related to sample collection. *Part Type:* groups *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set*

ID: NA Unit Set: naUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection date time

(collDT). For grab samples this is the date, time and timezone the sample was taken. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: datetime Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Collection date time end

(collDTEnd). For integrated time average samples this is the date, time and timezone the sample was finished being taken. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: datetime Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Collection date time start

(collDTStart). For integrated time averaged samples this is the date, time and timezone the sample was started being taken. *Part Type: at-*

Parts

tributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: datetime *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sample collection set

(collectSet). Methods for collection samples. Sample collection methods include water, air, and surface. See the attribute ‘Compartment set’ of the sample collection method. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: collectSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Collection number

(collNum). The number of subsamples that were combined to create the sample. Use NA for continuous, proportional or passive sampling. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Collection number and period

(collNumPer). Composite collection number.period. This a composite measure that combines ‘collection period’ (collPer) and ‘collection number’ (collNum). Collection number and period can be used to reduce the number of table headings during data collection. In data storage, ‘collNum’ and ‘collPeriod’ should be transformed to ‘collPer’ and ‘collNum’. For Grab sample, Surface swab, or Area proportional sample just enter 1

In case of Composite, Flow-proportional collection, the number of subsamples followed by a dot, followed by the period of the sampling, in hours. For example, for 4 composite subsamples covering a 8 hours period, the entry would be 4.2 For time-proportional 24 subsamples collected every hour, the entry would be 24.1 For volume-proportional, the entry for 24 subsamples collected over 24 hours should be 24.24, in this case the period part of the entry should be understood as the total collection time instead of the periodicity of the collection.

For a COSCA ball or Moore swab collecting for 24 hours, the entry would be 1.24 *Part Type:* attributes *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* 1:1 *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 6 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Collection period

(collPer). Collection period. The time period over which the sample was collected, in hours. Alternatively, use collectionStart and collectionEnd. Examples: 1 hour, 6 hours, 24 hours *Part Type:* attributes *Domain:*

Parts

naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample collection type

(collType). The type of collection. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: collectSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 12 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Co-located sample

(colocated). Second or multiples samples collected at same location but different time (water, air) or at a nearby location (soil, sediment). Sent blind to laboratory. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Column name depreciated

(columnNameDep). Name for the column Depreciated in ODM version 2. Look-up tables are now incorporated into parts and sets. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Composite sample - archival

(comp). A composite sample, usually generated by an autosampler. This is intended for use in updating and mapping archival data into the newest version of the ODM. Generally, autosamplers do time or flow-proportional sampling, and so those collection methods should be used instead. This is only for composite/autosampler samples where these additional details have not yet been specified. Use collectionPeriod to describe how many hours the sample was taken, and collectionNum to record the number of samples. *Part Type:* categories *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Composite grab sample of 3

(comp3). A grab-composite sample composed of 3 separate grab samples. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Composite grab sample of 3 depreciate

(comp3dep). A grab-composite sample composed of 3 separate grab samples. Deprecated in version 2, please do not use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Composite 3hr grab sample

(comp3h). A 3-hour composite with 3 grab samples each taken once per hour, generally manually performed. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:*

NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Composite 3hr grab sample depreciate

(comp3hdep). A 3-hour composite with 3 grab samples each taken once per hour, generally manually performed. Deprecated in version 2, please do not use. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Composite 8hr grab sample depreciated

(comp8h). An 8-hour composite with 8 grab samples each taken once per hour, generally manually performed. Deprecated in version 2, please do not use. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Composite 8hr grab sample deprecated

(comp8hDep). A 8-hour composite with 8 grab samples each taken once per hour, generally manually performed. *Part Type*: categories *Domain*:

Parts

naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Compartment

(compartment). The attribute identifying the substance from which where a sample was taken. For more information, see partID = compartment. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Compartments

(compartments). The substance from which a sample was taken. For example, wastewater has component of “water”. Compartments are attributes of measures, methods, units, and aggregations. *Part Type*: part-Type *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment sets

(compartmentSets). Sets of compartments. Compartment sets are used to identify when a measure can be recorded for more than one compartment. For example, SARS-CoV-2 can be measured in people (humans), water, surface, or air. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample concentrate

(conc). Sample concentrate. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Case report date

(conCase). Date that the numbers were reported publicly for a confirmed case. Typically, reported data and this measure is most commonly reported and used. *Part Type*: units *Domain*: bio *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale

Parts

Quality Set: NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Is concurrent to

(concurrent). Specifies that the object step or protocol occurs at the same time as the subject step or protocol in the overall protocol container. Use only for protocolRelationships table relationshipID. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Water conductivity

(cond). Measurement of conductivity of sample or site. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* siteFeat *Class:* conductivity *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* conductivityUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Conductivity class

(conductivity). Measures and methods related to conductivity. *Part Type:* classes *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* conductivity *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* conductivityUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Conductivity unit set

(conductivityUnitSet). Unit set related to conductivity measurements. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* conductivityUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Number of confirmed cases

(conf). Number of confirmed cases. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Confirmed cases episode date

(confEp). Episode date is the earliest of onset, test or reported date for a confirmed case. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Onset date

(confOn). Earliest that symptoms were reported for a confirmed case. This data is often not known and reported. In lieu, episode is used. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Contact ID

(contactID). A unique identifier for a given contact person. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Contact name

(contactName). Contact person or group, for the lab. Depreciated in ODM version 2. see 'name' *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 100 Part Status: deprecated.* First released: 1.0.0. Last updated: 2.0.0.

Contact table

(contacts). The table that contains information about a contact; a person who is the contact of a site or laboratory. adapt from version 1 documentation *Part Type: tables Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: programDescr Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Contacts table column order

(contactsOrder). Specifies the order of the columns in a Contacts table. *Part Type: tableSupport Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class:*

Parts

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Contact table required headers

(contactsRequired). Specifies the columns required in a Contacts table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Contact table set

(contactTableSet). Tables are where measures, methods and attributes are recorded. Tables represent the main entities of the environmental surveillance. *Part Type:* dictSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Confirmed cases test date

(conTest). Date that the covid-19 test was performed for a confirmed case. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* dateTime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Cphd ID

(corFcil). Correctional facility *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

COSCa ball

(cosca). COSCa passive sampling device. Use collectionPeriod to describe how many hours the ball was used to collect the sample *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* <https://pubs.rsc.org/en/content/articlelanding/2021/ew/d1ew00207d> *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Parts

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Country

(country). The country where a site or organization is located; part of the address. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Countries or sovereign states

(countryLevel). Countries or sovereign states *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: ISCU Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Districts, counties, regions

(countyLevel). Districts, counties, regions *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: ISCU Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar*

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Covid-19

(cov). Covid-19 infection. Any form of Covid-19 human infection - including testing for Covid-19, symptomatic Covid-19, asymptomatic Covid-19, hospitalized Covid-19, long-Covid-19, etc. If needed, define the specific form Covid-19 with units. See populationUnits. *Part Type:* measurements *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* humanCompartmentSet *Group:* sarsCov2 *Class:* disease *Nomenclature:* ICD *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* populationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

SARS-CoV-2 measure

(cov2Me). Measure the amount of SARS-CoV-2 virus. This measure should only be used if there is no other SARS-CoV-2 measure. See groupID = sarsCov2. Use a note to describe the specific measure used. *Part Type:* measurements *Domain:* bio *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

SARS-CoV-2-B.1.1.7

(covB117). Variant B.1.1.7 of the SARS-CoV-2 virus. *Part Type*: measurements *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: pangolin *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-B.1.351

(covB135). Variant B.1.351 of the SARS-CoV-2 virus. *Part Type*: measurements *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: pangolin *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-E

(covE). SARS-CoV-2 E gene. *Part Type*: measurements *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N1

(covN1). SARS-CoV-2 nucleocapsid gene, allele 1. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: any-CompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N2

(covN2). SARS-CoV-2 nucleocapsid gene, allele 2. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: any-CompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N3

(covN3). SARS-CoV-2 nucleocapsid gene, allele 3. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: any-CompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Parts

SARS-CoV-2-P.1

(covP1). Variant P.1 of the SARS-CoV-2 virus. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: pangolin *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-RdRp

(covRdrp). SARS-CoV-2 RdRp gene. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Covid-19 population measurement date

(cphDate). date of reporting for covid-19 measure. Depreciated in ODM version 2. The measure identification is now measureReportID with report dates as for any measure report. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0

Maximum Length: 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Cphd ID

(cphid). Unique identifier for the table. Depreciated in ODM version 2. The measure identification is now measureReportID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

crAssphage-N

(cra). crAssphage virus capsid protein gene region *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Crosswalk table set

(crosswalkTableSet). Tables used to translate to and from other environmental dictionaries and models. *Part Type:* dictSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet

Parts

Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Cycle threshold or quantification cycle (Ct or Cq)

(ct). Cycle thresholds in a PCR assay. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Cumulative count

(cuCo). Units for describing a population measure of a cumulative count of cases of a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Custody Contact ID

(custodyCont). A unique identifier for data custodians. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Data custodian ID

(custodyID). The data custodian of a database. Use Organization ID to populate this field, and the organizations table to describe contact information and other details for the data custodian. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d377y delta-variant gene target

(d377y). d377y delta-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference*

Parts

Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d3g omicron-variant gene target { #d3g }

(d3g). d3g omicron-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d63g delta-variant gene target

(d63g). d63g delta-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d796y omicron-variant gene target { #d796y }

(d796y). d796y omicron-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit*

Set: geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA
Reference Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal
Maximum Value: seeUnitVal *Minimum Length:* NA *Maximum Length:*
 NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

d950n delta-variant gene target

(d950n). d950n delta-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Daily count

(daiCo). Units for describing a population measure of a daily count of new cases of a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dataset creation date

(datasetDate). Specifies the date a given dataset was created. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dataset ID

(datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables. Where possible, datasetID name should correspond the original data custodian responsible for generating environmental data. (suggestion for a default name convention.) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dataset table

(datasets). A report table for capturing details about data's parental data set and data custodians. Supplying attribution for data collectors. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Datasets table column order

(datasetsOrder). Specifies the order of the columns in a Datasets table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset table required headers

(datasetsRequired). Specifies the columns required in a Datasets table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Data types

(dataTypes). The data type for a part. Data types used in the ODM include: varchar, boolean, float, category, date, time, datetime, url, email. dataType corresponds the the entry or cell within a data table, most

Parts

commonly within a report table. If the data entry has a unit, then the dataType corresponds to the unit and the dataType refers to ‘unit’. If the data entry is a category, then the dataType refers to ‘category’. All categories are varchar. Otherwise the dataType is identified within the part entry. TBA: dataType for dictionary. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Date

(date). Date date on which the assayMethod was created or updated (for version update). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Datetime data type

(datetime). The data type for date and time data. Only used for the dictionary entries of parts. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum*

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Days

(days). A unit for indicating a length of time in days. *Part Type:* units *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* time *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* timeUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ddcov_e sars-cov-2 gene target

(ddcovE). ddcov_e sars-cov-2 gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ddcov_n sars-cov-2 gene target { #ddcovN }

(ddcovN). ddcov_n sars-cov-2 gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference*

Parts

Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Deaths

(death). Units for describing a population measure of patients who have died from a given cause. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

del143/145

(del143). 143 or 145 deletion omicron-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

del 157/158

(del157). 157 or 158 deletion delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:*

anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

del2084/2084

(del2084). 2084 or 2084 deletion omicron-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

del212/212

(del212). 212 or 212 deletion omicron-variant gene target. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

del3674/3676

(del3674). 3674 or 3676 deletion omicron-variant gene target. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

del69/70

(del6970). 69 or 70 deletion omicron-variant gene target. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Delta

(delta). B.1.617.2 *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Deprecated

(deprecated). Indicator to say that a part is no longer in current use in the model. See `partID = status` *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Derived sample

(derived). Specifies a sample that is derrived or made from another sample or material. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Description

(descr). A detailed description as an attribute for describing results or program elements. A description of the part that serves a clear presentation of the part to a wide audience including technical and not technical staff. A description should allow any person who generates environment surveillance data to know how to identify how to record their data elements in the ODM. `partReference` can be used to further describe a part. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment*

Parts

Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 200 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Description of change

(descrChange). A description of change in a part from the previous version. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Desk or counter

(desk). Desk, table, countertop or other flat working surface. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dewatered solids

(deso). Dewatered solids. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Detected

(det). Substance detected or not detected. TRUE = detected, FALSE = not detected *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* 0 *Maximum Value:* 1 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Access to details deprecated

(detailsDep). More details on the existing confidentiality requirements of this measurement. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum*

Parts

Length: 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Development

(development). Indicator that a part is under development use. See *partID* = *status* *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dictionary reference table set

(dictionaryRefTableSet). Reference or look-up tables. For example, the parts table describe all elements of the ODM, including tables, table headers, measures, methods, categories, and units. sets are collections of parts. For example, units can be grouped together in a unitSet. languages and translations support translations. *Part Type:* dictSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Additional dictionary sheets

(dictionarySupport). Additional sheets for the Excel version of the ODM dictionary. Sheets or tabs for the ODM Excel dictionary include: all ERD tables (look-up tables and data entry tables) and additional supportig tables. *Part Type:* partType *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dictionary set type

(dictSet). Sets used to describe and group dictionary tables. *Part Type:* dictSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dilution factor

(dilFact). Specifies the extent to which a sample or aliquot was diluted prior to analysis. Dilution factor is reported as a unitless measure, where a value of 10 indicates a 10:1 dilution. *Part Type:* measurements *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* dilution *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* unitlessUnitSet

Parts

Aggregation Scale: seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Point dillutions

(dillute). Exact concentration or dillutions for generating a standard curve. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* dilution *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dilution Class

(dilution). Measures and methods relating to dilutions. *Part Type:* classes *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* dilution *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* unitlessUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Diseases (human) class

(disease). Measure and methods related to disease or infection in humans. *Part Type:* classes *Domain:* bio *Speciment ID:* poSpecimenSet

Compartment Set: humanCompartmentSet *Group:* naGroup *Class:* disease *Nomenclature:* ICD *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* populationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dissolved gas concentration unit set

(dissGasUnitSet). Unit set for carbon dioxide concentrations measurements in water or air. *Part Type:* unitSets *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* airWaterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* dissGasUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Depth of coverage

(doc). The sequencing read depth. Used to interpret the confidence in a presence or amount of a variant. *Part Type:* units *Domain:* bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Domain

(domain). Domain is the highest level of describing of a measure. The domain ID that corresponds to a given part. Mostly applicable for measures, methods, units, and aggregations. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Domains

(domains). There are three domain types: biologic (i.e. Covid-19, chemical (i.e. nitrogen), physical measure (temperature). *Part Type*: part-Type *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Higher education dormitory or residential building

(dorm). Higher education dormitory or residential building *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

e156g delta-variant gene target

(e156g). e156g delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

e484a omicron-variant gene target { #e484a }

(e484a). e484a omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Escherichia coli

(ecoli). Concentration of bacteria that are passed through the faecal excrement of humans, livestock and wildlife *Part Type:* measurements *Domain:*

Parts

bio Speciment ID: siSaSpecimenSet *Compartment Set:* surfaceWaterCompartmentSet *Group:* bactMisc *Class:* bacteria *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* bacteriaUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Education

(education). A measure or sample taken for education or training. Use this purpose, for example, when teaching how to use the PHES-ODM. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Efficiency

(efficient). The efficiency reported for a standard curve. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Contact email

(email). Contact e-mail address, for the lab. *Part Type:* attributes
Domain: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:*
 naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*
 naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*
 naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA
Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*
Length: 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0.
 Last updated: 2.0.0.

European Nucleotide Association

(ena). ENA header *Part Type:* crosswalkTableSet *Domain:* naDomain
Speciment ID: naSpecimenSet *Compartment Set:* naCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology*
Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar
Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maxi-*
imum Length: 30 *Part Status:* active. First released: 2.0.0. Last updated:
 2.0.0.

European Nucleotide Association (ENA) - notes

(enaNotes). ENA notes *Part Type:* crosswalkTableSet *Domain:* naDomain
Speciment ID: naSpecimenSet *Compartment Set:* naCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology*
Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar
Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maxi-*
imum Length: 30 *Part Status:* active. First released: 2.0.0. Last updated:
 2.0.0.

Parts

Rainfall

(envRnF). Rainfall, i.e. amount of precipitation in the form of rain. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* precipitation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* precipitationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Ground snow depth

(envSnwD). Total depth of snow on the ground. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* precipitation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* precipitationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Snowfall

(envSnwF). Snowfall, i.e. amount of precipitation in the form of snow. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* precipitation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* precipitationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA

Maximum Length: NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Episode date

(epiDate). Episode date is the earliest of onset, test or reported date. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

ERD table set

(erdTableSet). All tables listed in the Entity Relationship Diagram. The full ODM model is commonly referred to as “long” tables as it stores data with one measurement per row. *Part Type:* dictSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sarbecovirus-specific E sars-cov-2 gene target

(eSarbec). Sarbecovirus-specific E sars-cov-2 gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:*

Parts

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Estimated frequency of reads

(estFreqReads). Estimated frequency of reads in a sequencing assay *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: sequence *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Estuary, natural water body

(estuary). Estuary, natural water body *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Experiment Failed

(expFail). PCR experiment failed. No value reported. Value should be blank for a measure with this qualityFlag value. *Part Type:* qualityIndicators *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

4s method

(extract4s). Nucleic acid extraction performed using the 4s method. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* <https://www.protocols.io/view/v-4-direct-wastewater-rna-capture-and-purification-bpdfmi3n> *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nucleic acid extraction method

(extraction). Description of the nucleic acid extraction method. Description of the method used to extract the sample *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* extractSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data*

Parts

Type: categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Nucleic Acid Extraction set

(extractSet). set used for storing all the valid category values for the nucleic acid extraction method. *Part Type:* mmaSets *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* extractSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Extraction volume of sample

(exvol). Extraction volume of sample. Size of the sample that is analyzed. *Part Type:* measurements *Domain:* phy *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Fecal matter

(faeces). Fecal matter. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology*

Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

FALSE

(false). Boolean data type = FALSE Use these values and their set for any boolean measure or attribute. Use only “FALSE” (case sensitive) *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 5 *Maximum Length:* 5 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Field sample

(field). Specifies a sample taken from the field; directly collected from an area for testing. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Field sample replicate

(fieldReplicate). A sample divided into two or more homogeneous parts. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

File location of polygon

(fileLocation). The location of the file containing the geometry of the polygon. File path specified, or a URL *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* blob *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 65535 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Filtration

(filt). Describes solid separation from a wastewater sample via filtration. Proceeds further concentration or analysis of the liquid filtrate. *Part Type:* categories *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:*

NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

First Nation

(fiNa). Used to categorize a sample that is a First Nation, or on reserve lands. Likely for internal use only, Indigenous data can and should not be shared without explicit consent of the nation in question. *Part Type*: categories *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

First name of contact

(firstName). Specifies the first name of a given contact. *Part Type*: attributes *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

First released version

(firstReleased). The version in which a part was first released *Part Type*: partSupport *Domain*: naDomain *Specimen ID*: naSpecimenSet *Compartment*:

Parts

ment Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Foreign key

(fK). Foreign key for a table. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed

(flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution. The reported concentration estimate is the updated result after addressing inhibition. *Part Type:* qualityIndicators *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination

(flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. *Part Type:* qualityIndicators *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

FI - Inhibition present and unaddressed

(flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. *Part Type:* qualityIndicators *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

J - Weak signal extrapolation

(flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of the standard curve in the non-linear region. *Part Type:* qualityIndicators *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:*

Parts

naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

ND - Non-detect

(flagND). No amplification occurred in the reaction; non-detect. For the value of a non-detected measure, report the actual value even if it is below the limit of detection. The flag here ensures that it's recorded as a non-detect regardless. In instances where the original data doesn't record a value, but only has the flag, please populate the value field with either a 0 if dealing with variant percentages, and a 1 for all other purposes to further indicate a null result. *Part Type: qualityIndicators Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

UJ - Trace signal extrapolation

(flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., "trace" signal observed); the reported value is based on the extrapolation of the standard curve in the non-linear region. *Part Type: qualityIndicators Domain: bio Speciment ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set*

ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable

(flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. *Part Type: qualityIndicators Domain: bio Specimen ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Float data type

(float). The data type for float data. Only used for the dictionary entries of parts. *Part Type: dataTypes Domain: naDomain Specimen ID: naSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 15 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Flow-normalized mean

(floMean). Mean measure normalized to wastewater flow. Mostly used for reporting specific alleles (N1, N2, E, etc.) *Part Type: aggregations*

Parts

Domain: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Floor

(floor). Floor of a building or room. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Flow rate

(floRate). Wastewater volumetric flow rate at the sample collection location over the 24-hr period during which the sample was collected. If only an instantaneous flow measurement is available, it may be reported in units of million gallons per day. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* siteFeat *Class:* flow *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* flowUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Flow class

(flow). Measures and methods related to flow. *Part Type:* classes *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* flow *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* flowUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Flow proportional 24hr sample depreciated

(flow24hDep). A flow proportional 24-hour composite sample generally collected by an autosampler. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Flow proportional sample

(flowPr). A flow proportional composite sample generally collected by an autosampler. Use collectionPeriod to describe how many hours the sample was taken. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA

Parts

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Volume flow rate unit set

(flowUnitSet). Unit set for volume flow measurements. *Part Type:* unit-Sets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* flowUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Flow volume

(flowVol). Volume of influent. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* siteFeat *Class:* flow *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* flowUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus measure

(flu). General influenza virus measure This measure should only be used if there is no other influenza virus measure. See groupID = virusMisc. Use a note to describe the specific measure used. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature

Ontology Reference: http://purl.obolibrary.org/obo/MONDO_0005812
Category Set ID: NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar
Minimum Value: seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus A1

(fluA1). Influenza virus A1 type *Part Type:* measurements *Domain:* bio
Speciment ID: saSpecimenSet *Compartment Set:* anyCompartmentSet
Group: virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData
Minimum Value: seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

influenza virus A2

(fluA2). influenza virus A2 type *Part Type:* measurements *Domain:* bio
Speciment ID: saSpecimenSet *Compartment Set:* anyCompartmentSet
Group: virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData
Minimum Value: seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Influenza virus B

(fluB). Influenza virus B type *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Fluidigm digital PCR

(fluidpcr). Describes a PCR analysis done using FluidIGM's digital PCR technology. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Foggy

(foggy). Qualitative category for the weather measure, specifying a foggy or hazy day. *Part Type*: categories *Domain*: phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Fraction analyzed

(fraction). Fraction of the sample that is analyzed. *Part Type:* attributes
Domain: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:*
 naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*
 naNomenclature *Ontology Reference:* NA *Category Set ID:* fractionSet
Unit Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference*
Link: NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:*
 NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First
 released: 1.0.0. Last updated: 2.0.0.

Sample fraction set

(fractionSet). set for the fraction of the sample (solid, liquid, etc.).
Part Type: mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet
Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set*
ID: fractionSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:*
 NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA
Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part*
Status: active. First released: 2.0.0. Last updated: 2.0.0.

Freyja Script

(freyja). Freyja Script for sequencing. *Part Type:* categories *Domain:*
 bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology*
Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* <https://github.com/andersen-lab/Freyja> *Data Type:* varchar *Minimum Value:* NA *Maximum Value:*
 NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First
 released: 2.0.0. Last updated: 2.0.0.

Parts

F-Specific RNA bacteriophages

(frna). A measure for amount of F-Specific RNA bacteriophages. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

F-Specific RNA bacteriophages, G2

(frnaG2). A measure for amount of G2 F-Specific RNA bacteriophages. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample frozen

(frozen). Sample was frozen before analysis or processing *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*:

NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample temperature

(fst). Temperature that the sample is stored at while it is being sampled. This field is mainly relevant for composite samples which are either kept at ambient temperature or refrigerated while being sampled. *Part Type*: measurements *Domain*: phy *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: temperature *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: temperatureUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Full dictionary sheets set

(fullDictionarySheetSet). Worksheets in the full Excel dictionary. The full dictionary is 'ODM_full-dictionary.xlsx' *Part Type*: dictSets *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Funder Contact ID

(funderCont). A unique identifier for funders. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*:

Parts

naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Funding agency ID

(funderID). The funding agency of the dataset. Use Organization ID to populate this field, and the organizations table to describe contact information and other details for the funding agency. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

g215c delta-variant gene target

(g215c). g215c delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

g339d omicron-variant gene target { #g339d }

(g339d). g339d omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

g496s omicron-variant gene target { #g496s }

(g496s). g496s omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

g662s delta-variant gene target

(g662s). g662s delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Gamma

(gam). P.1 *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Gas class

(gas). Measures and methods relating to gas and gases. *Part Type*: classes *Domain*: che *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: gas *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: dissGasUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per copy of crAssphage

(gcCrA). Gene or variant copies per copy of crAssphage. *Part Type*: units *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

gene copies per day per 100,000

(gcD100). The unit for measures reflecting the gene copies per day per 100,000 people in the population. *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per gram solids

(gCGS). Gene or variant copies per gram solids. *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per L

(gcL). Gene or variant copies per litre. *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer

Parts

Minimum Value: 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per mL

(gcMl). Gene or variant copies per millilitre of solution. *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per PMMoV copy

(gcPpmov). Gene or variant copies per copy of PMMoV. *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Genetics unit set

(geneticUnitSet). Unit set for genetic-related measurements. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

geneticUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

General missingness set

(genMissingnessSet). The general set for missingness values. *Part Type*: missingnessSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Generic quality flag set

(genQualitySet). A quality set to specify any generic quality concerns about a measure or sample. *Part Type*: qualityIndSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

European Petroleum Survey Group Coordinates

(geoEPSG). The unique EPSG code specifying a given geospatial area. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet

Parts

Compartment Set: naCompartmentSet *Group:* miscAttr *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 2.0.0. Last updated: 2.0.0.

Latitude

(geoLat). Geographical location, latitude in decimal coordinates, ie.: (45.424721) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* miscAttr *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* -90 *Maximum Value:* 90 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Longitude

(geoLong). Geographical location, longitude in decimal coordinates, ie.: (-75.695000) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* miscAttr *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* -180 *Maximum Value:* 180 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Type of geography

(geoType). Type of geography that is represented by the polygon. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* geoTypeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Geographic set

(geoTypeSet). set for different type of geographic components. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* geoTypeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Well-known text

(geoWKT). Well-known text of the polygon *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 63 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Gram per cubic metre

(gm3). Density unit. Used for absolute humidity and other measures. *Part Type*: units *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: airCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Geometric mean

(gmn). Geometric mean. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Government agency

(govt). The category of organization type used for government agencies, programs, or crown-owned bodies. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0

Maximum Length: 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Grams

(grams). A unit of mass or weight. *Part Type:* units *Domain:* phy *Specimen ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Gravity settling

(graSet). Describes solid separation from a wastewater sample where the sample material is allowed to settle by gravity, and then separated. *Part Type:* categories *Domain:* allDo *Specimen ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Grab sample

(grb). A single large representative grab sample. If the sample was collected over a series of hours or is a composite sample, please use partID = comp *Part Type:* categories *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:*

Parts

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Gromstole 1.0 Script

(gromstole). Gromstole 1.0 Script for sequencing. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* <https://github.com/PoonLab/gromstole> *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Group

(group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations. A collection of related measures. For example, SARS-CoV-2 is a group. Within the SARS-CoV-2 group, there are measure classes tht include RNA alleles (N1, N2, E, etc.), mutations, (E484K0), an entire sequence, viral proteins, etc. Currently groups are used for measures only. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Groups

(groups). A collection of related measures. Used primary to group measurements and methods, this helps pare down the drop down list for a given measurement or *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Header

(header). Header for a table. Also known as a table variable or entiyp relationship ‘attribute’. Header is the top row or the variable name in a table. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Health administration or planning agency

(healthAdm). Health adminstrative or planning organization. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*:

Parts

NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Heat inactivated sars-cov-2 virus spike target

(heatInacSARS). Heat inactivated SARS-CoV-2 virus is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

hep g armored rna

(hepGRna). Measure of the amount Hepatitis G Armored RNA. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

hep g armored rna spike target

(hepGRnaMat). Hepatitis G armored RNA is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality*

Set: NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

High severity

(high). Indicates a very sever quality issue, likely meaning the data should not be reported. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sewer Network Health Region

(hlthReg). Health region served by the sewer network *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Hand Measurement

(hma). Handheld measurement analyzer. A handheld measurement analyzer. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:*

Parts

naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Hollow fiber dead end ultrafiltration

(hollowFiberUF). Hollow fiber dead end ultrafiltration *Part Type*: categories *Domain*: che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

General Hospital Admissions

(hosa). Hospital admissions or patients newly admitted to hospital. *Part Type*: units *Domain*: bio *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Hospital Census

(hosc). Hospital census or the number of people admitted with an ailment. *Part Type*: units *Domain*: bio *Speciment ID*: poSpecimenSet *Compartment*

Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Hospital

(hosptl). Hospital *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Holding tank wastewater depreciated

(hoTaWa). Wastewater from a holding tank, such as from an airplane or ship Depreciated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Hours

(hours). A unit for indicating a length of time in hours. *Part Type*: units *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: anyCompartmentSet *Group*: miscMeas *Class*: time *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Holding tank wastewater

(htSam). Wastewater sampled from a holding tank, such as from an airplane or ship *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Holding tank

(htSite). Holding tank *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Human compartment

(hum). A measure or observation made about a human. *Part Type:* compartments *Domain:* bio *Speciment ID:* anySpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Human compartment set

(humanCompartmentSet). A compartment set for measures and methods in the human compartment. *Part Type:* compartmentSets *Domain:* bio *Speciment ID:* anySpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Humidity class

(humid). Measures and methods related to humidity. *Part Type:* classes *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* airCompartmentSet *Group:* naGroup *Class:* humid *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Parts

Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

i1566v omicron-variant gene target { #i1566v }

(i1566v). i1566v omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

i3758v omicron-variant gene target { #i3758v }

(i3758v). i3758v omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

i82t delta-variant gene target

(i82t). i82t delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data*

Type: seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

International Classification of Diseases

(icd). Classification system for diseases in humans. *Part Type:* nomenclatures *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* disease *Nomenclature:* ICD *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Intensive care unit patients

(icu). Units for describing a population measure of patients who are in intensive care due to a given cause. Unit for a population or disease measure to explain population-level clinical effects. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Index

(index). Index number in case the measurement was taken multiple times. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet

Parts

Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 50 *Part Status:* active.
First released: 1.1.0. Last updated: 2.0.0.

Inhibition set

(inhibitionSet). Category set for inhibition methods. *Part Type:* mmaSets *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* inhibitionSet
Unit Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Inhibition measure

(inhibMe). Parameter to report whether or not inhibition was detected in the sample. Detected = TRUE, not detected = FALSE. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* booleanUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA
Reference Link: NA *Data Type:* categorical *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Inhibition method

(inhibMeth). Description of the method used to evaluate molecular inhibition. Description of the inhibition parameters. *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* inhibitionSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Innovaprep ultrafiltration

(innovaprepUF). Innovaprep ultrafiltration *Part Type:* categories *Domain:* che *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Input

(input). Input for a table. Indicates if a part can be used in a table. “Input” is used to indicate if a part can be used as an entry in a table (something with partType = table). For example, covN1 is a specific measure that can be entered in the measureID field of the measures table. Therefore, covN1 has a value “input” in the measures column in the parts list. ODM users can generate their own custom template by modifying using the value “input” for only the parts that are applicable for their program.

Parts

Part Type: categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Instrument ID

(instrumentID). A unique identifier for an instrument. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Instrument table

(instruments). The table that contains information about instruments. adapt from version 1 documentation *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Instruments table column order

(instrumentsOrder). Specifies the order of the columns in the Instruments table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Instrutment table required headers

(instrumentsRequired). Specifies the columns required in the Instruments table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other instrument

(instrumentTypeOther). Type of instrument other than those included in the PHES-ODM. An other type of measurement instrument. Add description to notes. See documentation for how to request new instruments added to the dictionary. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum*

Parts

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Instrument Type

(insType). Type of instrument used to perform the measurement. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* insTypeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Describe other instrument type, if applicable

(insTypeOth). Description of the instrument in case it is not listed in instrumentType. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Instrument set

(insTypeSet). List of instruments that are used for measures and methods *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* insTypeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:*

NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Integer data type

(integer). The data type for integers. Only used for the dictionary entries of parts. *Part Type*: dataTypes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Intercept

(inter). Intercept value of the calibration curve. Used for storing calibration curve information, potentially also applicable for other curve types. *Part Type*: measurements *Domain*: allDo *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: standard-Curve *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Combined ip2 and ip4 sars-cov-2 gene target

(ip2ip4). Combined ip2 and ip4 sars-cov-2 gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*:

Parts

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-1

(iso6391). The first part of the ISO 639 series of international standards for language codes. Part 1 covers the registration of two-letter codes. There are 183 two-letter codes registered as of June 2021. The registered codes cover the world's major languages. Part of the metadata for a languageID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2B

(iso6392B). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'B' specifies the bibliographic code (B code). Part of the metadata for a languageID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA

Data Type: categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2T

(iso6392T). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The ‘T’ specifies the terminological code (T code). Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-3

(iso6393). A set of international standards that lists short codes for language names. ISO 639-3 extends the ISO 639-2 alpha-3 codes with an aim to cover all known natural languages. Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-6

(iso6396). A set of international standards that lists short codes for language names. ISO 639-6 builds off ISO639-3 with the use of four-letter codes, and allowing users to differentiate between variants of languages and language families, such as historical vs. revived versions of languages. Part of the metadata for a languageID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

k856r omicron-variant gene target { #k856r }

(k856r). k856r omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Kilogram per second

(kgS). Kilograms per second. *Part Type*: units *Domain*: phy *Speciment ID*: siSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA

Category Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Kilolitres

(kl). Kilolitres of volume *Part Type:* units *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Litres

(l). Litres of volume *Part Type:* units *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

l452r delta-variant gene target

(l452r). l452r delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet

Parts

Group: sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

1981f omicron-variant gene target { #1981f }

(1981f). 1981f omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory

(lab). Laboratory for environmental testing. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Lab ID default deprecated

(labDefDep). Used as default when a new sample is created by this reporter. See ID in Lab table. Deprecated as of version 2, no defaults please specify. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Laboratory duplicate

(labDuplicate). Second (time or more) processing and analysis of sample. Usually for general chemistry or metals analyses. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Lab ID

(labID). Unique identifier for a laboratory. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0

Parts

Maximum Length: 10 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Lagoon system

(lagoon). Lagoon system for extensive wastewater treatment *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Lake, natural water body

(lake). Lake, natural water body *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Lambda

(lambda). C.37 *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* who *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum*

Value: seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA
Maximum Length: NA *Part Status:* active. First released: 2.0.0. Last
 updated: 2.0.0.

Language ID

(lang). Language code for translation purposes. Specifies the language for each translation, other than the default English. Follow the ISO639-3 codes for now. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Language family

(langFam). Specifies the language family of a given language for translation and language tracking purposes. Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Language name

(langName). Specifies the name of the language in roman alphabet characters for translation and language tracking purposes. Part of the metadata

Parts

for a languageID. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Language Look-up table

(languages). Look up table for all languages, used to give structure to the translation table. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: lookup *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Language table column order

(languagesOrder). Specifies the order of the columns in the Languages table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Language table required headers

(languagesRequired). Required headers in the Languages table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Last edited

(lastEdited). The date the entry was last updated. Use lastEdited if an entry is updated. Leave lastEdited blank or 'NA' for the first entry. Updates can include additions or revisions of the entry for any reason. For example, a wastewater measure was repeated later with an improved method. You can revise the entry by updating the value and then adding the date the new measure was performed to 'lastEdited'. To ensure data provenance, the best practice in this example is to generate a new entry with the same measureID as the original measureID. The original and new measures are kept in the database with a date in the 'lastEdited' field for the new entry but not the original one. Some databases may choose a delete-and-replace approach where the original entry is deleted and replaced by a new entry. In this approach, the 'lastEdited' field indicates the delete and replace occurred. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Last name of contact

(lastName). Specifies the last name of a given contact. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Last updated version

(lastUpdated). The version in which the part was last updated. Any change to the part of list will result in a change in the lastUpdated field to the dictionary version where the update occurred. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

LaTeX expression

(latExp). LaTeX expression used to generate formulas, symbols, etc. Mainly relevant for units. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 20 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage

(IBC). The percentage of the genome covered by reads (the breadth) is low. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory control sample duplicate

(lcsd). Known amounts of an analyte or representative compounds are added to a second “clean” matrix (lab water or clean sand) in laboratory. Duplicate of LCS *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage

(IDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

ment Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Leaked sample

(leaked). Sample leaked, some volume and material was lost. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Link

(li). Link to an external reference that describes the geometry of the polygon. Deprecated as of version 2, please avoid use. Use `referenceLink` instead. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

License

(license). The license of a dataset. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Linear scale aggregation set

(linearAggrSet). The aggregation set that contains all aggregations that exist on the linear scale. *Part Type:* aggregationSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Liquid fraction

(liq). Liquid fraction of a sample. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

Parts

0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

List worksheet

(lists). Parts lists for template dropdowns and other documentation. Sheets or tabs for the ODM Excel dictionary include: all ERD tables (look-up tables and data entry tables) and additional supportig tables. *Part Type:* dictionarySupport *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

List set

(listSet). PartTypes that contain categories, sets, or lists. Used to generate documentation. *Part Type:* miscSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* airCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* windSpeedUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Access to local ha

(localHA). If this is ‘no’, the, data will not be available to local health authorities. If missing, data will be available to local health authorities.

Part Type: attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Limit of detection (LOD)

(lod). Limit of detection. Limit of detection (LOD) for this method, if one exists. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Limit of detection (LOD) - sewage

(lodSewa). Limit of detection for sewage samples. Limit of detection (LOD) for sewage *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Logarithmic scale aggregation set

(logAggrSet). The aggregation set that contains all aggregations that exist on the logarithmic scale (rather than linear or qualitative). *Part Type:* aggregationSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Dictionary tables

(lookup). Tables to describe and support the ODM dictionary. These tables hold information on sets, all parts, and support multiple languages and translations. The tables are updated by the ODM development team when a new ODM version is created. *Part Type:* classes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* lookup *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Limit of quantification (LOQ)

(loq). Loq Limit of quantification (LOQ) for this method if one exists. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Low severity

(low). A marker for low severity *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval

(lowerCI95). Specifies the the lower limit of a 95% confidence interval. Should typically be linked to an upper limit, and some indication of data spread (ie. standard deviation), and a mean/median measure. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Low-volume sample

(lowVol). Sample is low-volume, but was run regardless. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Long-term care facility

(ltcf). A residential healthcare facility that provides 24-medical care. These are also called skilled nursing facilities. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Long-term care - assisted living or retirement home

(ltcfAl). A residential facility that provides assistance with daily care but generally does not provide skilled nursing care. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other long-term care

(ltcfO). Other residential facilities that provide daily and/or medical care, but are not defined as nursing home/skilled nursing facilities or assisted living facilities. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Life technologies digital PCR

(ltDpcr). Describes a PCR analysis done using Life Technologies' digital PCR technology. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Lysis buffer

(lysi). Lysis buffer. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale:*

Parts

NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Cubic metres per day

(m3D). Cubic metres per day. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per hour

(m3H). Cubic metres per hour. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per second

(m3S). Cubic metre per second. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum*

Length: NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Manufacturer

(manufacturer). Manufacturer of an instrument. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Mate pair layout

(mateLay). Specifies the mate pair layout for a sequencing method. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Maximum length

(maxLength). The maximum length of the value of a part or measure. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

Parts

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum value

(maxVal). Highest value in a range of values for a measure. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Maximum value part support

(maxValue). The maximum value of a part. *Part Type*: partSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Arithmetic mean

(me). Arithmetic mean. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology*

Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Measurement group

(measGrp). A group of measures that cannot be otherwise categorized. *Part Type:* groups *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other Measure

(measOth). Other measure, not otherwise specified in measures. Add description to categoryOther. *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Measures quality set

(measQualitySet). A quality set for any measure, includes all the quality flag options for any measure. Used as an “any quality set” for measures, but excludes the sample quality flags to avoid confusion and data entry errors. *Part Type*: qualityIndSets *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Measure

(measure). A measurement or observation of any substance including a biological, physical or chemical substance. An measurement is recorded from a specimen (i.e. a site, sample, person or population). Measures are organized into groups and classes (sub-groups). *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 12 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Measures

(measurements). The attribute to describe a part type of measures. All measures have ‘partType’ = “measure”. *Part Type*: partType *Domain*:

naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Report ID

(measureRepID). Unique identifier for a measurement. Report IDs cannot be repeated. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Measure report table

(measures). The table that contains information and details about a given measure adapt from version 1 documentation *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: results *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Measure set report table

(measureSets). The table that identifies sets of measures. Examples of measure sets include a set of replicates, dilutions (used to generate a Ct curve) or variants that are identified in a single sample. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* results *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measure sets table column order

(measureSetsOrder). Specifies the order of the columns in a Measure Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measure Set table required headers

(measureSetsRequired). Specifies the columns required in a Measure Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measures table column order

(measuresOrder). Specifies the order of the columns in a Measures table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Measures table required headers

(measuresRequired). Specifies the columns required in a Measures table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Report set ID

(measureSetRepID). Unique identifier that links together a group of related measures. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality*

Parts

Set: NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Median

(med). Median. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Normalized arithmetic mean

(menr). Arithmetic mean, normalized. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Other Measure Description

(meOthDe). Description of other measure (measOtherDesc). *Part Type:* measurements *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set*

ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 60 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Measure license

(measureLic). Specifies the access and use licensing for a given single measurement. Populated by partID = licenseID *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Method

(method). A procedure for collecting a sample or performing a measure. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Concentration method

(methodConc). Description of the method to concentrate a wastewater sample. Description of the method used to concentrate the sam-

Parts

ple *Part Type*: methods *Domain*: allDo *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: methodConcSet *Unit Set*: unitlessUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 100 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration method set

(methodConcSet). A set a concentration methods. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Methods

(methods). Procedures or steps for collecting samples or performing measures. For example, ‘methodExtract’ is a method that describes the how a sample was extracted. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2

(mfAcidmgcl2). Membrane filtration with acidification and mgcl2 *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2, membrane recombined with separated solids

(mfAcidmgcl2SS). Membrane filtration with acidification and mgcl2, membrane recombined with separated solids *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with addition of mgcl2

(mfMgcl2). Membrane filtration with addition of mgcl2 *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum*

Parts

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0.
Last updated: 2.0.0.

Membrane filtration with mgcl2 and separated solids

(mfMgcl2SS). Membrane filtration with addition of mgcl2, membrane recombined with separated solids *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment

(mfNoAmend). Membrane filtration with no amendment *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment, membrane recombined with separated solids

(mfNoAmendSS). Membrane filtration with no amendment, membrane recombined with separated solids *Part Type: categories Domain: naDomain*

Speciment ID: naSpecimenSet *Compartment Set:* naCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology*
Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation*
Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar
Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum*
Length: 30 *Part Status:* active. First released: 2.0.0. Last updated:
 2.0.0.

Membrane filtration with sample acidification

(mfSampleAcid). Membrane filtration with sample acidification *Part Type:*
 categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment*
Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomen-*
clature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA
Unit Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference*
Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:*
 NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First
 released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with sample acidification, membrane recombined with separated solids

(mfSampleAcidSS). Membrane filtration with sample acidification, mem-
 brane recombined with separated solids *Part Type:* categories *Domain:*
 naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet
Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature
Ontology Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Ag-*
gregation Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:*
 varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0
Maximum Length: 30 *Part Status:* active. First released: 2.0.0. Last
 updated: 2.0.0.

Parts

Millions of gallons per day (MG/D)

(mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to measure flow. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre

(mgL). Milligrams per litre. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Murine Hepatitis Virus

(mhv). A measure for amount of murine hepatitis viurs. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Multiple issues

(mI). Multiple issues have arisen in the sequencing process. When using the ‘multiple issues’ quality flag, please specify the issues and details in the notes section. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Mid-level severity

(mid). A marker for med-level severity, where there are some concerns. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Minimum length

(minLength). The maximum value of measure. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Parts

Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Minutes

(minutes). A unit for indicating a length of time in minutes. *Part Type:* units *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* time *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Minimum value

(minVal). Lowest value in a range of values for a measure. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Minimum value part support

(minValue). The minimum value of part. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:*

varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous attribute group

(miscAttr). A group of miscellaneous measurement-like attributes. Examples of these would be longitude and latitude, or EPSG coordinates. *Part Type:* groups *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous measure group

(miscMeas). A group of measures that cannot be otherwise categorized. *Part Type:* groups *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Missingness

(missingness). The part type for missingness values. Only used for the dictionary entries of parts. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:*

Parts

naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 15 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Missingness sets

(missingnessSets). Missingness sets for measures, methods or attributes. Inputted data into an ODM field may have missing data. The missingness set lists valid missing categories for the measure, method, or attribute. *Part Type: partType Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 15 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Mixed/homogenized sample

(mix). Mixed or homogenized sample or fraction analyzed. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Millilitres

(ml). Millilitres of volume *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Megalitres per day (ML/d)

(mld). Megalitres per day. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres

(mm). Unit part for the SI unit of millimetres. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Parts

Measure, method, or attribute

(mmaSet). The set for a measure, method, or attribute. Only applicable for categorical fields, and the categories that populate them. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 12 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Measure, method, or attribute sets

(mmaSets). A set of categories. For example, site type has a list of different sites, such as a wastewater treatment plant, a septic tank, or a pumping station. Only used for the dictionary entries of parts. *Part Type*: part-Type *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Model

(model). Model number or version of the instrument. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Months

(months). A unit for indicating a length of time in months. *Part Type:* units *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* time *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Moore swab passive sample

(moorSw). Moore swab passive sample. Use collectionPeriod to describe how many hours the Moore swab was used to collect the sample *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Metres per second

(ms). metres per second. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:*

Parts

naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage

(ms2Col). Measure of the amount of ms2 coliphage. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage spike target

(ms2ColMat). ms2 coliphage is used as the recovery efficiency control target. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Matrix spike duplicate

(msd). A known amounts of an analyte or representative compounds are added in the laboratory to a second aliquot of the sample used for matrix

spike. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Moore swab passive sample depreciated

(mspsDep). Moore swab passive sample. Deprecated in version 2, please do not use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Major sewer pipeline

(mSwrPpl). Major sewer pipeline *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Mu

(mu). B.1.621 *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

murine coronavirus

(muCo). Measure of the amount of murine coronavirus. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

murine coronavirus spike target

(muCoMat). Murine coronavirus is used as the recovery efficiency control target. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple fraction

(mul). Multiple fractions were analyzed separately and aggregated together post-analysis. This fractionID is intended for use for archival data and should not be used for newly added new moving forward (2022). Please also specify in the notes section which fractions were used for the multiple fractions (ex. solid and liquid) to avoid loss of information. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Multiple Purpose

(multiple). A measure or sample taken for multiple purposes, not easily captured by the other purpose categories. Useful for parent samples where various subsamples will have different purposes/ *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Municipalities or communes

(municipalLevel). Municipalities or communes *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* na-

Parts

CompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: ISCU *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Municipality

(municip). A complete municipality, this specifies an entire metropolitan area, either a city, town, etc. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Mutations class

(mutation). Measures and methods related to mutations. *Part Type*: classes *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

N sars-cov-2 gene target

(n). N sars-cov-2 gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Combined N1 and N2 sars-cov-2 gene target

(n1n2). Combined N1 and N2 sars-cov-2 gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n211i omicron-variant gene target { #n211i }

(n211i). n211i omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n679k omicron-variant gene target { #n679k }

(n679k). n679k omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n856k omicron-variant gene target { #n856k }

(n856k). n856k omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

n969k omicron-variant gene target { #n969k }

(n969k). n969k omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Not applicable

(NA). The field for which the expected value is not a property of the described object. Missing value indicator for missing data with no explanation. *Part Type:* missingness *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation not applicable

(naAggr). Not applicable for aggregations. Used for parts for which an aggregation value is not applicable or appropriate. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scale not applicable

(naAggrScale). Not application for aggregation sets. Used for parts for which an aggregation scale value is not applicable or appropriate. *Part Type:* aggregationScales *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality*

Parts

Set: NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation set not applicable

(naAggrSet). Not application for aggregation sets. Used for parts for which an aggregation set value is not applicable or appropriate. *Part Type:* aggregationSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Class not applicable

(naClass). Class is not applicable. Use this for all parts that don't have classes (most non-measure or method parts) *Part Type:* classes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Compartment not applicable

(naCompartment). Compartment not applicable. *Part Type:* compartments *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment*

Set: naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 13 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Compartment set not applicable

(naCompartmentSet). Compartment not applicable. *Part Type:* compartmentSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 13 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Domain not applicable

(naDomain). Not applicable. Use ‘not applicable’ when there is no domain for the part *Part Type:* domains *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Group not applicable

(naGroup). Used when there is no group; ie. group is not applicable. Use 'not applicable' when there is no group for the part *Part Type*: groups *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Name

(name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity. Name of any entity. Use a prefix to distinguish names from different report tables. See fully specified name (FSN) list. i.e. Si-name = Site name, Po-name = Polygon name, Or-name = Organization name, Me-name = Method-name. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Not a number

(nan). The outcome of a measurement is not a valid number (plus or minus infinity, error, ...) Missing values indicator for when a numeric value is

expected but not given. *Part Type*: missingness *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Nomenclature not applicable

(naNomenclature). Not applicable. Use ‘not applicable’ when there is no nomenclature for the part *Part Type*: nomenclatures *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen not applicable

(naSpecimen). Non applicable specimen. Use ‘not applicable’ when there is no specimen/set for the part *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Specimen set not applicable

(naSpecimenSet). A specimen set for when specimen/specimen set is not applicable. Use 'not applicable' when there is no specimen/set for the part
Part Type: compartmentSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Native Name

(natName). The native name of the language, i.e. what the language is called by its speakers. Part of the metadata for a languageID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Unit not applicable

(naUnit). Not appliable for units. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Unit set not applicable

(naUnitSet). Not applicable for unit sets. *Part Type*: unitSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

National Center for Biotechnology Information

(ncbi). NCBI header *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

NCBI - notes

(ncbiNotes). NCBI notes *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*:

Parts

varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Neighbourhood

(neigh). A municipal neighbourhood, this specifies a sub-section of a larger municipality. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Local administrative units or neighborhoods

(neighborLevel). Local administrative units or neighborhoods *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* ISCU *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population of newly vaccinated persons

(newVax). Population of newly vaccinated persons *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

NextClade nomenclature

(nextclade). Specifies variant or genetic nomenclature as set out by NextClade. *Part Type*: nomenclatures *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Normalized geometric mean

(ngmn). Geometric mean, normalized. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Ammonium Nitrogen

(nH4N). Ammonium nitrogen concentration, as N. *Part Type*: measurements *Domain*: che *Speciment ID*: siSaSpecimenSet *Compartment Set*:

Parts

waterCompartmentSet *Group*: miscMeas *Class*: standardConc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

niid_2019-ncov_n sars-cov-2 gene target

(niid19). niid_2019-ncov_n sars-cov-2 gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns

(noConcern). A flag to indicate there is are no quality concern about the measure or sample. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: 1 *Maximum Value*: 1 *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sample not labelled

(noLabel). Sample had no label *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

No liquid concentration, liquid recombined with separated solids

(noliquid). No liquid concentration, liquid recombined with separated solids *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclature

(nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases. Currently only used for class = variant, mutation, or disease. Valid options are currently restricted to 'NextClade', 'Pangolin', 'WHO', or 'ICD'. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*:

Parts

NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclatures

(nomenclatures). A classification system to report the measure class. See *partID* = *nomenclatureID*. *Part Type*: *partType* *Domain*: *naDomain* *Speciment ID*: *naSpecimenSet* *Compartment Set*: *naCompartmentSet* *Group*: *naGroup* *Class*: *naClass* *Nomenclature*: *naNomenclature* *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: *naUnitSet* *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

NORMAN

(norman). NORMAN header *Part Type*: *crosswalkTableSet* *Domain*: *naDomain* *Speciment ID*: *naSpecimenSet* *Compartment Set*: *naCompartmentSet* *Group*: *naGroup* *Class*: *naClass* *Nomenclature*: *naNomenclature* *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: *naUnitSet* *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: *varchar* *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

NORMAN - notes

(normanNote). Norman notes *Part Type*: *crosswalkTableSet* *Domain*: *naDomain* *Speciment ID*: *naSpecimenSet* *Compartment Set*: *naCompartmentSet* *Group*: *naGroup* *Class*: *naClass* *Nomenclature*: *naNomenclature*

Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Norovirus G1

(norog1). Norovirus genogroup 1 *Part Type: measurements Domain: bio Speciment ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: virusMisc Class: allele Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: geneticUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: seeUnitData Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Norovirus G2

(norog2). Norovirus genogroup 2 *Part Type: measurements Domain: bio Speciment ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: virusMisc Class: allele Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: geneticUnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: seeUnitData Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Notes

(notes). A note used to describe details that are not captured in other attributes There is no recommended structure of a note. Use notes as

Parts

needed. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 200 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample missing time stamp

(noTime). Sample is missing the autosampler time; incomplete metadata on collection. *Part Type*: qualityIndicators *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Not reported

(nr). A value could have been recorded, however, it was not. Missing value indicator for data that is absent because it was not reported. *Part Type*: missingness *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported, Not applicable missing set

(nrNAMissingnessSet). Not reported, Not applicable missing set. *Part Type:* missingnessSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sarbecovirus-specific N sars-cov-2 gene target

(nSarbec). Sarbecovirus-specific N sars-cov-2 gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nephelometric turbidity unit

(ntu). Nephelometric Turbidity Units, a unit used in the measurement of turbidity (see “turbidity” under measureIDs). Nephelometric turbidity units (NTU) are based on white light (400–680nm) and 90° incident angle. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference*

Parts

Link: <https://www.sciencedirect.com/topics/engineering/nephelometric-turbidity-unit> *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nuclisens automated magnetic bead extraction kit

(nucAuto). Nucleic acid extraction performed using the nuclisens automated magnetic bead extraction kit. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nuclisens manual magnetic bead extraction kit

(nucManu). Nucleic acid extraction performed using the nuclisens manual magnetic bead extraction kit. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Null

(null). A logical representation of a statement that is neither TRUE nor FALSE. Missing value indicator when the value is neither TRUE nor FALSE. *Part Type:* missingness *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

National Wastewater Surveillance System

(nwss). NWSS header *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

NWSS - notes

(nwssNotes). NWSS notes *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

Parts

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

NWSS - process

(nwssProcess). NWSS process *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Water

(nww). Non-wastewater, coming from any kind of water body *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

coronavirus OC43

(oc43). Measure of the amount of coronavirus OC43. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference*

Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

oc43 spike target

(oc43Mat). Human coronavirus OC43 is used as the recovery efficiency control target. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ocean, natural water body

(ocean). Ocean, natural water body *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Other Collection deprecated

(oColDep). Other type of collection method. Add description to collectionOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or

Parts

at phesd_odm@ohri.ca *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Lab Analysis

(ola). Offline laboratory analysis. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Omicron BA.1

(omicr1). Omicron B.1.1.529 *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* who *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2

(omicr2). Omicron BA.2 *Part Type*: measurements *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2.75

(omicr275). Omicron BA.2.75 *Part Type*: measurements *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.4

(omicr4). Omicron BA.4 *Part Type*: measurements *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Omicron BA.5

(omicr5). Omicron BA.5 *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: variant *Nomenclature*: who *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Onset date

(onsDate). Earliest that symptoms were reported for this case. This data is often not known and reported. In lieu, episode is used. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Online Sensor

(onse). Online sensor An online sensor *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0

Maximum Length: 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Ontology reference

(ontologyRef). Ontology reference for a part. This field contains a link to an existing ontology reference for the part. ODM content spans several existing ontologies. ODM does not strive to generate a new ontology, rather ODM seeks to harmonize and extend dictionaries for application in environmental and public health surveillance. *Part Type:* part-Support *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 200 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Orthophosphates

(ophos). Ortho-phosphate concentration. *Part Type:* measurements *Domain:* che *Speciment ID:* siSaSpecimenSet *Compartment Set:* water-CompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Other residential building

(orb). Individual residential buildings or institutions not captured in other categories *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

OR Boolean aggregation

(orBoo). “OR” aggregation. If any value in the aggregation is “TRUE” then the OR aggregation is also “TRUE”. If all values in the aggregation is “FALSE” then the OR aggregation is also “FALSE”. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ORF1a sars-cov-2 gene target

(orf1a). ORF1a sars-cov-2 gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal

Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ORF1ab sars-cov-2 gene target

(orf1ab). ORF1ab sars-cov-2 gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ORF1b sars-cov-2 gene target

(orf1b). ORF1b sars-cov-2 gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Organization ID

(organizationID). A unique identifier for the organization to which the reporter is affiliated. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:*

Parts

NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 100 *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Organization table

(organizations). The table that contains information about a laboratory. adapt from version 1 documentation *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: programDescr *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Organizations table column order

(organizationsOrder). Specifies the order of the columns in a Organizations table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization table required headers

(organizationsRequired). Specifies the columns required in a Organizations table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup

Class: naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Organization level

(orgLevel). The geographic level of an organization. There are six levels based on International Standard Classification of Administrative Units (ISCU). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* orgLevelSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Organization level set

(orgLevelSet). Categories of organization levels. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Organization sector

(orgSector). The sector of an organization *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: orgSectorSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization sector set

(orgSectorSet). Cateogries of organization sectors. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization Type

(orgType). Specifies the type or purpose of a given organization. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: orgTypeSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Organization type set

(orgTypeSet). The set for storing organization types. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample origin

(origin). An attribute of a sample specifying the origin. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample Origin set

(originSet). The set for storing the valid categorical values of sample origin. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Other Collection deprecated

(otco). Other type of collection method. Add description to collectionOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Other aggregation scale

(othAgg). Specifies an aggregation scale that can't be described as either "qualitative" or "quantitative". *Part Type*: aggregationScales *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Other

(other). Other Other category. A categorical response for several variables in version 1. Also accompanied by a free text field. Deprecated in version 2, with common text fields added as new categories in version 2.0.0

Part Type: categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation set

(otherAggrSet). Aggregation set used for unitless measures. *Part Type:* aggregationSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation deppricated

(otherDep). Other aggregation method. Add description to aggregationOther Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA

Parts

Maximum Length: NA *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Access to other prov deprecated

(otherProvDep). If this is ‘no’, this data will not be available to other data providers not listed before. If missing, data will be available to other data providers not listed before *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Other report table set

(otherReportTableSet). Additional tables that form the full data model. *Part Type:* dictSets *Domain:* naDomain *Specimen ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Variants Other

(otherV). Used for reporting variants detected for which a specific measureID doesn’t already exist. The ODM team will periodically migrate collected responses in this category into set measureIDs. Please write

in an official name for the new variant that was detected, and specify which nomenclature you're using with the nomenclatureID. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* variant *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other Site - sample depreciated

(otsisaDep). Other type of site. Add description to typeOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at ph-esd_odm@ohri.ca *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Other Site - Site depreciated

(otsisiDep). Other site type. Add description to typeOther *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA

Parts

Data Type: NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Other Unit deprecated

(otunDep). Other measurement of viral copies or wastewater treatment plant parameter. Add description to UnitOther. Deprecated as of version 2, please avoid use. Specify in type, or contact the PHES-ODM research team to add your specific item via a github issue (<https://github.com/Big-Life-Lab/PHES-ODM/issues>), or at phesd_odm@ohri.ca *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Outbreak

(outb). Measure to indicate outbreak status. Given that when outbreaks occur, full data on case counts may not be available or completely accurate, using this measure indicates that case counts may be approximate. The datetime for the measure helps indicate that start and end periods for an outbreak. *Part Type:* measurements *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* humanCompartmentSet *Group:* miscMeas *Class:* outbreak *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* outbreakSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak End

(outbEnd). Indicates an outbreak has ended. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Outbreak - on-going

(outbOngoing). Indicates an outbreak is on-going. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Outbreak class

(outbreak). Measures and methods relating to public health outbreaks. *Part Type: classes Domain: bio Speciment ID: poSpecimenSet Compartment Set: anyCompartmentSet Group: naGroup Class: outbreak Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: unitlessUnitSet Aggregation Scale: seeUnitAggScale Quality Set:*

Parts

NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak set

(outbreakSet). set for the valid values of the outbreak measure. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: outbreakSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak start

(outbStart). Indicates an outbreak began. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

p100l delta-variant gene target

(p100l). p100l delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet

Group: sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

p2046l delta-variant gene target

(p2046l). p2046l delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

p2287s delta-variant gene target

(p2287s). p2287s delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

p3395h omicron-variant gene target { #p3395h }

(p3395h). p3395h omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

p681r delta-variant gene target

(p681r). p681r delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Paired Layout

(pairLay). Specifies the paired layout for a sequencing method *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Pangolin nomenclature

(pangolin). Specifies variant or genetic nomenclature as set out by the Phylogenetic Assignment of Named Global Outbreak Lineages (Pangolin). *Part Type:* nomenclatures *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parent dataset ID

(parDatasetID). The datasetID that is a parent to another datasetID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parent sample ID

(parent). If this sample has been pooled into one big sample for analysis this indicates the sample of the larger pooled sample. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Parts

*Length: 0 Maximum Length: 30 Part Status: **deprecated**.* First released: 1.1.0. Last updated: 2.0.0.

Parent Site ID

(parSiteID). The siteID that is a parent to another siteID; usually the sub-site will be contained in the parentSite. An example would be surface testing in a LTC facility. The LTC facility would be the parentSiteID, and the various rooms/surfaces would be the sub-sites, or child-sites. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Part description

(partDesc). The description of the part. Provides description of the part, used for translation. *Part Type: partSupport Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active.* First released: 2.0.0. Last updated: 2.0.0.

Part identifier

(partID). The unique identify of any entity within the dictionary. Every entity in the ODM has a partID. The partID must be unique with no du-

plcates. partIDs are short. partIDs that are constructed by different part attributes can be up to 28 characters (check max), but each part attribute should be 1 to 5 characters. partIDs are machine readable, but a knowledgeable user can understand what they represent. partID are used, for example, as a header in an Excel input table. Numbers are allowed but there are no other special characters other than ‘:’. ‘:’ is a reserved character that is used to generate a partID using a component or table attributes. A partID for a subComponent is automatically generated using attributes of a component. For example, the partID for the component ‘SARS-CoV-2’ is ‘covid’. ‘Covid’ can have many subComponents including different groups (Covid-19-al, Covid-19-va, covid-pr (allele, variant or protein) or units and aggregations (i.e. covid-al-N1-mean-cpl representing covid allele N1, mean observation of viral copies per litre of effluent). A similar approach to naming subComponent is used for naming parts in a reprotTable. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 28 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Part instruction

(partInstr). Additional notes and instructions on how a part is used and/or defined. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 200 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Label

(partLabel). A human readable label of a part. Typically, a part label has no acrynomns (every word is spelled out). Equivalent to a LOINC common name. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts Look-up table

(parts). Look up table containing all parts in of the data model. Contains all parts, including self-referential parts. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* lookup *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts table column order

(partsOrder). Order of headers in the Parts table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Parts table required headers

(partsRequired). Required headers in the Parts table. *Part Type: table-Support Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Part types

(partType). Part types describe the purpose or use of the part. Enivornment data has three main part types: measure (i.e. covN1 viral measures, wasteawater flow rate, temperature), method (e.g. how the measure was taken), and attribute (such as a site name). See description of each of the part types within categorySetID = partTypeCatSets. *Part Type: part-Type Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: categorical Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 15 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

Parts

Primary Clarifier Effluent deprecated

(pceDep). Effluent obtained after primary clarifiers. Deprecated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Postal or Zip Code

(pCode). The zip code or postal code for a given address, specifying a specifc geographic area. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qPCR Class

(pcr). Measures and methods relating to qPCR. *Part Type:* classes *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

PCR method

(pcrmeth). Description of the PCR method. Description of the PCR method used *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* pcrSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

PCR quality set

(pcrQualitySet). Quality set for PCR measures. *Part Type:* qualityInd-Sets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

PCR sequencingsSelection method

(pcrSeq). Specifies the PCR selection method for sequencing, ie. that a pre-determined primer was used. Primer should be specified if this is selected. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

Parts

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

PCR Method set

(pcrSet). The set capturing containing all of the valid methods for the PCR methodID. *Part Type*: mmaSets *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: pcrSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Primary Clarifier Sludge depreciated

(pcs). Sludge produced by primary clarifiers. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Primary clarifier effluent

(pEfflu). Effluent obtained after primary clarifiers *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*:

naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Polyethyleneglycol (PEG) precipitation

(peg). Peg precipitation *Part Type*: categories *Domain*: che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Percent

(perc). Percentage. *Part Type*: units *Domain*: allDo *Speciment ID*: anySpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Percent recovery

(percRec). Percent of the surrogate recovery for a recovery efficiency control assay. Use this for reporting the results of any recovery control assay. *Part Type*: units *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass

Parts

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

pH

(ph). pH measurement *Part Type:* measurements *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* pHClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* unitlessUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 1 *Maximum Value:* 14 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Access to phac

(phac). If this is 'no', the data will not be available to employees of the Public Health Agency of Canada - PHAC. If missing, data will be available to employees of the Public Health Agency of Canada - PHAC. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

PHAGE

(phage). PHAGE header *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

PHAGE - notes

(phageNotes). PHAGE notes *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

pH class

(pHClass). Measures and methods relating to pH. *Part Type:* classes *Domain:* che *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* pHClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* unitlessUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Phenol chloroform

(phenCl). Nucleic acid extraction performed using phenol chloroform. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Pseudomonas virus phi6

(phi6). Measure of the amount of pseudomonas virus phi6. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

phi6 spike target

(phi6Mat). Pseudomonas virus phi6is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Contact phone

(phone). Contact phone number, for the lab. *Part Type:* attributes
Domain: naDomain *Speciment ID:* naSpecimenSet *Compartment Set:*
 naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:*
 naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*
 naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA
Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*
Length: 10 *Maximum Length:* 12 *Part Status:* active. First released: 1.0.0.
 Last updated: 2.0.0.

Total phosphorous

(phos). Total phosphorous *Part Type:* measurements *Domain:* phy *Spec-*
iment ID: siSaSpecimenSet *Compartment Set:* waterCompartmentSet
Group: miscMeas *Class:* standardConc *Nomenclature:* naNomenclature
Ontology Reference: NA *Category Set ID:* NA *Unit Set:* stdConcentratio-
 nUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference*
Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maxi-*
imum Value: seeUnitVal *Minimum Length:* NA *Maximum Length:* NA
Part Status: active. First released: 1.1.0. Last updated: 2.0.0.

Total Phosphates

(phostot). Total phosphates *Part Type:* measurements *Domain:* phy
Speciment ID: siSaSpecimenSet *Compartment Set:* waterCompartmentSet
Group: miscMeas *Class:* standardConc *Nomenclature:* naNomenclature
Ontology Reference: NA *Category Set ID:* NA *Unit Set:* stdConcentratio-
 nUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference*
Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maxi-*
imum Value: seeUnitVal *Minimum Length:* NA *Maximum Length:* NA
Part Status: active. First released: 1.1.0. Last updated: 2.0.0.

Parts

PHRED quality score

(phred). PHRED Quality Score For Sequencing *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: sequence *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Physical property

(phy). A physical property or object not characterized by life or chemistry. *Part Type*: domains *Domain*: phy *Speciment ID*: anySpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Physical class

(physical). Measures and methods related to generic physical properties. *Part Type*: classes *Domain*: phy *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: physical *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: volumeUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Primary key

(pK). Primary key for a table. All report tables have a primary key. Each row in a report table has unique identifier recorded in the primary key. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean

(pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Mostly used for reporting specific alleles (N1, N2, E, etc.) *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-normalized mean

(pmmovNorm). Mean measure normalized to amount of PMMoV. Mostly used for reporting specific alleles (N1, N2, E, etc.) *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*:

Parts

naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Population

(po). An measure or observation for a geographic area or population. Examples include human compartment measures such as case numbers, hospitalization rates for diseases or surface or air compartment measures such as snow coverage or ambient temperature. *Part Type*: specimens *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Post-Grit depreciated

(pogr). Raw wastewater after a treatment plant's headworks. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Polygon ID

(polygonID). Unique identifier for the polygon. A polygon is a geographic area representing the capment area of an environmental sample. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Polygon table

(polygons). The table that contains information about the geometry of a geographic area. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Polygons table column order

(polygonsOrder). Specifies the order of the columns in a Polygons table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum*

Parts

Value: NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Polygon table required headers

(polygonsRequired). Specifies the columns required in a Polygons table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Polygon Population

(polyPop). An attribute of a polygon, which specifies the population of that polygon. A rough estimate of the number of human residents. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Pooled

(pooled). Is this a pooled sample, and therefore composed of multiple child samples obtained at different sites *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature

Ontology Reference: NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Type of date for case reporting deprecated

(popDateTypeDep). Type of date used for confirmed cases. Typically, report or episode are reported. Onset and test date is not usually reported within aggregate data. *Part Type:* attributes *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* humanCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* methodConcSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Case reporting date set

(popDateTypeSet). set for the date of case reporting. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* popDateTypeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Population equivalents

(popEq). A unit for measuring of design capacity for wastewater treatment plants, represented as the approximate amount of water a single person would use. *Part Type:* units *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Population group

(popGrp). A group of measures/methods related to population-level factors. Examples might be case numbers, hospitalization rates, etc. *Part Type:* groups *Domain:* allDo *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* popGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population Served

(popServ). An attribute of a site, which specifies the population of/population served by a given site. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:*

categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Population unit set

(populationUnitSet). Unit set for hospital-related measurements. *Part Type*: unitSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: humanCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: populationUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Population or sample set

(poSaSpecimenSet). A specimen set that inculdes population or sample specimen. *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: poSaSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Population or site set

(poSiSpecimenSet). A specimen set that inculdes population or site specimen. *Part Type*: specimenSets *Domain*: naDomain *Speciment ID*: poSiSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA

Parts

Category Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population specimen set

(poSpecimenSet). A specimen set that includes only a population specimen. *Part Type:* specimenSets *Domain:* naDomain *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Percent positive

(pp). Percent positive of sample measures. Can use for Moore swabs, etc. *Part Type:* units *Domain:* bio *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

parts per million

(ppm). Parts per million. *Part Type:* units *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Refer-*

ence: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

PMMoV-CP

(ppmv). Pepper mild mottle virus capsid protein gene region *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Percent positivity rate

(pprt). Percent positivity rate of tests conducted within a day in a given region. *Part Type*: units *Domain*: bio *Speciment ID*: poSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Percent primary sludge

(pps). Percentage of total solids, for primary sludge. *Part Type*: units *Domain*: phy *Speciment ID*: siSaSpecimenSet *Compartment Set*: waterCom-

Parts

partmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: 100 *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Precipitation class

(precipitation). Measures and methods related to precipitation. *Part Type*: classes *Domain*: phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: precipitation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: precipitationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Precipitation unit set

(precipitationUnitSet). Unit set for percipitation measurements. *Part Type*: unitSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: precipitationUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Pretreatment

(pretreat). Was the sample chemically treated in anyway with the addition of stabilizers or other Pretreatment method can be describe using a method. See methodID. *Part Type:* methods *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Genetic primer

(primer). Method ID used for indicating the primer used for a PCR or sequencing assay. *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* primerSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Genetic primer set

(primerSet). The set for genetic primers used in sequencing or PCR assays. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* primerSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA

Parts

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Private sector

(priv). The category of organization type used for private sector or non-profit groups that are not otherwise captured by the academic category. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Processing group

(procGrp). A group of measures/methods related to processing samples for analysis. *Part Type:* groups *Domain:* allDo *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Program decription tables

(programDescr). Tables used to describe surveillance and testing programs. Program description tables record metadata on the organizations, locations, and protocols. Information include the name, location, and other contact information. These tables are usually completed once and they are

updated only when contact information changes. *Part Type:* classes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

promega automated tna kit

(promAuto). Nucleic acid extraction performed using the promega automated tna kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

promega ht tna kit

(promHt). Nucleic acid extraction performed using the promega ht tna kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

promega manual tna kit

(promManu). Nucleic acid extraction performed using the promega manual tna kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

promega wastewater large volume tna capture kit

(promWW). Nucleic acid extraction performed using the promega wastewater large volume tna capture kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Proportion of total

(prop). Proportion as a percent of total. *Part Type:* units *Domain:* naDomain *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of variant in sample

(propV). Proportion of a variant as percent of total variants. *Part Type:* units *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Protocol ID

(protocolID). A unique identifier for a given protocol. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID container

(protocolIDContainer). Unique identifier for the protocol and the steps and other protocols that make it up. Protocol ID Container in populated by a protocolID (partID = protocolID) *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0

Parts

Maximum Length: 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID object

(protocolIDObj). The object of the relationship between one protocol and another protocol or step, within a given protocol ID container. Populate this field using Protocol ID (partID = protocolID). Together with 'protocolIDSub' or 'stepIDSub' and 'relationshipID', the protocol relationship is defined as: [protocolIDSubject or stepIDSubject] [relationshipID] [protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID subject

(protocolIDSub). The subject of the relationship between one protocol and another protocol or step, within a given protocol ID container. Populate this field using Protocol ID (partID = protocolID). Together with 'protocolIDObj' or 'stepIDObj' and 'relationshipID' the row of 'protocolRelationships' should say: [protocolIDSubject or stepIDSubject] [relationshipID] [protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol relationships table

(protocolRelationships). The table that contains the organizational information and details (such as order) about a given protocol, recorded as a series of protocol steps (protocolSteps). Protocols are a group of related protocol steps. For example a protocol step can describe a sample concentration, extraction, inhibition or a measurement. Protocols can include both methodID (a procedure) and measureIDs (a measure). An example of a method is centrifugation; where the centrifugation rotation speed is described as a measure. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol relationships table column order

(protocolRelationshipsOrder). Specifies the order of the columns in the Protocol Organization table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol Relationship table required headers

(protocolRelationshipsRequired). Specifies the columns required in the Protocol Organization table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol Relationships set

(protocolRelSet). set for valid values of relationshipID in the protocolRelationships table. Use only for protocolRelationships table relationshipID. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* protocolRelSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocols table

(protocols). The table for protocols. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum*

Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocols table column order

(protocolsOrder). Specifies the order of the columns in the Protocols table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocols table required headers

(protocolsRequired). Specifies the columns required in the Protocols table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol steps table

(protocolSteps). The table for collecting metadata on individual steps in a protocol, methodological process, or assay. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* programDescr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet

Parts

Aggregation Scale: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Protocol steps table column order

(protocolStepsOrder). Specifies the order of the columns in a Protocol Steps table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol steps table required headers

(protocolStepsRequired). Specifies the columns required in a Protocol Steps table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol table set

(protocolTableSet). Tables holding information on the methods used for sample collection or measurement. *Part Type*: dictSets *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: naCompartmentSet

Group: naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol version

(protocolVersion). Specifies the version of a method set. Version of the method set. Semantic versioning is recommended. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Access to prov ha

(provHA). If this is 'no', this data will not be available to provincial health authorities. If missing, data will be available to provincial health authorities. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Provisional report

(provisional). A provisional or interm report. Use for measurement that are not yet finalized. Typically this purpose is not publicly reported. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Primary clarifier sludge

(pSludge). Sludge produced by primary clarifiers *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Pumping station

(pStat). Pumping station *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar

Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.

Post-grit

(pstGrit). Raw wastewater after a treatment plant's headworks (post grit or removal of large solids at a treatment plant but prior to a primary clarifier) *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: waterCompartmentSet Group: colGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.*

Pretreatment description - depreciated

(ptDescDep). If preTreatment then describe the treatment that was performed. Depreciated in version 2, please do not use. Use pretreatment methodID and specify further in summary and/or notes, if necessary. *Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: procGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 255 Part Status: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.*

Parts

Number of positive tests

(ptot). Number of positive tests conducted within a day in a given region. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Public healith agency

(pubHealth). Public health agency. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Access to public

(public). If this is 'no', this data will not be available to the public. If missing, data will be available to the public. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

0 *Maximum Length*: 1 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Puro virus

(puro). Measure of the amount of puro virus. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Puro virus spike target

(puroMat). Puro Virus is used as the recovery efficiency control target. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Purpose

(purpose). The reason the measure or sample was taken. See allowed categories: report, quality control, validation study, test *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: purposeSet

Parts

Unit Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Purpose set

(purposeSet). Purpose set for a sample or measure. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* purposeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

q19e omicron-variant gene target { #q19e }

(q19e). q19e omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

q493r omicron-variant gene target { #q493r }

(q493r). q493r omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature

ture *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-
UnitSet Aggregation Scale: seeUnitAggScale *Quality Set*: NA *Reference*
Link: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum*
Value: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

q498r omicron-variant gene target { #q498r }

(q498r). q498r omicron-variant gene target *Part Type*: measurements *Do-*
main: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet
Group: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature
Ontology Reference: NA *Category Set ID*: NA *Unit Set*: genetic-
UnitSet Aggregation Scale: seeUnitAggScale *Quality Set*: NA *Reference*
Link: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum*
Value: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

q954h omicron-variant gene target { #q954h }

(q954h). q954h omicron-variant gene target *Part Type*: measurements *Do-*
main: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet
Group: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature
Ontology Reference: NA *Category Set ID*: NA *Unit Set*: genetic-
UnitSet Aggregation Scale: seeUnitAggScale *Quality Set*: NA *Reference*
Link: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum*
Value: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quaility assurance method

(qaqc). Quality assurance and quality control method. Description of the
 quality control steps taken *Part Type*: methods *Domain*: allDo *Speciment*

Parts

ID: anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Quality concerns

(qf1). A flag to indicate there is a quality concern, not otherwise specified. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* 0 *Maximum Value:* 0 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

qiagen allprep dna/rna kit

(qgDNARNA). Nucleic acid extraction performed using the qiagen allprep dna/rna kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen digital PCR

(qgDpcr). Describes a PCR analysis done using Qiagen's digital PCR technology. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen ez1 virus mini kit v2.0

(qgEz1). Nucleic acid extraction performed using the qiagen ez1 virus mini kit v2.0. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep powerfecal dna/rna kit

(qgPwrFecal). Nucleic acid extraction performed using the qiagen allprep powerfecal dna/rna kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

qiagen allprep powerviral dna/rna kit

(qgPwrViral). Nucleic acid extraction performed using the qiagen allprep powerviral dna/rna kit. *Part Type*: categories *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen powerwater kit

(qgPwrWtr). Nucleic acid extraction performed using the qiagen powerwater kit. *Part Type*: categories *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen qiaamp buffers with epoch columns

(qgQiAmp). Nucleic acid extraction performed using the qiagen qiaamp buffers with epoch columns. *Part Type*: categories *Domain*: bio *Specimen ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy kit

(qgRneasy). Nucleic acid extraction performed using the qiagen rneasy kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy powermicrobiome kit

(qgRneasyPwr). Nucleic acid extraction performed using the qiagen rneasy powermicrobiome kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quantitative PCR

(qpcr). Real-time PCR, also called ‘quantAggScale’ PCR *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Qualitative

(qualAggScale). The qualitative aggregation scale. *Part Type:* aggregationScales *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* qualAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality control

(qualityControl). A measure or sample taken for the purpose of quality control. Measures with this category are take to assess quality control. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality flag

(qualityFlag). A field for reporting any quality concerns - of lack thereof - for a sample or measure. Types of quality flag vary on context; please see the qulaityFlagSet to be confirm. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:*

categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Quality report ID

(qualityID). A unique identifier for a given quality report. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicators

(qualityIndicators). A measure of the quality of a reported value or sample. Only used for the dictionary entries of parts. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicator sets

(qualityIndSets). Sets of quality indicators or measures. For example, PCR have a quality measure set. *Part Type*: partType *Domain*:

Parts

naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Quality reports table

(qualityReports). The table for recording the various quality metrics and indicators for samples and measures. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: results *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Quality repots table column order

(qualityReportsOrder). Specifies the order of the columns in the Quality Reports table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Quality report table required headers

(qualityReportsRequired). Specifies the columns required in the Quality Reports table. *Part Type:* tableSupport *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quality set

(qualitySetID). The quality set that corresponds to a given part. Only applicable for samples and measures. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Quantitative

(quantAggScale). The “quantAggScale” aggregation scale. *Part Type:* aggregationScales *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum*

Parts

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 2.0.0. Last updated: 2.0.0.

R squared

(r2). R-squared value of the calibration curve. Used for storing calibration curve information, potentially also applicable for other curve types. *Part Type:* measurements *Domain:* allDo *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standard-Curve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

r203m delta-variant gene target

(r203m). r203m delta-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Raindance digital PCR

(rainDpcr). Describes a PCR analysis done using Raindance's digital PCR technology. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:*

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Rainy

(rainy). Qualitative category for the weather measure, specifying a rainy day. *Part Type:* categories *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Random sequencing selection method

(ranSeq). Specifies the random selection method for sequencing No primer set should be used if this is selected. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Ratio

(ratio). Ratio (unitless) Report as a real number or fraction. i.e 10:2 ratio is reported as 5. The specifics of the ratio (the numerator and the demoninator) are described in the measure. *Part Type*: units *Domain*: naDomain *Speciment ID*: anySpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Raw wastewater

(raws). Raw wastes water sample *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage at site

(rawWW). Wastewater without any form of treatment *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum*

Length: 0 Maximum Length: 30 Part Status: active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage downstream from a site

(rawWWdown). Downstream from a site. See partType = rawWWup.
*Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet
 Compartment Set: waterCompartmentSet Group: naGroup Class: naClass
 Nomenclature: naNomenclature Ontology Reference: NA Category Set ID:
 NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference
 Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA
 Minimum Length: 0 Maximum Length: 30 Part Status: active. First
 released: 2.0.0. Last updated: 2.0.0.*

Raw sewage upstream from a site

(rawWWup). Upstream from a site. Used when there is not direct access to a site's wastewater inconjuction with downstream sample. For example, if the site is a long-term care home, but ther is not access to the building or property cleanout. Use two sameple, one sample upstream and one sample downstream. *Part Type: categories Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: waterCompartmentSet Group: colGrp Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

IP2 rdrp sars-cov-2 gene target

(rdrpIP2). IP2 rdrp sars-cov-2 gene target *Part Type: measurements Do-
 main: bio Speciment ID: saSpecimenSet Compartment Set: anyCompartment*

Parts

mentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

IP4 rdrp sars-cov-2 gene target

(rdrpIP4). IP4 rdrp sars-cov-2 gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: genetic-UnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Date sample recieved

(recDate). The date the sample was received at the laboratory for analysis. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: datetime *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Recovered patients

(recov). Units for describing a population measure of patients who have recovered from a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type*: units *Domain*: bio *Speciment ID*: poSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Reference link

(refLink). Link to the reference material for a part. May link to literature on a method, measure, etc. Link to standard operating procedure for a measure or part. Part reference is a reference to a technical document that describes the part in detail. The reference could be a URL, DOI, or reference to a technical report. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 255 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Regular

(regular). A measure or sample taken for surveillance or epidemiology. A 'regular' measure is the default. Missing data will be recoded to this pur-

Parts

pose. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Regular reports table set

(regularReportTableSet). Tables used for daily reporting of new measurements and information on sample collection. *Part Type*: dictSets *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Relative humidity

(relHum). The unit of relative humidity, or the air-water mixture. The ratio of the partial pressure of water vapor in the mixture to the equilibrium vapor pressure of water over a flat surface. of pure water at a given temperature. *Part Type*: measurements *Domain*: phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: humid *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: relHumidUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum*

Length: NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Relative humidity unit set

(relHumidUnitSet). Unit set for relative humidity measurements. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* relHumidUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Representative grab sample

(repGrab). A single large representative grab sample. Deprecated in version 2, please do not use. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Replicate Type set

(replicateSet). The set for storing valid categorical values of replicate type. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set*

Parts

ID: replicateSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

replicate number

(repNum). The replicate number for a Ct or Cq value - used for specifying out a protocol standard curve. Only applies and used for aggregated standard curves, or curves that are used across multiple plates. Otherwise, these details can be captured in a measure report and specified by measure index. *Part Type:* measurements *Domain:* bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* see-UnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Primary reporting authority

(repOrg1). The primary or most responsible authority for routine surveillance reports or findings, or where the site is located. Use the 'organizations' table to describe attributes such as the name and address of the organization. Use 'orgType', 'orgLevel' and 'orgSector' to describe additional attributes. For example, whether the organization is a government agency ('orgType'), local ('orgLevel'), or public health ('orgSector'). If the attributes are not specified, the default is version 1.1.0 'publicHealthDepartment' defined as 'The public health department or region where the site is located. See also 'healthReg' if there is a separate regional health

care delivery authority.' *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Secondary reporting authority

(repOrg2). The secondary, additional or alternative authority for routine surveillance reports or findings, or where the site is located. Use the 'organizations' table to describe attributes such as the name and address of the organization. Use 'orgType', 'orgLevel' and 'orgSector' to describe additional attributes. For example, whether the organization is a government agency ('orgType'), local ('orgLevel'), or public health ('orgSector'). If the attributes are not specified, the default is version 1.1.0 'healthRegion', The health planning authority where site is located. defined as 'The public health department or region where the site is located.' *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Reportable

(reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:*

Parts

naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 10 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Report date

(reportDate). The date a measure was reported. This date is the first date the laboratory or entity reported the measure findings to the public health or other entity who ordered the measure or is responsible for action. Use the analyses end date if the measure is part of a research study or other purpose that does not have a reporting entity. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Reporter ID deprecated

(reporterIDDep). Unique identifier for the person or organization that is reporting the data. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.1.0. Last updated: 2.0.0.

Reporter table deprecated

(reportersDep). The table that contains information about a reporter of a sample, method, measure, or attribute. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Replicate Type

(repType). Attribute of a sample, specifying whether the sample is unique, or a replicate. And if it is a replicate, what type of replicate is it? *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* replicateSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Resuspend COSCa filter collection

(resCosca). Nucleic acid extraction via resuspension of a sample collected using the COSCa ball method. See `partID = cosca` for additional details. *Part Type:* categories *Domain:* bio *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:*

Parts

NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 10 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Results tables

(results). Tables used to record samples, measures and quality reports. These tables are updated whenever a new sample is taken or a measure is reported. Quality reports are used to report quality assurance and quality control for samples and measures. *Part Type*: classes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: results *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Retention pond

(retPond). Retention pond *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

River, natural water body

(river). River, natural water body *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Role of contact

(role). Specifies the organizational role of a given contact. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Rate per 100,000

(rP100). Units for describing a population measure of a rate of case incidence per 100,000 people of a given disease. Unit for a population or disease measure to explain population-level clinical effects. *Part Type*: units *Domain*: bio *Speciment ID*: poSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: 0 *Maximum Value*: NA

Parts

Minimum Length: NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Raw post-pasteurized wastewater

(rppw). Raw wastewater sample post-pasteurization *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Respiratory syncytial virus

(rsv). Respiratory syncytial virus. *Part Type:* measurements *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* virusMisc *Class:* allele *Nomenclature:* naNomenclature *Ontology Reference:* http://purl.obolibrary.org/obo/NCIT_C14267 *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

s2083i omicron-variant gene target { #s2083i }

(s2083i). s2083i omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA

Reference Link: NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal
Maximum Value: seeUnitVal *Minimum Length:* NA *Maximum Length:* NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

s371l omicron-variant gene target { #s371l }

(s371l). s371l omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

s373p omicron-variant gene target { #s373p }

(s373p). s373p omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

s375f omicron-variant gene target { #s375f }

(s375f). s375f omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet

Parts

UnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: seeUnitData Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

s477n omicron-variant gene target { #s477n }

(s477n). s477n omicron-variant gene target Part Type: measurements Domain: bio Speciment ID: saSpecimenSet Compartment Set: anyCompartmentSet Group: sarsCov2 Class: mutation Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: genetic-UnitSet Aggregation Scale: seeUnitAggScale Quality Set: NA Reference Link: NA Data Type: seeUnitData Minimum Value: seeUnitVal Maximum Value: seeUnitVal Minimum Length: NA Maximum Length: NA Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample

(sa). A measure made on a compartment or property from a sample of a substance. Part Type: specimens Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample material

(saMaterial). Type of sample. Part Type: attributes Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: colGrp Class: naClass Nomenclature: naNomenclature Ontology

Reference: NA *Category Set ID:* sampleMatSet *Unit Set:* naUnitSet
Aggregation Scale: NA *Quality Set:* NA *Reference Link:* NA *Data Type:*
 categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*
 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last
 updated: 2.0.0.

Sample ID

(sampleID). Unique identifier for a sample. Suggestion:siteID-date-index.
Part Type: attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet
Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set*
ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum*
Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
 First released: 1.1.0. Last updated: 2.0.0.

Sample ID object

(sampleIDObject). Populated by ‘sampleID’ - this specifies the object
 of a sample relationship. Use ‘sampleID’ for this field. Together with
 ‘sampleIDSub’ and ‘relationshipID’ the row of the ‘sampleRelationships’
 should say: [sampleIDSubject] is a[relationshipID] of [sampleIDObject]
Part Type: attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet
Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set*
ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* varchar *Minimum Value:* NA *Maximum*
Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
 First released: 2.0.0. Last updated: 2.0.0.

Sample ID subject

(sampleIDSubject). Populated by 'sampleID' - this specifies the subject of a sample relationship. Use 'sampleID' for this field. Together with 'sampleIDSub' and 'relationshipID' the row of the 'sampleRelationships' should say: [sampleIDSubject] is a[relationshipID] of [sampleIDObject]
Part Type: attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample material set

(sampleMatSet). set for the types of material that can be found in a sample.
Part Type: mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sampleMatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample quality set

(sampleQualitySet). Quality set for a sample. *Part Type:* qualityIndSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:*

NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA
Minimum Length: 0 *Maximum Length*: 30 *Part Status*: active. First
 released: 2.0.0. Last updated: 2.0.0.

Sample relationships table

(sampleRelationships). Table for recording the relationships between samples. Samples can be pooled or split. The sample relationships table holds information on parent-child relationships between samples, and allow for tracking sample lineage for single and pooled samples. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: results *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationships table column order

(sampleRelationshipsOrder). Specifies the order of the columns in a Sample Relationships table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationships table required headers

(sampleRelationshipsRequired). Specifies the columns required in a Sample Relationships table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationship

(sampleRelID). Attribute for specifying the relationship between a sample subject and object. Together with 'sampleIDSubject' and 'sampleRelID', a row of the 'sampleRelationships' should say: [sampleIDSubject] is a [sampleID] of [sampleIDObject] *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sampleRelSet, protocolRelSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationships set

(sampleRelSet). set for valid values of relationshipID in the sampleRelationships table. Use only for sampleRelationships table relationshipID. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* sampleRelSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample report table

(samples). The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. Samples can be combined, split, stored and reused. The ‘Sample relationships’ (sampleRelationships) is used to describe the relationships between samples. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* results *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sample shed

(sampleShed). A geographic area, physical space, or structure. A sample is taken from a sampleshed for a representative measurement of a substance(s). Examples: Airplane, Long-term care home, Neighbourhood, University campus, Municipality *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* shedSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:*

Parts

0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Samples table column order

(samplesOrder). Specifies the order of the columns in a Samples table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sample table required headers

(samplesRequired). Specifies the columns required in a Samples table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Collection other deprecated

(sampleTypeOther). Description for other type of method not listed in collection. deprecated. *Part Type*: methods *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA

Quality Set: NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2

(sarsCov2). A group of measures/methods related to the SARS-CoV-2 virus. SARS-CoV-2 measures have the following classes: proteins, alleles, variants, mutations, diseases. *Part Type:* groups *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* http://purl.obolibrary.org/obo/NCIT_C169076 *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample specimen set

(saSpecimenSet). A specimen set that inculdes only a sample specimen. *Part Type:* specimenSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sparse coverage

(sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. If

Parts

there are multiple issues to flag, use the multiple issues flag and explain the details in the notes section. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier effluent depreciated

(sce). Effluent obtained after secondary clarifiers. Deprecated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Standard curve frequency

(scf). A method for specifying the frequency over which a standard curve (or 'Master Curve') is used. *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* timeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:*

NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

School

(school). A school serving students in the kindergarten to 12th grade range
Part Type: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet
Compartment Set: anyCompartmentSet *Group*: colGrp *Class*: naClass
Nomenclature: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA
Unit Set: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA
Reference Link: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA
Minimum Length: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sequencing script version set

(scriptSet). The set for script versions *Part Type*: mmaSets *Domain*: bio
Speciment ID: saSpecimenSet *Compartment Set*: anyCompartmentSet
Group: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA
Category Set ID: scriptSet *Unit Set*: naUnitSet *Aggregation Scale*: NA
Quality Set: NA *Reference Link*: NA *Data Type*: seeUnitData
Minimum Value: seeUnitVal *Maximum Value*: seeUnitVal
Minimum Length: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Genetic sequencing script version

(scriptVersion). Used to specify the script or script version used to extract summary information (i.e. coverage and mutation frequency statistics, etc.) from the sequencing output. *Part Type*: methods *Domain*: bio
Speciment ID: saSpecimenSet *Compartment Set*: anyCompartmentSet

Parts

Group: procGrp *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* scriptSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier sludge deprecated

(scsDep). Sludge produced by secondary clarifiers. Deprecated in version 2, please do not use. Duplicate of the sampleMat categories - refer to those. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Standard deviation

(sd). Standard deviation. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Normalized standard deviation

(sdn). Standard deviation, normalized. *Part Type*: aggregations *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sea, natural water body

(sea). Sea, natural water body *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

See unit (aggregation scales)

(seeUnitAggScale). A value for aggregation scale to be used when the aggregation scale depends on that specified in the parts list entry for the unit. When you encounter this value, see the unit being used for that entry and follow that aggregation scale entry. *Part Type*: aggregation-Scales *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA

Parts

Data Type: varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

See unit (data type)

(seeUnitData). The data type entry when the data type for a given entry is specified by the unit used. Only used for measures, which is where the data type entry is used to validate the entry in the value field. For these cases, however, the data type is more dependant on the unit than the measure , and so this entry directs users to consult the data type of the unit. *Part Type:* dataTypes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

See unit (value)

(seeUnitVal). A value for maximum or minimum value in the dictionary, use when the maximum or minimum value is determined by the unit. Used for measures when the constraints on maximum and minimum value are specified by the unit, not the measure. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier effluent

(sEfflu). Effluent obtained after secondary clarifiers *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Access to self

(self). If this is 'no', this data will not be shown on the portal when this reporter logs in. If missing, data will be available to this reporter. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Date sample was sent

(sentDate). The date the sample was sent for analyses at a laboratory. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* datetime *Minimum Value:* NA *Maximum*

Parts

Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active.
First released: 2.0.0. Last updated: 2.0.0.

Septic tank wastewater

(septage). Wastewater from within a septic tank *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Septic tank

(septTnk). Septic tank *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sequencing Layout

(seqLay). The layout of genetic material used in sequencing. *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* seqLaySet

Unit Set: unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Layout set

(seqLaySet). The set for the layout of genetic material used in sequencing *Part Type:* mmaSets *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* seqLaySet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing quality set

(seqQualitySet). Quality set for sequencing measures *Part Type:* qualityIndSets *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method

(seqSel). The primer sequence selection method used in sequencing. *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* sequence *Nomenclature:*

Parts

naNomenclature *Ontology Reference*: NA *Category Set ID*: seqSelSet *Unit Set*: unitlessUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method set

(seqSelSet). The set for the selection method used in sequencing *Part Type*: mmaSets *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: seqSelSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy

(seqStrat). The sequencing strategy used for an analysis. *Part Type*: methods *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: sequence *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: seqStratSet *Unit Set*: unitlessUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy set

(seqStratSet). The set for the sequencing strategy used for an analysis *Part Type*: mmaSets *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment*

Set: anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* seqStratSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Genetic sequences class

(sequence). Measures and methods related to genetic sequencing. *Part Type:* classes *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* sequence *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Serial dilution

(serialDilution). The serial dilution method for assessing inhibition. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* dilution *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Sewer Sediments Depreciated

(seseDep). Sediments obtained in sewer. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Set ID

(setID). The unique identifier of a value set. A set is a group of units, attributes, categories, specimens, or compartments that can be used for a measure, method or attribute. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 28 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sets look-up table

(sets). Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together. This table manages many-to-many relationships and category recycling within the ODM. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* lookup *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit*

Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sets table column order

(setsOrder). Specifies the order of the columns in the Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 1 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sets table required headers

(setsRequired). Required headers in the Sets table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Settled solids

(settsol). Amount of settled solids from a wastewater sample. May be a volume or a weight, see the unitID for a given measure row for details on a given sample. *Part Type:* measurements *Domain:* phy *Speciment ID:*

Parts

saSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: proc-Grp *Class*: physical *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: volumeUnitSet *Aggregation Scale*: see-UnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Set type

(setType). The type of set. i.e. quality set, aggregation set, unit set *Part Type*: partSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Set Value

(setValue). The partID for the set value or category. Populated by any part where partType = aggregationSet, compartmentSet, qualitySet, specimenSet, or unitSet *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Severity indicator

(severity). An indicator of the severity or seriousness of a quality flag. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sevSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Severity set

(sevSet). A set for severity indicators. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* sevSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sewer network file link depreciated

(sewerNetworkFileBLOB). Link to a file that has any detailed information about the sewer network associated with the site (any format). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:*

Parts

NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**.
First released: 1.0.0. Last updated: 2.0.0.

S sars-cov-2 gene target

(sGene). S sars-cov-2 gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: allele *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active.
First released: 2.0.0. Last updated: 2.0.0.

Sampled set

(shedSet). The set for all valid values of sampled. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: siSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Ship

(ship). A cruise ship or other ship. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*:

varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Shipped on ice

(shipOnIce). Was the sample kept cool while being shipped to the lab *Part Type:* methods *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* booleanSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Short names

(shortName). Shortened names for tables and other important parts for use in wide names. *Part Type:* partSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Site

(si). A measure made on a compartment or property at a site. An example is the wastewater flow rate, taken at a wasteawater treatment plant. *Part Type:* specimens *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass

Parts

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Single

(sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. *Part Type:* aggregations *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Single Layout

(sinLay). Specifies the single layout for a sequencing method *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Site or sample specimen set

(siSaSpecimenSet). A specimen set that includes site or sample specimen. *Part Type:* specimenSets *Domain:* naDomain *Speciment ID:* siSaSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Site specimen set

(siSpecimenSet). A specimen set that includes only a site specimen. *Part Type:* specimenSets *Domain:* naDomain *Speciment ID:* siSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Site ID default

(siteDef). Used as default when a new sample is created by this reporter. See ID in Site table. Deprecated in version 2, please do not use. Just populate standard SiteID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:*

Parts

NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Site features group

(siteFeat). A group of environmental measures/methods and those pertaining to the features of a site. *Part Type*: groups *Domain*: allDo *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Site ID

(siteID). Unique identifier for the location where a sample was taken. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Site measure ID

(siteMeasureID). Unique identifier for wide table only. Use when all measures are performed on a single sample. Depreciate in version 2. *siteMeasures* *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*:

naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Sites table

(sites). The table that contains information about a site; the location where an environmental sample was taken. The site of an environmental sample. Information in the site table does not regularly change. Consider using the MeasureReport table if the information changes often. *Part Type*: tables *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: programDescr *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sites table column order

(sitesOrder). Specifies the order of the columns in a Sites table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Sites table required headers

(sitesRequired). Specifies the columns required in a Sites table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Site Type

(siteType). Type of site or institution where sample was taken. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* siteTypeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Site set

(siteTypeSet). set for the type of sampling site. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* siteTypeSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Skimmed milk flocculation

(skimMilkFloc). Skimmed milk flocculation *Part Type*: categories *Domain*: che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

slope

(slope). Slope value of the calibration curve. Used for storing calibration curve information, potentially also applicable for other curve types. *Part Type*: measurements *Domain*: allDo *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: standardCurve *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sewer network file link depreciated 2

(snfl). Link to a file that has any detailed information about the sewer network associated with the site (any format). Depreciate in version 2. This attribute has been changed to referenceLink. Use po_referenceLink to specify polygon link to sewer network file. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature

Parts

Ontology Reference: NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Snowy

(snowy). Qualitative category for the weather measure, specifying a snowy day. *Part Type:* categories *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Solid fraction

(sol). Solid fraction of a sample. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Solid seperation

(solidSep). Process used to separate solid and liquid phases of the sample. Solid seperation is used either prior to or in the absence of the concentra-

tion method specified in ‘methodConc’. *Part Type:* methods *Domain:* allDo *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* solidSeparationSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Solid separation set

(solidSeparationSet). set for the separation of solids. *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* solidSeparationSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Source Protocol ID

(sourceProtocol). A protocol that served as a basis for another protocol. This may be due to version changes, or a referencing a protocol that was later adapted. Populate with protocolID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step source ID

(sourceStep). Specifies the protocol step which serves as a basis for a given protocol step. Populate with a protocol step ID (stepID). *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity

(specHum). Measure for specific humidity, the unit is the ratio of the mass of water vapour and the mass of total air. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* humid *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* specHumidUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity unit set

(specHumidUnitSet). Unit set for specific humidity. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* specHumidUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA

Maximum Value: NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specifies

(specifies). Specifies that the object step or protocol occurs specifies details for the subject step or protocol in the overall protocol container. Use only for protocolRelationships table relationshipID. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specimen

(specimen). The substance or thing upon which the observation was made. Specimens include: site, sample, person and population. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 12 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Specimens

(specimens). Measures or observations are taken from three types of substances: populations (humans or geographic areas), samples, sites. See

Parts

partID = specimenID *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen sets

(specimenSets). Sets of specimens. Specimen sets are used when a measure or attribute can be recorded for more than one specimen. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 15 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Speed class

(speed). Measures and methods relating to speed. *Part Type*: classes *Domain*: phy *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: speed *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: windSpeedUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Spike matrix and recovery

(spike). The spike matrix and recovery method for assessing inhibition. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* pcr *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Spike material

(spikeMat). Material into which the recovery efficiency control target is spiked. *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeMatSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 4 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Spike material set

(spikeMatSet). set for spikeMat (aterial into which the recovery efficiency control target is spiked). *Part Type:* mmaSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeMatSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical

Parts

Minimum Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Recovery efficiency spike target

(spikeTarget). Method specifying the recovery efficiency control target. This matches on to the the spike material method ID. Should typically be indexed as a step proceeding the spikeMaterial methodID step. *Part Type:* methods *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeTargetSet *Unit Set:* unitlessUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Recovery efficiency spike target set

(spikeTargetSet). The set capturing containing all of the valid categories for the recovery efficiency control (or spike) target. *Part Type:* mmaSets *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* spikeTargetSet *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample spilled

(spill). Sample contents spilled from container. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier sludge

(sSludge). Sludge produced by secondary clarifiers *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Social services shelter

(sss). Other type of social services shelter *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Stabilization pond

(stabPnd). Specifies a site type which is a pond designed and built for wastewater treatment to reduce the organic content and remove pathogens from wastewater. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Standard concentrations class

(standardConc). Measures and methods relating to standard concentrations. *Part Type:* classes *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Standard curve class

(standardCurve). Measures and methods relating to generating or recording a standard curve. *Part Type:* classes *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* standardCurve *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:*

varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0
Maximum Length: 30 *Part Status*: active. First released: 2.0.0. Last
updated: 2.0.0.

Departments, states, or provinces

(stateProvLevel). Departments, states, or provinces *Part Type*: attributes
Domain: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: na-
CompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: ISCU
Ontology Reference: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Ag-*
gregation Scale: NA *Quality Set*: NA *Reference Link*: NA *Data Type*:
varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0
Maximum Length: 30 *Part Status*: active. First released: 2.0.0. Last
updated: 2.0.0.

State, Province, or Region

(stateProvReg). The state, province, or region where a site or organi-
zation is located; part of the address. *Part Type*: attributes *Domain*:
naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet
Group: naGroup *Class*: naClass *Nomenclature*: naNomenclature
Ontology Reference: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Ag-*
gregation Scale: NA *Quality Set*: NA *Reference Link*: NA *Data Type*:
varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0
Maximum Length: 30 *Part Status*: active. First released: 1.0.0. Last
updated: 2.0.0.

Status

(status). Whether the part is still active and can be used in the most
current ODM version. Values are ‘active’ or ‘inactive’. The status of

Parts

the part or list. ‘active’ means the part is used in the dictionary version. ‘inactive’ means the part has been retired and not longer used in the version. *Part Type*: partSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: statusSet *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Status set

(statusSet). A set for partID = Status to indicate whether a part is in current use or not. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Standard concentration unit set

(stdConcentrationUnitSet). Unit set for concentration measurements. *Part Type*: unitSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Storage temp

(sTemp). Temperature that the sample is stored at in Celsius. *Part Type:* measurements *Domain:* phy *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* temperature *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Protocol Step ID

(stepID). The unique identifier for a specific protocol step. Protocol Steps are the component parts of a larger protocols, the latter being linked to specific measure reports. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Step ID Object

(stepIDObj). The object of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Populate this field using Protocol Step ID (partID = stepID). Together with 'protocolIDSub' or 'stepIDSub' and 'relationshipID' the row of 'protocolRelationships' should say: [protocolIDSubject or stepIDSubject] [relationshipID] [protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type:* attributes *Domain:* naDomain *Specimen ID:*

Parts

naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup
Class: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA
Category Set ID: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA
Reference Link: NA *Data Type*: integer *Minimum Value*: 1
Maximum Value: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active.
First released: 2.0.0. Last updated: 2.0.0.

Step ID Subject

(stepIDSub). The subject of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Populate this field using Protocol ID (partID = stepID). Together with ‘protocolIDSub’ or ‘stepIDSub’ and ‘relationshipID’ the row of ‘protocolRelationships’ should say: [protocolIDSubject or stepIDSubject] [relationshipID] [protocolIDObject or stepIDObject]. Ex: stepIDSubject precedes protocolIDObject. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup
Class: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA
Category Set ID: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA
Reference Link: NA *Data Type*: varchar *Minimum Value*: NA
Maximum Value: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active.
First released: 2.0.0. Last updated: 2.0.0.

Method step parent ID

(stepProvenanceID). A method step that served as a basis for another method step. This may be due to version changes, or a referencing a method set that was later adapted. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup
Class: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA
Category Set ID: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA
Reference Link: NA *Data Type*:

varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step version

(stepVer). Specifies the version of a given protocol step. Version of the method step. Semantic versioning is recommended. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 50 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Sewer network file blob

(storTempDef). A file blob that has any detailed information about the sewer network associated with the site (any format). Deprecated in version 2, please do not use. Refer to WKT and EPSG instead. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* blob *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 65535 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

Storage time

(stoTim). Length of time that a sample was in storage. *Part Type*: measurements *Domain*: allDo *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: time *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: timeUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Sub-variant or lineage

(subVar). A unit used to report a specific genetic lineage, or to specify that a reported variant is a sub-variant of another. Combine in a measure report if necessary. *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: variant *Nomenclature*: pangolin *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Summary

(summ). Short description of the assay and how it is different from the other assay methods. *Part Type*: attributes *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum*

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 255 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Summary worksheet

(summary). Summary sheet of the dictionary.xlsx file. The summary sheet contains versioning and changelog. *Part Type:* dictionarySupport *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sunny

(sunny). Qualitative category for the weather measure, specifying a clear and sunny day. *Part Type:* categories *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment

(surf). A measure or observation made from a substance on a surface. *Part Type:* compartments *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category*

Parts

Set ID: NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Other surface

(surface). Surface other than floor or desk. See also categories for floor and desk. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment set

(surfaceCompartmentSet). A compartment set for measures and methods in the surface compartment. *Part Type:* compartmentSets *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Surface and water compartment set

(surfaceWaterCompartmentSet). A compartment set for measures and methods in the surface or water compartments. *Part Type:* compart-

mentSets *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* surfaceWaterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Surface swab

(surfSw). Surface swab. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* surfaceCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sewage truck

(swgTrck). Sewage truck *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Sewer sediment

(swrSed). Sediments obtained in sewer *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer catchment area

(swrSet). Sewer catchment area. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Synthetic sample

(synthetic). Specifies a synthetic, or artificially derived sample. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

t19r delta-variant gene target

(t19r). t19r delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

t3646a delta-variant gene target

(t3646a). t3646a delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

t547k omicron-variant gene target { #t547k }

(t547k). t547k omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

t9i omicron-variant gene target { #t9i }

(t9i). t9i omicron-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Tables

(tables). Tables are where measures, methods and attributes are recorded. Tables represent the main entities of the environmental surveillance. There An example of a table is ‘site’ a collection of attributes such as site name and address to describe where environment samples are taken. PHES-ODM has report tables, but users can create their own tables as well. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Table column set

(tableSet). The set for valid inputs in table columns. These categories are used in the parts list for columns named after tables to indicate the position and function of a given part in the entity relationship diagram of the data model. *Part Type*: mmaSets *Domain*: naDomain *Speciment ID*:

naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup
Class: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA
Category Set ID: tableSet *Unit Set*: naUnitSet *Aggregation Scale*: NA
Quality Set: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

TaqPath N sars-cov-2 gene target

(taqpatN). TaqPath N sars-cov-2 gene target *Part Type*: measurements
Domain: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet
Group: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA
Category Set ID: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale
Quality Set: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal
Maximum Value: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

TaqPath S sars-cov-2 gene target

(taqpatS). TaqPath N sars-cov-2 gene target *Part Type*: measurements
Domain: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet
Group: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA
Category Set ID: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale
Quality Set: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal
Maximum Value: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA
Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Sample temperature

(temp). Temperature of the sample measured in degrees Celcius *Part Type*: measurements *Domain*: phy *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: miscMeas *Class*: temperature *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: temperatureUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: float *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Temperature class

(temperature). Measures and methods related to temperatures. *Part Type*: classes *Domain*: phy *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: temperature *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: temperatureUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Temperature unit set

(temperatureUnitSet). Unit set for temperature measurements. *Part Type*: unitSets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: temperatureUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA

Reference Link: NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Nucleic acid template volume

(tempVol). The volume of DNA or RNA template used for PCR. *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Airport terminal

(terminal). Airport terminal sample shed category type *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Test date

(tesDate). Date that the covid-19 test was performed. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:*

Parts

naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: NA *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Number of tests performed

(test). Number of tests performed. *Part Type*: units *Domain*: bio *Speciment ID*: poSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Testing

(testing). A measure or sample taken to test a new method. These measures are typically used for internal lab uses and not reported to external partners. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

thermo magmax microbiome ultra nucleic acid isolation kit

(thermMag). Nucleic acid extraction performed using the thermo magmax microbiome ultra nucleic acid isolation kit. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet

Group: procGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Time class

(time). Measures and methods relating to time. *Part Type:* classes *Domain:* allDo *Speciment ID:* siSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* time *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* timeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Time proportional sample

(timePr). A time proportionalcomposite sample generally collected by an autosampler. Use collectionPeriod to describe how many hours the sample was taken. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Time unit set

(timeUnitSet). The unit set for measures of time. *Part Type:* unitSets *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* time *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* timeUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Total Kjeldahl Nitrogen

(tkn). A measure of Total Kjeldahl Nitrogen; ie. the sum of nitrogen bound in organic substances, nitrogen in ammonia (NH₃-N) and in ammonium (NH₄⁺-N) in the chemical analysis of soil, water, or waste water. *Part Type:* measurements *Domain:* che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Total Nitrogen

(tn). Total nitrogen concentration, as N. *Part Type:* measurements *Domain:* che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality*

Set: NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Time Proportional 24hr sample

(tp24s). A time proportional 24-hour composite sample generally collected by an autosampler. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* col-Grp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **deprecated**. First released: 1.0.0. Last updated: 2.0.0.

Translation look-up table

(translations). Look up table for translations of the description, label, and instruction for all parts. The default language if a translation is not specified is English. *Part Type:* tables *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* lookup *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Translation table column order

(translationsOrder). Specifies the order of the columns in the Translation table. *Part Type:* tableSupport *Domain:* naDomain *Speciment ID:*

Parts

naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 1 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Translation table required headers

(translationsRequired). Required headers in the Translations table. *Part Type*: tableSupport *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

trizol, zymo mag beads w/ zymo clean and concentrator

(trizol). Nucleic acid extraction performed using the trizol, zymo mag beads w/ zymo clean and concentrator method. *Part Type*: categories *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

TRUE

(true). Boolean data type = TRUE Use only “TRUE” (case sensitive)
Part Type: categories *Domain:* naDomain *Speciment ID:* naSpecimenSet
Compartment Set: naCompartmentSet *Group:* naGroup *Class:* naClass
Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA
Unit Set: naUnitSet *Aggregation Scale:* NA *Quality Set:* NA
Reference Link: NA *Data Type:* boolean *Minimum Value:* NA *Maximum Value:* NA
Minimum Length: 4 *Maximum Length:* 4 *Part Status:* active.
 First released: 1.1.0. Last updated: 2.0.0.

Total solids concentration

(ts). Total solids concentration *Part Type:* measurements *Domain:* che
Speciment ID: siSaSpecimenSet *Compartment Set:* waterCompartmentSet
Group: miscMeas *Class:* standardConc *Nomenclature:* naNomenclature
Ontology Reference: NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet
Aggregation Scale: seeUnitAggScale *Quality Set:* NA *Reference Link:* NA
Data Type: seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal
Minimum Length: NA *Maximum Length:* NA
Part Status: active. First released: 1.0.0. Last updated: 2.0.0.

Total suspended solids

(tss). Total suspended solids *Part Type:* measurements *Domain:* che
Speciment ID: siSaSpecimenSet *Compartment Set:* waterCompartmentSet
Group: miscMeas *Class:* standardConc *Nomenclature:* naNomenclature
Ontology Reference: NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet
Aggregation Scale: seeUnitAggScale *Quality Set:* NA *Reference Link:* NA
Data Type: seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal
Minimum Length: NA *Maximum Length:* NA
Part Status: active. First released: 1.0.0. Last updated: 2.0.0.

Parts

Concentration of total suspended solids

(tssConc). Total suspended solids concentration of the wastewater. *Part Type*: measurements *Domain*: phy *Specimen ID*: siSaSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: standard-Conc *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: stdConcentrationUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Turbidity

(turb). A measure for the turbidity of water, or a liquid sample, quantifying the relative clarity or opacity of a liquid. Highly indicative of water quality. *Part Type*: measurements *Domain*: phy *Specimen ID*: siSaSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: miscMeas *Class*: turbidity *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: turbidityUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Turbidity class

(turbidity). Measures and methods relating to water and wastewater turbidity. *Part Type*: classes *Domain*: che *Specimen ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: turbidity *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*:

NA *Unit Set*: turbidityUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Turbidity unit set

(turbidityUnitSet). Unit set for turbidity measurements. *Part Type*: unit-Sets *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: turbidityUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Tube weight empty

(tWeighE). The weight of the tube used for analysis while empty. *Part Type*: measurements *Domain*: phy *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: physical *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: volumeUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: integer *Minimum Value*: 0 *Maximum Value*: NA *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Tube weight full

(tWeighF). The weight of the tube used for analysis while full of analyte. *Part Type*: measurements *Domain*: phy *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: procGrp *Class*: physical

Parts

Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Type other depreciate

(tyOtDep). Description for other type of sample not listed in *Part Type*: methods *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Type other depreciate2

(tyOtDep2). Description of the site when the site is not listed. See site-Type. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 100 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

University campus

(uCampus). Universityor college campus - comprising an entire campus or part of a campus. See also 'Higher Education Domitory'. *Part Type:* cat-

egories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.1.0. Last updated: 2.0.0.

Undisclosed

(undisc). A value has been recorded, however it was not included in the dataset. Missing value indicator for data that has not been shared or disclosed with outside parties. *Part Type:* missingness *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Unique sample

(unique). A unique sample, not duplicated. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Unit

(unit). The units of a measurement. Different units that are used to describe measurement values. Units are combined into UnitSetIDs. Please requests new units by creating an issue in the ODM repository. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Unitless measure

(unitless). A unit for unitless measures, like pH. *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Unitless unit set

(unitlessUnitSet). Unit set for measurements that does not have units. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* unitlessUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:*

NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Units

(units). The unit of the measurement. Every measurement must have a unit. Meaning, ‘units’ are a mandatory field whenever a measurement is recorded. *Part Type*: partType *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Unit set

(unitSet). An identification of a set of units that can be used for a measure or method. *Part Type*: attributes *Domain*: che *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: procGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Unit sets

(unitSets). Sets of units. Contains units that are associated with part types such as measures. Sets are lists of input values. Each input value within a set is a part (has a partID) that can be reused between sets. For example, a set A could contain the values of ‘yes’ and ‘no’; and set B

Parts

could contain the values of ‘yes’, ‘no’, and ‘maybe.’ partTypes with sets include aggregations, compartments, missingness, quality flags, specimens, and units. Measures with categories and dataTypes have their own specific sets (catSetID) because their input values are unique and are not reused. *Part Type:* partType *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* categorical *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Ultracentrifugation

(untracent). Ultracentrifugation Should typically be linked to a lower limit, and some indication of data spread (ie. standard deviation), and a mean/median measure. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval

(upperCI95). Specifies the the upper limit of a 95% confidence interval. A catch-all category for upstream wastewater sampling sites, including major sewer pipelines and pumping stations. *Part Type:* aggregations *Domain:*

naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: colGrp *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.0.0. Last updated: 2.0.0.

Upstream sites

(upstream). A general site type for upstream wastewater sampling sites. *Part Type*: categories *Domain*: phy *Speciment ID*: siSaSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: quantAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 1.1.0. Last updated: 2.0.0.

Micro-Siemens per centimetre

(uSCm). Micro-siemens per centimetre. Deprecated in version 2, please do not use. MeasureID will be associated with site, with site specified as the specimen as well. *Part Type*: units *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: waterCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 12 *Part Status*: **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Parts

U site measure ID

(uSiteMeasureID). Unique identifier for each measurement for a site. *Part Type*: attributes *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

v2930l delta-variant gene target

(v2930l). v2930l delta-variant gene target *Part Type*: measurements *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: sarsCov2 *Class*: mutation *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: geneticUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Validation study

(validationStudy). A measure or sample taken for a validation study. *Part Type*: categories *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Value

(value). Value of a measure, observation or attribute. Only used for the dictionary entries of parts. *Part Type:* attributes *Domain:* naDomain *Specimen ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Variable character data type

(varchar). The data type for variable character data. Measured by sequencing. *Part Type:* dataTypes *Domain:* bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* 100 *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Frequency of variants detected

(varFreq). A description of the frequency of a variants detected. *Part Type:* units *Domain:* bio *Specimen ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* variant *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

Variants class

(variant). Measures and methods relating to variants. *Part Type:* classes *Domain:* bio *Speciment ID:* poSaSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* variant *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population with 1 dose of vaccine

(vax1). A measure of the population with a single dose of vaccine. *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population with 2 doses of vaccine

(vax2). Population with 2 doses of vaccine *Part Type:* units *Domain:* bio *Speciment ID:* poSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Population with 3 doses of vaccine

(vax3). Population with 3 doses of vaccine *Part Type:* units *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ODM V1 category

(version1Category). ODM V1 category *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

ODM V1 location

(version1Location). ODM V1 location *Part Type:* crosswalkTableSet *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Parts

ODM V1 table

(version1Table). ODM V1 table *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ODM version 1 to 2 changes

(version1to2Changes). Changes for part between ODM v1 and V2 *Part Type*: crosswalkTableSet *Domain*: naDomain *Speciment ID*: naSpecimenSet *Compartment Set*: naCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

ODM V1 variable

(version1Variable). ODM V1 variable *Part Type*: crosswalkTableSet *Domain*: bio *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: virusMisc *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous viruses group

(virusMisc). A group of measures/methods related to miscellaneous viruses. The miscallenous viruses are often measured for normalization or standardization and have only one measure or method. When a virus has many measures, they will be given their own groupID in updated dictionary versions. *Part Type:* groups *Domain:* phy *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* measGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Size in volume

(vol). Total volume of water or sludge sampled. Used for water or air testing *Part Type:* measurements *Domain:* phy *Speciment ID:* saSpecimenSet *Compartment Set:* airWaterCompartmentSet *Group:* colGrp *Class:* physical *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* volumeUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Volume proportional sample

(volPr). A volume proportional sample generally collected by an autosampler. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:*

Parts

naClass Nomenclature: naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Volume unit set

(volumeUnitSet). Unit set of volume measurements. *Part Type:* unitSets *Domain:* che *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Volatile suspended solids

(vss). Volatile suspended solids *Part Type:* measurements *Domain:* phy *Speciment ID:* saSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* measGrp *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Volatile suspended solids - ignition

(vssIg). A water quality measure, captured by via the loss on ignition of the mass of measured total suspended solids. This ignition generally takes place in an oven at a temperature of 550 °C to 600 °C. It represents the amount of volatile matter present in the solid fraction of the measured solution. *Part Type:* measurements *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* waterCompartmentSet *Group:* meas-Grp *Class:* standardConc *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* stdConcentrationUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 1 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment

(wat). A measure or observation made from a substance in the water, including wastewater. *Part Type:* compartments *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* NA *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Water depreciated

(water). Non-wastewater, coming from any kind of water body. *Part Type:* categories *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment*

Parts

Set: waterCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 15 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment set

(waterCompartmentSet). A compartment set for measures and methods in the water compartment. *Part Type:* compartmentSets *Domain:* phy *Speciment ID:* siSaSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* miscMeas *Class:* temperature *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Wastewater temperature

(watTemp). Temperature of the wastewater. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* weather *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* weathSet *Unit Set:* temperatureUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* float *Minimum Value:* 0 *Maximum Value:* 40 *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Weather

(weath). A measure for weather captured through qualitative categories; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources. *Part Type*: measurements *Domain*: phy *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: siteFeat *Class*: weather *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: weathSet *Unit Set*: unitlessUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: categorical *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Weather class

(weather). Measures and methods relating to weather. *Part Type*: classes *Domain*: allDo *Speciment ID*: siSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: weather *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: unitlessUnitSet *Aggregation Scale*: seeUnitAggScale *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Weather set

(weathSet). A set of the valid qualitative categories for the qualitative weather measure. *Part Type*: mmaSets *Domain*: bio *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*:

Parts

NA *Reference Link*: NA *Data Type*: seeUnitData *Minimum Value*: seeUnitVal *Maximum Value*: seeUnitVal *Minimum Length*: NA *Maximum Length*: NA *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Whole genome sequencing

(wgs). Specifies the whole genome sequencing (WGS) strategy for genetic sequencing *Part Type*: categories *Domain*: bio *Speciment ID*: naSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

WHO nomenclature

(who). Specifies variant or genetic nomenclature as set out by the World Health Organization. If there are multiple issues to flage, use the multiple issues flag and explain the details in the notes section. *Part Type*: nomenclatures *Domain*: naDomain *Speciment ID*: saSpecimenSet *Compartment Set*: anyCompartmentSet *Group*: naGroup *Class*: naClass *Nomenclature*: naNomenclature *Ontology Reference*: NA *Category Set ID*: NA *Unit Set*: naUnitSet *Aggregation Scale*: NA *Quality Set*: NA *Reference Link*: NA *Data Type*: varchar *Minimum Value*: NA *Maximum Value*: NA *Minimum Length*: 0 *Maximum Length*: 30 *Part Status*: active. First released: 2.0.0. Last updated: 2.0.0.

Wide 95% interval

(wI). The 95% interval is too wide, low confidence in results. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Wide name table

(wideNames). The table for wide names. *Part Type:* tables *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* siteFeat *Class:* speed *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* windSpeedUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* seeUnitData *Minimum Value:* seeUnitVal *Maximum Value:* seeUnitVal *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Wind Speed

(wind). A measure for wind speed; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* airCompartmentSet *Group:* siteFeat *Class:* speed *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* windSpeedUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum*

Parts

Value: NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Wind speed unit set

(windSpeedUnitSet). Unit set for wind speed measurements. *Part Type:* unitSets *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* naCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* **depreciated**. First released: 1.0.0. Last updated: 2.0.0.

Measure identification (wide table)

(wMeasureID). Unique identifier for wide table only. Use when all measures are performed on a single sample at the same time and same laboratory. Suggestion: siteID_sampleID_LabID_reportDate_ID. *Part Type:* attributes *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored incorrectly

(wrongStorage). Sample was stored inappropriately. *Part Type:* qualityIndicators *Domain:* naDomain *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* naGroup *Class:* naClass *Nomen-*

clature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored at wrong temperature

(wrongTemp). Sample was stored at an inappropriate temperature. *Part Type: qualityIndicators Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

World Sphere (W-SPHERE)

(wSphere). World Sphere header *Part Type: crosswalkTableSet Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: naCompartmentSet Group: naGroup Class: naClass Nomenclature: naNomenclature Ontology Reference: NA Category Set ID: NA Unit Set: naUnitSet Aggregation Scale: NA Quality Set: NA Reference Link: NA Data Type: varchar Minimum Value: NA Maximum Value: NA Minimum Length: 0 Maximum Length: 30 Part Status: active. First released: 2.0.0. Last updated: 2.0.0.*

World Sphere (W-SPHERE) notes

(wSphereNotes). W-Sphere notes *Part Type: crosswalkTableSet Domain: naDomain Speciment ID: naSpecimenSet Compartment Set: waterCom-*

Parts

partmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Wastewater treatment plant

(wwtp). A general site type for wastewater treatment plants. Used as a catch-all category for industrial or municipal wastewater treatment plants that carry combined or sanitary sewage, or lagoons. *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Back-up wastewater treatment site

(wwtpBack). Indicates a wastewater treatment plant where multiple labs are or were sampling, and where secondary samples are being taken taken. Maps to units of either million gallons per day, or population equivalents. *Part Type:* categories *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* siteFeat *Class:* capacity *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Wastewater treatment plant designed capacity

(wwtpCap). A measure for the designed capacity of a wastewater treatment plant, in terms of how much water it was designed to be able to hold and process. *Part Type:* measurements *Domain:* phy *Speciment ID:* siSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* capacity *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* capacityUnitSet *Aggregation Scale:* see-UnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Industrial wastewater treatment plant

(wwtpInd). Industrial wastewater treatment plant *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for combined sewage

(wwtpMuC). Municipal wastewater treatment plant for combined sewage *Part Type:* categories *Domain:* naDomain *Speciment ID:* naSpecimenSet *Compartment Set:* waterCompartmentSet *Group:* colGrp *Class:* naClass *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* NA *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA

Parts

Minimum Length: 0 *Maximum Length:* 30 *Part Status:* active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for sanitary dewage only

(wwtpMuS). Municipal wastewater treatment plant for sanitary sewage only *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* sarsCov2 *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

y505h omicron-variant gene target { #y505h }

(y505h). y505h omicron-variant gene target *Part Type:* measurements *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* geneticUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Years

(years). A unit for indicating a length of time in years. *Part Type:* units *Domain:* allDo *Speciment ID:* anySpecimenSet *Compartment Set:* anyCompartmentSet *Group:* miscMeas *Class:* time *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* quantAggScale *Quality Set:* NA *Reference Link:* NA

Data Type: integer *Minimum Value:* 0 *Maximum Value:* NA *Minimum Length:* NA *Maximum Length:* NA *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

zymo environ water rna kit/ zymo environ water rna kit (cat. r2042)

(zymoEnv). Nucleic acid extraction performed using the zymo environ water rna kit/ zymo environ water rna kit (cat. r2042). *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

zymo quick-rna fungal/bacterial miniprep #r2014

(zymoQuick). Nucleic acid extraction performed using the zymo quick-rna fungal/bacterial miniprep #r2014 kit. *Part Type:* categories *Domain:* bio *Speciment ID:* saSpecimenSet *Compartment Set:* anyCompartmentSet *Group:* procGrp *Class:* mutation *Nomenclature:* naNomenclature *Ontology Reference:* NA *Category Set ID:* NA *Unit Set:* naUnitSet *Aggregation Scale:* seeUnitAggScale *Quality Set:* NA *Reference Link:* NA *Data Type:* varchar *Minimum Value:* NA *Maximum Value:* NA *Minimum Length:* 0 *Maximum Length:* 30 *Part Status:* active. First released: 2.0.0. Last updated: 2.0.0.

Sets

ODM version v2.0.0.

Sets are lists of categories and other parts. For example, `Collection category set` ([collectionCatSet](#)) contains a list of the different collection types including: `Moore swab`, `Grab sample`, `Flow proportional sample`, etc.

Sets can have parts that are reused. For example, there are different quality measures for sample collection (e.g. `Sample not collected`), PCR measures (`Inhibition preset`), and sequencing measures (`Low depth of coverage`), but all quality sets share the `No quality concerns` indicator.

Aggregations

(aggregations). Statistical measures used to report a measure. Each measure reported as a number should be reported with an aggregation. The default aggregation is `single`, which corresponds to a single measurement.

AND Boolean aggregation (`andBoo`). “AND” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (`floMean`). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Geometric mean (gmn). Geometric mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum value (maxVal). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arithmetic mean (me). Arithmetic mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Median (med). Median. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized arithmetic mean (menr). Arithmetic mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum value (minVal). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation not applicable (naAggr). Not applicable for aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Normalized geometric mean (ngmn). Geometric mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

OR Boolean aggregation (orBoo). “OR” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation deprecated (otherDep). Other aggregation method. Add description to aggregationOther Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean (pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-normalized mean (pmmovNorm). Mean measure normalized to amount of PMMoV. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard deviation (sd). Standard deviation. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized standard deviation (sdn). Standard deviation, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Specifies the the upper limit of a 95% confidence interval. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation scales

(aggregationScales). The scale of an aggregation set. Aggregation scales include quantitative and qualitative.

Aggregation scale not applicable (naAggrScale). Not application for aggregation sets. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation scale (othAgg). Specifies an aggregation scale that can't be described as either "qualitative" or "quantitative". Status: active. First released: 2.0.0. Last updated: 2.0.0.

Qualitative (qualAggScale). The qualitative aggregation scale. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Quantitative (quantAggScale). The “quantAggScale” aggregation scale. Status: active. First released: 2.0.0. Last updated: 2.0.0.

See unit (aggregation scales) (seeUnitAggScale). A value for aggregation scale to be used when the aggregation scale depends on that specified in the parts list entry for the unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation sets

(aggregationSets). Sets of aggregations. Examples of aggregation sets include logarithm, linear and boolean.

Boolean aggregation set

(booleanAggrSet)Aggreagation set for boolean. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

AND Boolean aggregation (andBoo). “AND” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

OR Boolean aggregation (orBoo). “OR” aggregation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Linear scale aggregation set

(linearAggrSet)The aggregation set that contains all aggregations that exist on the linear scale. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (floMean). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum value (maxVal). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arithmetic mean (me). Arithmetic mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Median (med). Median. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized arithmetic mean (menr). Arithmetic mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum value (minVal). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean (pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-normalized mean (pmmovNorm). Mean measure normalized to amount of PMMoV. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard deviation (sd). Standard deviation. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized standard deviation (sdn). Standard deviation, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Specifies the the upper limit of a 95% confidence interval. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Logarithmic scale aggregation set

(logAggrSet)The aggregation set that contains all aggregations that exist on the logarithmic scale (rather than linear or qualitative). Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (floMean). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Geometric mean (gm). Geometric mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum value (maxVal). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum value (minVal). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized geometric mean (ngm). Geometric mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

PMMoV- and flow-normalized mean (pmFloMean). Mean measure normalized to amount of PMMoV and wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-normalized mean (pmmovNorm). Mean measure normalized to amount of PMMoV. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Specifies the the upper limit of a 95% confidence interval. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation set not applicable

(naAggrSet)Not application for aggregation sets. Used for parts for which an aggregation set value is not applicable or appropriate. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Aggregation not applicable (naAggr). Not applicable for aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other aggregation set

(otherAggrSet)Aggregation set used for unitless measures. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Flow-normalized mean (floMean). Mean measure normalized to wastewater flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Geometric mean (gmn). Geometric mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Lower limit of a 95% confidence interval (lowerCI95). Specifies the the lower limit of a 95% confidence interval. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Maximum value (maxVal). Highest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arithmetic mean (me). Arithmetic mean. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Median (med). Median. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized arithmetic mean (menr). Arithmetic mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Minimum value (minVal). Lowest value in a range of values for a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Aggregation not applicable (naAggr). Not applicable for aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Normalized geometric mean (ngmn). Geometric mean, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Standard deviation (sd). Standard deviation. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Normalized standard deviation (sdn). Standard deviation, normalized. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Single (sin). A value that is not an aggregate measurement (ie. not a mean, median, max or any other) and can be a replicate value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upper limit of a 95% confidence interval (upperCI95). Specifies the the upper limit of a 95% confidence interval. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Attributes

(attributes). Attributes describe the who, where, when, and why of environmental surveillance. There are attributes for domain, specimen, compartment, or table. Attributes are one of the three main components or ways to report environment surveillance data. The other two components are measures and methods.

Analysis date (aDate). Date the measurement was performed in the lab. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date end (aDateEnd). Date the measurement or analysis was completed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Analysis date start (aDateStart). Date the measurement or analysis was started. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Address Line 1 (addL1). Line 1 (the street name, number and direction) for a given address. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Address Line 2 (addL2). Line 2 (the unit number) for a given address. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Address ID (addressID). A unique identifier for an address. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Administrative regions (admRegLevel). Administrative regions Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation scale (aggregationScale). A scale used for an aggregation. Only applicable for measures and units. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Aggregation (aggregation). Statistical measures used to report a measure. Each aggregation has a corresponding value. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Aggregation set (aggregationSet). The aggregation set for a unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Alias depreciated (aliasDep). ID of an assay that is the same or similar. A comma separated list. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Alias ID depreciated (aliasIDDep). Alias id Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Access to all org (allOrgs). If this is 'no', this data will not be available to any partner organization. If missing, data will be available to the all organizations. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Assay method id default depreciated (amDefDep). Used as default when a new measurement is created for this lab. See ID in AssayMethod table. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Assay ID (asID). Links with the AssayMethod used to perform the analysis. Use instrument.ID for measures that are not viral measures. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

City (city). The city where a site or organization is located; part of the address. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Class (class). A unique identifier for a class, which is akin to a sub-group; it's a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection date time (collDT). For grab samples this is the date, time and timezone the sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Collection date time end (collDTend). For integrated time average samples this is the date, time and timezone the sample was finished being taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Collection date time start (collDTStart). For integrated time averaged samples this is the date, time and timezone the sample was started being taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Collection number (collNum). The number of subsamples that were combined to create the sample. Use NA for continuous, proportional or passive sampling. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection number and period (collNumPer). Composite collection number.period. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Collection period (collPer). Collection period. The time period over which the sample was collected, in hours. Alternatively, use collectionStart and collectionEnd. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample collection type (collType). The type of collection. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Column name depreciated (columnNameDep). Name for the column Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Compartment (compartment). The attribute identifying the substance from which where a sample was taken. For more information, see partID = compartment. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Contact ID (contactID). A unique identifier for a given contact person. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Contact name (contactName). Contact person or group, for the lab. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Country (country). The country where a site or organization is located; part of the address. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Countries or sovereign states (countryLevel). Countries or sovereign states Status: active. First released: 2.0.0. Last updated: 2.0.0.

Districts, counties, regions (countyLevel). Districts, counties, regions Status: active. First released: 2.0.0. Last updated: 2.0.0.

Covid-19 population measurement date (cphDate). date of reporting for covid-19 measure. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Cphd ID (cphid). Unique identifier for the table. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Custody Contact ID (custodyCont). A unique identifier for data custodians. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Data custodian ID (custodyID). The data custodian of a database. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset creation date (datasetDate). Specifies the date a given dataset was created. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Date (date). Date Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Description (descr). A detailed description as an attribute for describing results or program elements. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to details deprecated (detailsDep). More details on the existing confidentiality requirements of this measurement. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Domain (domain). Domain is the highest level of describing of a measure. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Attributes

Contact email (email). Contact e-mail address, for the lab. Status: active. First released: 1.0.0. Last updated: 2.0.0.

File location of polygon (fileLocation). The location of the file containing the geometry of the polygon. Status: active. First released: 1.0.0. Last updated: 2.0.0.

First name of contact (firstName). Specifies the first name of a given contact. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Fraction analyzed (fraction). Fraction of the sample that is analyzed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Funder Contact ID (funderCont). A unique identifier for funders. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Funding agency ID (funderID). The funding agency of the dataset. Status: active. First released: 2.0.0. Last updated: 2.0.0.

European Petroleum Survey Group Coordinates (geoEPSG). The unique EPSG code specifying a given geospatial area. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Latitude (geoLat). Geographical location, latitude in decimal coordinates, ie.: (45.424721) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Longitude (geoLong). Geographical location, longitude in decimal coordinates, ie.: (-75.695000) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Type of geography (geoType). Type of geography that is represented by the polygon. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Well-known text (geoWKT). Well-known text of the polygon Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Group (group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Index (index). Index number in case the measurement was taken multiple times. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Instrument ID (instrumentID). A unique identifier for an instrument. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Instrument Type (insType). Type of instrument used to perform the measurement. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Describe other instrument type, if applicable (insTypeOth). Description of the instrument in case it is not listed in instrumentType. Status: active. First released: 1.0.0. Last updated: 2.0.0.

ISO639-1 (iso6391). The first part of the ISO 639 series of international standards for language codes. Part 1 covers the registration of two-letter codes. There are 183 two-letter codes registered as of June 2021. The registered codes cover the world's major languages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2B (iso6392B). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'B' specifies the bibliographic code (B code). Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-2T (iso6392T). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'T' specifies the terminological code (T code). Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-3 (iso6393). A set of international standards that lists short codes for language names. ISO 639-3 extends the ISO 639-2 alpha-3 codes with an aim to cover all known natural languages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ISO639-6 (iso6396). A set of international standards that lists short codes for language names. ISO 639-6 builds off ISO639-3 with the use of four-letter codes, and allowing users to differentiate between variants of languages and language families, such as historical vs. revived versions of languages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Lab ID default deprecated (labDefDep). Used as default when a new sample is created by this reporter. See ID in Lab table. Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

Lab ID (labID). Unique identifier for a laboratory. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Language family (langFam). Specifies the language family of a given language for translation and language tracking purposes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Language name (langName). Specifies the name of the language in roman alphabet characters for translation and language tracking purposes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Last edited (lastEdited). The date the entry was last updated. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Last name of contact (lastName). Specifies the last name of a given contact. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Link (li). Link to an external reference that describes the geometry of the polygon. Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

License (license). The license of a dataset. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to local ha (localHA). If this is ‘no’, the data will not be available to local health authorities. If missing, data will be available to local health authorities. Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

Manufacturer (manufacturer). Manufacturer of an instrument. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure (measure). A measurement or observation of any substance including a biological, physical or chemical substance. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Report ID (measureRepID). Unique identifier for a measurement. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Report set ID (measureSetRepID). Unique identifier that links together a group of related measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure license (measureLic). Specifies the access and use licensing for a given single measurement. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Method (method). A procedure for collecting a sample or performing a measure. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measure, method, or attribute (mmaSet). The set for a measure, method, or attribute. Only applicable for categorical fields, and the categories that populate them. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Model (model). Model number or version of the instrument. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Municipalities or communes (municipalLevel). Municipalities or communes Status: active. First released: 2.0.0. Last updated: 2.0.0.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Native Name (natName). The native name of the language, i.e. what the language is called by its speakers. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Local administrative units or neighborhoods (neighborLevel). Local administrative units or neighborhoods Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclature (nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Notes (notes). A note used to describe details that are not captured in other attributes Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Organization level (orgLevel). The geographic level of an organization. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization sector (orgSector). The sector of an organization Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization Type (orgType). Specifies the type or purpose of a given organization. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Sample origin (origin). An attribute of a sample specifying the origin. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to other prov deprecated (otherProvDep). If this is 'no', this data will not be available to other data providers not listed before. If missing, data will be available to other data providers not listed before. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Parent dataset ID (parDatasetID). The datasetID that is a parent to another datasetID. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Parent sample ID (parent). If this sample has been pooled into one big sample for analysis this indicates the sample of the larger pooled sample. Status: depreciated. First released: 1.1.0. Last updated: 2.0.0.

Parent Site ID (parSiteID). The siteID that is a parent to another siteID; usually the sub-site will be contained in the parentSite. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Postal or Zip Code (pCode). The zip code or postal code for a given address, specifying a specific geographic area. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to phac (phac). If this is 'no', the data will not be available to employees of the Public Health Agency of Canada - PHAC. If missing, data will be available to employees of the Public Health Agency of Canada - PHAC. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Contact phone (phone). Contact phone number, for the lab. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Polygon ID (polygonID). Unique identifier for the polygon. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Polygon Population (polyPop). An attribute of a polygon, which specifies the population of that polygon. A rough estimate of the number of

human residents. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Pooled (pooled). Is this a pooled sample, and therefore composed of multiple child samples obtained at different sites Status: active. First released: 1.0.0. Last updated: 2.0.0.

Type of date for case reporting depreciated (popDateTypeDep). Type of date used for confirmed cases. Typically, report or episode are reported. Onset and test date is not usually reported within aggregate data. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Population Served (popServ). An attribute of a site, which specifies the population of/population served by a given site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID (protocolID). A unique identifier for a given protocol. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID container (protocolIDContainer). Unique identifier for the protocol and the steps and other protocols that make it up. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID object (protocolIDObj). The object of the relationship between one protocol and another protocol or step, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol ID subject (protocolIDSub). The subject of the relationship between one protocol and another protocol or step, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol version (protocolVersion). Specifies the version of a method set. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to prov ha (provHA). If this is 'no', this data will not be available to provincial health authorities. If missing, data will be available to provincial health authorities. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Sets

Pretreatment description - depreciated (ptDescDep). If preTreatment then describe the treatment that was performed. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Access to public (public). If this is 'no', this data will not be available to the public. If missing, data will be available to the public. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Purpose (purpose). The reason the measure or sample was taken. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality flag (qualityFlag). A field for reporting any quality concerns - of lack thereof - for a sample or measure. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality report ID (qualityID). A unique identifier for a given quality report. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality set (qualitySetID). The quality set that corresponds to a given part. Only applicable for samples and measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Date sample recieved (recDate). The date the sample was received at the laboratory for analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Primary reporting authority (repOrg1). The primary or most responsible authority for routine surveillance reports or findings, or where the site is located. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Secondary reporting authority (repOrg2). The secondary, additional or alternative authority for routine surveillance reports or findings, or where the site is located. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Reportable (reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Report date (reportDate). The date a measure was reported. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Reporter ID depreciated (reporterIDDep). Unique identifier for the person or organization that is reporting the data. Status: depreciated. First released: 1.1.0. Last updated: 2.0.0.

Replicate Type (repType). Attribute of a sample, specifying whether the sample is unique, or a replicate. And if it is a replicate, what type of replicate is it? Status: active. First released: 2.0.0. Last updated: 2.0.0.

Role of contact (role). Specifies the organizational role of a given contact. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample material (saMaterial). Type of sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample ID (sampleID). Unique identifier for a sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sample ID object (sampleIDObject). Populated by **sampleID** - this specifies the object of a sample relationship. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample ID subject (sampleIDSubject). Populated by **sampleID** - this specifies the subject of a sample relationship. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample relationship (sampleRelID). Attribute for specifying the relationship between a sample subject and object. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Sample shed (sampleShed). A geographic area, physical space, or structure. A sample is taken from a sampleShed for a representative measurement of a substance(s). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Access to self (self). If this is 'no', this data will not be shown on the portal when this reporter logs in. If missing, data will be available to this reporter. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Date sample was sent (sentDate). The date the sample was sent for analyses at a laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Set ID (setID). The unique identifier of a value set. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Set Value (setValue). The partID for the set value or category. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Severity indicator (severity). An indicator of the severity or seriousness of a quality flag. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sewer network file link depreciated (sewerNetworkFileBLOB). Link to a file that has any detailed information about the sewer network associated with the site (any format). Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Site ID default (siteDef). Used as default when a new sample is created by this reporter. See ID in Site table. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Site ID (siteID). Unique identifier for the location where a sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Attributes

Site measure ID (siteMeasureID). Unique identifier for wide table only. Use when all measures are performed on a single sample. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Site Type (siteType). Type of site or institution where sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer network file link depreciated 2 (snfl). Link to a file that has any detailed information about the sewer network associated with the site (any format). Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Source Protocol ID (sourceProtocol). A protocol that served as a basis for another protocol. This may be due to version changes, or a referencing a protocol that was later adapted. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step source ID (sourceStep). Specifies the protocol step which serves as a basis for a given protocol step. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen (specimen). The substance or thing upon which the observation was made. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Departments, states, or provinces (stateProvLevel). Departments, states, or provinces Status: active. First released: 2.0.0. Last updated: 2.0.0.

State, Province, or Region (stateProvReg). The state, province, or region where a site or organization is located; part of the address. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Protocol Step ID (stepID). The unique identifier for a specific protocol step. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Step ID Object (stepIDObj). The object of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Step ID Subject (stepIDSub). The subject of the relationship between one protocol step and another step or protocol, within a given protocol ID container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Method step parent ID (stepProvenanceID). A method step that served as a basis for another method step. This may be due to version changes, or a referencing a method set that was later adapted. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol step version (stepVer). Specifies the version of a given protocol step. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer network file blob (storTempDef). A file blob that has any detailed information about the sewer network associated with the site (any format). Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Summary (summ). Short description of the assay and how it is different from the other assay methods. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Type other depreciate2 (tyOtDep2). Description of the site when the site is not listed. See siteType. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Unit (unit). The units of a measurement. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Unit set (unitSet). An identification of a set of units that can be used for a measure or method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

U site measure ID (uSiteMeasureID). Unique identifier for each measurement for a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Value (value). Value of a measure, observation or attribute. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure identification (wide table) (wMeasureID). Unique identifier for wide table only. Use when all measures are performed on a single sample at the same time and same laboratory. Suggestion: siteID_sampleID_LabID_reportDate_ID. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Classes

(classes). A class is a collection of one or more related measures or methods within a group. A group can have one more more classes to describe different parts or the class.

Alleles class (allele). Measures and methods related to alleles. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Bacteria Class (bacteria). Measures and methods relating to bacteria. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Capacity class (capacity). Measures and methods relating to capacity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Conductivity class (conductivity). Measures and methods related to conductivity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dilution Class (dilution). Measures and methods relating to dilutions. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Diseases (human) class (disease). Measure and methods related to disease or infection in humans. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Flow class (flow). Measures and methods related to flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gas class (gas). Measures and methods relating to gas and gases. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Humidity class (humid). Measures and methods related to humidity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dictionary tables (lookup). Tables to describe and support the ODM dictionary. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Mutations class (mutation). Measures and methods related to mutations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Class not applicable (naClass). Class is not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak class (outbreak). Measures and methods relating to public health outbreaks. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qPCR Class (pcr). Measures and methods relating to qPCR. Status: active. First released: 2.0.0. Last updated: 2.0.0.

pH class (pHClass). Measures and methods relating to pH. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Physical class (physical). Measures and methods related to generic physical properties. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Precipitation class (precipitation). Measures and methods related to precipitation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Program decription tables (programDescr). Tables used to describe surveillance and testing programs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Results tables (results). Tables used to record samples, measures and quality reports. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Genetic sequences class (sequence). Measures and methods related to genetic sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Speed class (speed). Measures and methods relating to speed. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard concentrations class (standardConc). Measures and methods relating to standard concentrations. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard curve class (standardCurve). Measures and methods relating to generating or recording a standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Temperature class (temperature). Measures and methods related to temperatures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Time class (time). Measures and methods relating to time. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Turbidity class (turbidity). Measures and methods relating to water and wastewater turbidity. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Variants class (variant). Measures and methods relating to variants. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Weather class (weather). Measures and methods relating to weather. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Compartments

(compartments). The substance from which a sample was taken. For example, wastewater has component of “water”. Compartments are attributes of measures, methods, units, and aggregations.

Human compartment (hum). A measure or observation made about a human. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment not applicable (naCompartment). Compartment not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment (wat). A measure or observation made from a substance in the water, including wastewater. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Compartment sets

(compartmentSets). Sets of compartments. Compartment sets are used to identify when a measure can be recorded for more than one compartment. For example, SARS-CoV-2 can be measured in people (humans), water, surface, or air.

Air compartment set

(airCompartmentSet)A compartment set for measures and methods in the air compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Air and surface compartment set

(airSurfaceCompartmentSet)A compartment set for measures and methods in the air or surface compartments. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Air and water compartment set

(airWaterCompartmentSet)A compartment set for measures and methods in the air or water compartments. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Any compartment set

(anyCompartmentSet)A compartment set for measures and methods in any compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sets

Air compartment (air). A measure or observation made from a substance in the air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment (wat). A measure or observation made from a substance in the water, including wastewater. Status: deprecated. First released: 1.0.0. Last updated: 2.0.0.

Human compartment set

(humanCompartmentSet)A compartment set for measures and methods in the human compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Human compartment (hum). A measure or observation made about a human. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Compartment set not applicable

(naCompartmentSet)Compartment not applicable. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Compartment not applicable (naCompartment). Compartment not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface compartment set

(surfaceCompartmentSet)A compartment set for measures and methods in the surface compartment. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Surface and water compartment set

(surfaceWaterCompartmentSet)A compartment set for measures and methods in the surface or water compartments. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Surface compartment (surf). A measure or observation made from a substance on a surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water compartment (wat). A measure or observation made from a substance in the water, including wastewater. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Water compartment set

(waterCompartmentSet)A compartment set for measures and methods in the water compartment. Status: active. First Released: 1.0.0. Last updated: 2.0.0.

Water compartment (wat). A measure or observation made from a substance in the water, including wastewater. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Data types

(dataTypes). The data type for a part. Data types used in the ODM include: varchar, boolean, float, category, date, time, datetime, url, email. dataType corresponds to the entry or cell within a data table, most commonly within a report table. If the data entry has a unit, then the dataType corresponds to the unit and the dataType refers to 'unit'. If the data entry is a category, then the dataType refers to 'category'. All categories are varchar. Otherwise the dataType is identified within the part entry. TBA: dataType for dictionary.

Binary Large Object (BLOB) data type (blob). The data type for Binary Large Object (BLOB) data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Boolean data type (boolean). The data type for boolean/binary data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Categorical data type (categorical). The data type for categorical data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Datetime data type (datetime). The data type for date and time data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Float data type (float). The data type for float data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Integer data type (integer). The data type for integers. Status: active. First released: 2.0.0. Last updated: 2.0.0.

See unit (data type) (seeUnitData). The data type entry when the data type for a given entry is specified by the unit used. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Variable character data type (varchar). The data type for variable character data. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Domains

(domains). There are three domain types: biologic (i.e. Covid-19, chemical (i.e. nitrogen), physical measure (temperature).

All domains (allDo). Domain that specifies that it could apply to all domains; biological, chemical, and physical. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Biologic (bio). A living organism or biological substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Chemical (che). A chemical compound. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Domain not applicable (naDomain). Not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Physical property (phy). A physical property or object not characterized by life or chemistry. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Groups

(groups). A collection of related measures. Used primarily to group measurements and methods, this helps pare down the drop down list for a given measurement or

Miscellaneous bacteria group (bactMisc). A group of measures/methods related to miscellaneous bacteria. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Collection group (colGrp). A group of measurement-like attributes related to sample collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measurement group (measGrp). A group of measures that cannot be otherwise categorized. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous attribute group (miscAttr). A group of miscellaneous measurement-like attributes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous measure group (miscMeas). A group of measures that cannot be otherwise categorized. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Group not applicable (naGroup). Used when there is no group; ie. group is not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population group (popGrp). A group of measures/methods related to population-level factors. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Processing group (procGrp). A group of measures/methods related to processing samples for analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

SARS-CoV-2 (sarsCov2). A group of measures/methods related to the SARS-CoV-2 virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site features group (siteFeat). A group of environmental measures/methods and those pertaining to the features of a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Miscellaneous viruses group (virusMisc). A group of measures/methods related to miscellaneous viruses. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measures

(measurements). The attribute to describe a part type of measures. All measures have `partType` = “measure”.

a1306s delta-variant gene target (a1306s). a1306s delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a1918v delta-variant gene target (a1918v). a1918v delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a2710t omicron-variant gene target (a2710t). a2710t omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a63t omicron-variant gene target (a63t). a63t omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

a67v omicron-variant gene target (a67v). a67v omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Absolute humidity (absHum). A measure of the total mass of water vapour present in the air per volume of air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Environmental temperature (airTemp). Environmental temperature. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Arrival temperature (arTemp). The temperature of a sample upon arrival at the laboratory for analysis. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

bcov (bcov). Measure of the amount of bovine coronavirus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Beta (beta). B.1.351 Status: active. First released: 2.0.0. Last updated: 2.0.0.

5-Day carbonaceous biochemical oxygen demand (bod5c). The quantity of oxygen utilized for the biochemical degradation of organic matter under standard laboratory procedures in five (5) days in the presence of a nitrification inhibitor, expressed in milligrams per litre (mg/l). Status: active. First released: 1.1.0. Last updated: 2.0.0.

5-day total biochemical oxygen demand (bod5t). 5 day total biochemical oxygen demand. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Bovine respiratory syncytial virus group (brsv). Bovine respiratory syncytial virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

BRSV-N (brsvN). bovine respiratory syncytial virus capsid protein gene region Status: active. First released: 1.0.0. Last updated: 2.0.0.

Carbon Dioxide (co2). A measure of an amount or concentraion of carbon dioxide. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Chemical Oxygen Demand (cod). Chemical oxygen demand. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Water conductivity (cond). Measurement of conductivity of sample or site. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Covid-19 (cov). Covid-19 infection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

SARS-CoV-2 measure (cov2Me). Measure the amount of SARS-CoV-2 virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

SARS-CoV-2-B.1.1.7 (covB117). Variant B.1.1.7 of the SARS-CoV-2 virus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-B.1.351 (covB135). Variant B.1.351 of the SARS-CoV-2 virus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-E (covE). SARS-CoV-2 E gene. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N1 (covN1). SARS-CoV-2 nucleocapsid gene, allele 1. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N2 (covN2). SARS-CoV-2 nucleocapsid gene, allele 2. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-N3 (covN3). SARS-CoV-2 nucleocapsid gene, allele 3. Status: active. First released: 1.0.0. Last updated: 2.0.0.

SARS-CoV-2-P.1 (covP1). Variant P.1 of the SARS-CoV-2 virus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

SARS-CoV-2-RdRp (covRdRp). SARS-CoV-2 RdRp gene. Status: active. First released: 1.0.0. Last updated: 2.0.0.

crAssphage-N (cra). crAssphage virus capsid protein gene region Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cycle threshold or quantification cycle (Ct or Cq) (ct). Cycle thresholds in a PCR assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

d377y delta-variant gene target (d377y). d377y delta-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

d3g omicron-variant gene target (d3g). d3g omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

d63g delta-variant gene target (d63g). d63g delta-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

d796y omicron-variant gene target (d796y). d796y omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

d950n delta-variant gene target (d950n). d950n delta-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ddcov_e sars-cov-2 gene target (ddcovE). ddcov_e sars-cov-2 gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ddcov_n sars-cov-2 gene target (ddcovN). ddcov_n sars-cov-2 gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del143/145 (del143). 143 or 145 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del 157/158 (del157). 157 or 158 deletion delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

del2084/2084 (del2084). 2084 or 2084 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del212/212 (del212). 212 or 212 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del3674/3676 (del3674). 3674 or 3676 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

del69/70 (del6970). 69 or 70 deletion omicron-variant gene target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Delta (delta). B.1.617.2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dilution factor (dilFact). Specifies the extent to which a sample or aliquot was diluted prior to analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Point dillutions (dillute). Exact concentration or dillutions for generating a standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

e156g delta-variant gene target (e156g). e156g delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

e484a omicron-variant gene target (e484a). e484a omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Escherichia coli (ecoli). Concentration of bacteria that are passed through the faecal excrement of humans, livestock and wildlife Status: active. First released: 2.0.0. Last updated: 2.0.0.

Efficiency (efficient). The efficiency reported for a standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rainfall (envRnF). Rainfall, i.e. amount of precipitation in the form of rain. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ground snow depth (envSnwD). Total depth of snow on the ground. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Snowfall (envSnwF). Snowfall, i.e. amount of precipitation in the form of snow. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sarbecovirus-specific E sars-cov-2 gene target (eSarbec). Sarbecovirus-specific E sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Estimated frequency of reads (estFreqReads). Estimated frequency of reads in a sequencing assay Status: active. First released: 2.0.0. Last updated: 2.0.0.

Extraction volume of sample (exvol). Extraction volume of sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Flow rate (floRate). Wastewater volumetric flow rate at the sample collection location over the 24-hr period during which the sample was collected. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Flow volume (flowVol). Volume of influent. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Influenza virus measure (flu). General influenza virus measure Status: active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus A1 (fluA1). Influenza virus A1 type Status: active. First released: 2.0.0. Last updated: 2.0.0.

influenza virus A2 (fluA2). influenza virus A2 type Status: active. First released: 2.0.0. Last updated: 2.0.0.

Influenza virus B (fluB). Influenza virus B type Status: active. First released: 2.0.0. Last updated: 2.0.0.

F-Specific RNA bacteriophages (frna). A measure for amount of F-Specific RNA bacteriophages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

F-Specific RNA bacteriophages, G2 (frnaG2). A measure for amount of G2 F-Specific RNA bacteriophages. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample temperature (fst). Temperature that the sample is stored at while it is being sampled. This field is mainly relevant for composite samples which are either kept at ambient temperature or refrigerated while being sampled. Status: active. First released: 1.0.0. Last updated: 2.0.0.

g215c delta-variant gene target (g215c). g215c delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

g339d omicron-variant gene target (g339d). g339d omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

g496s omicron-variant gene target (g496s). g496s omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

g662s delta-variant gene target (g662s). g662s delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gamma (gam). P.1 Status: active. First released: 2.0.0. Last updated: 2.0.0.

hep g armored rna (hepGRna). Measure of the amount Hepatitis G Armored RNA. Status: active. First released: 2.0.0. Last updated: 2.0.0.

i1566v omicron-variant gene target (i1566v). i1566v omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

i3758v omicron-variant gene target (i3758v). i3758v omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

i82t delta-variant gene target (i82t). i82t delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Inhibition measure (inhibMe). Parameter to report whether or not inhibition was detected in the sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Intercept (inter). Intercept value of the calibration curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Combined ip2 and ip4 sars-cov-2 gene target (ip2ip4). Combined ip2 and ip4 sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

k856r omicron-variant gene target (k856r). k856r omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

l452r delta-variant gene target (l452r). l452r delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

l981f omicron-variant gene target (l981f). l981f omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Lambda (lambda). C.37 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Limit of detection (LOD) (lod). Limit of detection. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sets

Limit of detection (LOD) - sewage (lodSewa). Limit of detection for sewage samples. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Limit of quantification (LOQ) (loq). Loq Status: active. First released: 1.1.0. Last updated: 2.0.0.

Other Measure (measOth). Other measure, not otherwise specified in measures. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Other Measure Description (meOthDe). Description of other measure (measOtherDesc). Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Murine Hepatitis Virus (mhv). A measure for amount of murine hepatitis viurs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage (ms2Col). Measure of the amount of ms2 coliphage. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Mu (mu). B.1.621 Status: active. First released: 2.0.0. Last updated: 2.0.0.

murine coronavirus (muCo). Measure of the amount of murine coronavirus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

N sars-cov-2 gene target (n). N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Combined N1 and N2 sars-cov-2 gene target (n1n2). Combined N1 and N2 sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n211i omicron-variant gene target (n211i). n211i omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n679k omicron-variant gene target (n679k). n679k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n856k omicron-variant gene target (n856k). n856k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

n969k omicron-variant gene target (n969k). n969k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Ammonium Nitrogen (nH4N). Ammonium nitrogen concentration, as N. Status: active. First released: 1.0.0. Last updated: 2.0.0.

niid_2019-ncov_n sars-cov-2 gene target (niid19). niid_2019-ncov_n sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Norovirus G1 (noroG1). Norovirus genogroup 1 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Norovirus G2 (noroG2). Norovirus genogroup 2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sarbecovirus-specific N sars-cov-2 gene target (nSarbec). Sarbecovirus-specific N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

coronavirus OC43 (oc43). Measure of the amount of coronavirus OC43. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.1 (omicr1). Omicron B.1.1.529 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2 (omicr2). Omicron BA.2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.2.75 (omicr275). Omicron BA.2.75 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.4 (omicr4). Omicron BA.4 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Omicron BA.5 (omicr5). Omicron BA.5 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Orthophosphates (ophos). Ortho-phosphate concentration. Status: active. First released: 1.0.0. Last updated: 2.0.0.

ORF1a sars-cov-2 gene target (orf1a). ORF1a sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

ORF1ab sars-cov-2 gene target (orf1ab). ORF1ab sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

ORF1b sars-cov-2 gene target (orf1b). ORF1b sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Variants Other (otherV). Used for reporting variants detected for which a specific measureID doesn't already exist. The ODM team will periodically migrate collected responses in this category into set measureIDs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak (outb). Measure to indicate outbreak status. Given that when outbreaks occur, full data on case counts may not be available or completely accurate, using this measure indicates that case counts may be approximate. Status: active. First released: 2.0.0. Last updated: 2.0.0.

p100l delta-variant gene target (p100l). p100l delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p2046l delta-variant gene target (p2046l). p2046l delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p2287s delta-variant gene target (p2287s). p2287s delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p3395h omicron-variant gene target (p3395h). p3395h omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

p681r delta-variant gene target (p681r). p681r delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

pH (ph). pH measurement Status: active. First released: 1.1.0. Last updated: 2.0.0.

Pseudomonas virus phi6 (phi6). Measure of the amount of pseudomonas virus phi6. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Total phosphorous (phos). Total phosphorous Status: active. First released: 1.1.0. Last updated: 2.0.0.

Total Phosphates (phostot). Total phosphates Status: active. First released: 1.1.0. Last updated: 2.0.0.

PHRED quality score (phred). PHRED Quality Score For Sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

PMMoV-CP (ppmv). Pepper mild mottle virus capsid protein gene region Status: active. First released: 1.0.0. Last updated: 2.0.0.

Puro virus (puro). Measure of the amount of puro virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

q19e omicron-variant gene target (q19e). q19e omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

q493r omicron-variant gene target (q493r). q493r omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

q498r omicron-variant gene target (q498r). q498r omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

q954h omicron-variant gene target (q954h). q954h omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

R squared (r2). R-squared value of the calibration curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

r203m delta-variant gene target (r203m). r203m delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

IP2 rdrp sars-cov-2 gene target (rdrpIP2). IP2 rdrp sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

IP4 rdrp sars-cov-2 gene target (rdrpIP4). IP4 rdrp sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Relative humidity (relHum). The unit of relative humidity, or the air-water mixture. The ratio of the partial pressure of water vapor in the mixture to the equilibrium vapor pressure of water over a flat surface. of pure water at a given temperature. Status: active. First released: 2.0.0. Last updated: 2.0.0.

replicate number (repNum). The replicate number for a Ct or Cq value - used for specifying out a protocol standard curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Respiratory syncytial virus (rsv). Respiratory syncytial virus. Status: active. First released: 2.0.0. Last updated: 2.0.0.

s2083i omicron-variant gene target (s2083i). s2083i omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s371l omicron-variant gene target (s371l). s371l omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s373p omicron-variant gene target (s373p). s373p omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s375f omicron-variant gene target (s375f). s375f omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

s477n omicron-variant gene target (s477n). s477n omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Settled solids (settsol). Amount of settled solids from a wastewater sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

S sars-cov-2 gene target (sGene). S sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

slope (slope). Slope value of the calibration curve. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity (specHum). Measure for specific humidity, the unit is the ratio of the mass of water vapour and the mass of total air. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Storage temp (sTemp). Temperature that the sample is stored at in Celsius. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Storage time (stoTim). Length of time that a sample was in storage. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sub-variant or lineage (subVar). A unit used to report a specific genetic lineage, or to specify that a reported variant is a sub-variant of another. Status: active. First released: 2.0.0. Last updated: 2.0.0.

t19r delta-variant gene target (t19r). t19r delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

t3646a delta-variant gene target (t3646a). t3646a delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

t547k omicron-variant gene target (t547k). t547k omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

t9i omicron-variant gene target (t9i). t9i omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

TaqPath N sars-cov-2 gene target (taqpatN). TaqPath N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

TaqPath S sars-cov-2 gene target (taqpatS). TaqPath N sars-cov-2 gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample temperature (temp). Temperature of the sample measured in degrees Celcius Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Nucleic acid template volume (tempVol). The volume of DNA or RNA template used for PCR. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Total Kjeldahl Nitrogen (tkn). A measure of Total Kjeldahl Nitrogen; ie. the sum of nitrogen bound in organic substances, nitrogen in ammonia (NH₃-N) and in ammonium (NH₄⁺-N) in the chemical analysis of soil, water, or waste water. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Total Nitrogen (tn). Total nitrogen concentration, as N. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Total solids concentration (ts). Total solids concentration Status: active. First released: 1.0.0. Last updated: 2.0.0.

Total suspended solids (tss). Total suspended solids Status: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration of total suspended solids (tssConc). Total suspended solids concentration of the wastewater. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Turbidity (turb). A measure for the turbidity of water, or a liquid sample, quantifying the relative clarity or opacity of a liquid. Highly indicative of water quality. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Tube weight empty (tWeighE). The weight of the tube used for analysis while empty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Tube weight full (tWeighF). The weight of the tube used for analysis while full of analyte. Status: active. First released: 2.0.0. Last updated: 2.0.0.

v2930l delta-variant gene target (v2930l). v2930l delta-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Size in volume (vol). Total volume of water or sludge sampled. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Volatile suspended solids (vss). Volatile suspended solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Volatile suspended solids - ignition (vssIg). A water quality measure, captured by via the loss on ignition of the mass of measured total suspended solids. This ignition generally takes place in an oven at a temperature of 550 °C to 600 °C. It represents the amount of volatile matter present in the solid fraction of the measured solution. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wastewater temperature (watTemp). Temperature of the wastewater. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Weather (weath). A measure for weather captured through qualitative categories; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wind Speed (wind). A measure for wind speed; part of a larger suite of weather measures that may be useful when collecting samples from outdoor sources. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wastewater treatment plant designed capacity (wwtpCap). A measure for the designed capacity of a wastewater treatment plant, in terms of how much water it was designed to be able to hold and process. Status: active. First released: 1.0.0. Last updated: 2.0.0.

y505h omicron-variant gene target (y505h). y505h omicron-variant gene target Status: active. First released: 2.0.0. Last updated: 2.0.0.

Methods

(methods). Procedures or steps for collecting samples or performing measures. For example, `methodExtract` is a method that describes the how a sample was extracted.

Sets

Nucleic acid extraction method (extraction). Description of the nucleic acid extraction method. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Inhibition method (inhibMeth). Description of the method used to evaluate molecular inhibition. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration method (methodConc). Description of the method to concentrate a wastewater sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

PCR method (permeth). Description of the PCR method. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Pretreatment (pretreat). Was the sample chemically treated in anyway with the addition of stabilizers or other Status: active. First released: 1.0.0. Last updated: 2.0.0.

Genetic primer (primer). Method ID used for indicating the primer used for a PCR or sequencing assay. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quaility assurance method (qaqc). Quality assurance and quality control method. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Collection other depreciated (sampleTypeOther). Description for other type of method not listed in collection. depreciated. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Standard curve frequency (scf). A method for specifying the frequency over which a standard curve (or 'Master Curve') is used. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Genetic sequencing script version (scriptVersion). Used to specify the script or script version used to extract summary information (i.e. coverage and mutation frequency statistics, etc.) from the sequencing output. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Layout (seqLay). The layout of genetic material used in sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method (seqSel). The primer sequence selection method used in sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy (seqStrat). The sequencing strategy used for an analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Shipped on ice (shipOnIce). Was the sample kept cool while being shipped to the lab Status: active. First released: 1.0.0. Last updated: 2.0.0.

Solid seperation (solidSep). Process used to separate solid and liquid phases of the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Spike material (spikeMat). Material into which the recovery efficiency control target is spiked. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Recovery efficiency spike target (spikeTarget). Method specifying the recovery efficiency control target. This matches on to the the spike material method ID. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Type other depreciate (tyOtDep). Description for other type of sample not listed in Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Missingness

(missingness). The part type for missingness values. Only used for the dictionary entries of parts.

Not applicable (NA). The field for which the expected value is not a property of the described object. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not a number (nan). The outcome of a measurement is not a valid number (plus or minus infinity, error, ...) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported (nr). A value could have been recorded, however, it was not. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Null (null). A logical representation of a statement that is neither TRUE nor FALSE. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Undisclosed (undisc). A value has been recorded, however it was not included in the dataset. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Missingness sets

(missingnessSets). Missingness sets for measures, methods or attributes. Inputted data into an ODM field may have missing data. The missingness set lists valid missing categories for the measure, method, or attribute.

General missingness set

(genMissingnessSet)The general set for missingness values. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Measure, method, or attribute sets

Not applicable (NA). The field for which the expected value is not a property of the described object. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not a number (nan). The outcome of a measurement is not a valid number (plus or minus infinity, error, ...) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported (nr). A value could have been recorded, however, it was not. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Null (null). A logical representation of a statement that is neither TRUE nor FALSE. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Undisclosed (undisc). A value has been recorded, however it was not included in the dataset. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported, Not applicable missing set

(nrNAMissingnessSet)Not reported, Not applicable missing set. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Not applicable (NA). The field for which the expected value is not a property of the described object. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Not reported (nr). A value could have been recorded, however, it was not. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measure, method, or attribute sets

(mmaSets). A set of categories. For example, site type has a list of different sites, such as a wastewater treatment plant, a septic tank, or a pumping station. Only used for the dictionary entries of parts.

Boolean value set

(booleanSet)Set that contains the valid possible values for a boolean measure (TRUE or FALSE). Status: active. First Released: 2.0.0. Last updated: 2.0.0.

FALSE (false). Boolean data type = FALSE Status: active. First released: 1.1.0. Last updated: 2.0.0.

TRUE (true). Boolean data type = TRUE Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sample collection set

(collectSet)Methods for collection samples. Sample collection methods include water, air, and surface. See the attribute **Compartment set** of the sample collection method. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Area proportional sample (areaPr). An area proportional sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Composite sample - archival (comp). A composite sample, usually generated by an autosampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

COSCa ball (cosca). COSCa passive sampling device. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Flow proportional sample (flowPr). A flow proportional composite sample generally collected by an autosampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Grab sample (grb). A single large representative grab sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Moore swab passive sample (moorSw). Moore swab passive sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Surface swab (surfSw). Surface swab. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Time proportional sample (timePr). A time proportional composite sample generally collected by an autosampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Volume proportional sample (volPr). A volume proportional sample generally collected by an autosampler. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nucleic Acid Extraction set

(extractSet)set used for storing all the valid category values for the nucleic acid extraction method. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Amicon filter, extract with MP96 (amiMP96). Nucleic acid extraction, usually used for the liquid fraction of a wastewater sample, using amicon filtration and using an MP96 for final extraction. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Bead beating, extract with MP96 (bbMP96). Nucleic acid extraction, usually used for solid fraction of a wastewater sample, using bead beating and using an MP96 for final extraction. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Chemagic viral dna/rna 300 kit (chemVir). Nucleic acid extraction performed using the chemagic viral dna/rna 300 kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

4s method (extract4s). Nucleic acid extraction performed using the 4s method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Nuclisens automated magnetic bead extraction kit (nucAuto). Nucleic acid extraction performed using the nuclisens automated magnetic bead extraction kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nuclisens manual magnetic bead extraction kit (nucManu). Nucleic acid extraction performed using the nuclisens manual magnetic bead extraction kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Phenol chloroform (phenCl). Nucleic acid extraction performed using phenol chloroform. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega automated tna kit (promAuto). Nucleic acid extraction performed using the promega automated tna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega ht tna kit (promHt). Nucleic acid extraction performed using the promega ht tna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega manual tna kit (promManu). Nucleic acid extraction performed using the promega manual tna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

promega wastewater large volume tna capture kit (promWW). Nucleic acid extraction performed using the promega wastewater large volume tna capture kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep dna/rna kit (qgDNARNA). Nucleic acid extraction performed using the qiagen allprep dna/rna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen ez1 virus mini kit v2.0 (qgEz1). Nucleic acid extraction performed using the qiagen ez1 virus mini kit v2.0. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep powerfecal dna/rna kit (qgPwrFecal). Nucleic acid extraction performed using the qiagen allprep powerfecal dna/rna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen allprep powerviral dna/rna kit (qgPwrViral). Nucleic acid extraction performed using the qiagen allprep powerviral dna/rna kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen powerwater kit (qgPwrWtr). Nucleic acid extraction performed using the qiagen powerwater kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen qiaamp buffers with epoch columns (qgQiAmp). Nucleic acid extraction performed using the qiagen qiaamp buffers with epoch columns. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy kit (qgRneasy). Nucleic acid extraction performed using the qiagen rneasy kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen rneasy powermicrobiome kit (qgRneasyPwr). Nucleic acid extraction performed using the qiagen rneasy powermicrobiome kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Resuspend COSCa filter collection (resCosca). Nucleic acid extraction via resuspension of a sample collected using the COSCa ball method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

thermo magmax microbiome ultra nucleic acid isolation kit (thermMag). Nucleic acid extraction performed using the thermo magmax microbiome ultra nucleic acid isolation kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

trizol, zymo mag beads w/ zymo clean and concentrator (trizol). Nucleic acid extraction performed using the trizol, zymo mag beads w/ zymo clean and concentrator method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

zymo environ water rna kit/ zymo environ water rna kit (cat. r2042) (zymoEnv). Nucleic acid extraction performed using the zymo environ water rna kit/ zymo environ water rna kit (cat. r2042). Status: active. First released: 2.0.0. Last updated: 2.0.0.

zymo quick-rna fungal/bacterial miniprep #r2014 (zymoQuick). Nucleic acid extraction performed using the zymo quick-rna fungal/bacterial miniprep #r2014 kit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample fraction set

(fractionSet)set for the fraction of the sample (solid, liquid, etc.). Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Liquid fraction (liq). Liquid fraction of a sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Mixed/homogenized sample (mix). Mixed or homogenized sample or fraction analyzed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Solid fraction (sol). Solid fraction of a sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Geographic set

(geoTypeSet)set for different type of geographic components. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sewer Network Health Region (hlthReg). Health region served by the sewer network Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewer catchment area (swrSet). Sewer catchment area. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Inhibition set

(inhibitionSet)Category set for inhibition methods. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Serial dilution (serialDilution). The serial dilution method for assessing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Spike matrix and recovery (spike). The spike matrix and recovery method for assessing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Instrument set

(insTypeSet)List of instruments that are used for measures and methods Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Atline Analyzer (aas). Atline analyzer with sampler. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Hand Measurement (hma). Handheld measurement analyzer. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Other instrument (instrumentTypeOther). Type of instrument other than those included in the PHES-ODM. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lab Analysis (ola). Offline laboratory analysis. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Online Sensor (onse). Online sensor Status: active. First released: 1.0.0. Last updated: 2.0.0.

Concentration method set

(methodConcSet)A set a concentration methods. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Aloh3 precipitation (aloh3). Aloh3 precipitation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Amicon ultrafiltration (amiconUf). Amicon ultrafiltration Status: active. First released: 2.0.0. Last updated: 2.0.0.

Beef extract flocculation (beeExtractFloc). Beef extract flocculation. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Centricon ultrafiltration (centriconUf). Centricon ultrafiltration. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Ceres nanotrap (ceres). Ceres nanotrap. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Hollow fiber dead end ultrafiltration (hollowFiberUF). Hollow fiber dead end ultrafiltration Status: active. First released: 2.0.0. Last updated: 2.0.0.

Innovaprep ultrafiltration (innovaprepUF). Innovaprep ultrafiltration Status: active. First released: 2.0.0. Last updated: 2.0.0.

Polyethyleneglycol (PEG) precipitation (peg). Peg precipitation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Skimmed milk flocculation (skimMilkFloc). Skimmed milk flocculation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Ultracentrifugation (untracent). Ultracentrifugation Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization level set

(orgLevelSet)Categories of organization levels. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Administrative regions (admRegLevel). Administrative regions Status: active. First released: 2.0.0. Last updated: 2.0.0.

Countries or sovereign states (countryLevel). Countries or sovereign states Status: active. First released: 2.0.0. Last updated: 2.0.0.

Districts, counties, regions (countyLevel). Districts, counties, regions Status: active. First released: 2.0.0. Last updated: 2.0.0.

Municipalities or communes (municipalLevel). Municipalities or communes Status: active. First released: 2.0.0. Last updated: 2.0.0.

Local administrative units or neighborhoods (neighborLevel). Local administrative units or neighborhoods Status: active. First released: 2.0.0. Last updated: 2.0.0.

Departments, states, or provinces (stateProvLevel). Departments, states, or provinces Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization sector set

(orgSectorSet)Categories of organization sectors. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Airport (airport). Airport sample shed category type Status: active. First released: 2.0.0. Last updated: 2.0.0.

Long-term acute care hospital (ccc). Acute care hospitals, or complex continuing care, that provide care for patients with average length of stay longer than 25 days. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Higher education dormitory or residential building (dorm). Higher education dormitory or residential building Status: active. First released: 2.0.0. Last updated: 2.0.0.

First Nation (fiNa). Used to categorize a sampled that is a First Nation, or on reserve lands. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Health administration or planning agency (healthAdm). Health administrative or planning organization. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory (lab). Laboratory for environmental testing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Long-term care facility (ltcf). A residential healthcare facility that provides 24-medical care. These are also called skilled nursing facilities. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Long-term care - assisted living or retirement home (ltcfAl). A residential facility that provides assistance with daily care but generally does not provide skilled nursing care. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Public health agency (pubHealth). Public health agency. Status: active. First released: 2.0.0. Last updated: 2.0.0.

School (school). A school serving students in the kindergarten to 12th grade range Status: active. First released: 1.0.0. Last updated: 2.0.0.

Social services shelter (sss). Other type of social services shelter Status: active. First released: 2.0.0. Last updated: 2.0.0.

University campus (uCampus). University or college campus - comprising an entire campus or part of a campus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Organization type set

(orgTypeSet)The set for storing organization types. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Academic institution (academ). The category of organization type used for academic institutions or research groups. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Government agency (govt). The category of organization type used for government agencies, programs, or crown-owned bodies. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Private sector (priv). The category of organization type used for private sector or non-profit groups that are not otherwise captured by the academic category. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample Origin set

(originSet)The set for storing the valid categorical values of sample origin. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Derived sample (derived). Specifies a sample that is derrived or made from another sample or material. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample (field). Specifies a sample taken from the field; directly collected from an area for testing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Synthetic sample (synthetic). Specifies a synthetic, or artificially derived sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Outbreak set

(outbreakSet)set for the valid values of the outbreak measure. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Outbreak End (outbEnd). Indicates an outbreak has ended. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak - on-going (outbOngoing). Indicates an outbreak is on-going. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Outbreak start (outbStart). Indicates an outbreak began. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PCR Method set

(pcrSet)The set capturing containing all of the valid methods for the PCR methodID. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Digital droplet emulsification PCR (bioRadDdpcr). Describes a PCR analysis done using BioRad's digital droplet emulsification PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Fluidigm digital PCR (fluiDpcr). Describes a PCR analysis done using Fluidigm's digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Life technologies digital PCR (ltDpcr). Describes a PCR analysis done using Life Technologies' digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

qiagen digital PCR (qgDpcr). Describes a PCR analysis done using Qiagen's digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quantitative PCR (qpcr). Real-time PCR, also called 'quantAggScale' PCR Status: active. First released: 2.0.0. Last updated: 2.0.0.

Raindance digital PCR (rainDpcr). Describes a PCR analysis done using Raindance's digital PCR technology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Genetic primer set

(primerSet)The set for genetic primers used in sequencing or PCR assays. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Use of the articV3 primer (articV3). Initial implementation of an ARTIC bioinformatics platform for nanopore sequencing of nCoV2019 novel coronavirus; Artic Foundation v3 primer. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Use of the articV4 primer (articV4). Artic V4 sequencing primer for VOCs and VOIs. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol Relationships set

(protocolRelSet)set for valid values of relationshipID in the protocolRelationships table. Use only for protocolRelationships table relationshipID. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Is Before (before). Specifies that the object step or protocol occurs before the subject step or protocol in the overall protocol container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Is concurrent to (concurrent). Specifies that the object step or protocol occurs at the same time as the subject step or protocol in the overall protocol container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Specifies (specifies). Specifies that the object step or protocol occurs specifies details for the subject step or protocol in the overall protocol container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Purpose set

(purposeSet)Purpose set for a sample or measure. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Education (education). A measure or sample taken for education or training. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple Purpose (multiple). A measure or sample taken for multiple purposes, not easily captured by the other purpose categories. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Provisional report (provisional). A provisional or interm report. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality control (qualityControl). A measure or sample taken for the purpose of quality control. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Regular (regular). A measure or sample taken for surveillance or epidemiology. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Testing (testing). A measure or sample taken to test a new method. These measures are typically used for internal lab uses and not reported to external partners. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Validation study (validationStudy). A measure or sample taken for a validation study. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Replicate Type set

(replicateSet)The set for storing valid categorical values of replicate type. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Co-located sample (colocated). Second or multiples samples collected at same location but different time (water, air) or at a nearby location (soil, sediment). Sent blind to laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Field sample replicate (fieldReplicate). A sample divided into two or more homogeneous parts. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory duplicate (labDuplicate). Second (time or more) processing and analysis of sample. Usually for general chemistry or metals analyses. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory control sample duplicate (lcsd). Known amounts of an analyte or representative compounds are added to a second “clean” matrix (lab water or clean sand) in laboratory. Duplicate of LCS Status: active. First released: 2.0.0. Last updated: 2.0.0.

Matrix spike duplicate (msd). A known amounts of an analyte or representative compounds are added in the laboratory to a second aliquot of the sample used for matrix spike. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unique sample (unique). A unique sample, not duplicated. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample material set

(sampleMatSet)set for the types of material that can be found in a sample. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sets

Air filter (afu). Air filter as part of filtration or circulating unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Desk or counter (desk). Desk, table, countertop or other flat working surface. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Fecal matter (faeces). Fecal matter. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Floor (floor). Floor of a building or room. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Holding tank wastewater (htSam). Wastewater sampled from a holding tank, such as from an airplane or ship Status: active. First released: 2.0.0. Last updated: 2.0.0.

Water (nww). Non-wastewater, coming from any kind of water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Primary clarifier effluent (pEfflu). Effluent obtained after primary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Primary clarifier sludge (pSludge). Sludge produced by primary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Post-grit (pstGrit). Raw wastewater after a treatment plant's headworks (post grit or removal of large solids at a treatment plant but prior to a primary clarifier) Status: active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage at site (rawWW). Wastewater without any form of treatment Status: active. First released: 1.0.0. Last updated: 2.0.0.

Raw sewage downstream from a site (rawWWdown). Downstream from a site. See partType = rawWWup. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Raw sewage upstream from a site (rawWWup). Upstream from a site. Used when there is not direct access to a site's wastewater conjunction

with downstream sample. For example, if the site is a long-term care home, but there is not access to the building or property cleanout. Use two samples, one sample upstream and one sample downstream. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Secondary clarifier effluent (sEfflu). Effluent obtained after secondary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Septic tank wastewater (septage). Wastewater from within a septic tank Status: active. First released: 1.0.0. Last updated: 2.0.0.

Secondary clarifier sludge (sSludge). Sludge produced by secondary clarifiers Status: active. First released: 1.0.0. Last updated: 2.0.0.

Other surface (surface). Surface other than floor or desk. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sewer sediment (swrSed). Sediments obtained in sewer Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample relationships set

(sampleRelSet)set for valid values of relationshipID in the sampleRelationships table. Use only for sampleRelationships table relationshipID. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Child relationship (child). Indicated that this is a sample generated from (an)other sample(s) either because of pooling or sub-sampling. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Co-located sample (colocated). Second or multiples samples collected at same location but different time (water, air) or at a nearby location (soil, sediment). Sent blind to laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Field sample replicate (fieldReplicate). A sample divided into two or more homogeneous parts. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory duplicate (labDuplicate). Second (time or more) processing and analysis of sample. Usually for general chemistry or metals analyses. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Laboratory control sample duplicate (lcsd). Known amounts of an analyte or representative compounds are added to a second “clean” matrix (lab water or clean sand) in laboratory. Duplicate of LCS Status: active. First released: 2.0.0. Last updated: 2.0.0.

Matrix spike duplicate (msd). A known amounts of an analyte or representative compounds are added in the laboratory to a second aliquot of the sample used for matrix spike. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing script version set

(scriptSet)The set for script versions Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Freyja Script (freyja). Freyja Script for sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gromstole 1.0 Script (gromstole). Gromstole 1.0 Script for sequencing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Layout set

(seqLaySet)The set for the layout of genetic material used in sequencing Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Mate pair layout (mateLay). Specifies the mate pair layout for a sequencing method. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Paired Layout (pairLay). Specifies the paired layout for a sequencing method Status: active. First released: 2.0.0. Last updated: 2.0.0.

Single Layout (sinLay). Specifies the single layout for a sequencing method Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing selection method set

(seqSelSet)The set for the selection method used in sequencing Status: active. First Released: 2.0.0. Last updated: 2.0.0.

PCR sequencingsSelection method (pcrSeq). Specifies the PCR selection method for sequencing, ie. that a pre-determined primer was used. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Random sequencing selection method (ranSeq). Specifies the random selection method for sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing Strategy set

(seqStratSet)The set for the sequencing strategy used for an analysis Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Amplicon sequencing (amp). Specifies the amplicon strategy for genetic sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Whole genome sequencing (wgs). Specifies the whole genome sequencing (WGS) strategy for genetic sequencing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Severity set

(sevSet)A set for severity indicators. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

High severity (high). Indicates a very severe quality issue, likely meaning the data should not be reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low severity (low). A marker for low severity Status: active. First released: 2.0.0. Last updated: 2.0.0.

Mid-level severity (mid). A marker for med-level severity, where there are some concerns. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sampled set

(shedSet)The set for all valid values of sampledshed. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Airplane (airpln). Airplane sample shed category type Status: active. First released: 1.0.0. Last updated: 2.0.0.

Airport (airport). Airport sample shed category type Status: active. First released: 2.0.0. Last updated: 2.0.0.

Long-term acute care hospital (ccc). Acute care hospitals, or complex continuing care, that provide care for patients with average length of stay longer than 25 days. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Child day care (cdc). Child day care facility. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Cphd ID (corFcil). Correctional facility Status: active. First released: 1.0.0. Last updated: 2.0.0.

Higher education dormitory or residential building (dorm). Higher education dormitory or residential building Status: active. First released: 2.0.0. Last updated: 2.0.0.

First Nation (fiNa). Used to categorize a sampled that is a First Nation, or on reserve lands. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Hospital (hosptl). Hospital Status: active. First released: 1.0.0. Last updated: 2.0.0.

Long-term care facility (ltcf). A residential healthcare facility that provides 24-medical care. These are also called skilled nursing facilities. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Long-term care - assisted living or retirement home (ltcfAl). A residential facility that provides assistance with daily care but generally does not provide skilled nursing care. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other long-term care (ltcfO). Other residential facilities that provide daily and/or medical care, but are not defined as nursing home/skilled nursing facilities or assisted living facilities. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Municipality (municp). A complete municipality, this specifies an entire metropolitan area, either a city, town, etc. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Neighbourhood (neigh). A municipal neighbourhood, this specifies a sub-section of a larger municipality. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other residential building (orb). Individual residential buildings or institutions not captured in other categories Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

School (school). A school serving students in the kindergarten to 12th grade range Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ship (ship). A cruise ship or other ship. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Social services shelter (sss). Other type of social services shelter Status: active. First released: 2.0.0. Last updated: 2.0.0.

Airport terminal (terminal). Airport terminal sample shed category type Status: active. First released: 2.0.0. Last updated: 2.0.0.

University campus (uCampus). University or college campus - comprising an entire campus or part of a campus. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Site set

(siteTypeSet)set for the type of sampling site. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Branch sewer pipeline (bSwrPpl). Specifies a site type that is collection pipe that run lateral to other municipal sewer lines, allowing drainage into the main sewer. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Building cleanout (buildCO). Specifies a site type that is a capped pipe that connects to a building's main sewer line. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Estuary, natural water body (estuary). Estuary, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Holding tank (htSite). Holding tank Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lagoon system (lagoon). Lagoon system for extensive wastewater treatment Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lake, natural water body (lake). Lake, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Major sewer pipeline (mSwrPpl). Major sewer pipeline Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ocean, natural water body (ocean). Ocean, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Pumping station (pStat). Pumping station Status: active. First released: 1.0.0. Last updated: 2.0.0.

Retention pond (retPond). Retention pond Status: active. First released: 1.0.0. Last updated: 2.0.0.

River, natural water body (river). River, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sea, natural water body (sea). Sea, natural water body Status: active. First released: 1.0.0. Last updated: 2.0.0.

Septic tank (septTnk). Septic tank Status: active. First released: 1.0.0. Last updated: 2.0.0.

Stabilization pond (stabPnd). Specifies a site type which is a pond designed and built for wastewater treatment to reduce the organic content and remove pathogens from wastewater. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sewage truck (swgTrck). Sewage truck Status: active. First released: 1.0.0. Last updated: 2.0.0.

Upstream sites (upstream). A general site type for upstream wastewater sampling sites. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sets

Wastewater treatment plant (wwtp). A general site type for wastewater treatment plants. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Back-up wastewater treatment site (wwtpBack). Indicates a wastewater treatment plant where multiple labs are or were sampling, and where secondary samples are being taken. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Industrial wastewater treatment plant (wwtpInd). Industrial wastewater treatment plant. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for combined sewage (wwtpMuC). Municipal wastewater treatment plant for combined sewage. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Municipal wastewater treatment plant for sanitary sewage only (wwtpMuS). Municipal wastewater treatment plant for sanitary sewage only. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Solid separation set

(solidSeparationSet)set for the separation of solids. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Centrifugation (cent). Describes solid separation from a wastewater sample via centrifugation. Likely connected to other method steps prior to analysis. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Filtration (filt). Describes solid separation from a wastewater sample via filtration. Proceeds further concentration or analysis of the liquid filtrate. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gravity settling (graSet). Describes solid separation from a wastewater sample where the sample material is allowed to settle by gravity, and then separated. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2 (mfAcidmgcl2). Membrane filtration with acidification and mgcl2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with acidification and mgcl2, membrane recombined with separated solids (mfAcidmgcl2SS). Membrane filtration with acidification and mgcl2, membrane recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with addition of mgcl2 (mfMgcl2). Membrane filtration with addition of mgcl2 Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with mgcl2 and separated solids (mfMgcl2SS). Membrane filtration with addition of mgcl2, membrane recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment (mfNoAmend). Membrane filtration with no amendment Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with no amendment, membrane recombined with separated solids (mfNoAmendSS). Membrane filtration with no amendment, membrane recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Membrane filtration with sample acidification (mfSampleAcid). Membrane filtration with sample acidification Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Membrane filtration with sample acidification, membrane recombined with separated solids (mfSampleAcidSS). Membrane filtration with sample acidification, membrane recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

No liquid concentration, liquid recombined with separated solids (noliquid). No liquid concentration, liquid recombined with separated solids Status: active. First released: 2.0.0. Last updated: 2.0.0.

Spike material set

(spikeMatSet)set for spikeMat (aterial into which the recovery efficiency control target is spiked). Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Clarified sample (clari). Clarified sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample concentrate (conc). Sample concentrate. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Dewatered solids (deso). Dewatered solids. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Lysis buffer (lysi). Lysis buffer. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Raw wastewater (raws). Raw wastes water sample Status: active. First released: 1.0.0. Last updated: 2.0.0.

Raw post-pasteurized wastewater (rppw). Raw wastewater sample post-pasteurization Status: active. First released: 1.0.0. Last updated: 2.0.0.

Recovery efficiency spike target set

(spikeTargetSet) The set capturing containing all of the valid categories for the recovery efficiency control (or spike) target. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

bcov culture spike target (bcovCul). Cultured bovine coronavirus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

bcov vaccine spike target (bcovVac). The bovine coronavirus vaccine is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

brsv culture spike target (brsvCul). Cultured bovine respiratory syncytial virus (BRSV) is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

brsv vaccine spike target (brsvVac). The bovine respiratory syncytial virus (BRSV) vaccine is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Heat inactivated sars-cov-2 virus spike target (heatInacSARS). Heat inactivated SARS-CoV-2 virus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

hep g armored rna spike target (hepGRnaMat). Hepatitis G armored RNA is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ms2 coliphage spike target (ms2ColMat). ms2 coliphage is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

murine coronavirus spike target (muCoMat). Murine coronavirus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

oc43 spike target (oc43Mat). Human coronavirus OC43 is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

phi6 spike target (phi6Mat). Pseudomonas virus phi6 is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Puro virus spike target (puroMat). Puro Virus is used as the recovery efficiency control target. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Status set

(statusSet) A set for partID = Status to indicate whether a part is in current use or not. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Active (active). Indicator that a part is in current use. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Deprecated (deprecated). Indicator to say that a part is no longer in current use in the model. See partID = status Status: active. First released: 2.0.0. Last updated: 2.0.0.

Development (development). Indicator that a part is under development use. See partID = status Status: active. First released: 2.0.0. Last updated: 2.0.0.

Table column set

(tableSet)The set for valid inputs in table columns. These categories are used in the parts list for columns named after tables to indicate the position and function of a given part in the entity relationship diagram of the data model. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Foreign key (fK). Foreign key for a table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Header (header). Header for a table. Also known as a table variable or entity relationship ‘attribute’. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Input (input). Input for a table. Indicates if a part can be used in a table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Primary key (pK). Primary key for a table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Weather set

(weathSet)A set of the valid qualitative categories for the qualitative weather measure. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Cloudy (cloudy). Qualitative category for the weather measure, specifying an overcast day with no rain. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Foggy (foggy). Qualitative category for the weather measure, specifying a foggy or hazy day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rainy (rainy). Qualitative category for the weather measure, specifying a rainy day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Snowy (snowy). Qualitative category for the weather measure, specifying a snowy day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sunny (sunny). Qualitative category for the weather measure, specifying a clear and sunny day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclatures

(nomenclatures). A classification system to report the measure class. See `partID = nomenclatureID`.

International Classification of Diseases (icd). Classification system for diseases in humans. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nomenclature not applicable (naNomenclature). Not applicable. Status: active. First released: 2.0.0. Last updated: 2.0.0.

NextClade nomenclature (nextclade). Specifies variant or genetic nomenclature as set out by NextClade. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Pangolin nomenclature (pangolin). Specifies variant or genetic nomenclature as set out by the Phylogenetic Assignment of Named Global Outbreak Lineages (Pangolin). Status: active. First released: 2.0.0. Last updated: 2.0.0.

WHO nomenclature (who). Specifies variant or genetic nomenclature as set out by the World Health Organization. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicators

(qualityIndicators). A measure of the quality of a reported value or sample. Only used for the dictionary entries of parts.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed (flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution. The reported concentration estimate is the updated result after addressing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination (flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

FI - Inhibition present and unaddressed (flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

J - Weak signal extrapolation (flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ND - Non-detect (flagND). No amplification occurred in the reaction; non-detect. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UJ - Trace signal extrapolation (flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., “trace” signal observed); the reported value is based on the extrapolation of the standard curve in

Sets

the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable (flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample frozen (frozen). Sample was frozen before analysis or processing. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage (IBC). The percentage of the genome covered by reads (the breadth) is low. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage (IDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Leaked sample (leaked). Sample leaked, some volume and material was lost. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low-volume sample (lowVol). Sample is low-volume, but was run regardless. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple issues (mI). Multiple issues have arisen in the sequencing process. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there are no quality concerns about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample not labelled (noLabel). Sample had no label. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicator sets

Quality concerns (qf1). A flag to indicate there is a quality concern, not otherwise specified. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sparse coverage (sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample spilled (spill). Sample contents spilled from container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wide 95% interval (wI). The 95% interval is too wide, low confidence in results. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored incorrectly (wrongStorage). Sample was stored inappropriately. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored at wrong temperature (wrongTemp). Sample was stored at an inappropriate temperature. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Quality indicator sets

(qualityIndSets). Sets of quality indicators or measures. For example, PCR have a quality measure set.

Generic quality flag set

(genQualitySet) A quality set to specify any generic quality concerns about a measure or sample. Status: active. First Released: 1.0.0. Last updated: 2.0.0.

Multiple issues (mI). Multiple issues have arisen in the sequencing process. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

No quality concerns (noConcern). A flag to indicate there is no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Quality concerns (qf1). A flag to indicate there is a quality concern, not otherwise specified. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measures quality set

(measQualitySet) A quality set for any measure, includes all the quality flag options for any measure. Used as an “any quality set” for measures, but excludes the sample quality flags to avoid confusion and data entry errors. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed (flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution. The reported concentration estimate is the updated result after addressing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination (flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

FI - Inhibition present and unaddressed (flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

J - Weak signal extrapolation (flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of

the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ND - Non-detect (flagND). No amplification occurred in the reaction; non-detect. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UJ - Trace signal extrapolation (flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., “trace” signal observed); the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable (flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage (lBC). The percentage of the genome covered by reads (the breadth) is low. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage (lDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Multiple issues (mI). Multiple issues have arisen in the sequencing process. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there is are no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Quality concerns (qf1). A flag to indicate there is a quality concern, not otherwise specified. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sparse coverage (sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. Status: active. First released: 2.0.0. Last updated: 2.0.0.

PCR quality set

(pcrQualitySet)Quality set for PCR measures. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

AI - Inhibition present but addressed (flagAI). The original sample was inhibited; however, the inhibition has been addressed through dilution. The reported concentration estimate is the updated result after addressing inhibition. Status: active. First released: 2.0.0. Last updated: 2.0.0.

B - Trace levels of contamination (flagB). Analytical result may be subject to “trace” levels of contamination; the target analyte was also detected in negative controls on the same run as the sample. Status: active. First released: 2.0.0. Last updated: 2.0.0.

FI - Inhibition present and unaddressed (flagFI). The original sample was inhibited; however, the inhibition has not been successfully addressed. Therefore, no concentration estimate has been reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

J - Weak signal extrapolation (flagJ). Analytical result falls below the lowest concentration of the experiment-specific standard curve but above the y-intercept value; the reported value is based on the extrapolation of

the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

ND - Non-detect (flagND). No amplification occurred in the reaction; non-detect. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UJ - Trace signal extrapolation (flagUJ). Observed quantitation cycle is greater than the experiment-specific standard curve intercept value but evidence of clear amplification was present (i.e., “trace” signal observed); the reported value is based on the extrapolation of the standard curve in the non-linear region. Status: active. First released: 2.0.0. Last updated: 2.0.0.

UQ - Unquantifiable (flagUQ). Unquantifiable, Ct value exceeds the maximum value of the standard curve. There was a detect, but we cannot quantify it with certainty. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there is are no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample quality set

(sampleQualitySet)Quality set for a sample. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sample frozen (frozen). Sample was frozen before analysis or processing Status: active. First released: 2.0.0. Last updated: 2.0.0.

Leaked sample (leaked). Sample leaked, some volume and material was lost. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low-volume sample (lowVol). Sample is low-volume, but was run regardless. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

No quality concerns (noConcern). A flag to indicate there is are no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample not labelled (noLabel). Sample had no label Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample spilled (spill). Sample contents spilled from container. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored incorrectly (wrongStorage). Sample was stored inappropriately. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample stored at wrong temperature (wrongTemp). Sample was stored at an inappropriate temperature. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sequencing quality set

(seqQualitySet)Quality set for sequencing measures Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Experiment Failed (expFail). PCR experiment failed. No value reported. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low breadth of coverage (IBC). The percentage of the genome covered by reads (the breadth) is low. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Low depth of coverage (IDC). Poor coverage, specifically an insufficient number of reads or too many sequencing reads that are mapped incorrectly. Status: active. First released: 2.0.0. Last updated: 2.0.0.

No quality concerns (noConcern). A flag to indicate there is no quality concern about the measure or sample. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample missing time stamp (noTime). Sample is missing the autosampler time; incomplete metadata on collection. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sparse coverage (sC). Sequencing shows sparse coverage, ie. numerical breadth of coverage shows as adequate but coverage plots reveal missing genome pieces. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimens

(specimens). Measures or observations are taken from three types of substances: populations (humans or geographic areas), samples, sites. See `partID = specimenID`

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen sets

(specimenSets). Sets of specimens. Specimen sets are used when a measure or attribute can be recorded for more than one specimen.

Sets

Any specimen set

(anySpecimenSet) A specimen set that includes any specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specimen set not applicable

(naSpecimenSet) A specimen set for when specimen/specimen set is not applicable. Use 'not applicable' when there is no specimen/set for the part Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Specimen not applicable (naSpecimen). Non applicable specimen. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population or sample set

(poSaSpecimenSet) A specimen set that includes population or sample specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population or site set

(poSiSpecimenSet) A specimen set that includes population or site specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population specimen set

(poSpecimenSet) A specimen set that includes only a population specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Population (po). An measure or observation for a geographic area or population. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sample specimen set

(saSpecimenSet) A specimen set that includes only a sample specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Site or sample specimen set

(siSaSpecimenSet) A specimen set that includes site or sample specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sample (sa). A measure made on a compartment or property from a sample of a substance. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Site specimen set

(siSpecimenSet) A specimen set that includes only a site specimen. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Site (si). A measure made on a compartment or property at a site. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Tables

(tables). Tables are where measures, methods and attributes are recorded. Tables represent the main entities of the environmental surveillance. There An example of a table is **site** a collection of attributes such as site name and address to describe where environment samples are taken. PHES-ODM has report tables, but users can create their own tables as well.

Address table (addresses). The table that contains information about addresses. Addresses can recorded for sites, organizations or contacts (individuals). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Contact table (contacts). The table that contains information about a contact; a person who is the contact of a site or laboratory. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Dataset table (datasets). A report table for capturing details about data's parental data set and data custodians. Supplying attribution for data collectors. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Instrument table (instruments). The table that contains information about instruments. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Language Look-up table (languages). Look up table for all languages, used to give structure to the translation table. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Measure report table (measures). The table that contains information and details about a given measure Status: active. First released: 1.0.0. Last updated: 2.0.0.

Measure set report table (measureSets). The table that identifies sets of measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Organization table (organizations). The table that contains information about a laboratory. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Parts Look-up table (parts). Look up table containing all parts in of the data model. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Polygon table (polygons). The table that contains information about the geometry of a geographic area. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Protocol relationships table (protocolRelationships). The table that contains the organizational information and details (such as order) about

Sets

a given protocol, recorded as a series of protocol steps (protocolSteps). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocols table (protocols). The table for protocols. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Protocol steps table (protocolSteps). The table for collecting metadata on individual steps in a protocol, methodological process, or assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Quality reports table (qualityReports). The table for recording the various quality metrics and indicators for samples and measures. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Reporter table depreciated (reportersDep). The table that contains information about a reporter of a sample, method, measure, or attribute. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Sample relationships table (sampleRelationships). Table for recording the relationships between samples. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sample report table (samples). The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets look-up table (sets). Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sites table (sites). The table that contains information about a site; the location where an environmental sample was taken. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Translation look-up table (translations). Look up table for translations of the description, label, and instruction for all parts. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wide name table (wideNames). The table for wide names. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Units

(units). The unit of the measurement. Every measurement must have a unit. Meaning, **units** are a mandatory field whenever a measurement is recorded.

Number of active cases (acti). Number of active cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Breadth of coverage ($\geq 5\times$ depth) (boc). Positions with read depth greater or equal to 5. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Degrees Celcius (cel). Degrees Celsius. Status: active. First released: 1.1.0. Last updated: 2.0.0.

CFU per 100 ml (cfu). Colony forming units per 100 ml of filtered sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Centimetres (cm). Unit part for the SI unit of centimetres. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Confirmed Cases Case report date (conCase). Date that the numbers were reported publicly for a confirmed case. Typically, reported data and this measure is most commonly reported and used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Number of confirmed cases (conf). Number of confirmed cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases episode date (confEp). Episode date is the earliest of onset, test or reported date for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Confirmed Cases Onset date (confOn). Earliest that symptoms were reported for a confirmed case. This data is often not known and reported. In lieu, episode is used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases test date (conTest). Date that the covid-19 test was performed for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cumulative count (cuCo). Units for describing a population measure of a cumulative count of cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Daily count (daiCo). Units for describing a population measure of a daily count of new cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Days (days). A unit for indicating a length of time in days. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Deaths (death). Units for describing a population measure of patients who have died from a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Detected (det). Substance detected or not detected. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Depth of coverage (doc). The sequencing read depth. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per copy of crAssphage (gcCrA). Gene or variant copies per copy of crAssphage. Status: active. First released: 1.0.0. Last updated: 2.0.0.

gene copies per day per 100,000 (gcD100). The unit for measures reflecting the gene copies per day per 100,000 people in the population. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per gram solids (gCGS). Gene or variant copies per gram solids. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per L (gcL). Gene or variant copies per litre. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per mL (gcMl). Gene or variant copies per millilitre of solution. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per PMMoV copy (gcPpmov). Gene or variant copies per copy of PMMoV. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gram per cubic metre (gm3). Density unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Grams (grams). A unit of mass or weight. Status: active. First released: 2.0.0. Last updated: 2.0.0.

General Hospital Admissions (hosa). Hospital admissions or patients newly admitted to hospital. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Hospital Census (hosc). Hospital census or the number of people admitted with an ailment. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Hours (hours). A unit for indicating a length of time in hours. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Intensive care unit patients (icu). Units for describing a population measure of patients who are in intensive care due to a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Kilogram per second (kgS). Kilograms per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Kilolitres (kl). Kilolitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Litres (l). Litres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Cubic metres per day (m3D). Cubic metres per day. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per hour (m3H). Cubic metres per hour. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per second (m3S). Cubic metre per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Millions of gallons per day (MG/D) (mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to measure flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre (mgL). Milligrams per litre. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Minutes (minutes). A unit for indicating a length of time in minutes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millilitres (ml). Millilitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Megalitres per day (ML/d) (mld). Megalitres per day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres (mm). Unit part for the SI unit of millimetres. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Months (months). A unit for indicating a length of time in months. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Metres per second (ms). metres per second. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unit not applicable (naUnit). Not applicable for units. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population of newly vaccinated persons (newVax). Population of newly vaccinated persons Status: active. First released: 2.0.0. Last updated: 2.0.0.

Nephelometric turbidity unit (ntu). Nephelometric Turbidity Units, a unit used in the measurement of turbidity (see “turbidity” under measureIDs). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Other Unit deprecated (otunDep). Other measurement of viral copies or wastewater treatment plant parameter. Add description to UnitOther. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Percent (perc). Percentage. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Percent recovery (percRec). Percent of the surrogate recovery for a recovery efficiency control assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Population equivalents (popEq). A unit for measuring of design capacity for wastewater treatment plants, represented as the approximate amount of water a single person would use. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent positive (pp). Percent positive of sample measures. Status: active. First released: 1.0.0. Last updated: 2.0.0.

parts per million (ppm). Parts per million. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent positivity rate (pprt). Percent positivity rate of tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent primary sludge (pps). Percentage of total solids, for primary sludge. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of total (prop). Proportion as a percent of total. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Proportion of variant in sample (propV). Proportion of a variant as percent of total variants. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Number of positive tests (ptot). Number of positive tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Ratio (ratio). Ratio (unitless) Status: active. First released: 2.0.0. Last updated: 2.0.0.

Recovered patients (recov). Units for describing a population measure of patients who have recovered from a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rate per 100,000 (rP100). Units for describing a population measure of a rate of case incidence per 100,000 people of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Number of tests performed (test). Number of tests performed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Unitless measure (unitless). A unit for unitless measures, like pH. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Micro-Siemens per centimetre (uSCm). Micro-siemens per centimetre. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Frequency of variants detected (varFreq). A description of the frequency of a variants detected. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 1 dose of vaccine (vax1). A measure of the population with a single dose of vaccine. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 2 doses of vaccine (vax2). Population with 2 doses of vaccine Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 3 doses of vaccine (vax3). Population with 3 doses of vaccine Status: active. First released: 2.0.0. Last updated: 2.0.0.

Years (years). A unit for indicating a length of time in years. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unit sets

(unitSets). Sets of units. Contains units that are associated with part types such as measures. Sets are lists of input values. Each input value within a set is a part (has a partID) that can be reused between sets. For example, a set A could contain the values of ‘yes’ and ‘no’; and set B could contain the values of ‘yes’, ‘no’, and ‘maybe.’ partTypes with sets include aggregations, compartments, missingness, quality flags, specimens, and units. Measures with categories and dataTypes have their own specific sets (catSetID) because their input values are unique and are not reused.

Absolute humidity unit set

(absHumidUnitSet)Unit set for absolute humidity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Gram per cubic metre (gm3). Density unit. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Allele unit set

(alleleUnitSet)Unit set for alleles. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Sets

Cycle threshold or quantification cycle (Ct or Cq) (ct). Cycle thresholds in a PCR assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Detected (det). Substance detected or not detected. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Depth of coverage (doc). The sequencing read depth. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per copy of crAssphage (gcCrA). Gene or variant copies per copy of crAssphage. Status: active. First released: 1.0.0. Last updated: 2.0.0.

gene copies per day per 100,000 (gcD100). The unit for measures reflecting the gene copies per day per 100,000 people in the population. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per gram solids (gCGS). Gene or variant copies per gram solids. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per mL (gcMl). Gene or variant copies per millilitre of solution. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Gene copies per PMMoV copy (gcPpmov). Gene or variant copies per copy of PMMoV. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Bacteria unit set

(bacteriaUnitSet)Unit set for bacteria-related measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

CFU per 100 ml (cfu). Colony forming units per 100 ml of filtered sample. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Capacity unit set

(capacityUnitSet)Unit set for capacity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Millions of gallons per day (MG/D) (mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to measure flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Megalitres per day (ML/d) (mld). Megalitres per day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population equivalents (popEq). A unit for measuring of design capacity for wastewater treatment plants, represented as the approximate amount of water a single person would use. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Conductivity unit set

(conductivityUnitSet)Unit set related to conductivity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Micro-Siemens per centimetre (uSCm). Micro-siemens per centimetre. Status: depreciated. First released: 1.0.0. Last updated: 2.0.0.

Dissolved gas concentration unit set

(dissGasUnitSet)Unit set for carbon dioxide concentrations measurements in water or air. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre (mgL). Milligrams per litre. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Sets

parts per million (ppm). Parts per million. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Volume flow rate unit set

(flowUnitSet)Unit set for volume flow measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Kilogram per second (kgS). Kilograms per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cubic metres per day (m3D). Cubic metres per day. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per hour (m3H). Cubic metres per hour. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Cubic metres per second (m3S). Cubic metre per second. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Millions of gallons per day (MG/D) (mgd). A unit for measuring of design capacity for wastewater treatment plants, represented as millions of gallons per day. Also can be used to measure flow. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Megalitres per day (ML/d) (mld). Megalitres per day. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Genetics unit set

(geneticUnitSet)Unit set for genetic-related measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Breadth of coverage ($\geq 5\times$ depth) (boc). Positions with read depth greater or equal to 5. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Detected (det). Substance detected or not detected. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Depth of coverage (doc). The sequencing read depth. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Gene copies per L (gcL). Gene or variant copies per litre. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent recovery (percRec). Percent of the surrogate recovery for a recovery efficiency control assay. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of total (prop). Proportion as a percent of total. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Proportion of variant in sample (propV). Proportion of a variant as percent of total variants. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Unit set not applicable

(naUnitSet)Not applicable for unit sets. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Unit not applicable (naUnit). Not applicable for units. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population unit set

(populationUnitSet)Unit set for hospital-related measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Number of active cases (acti). Number of active cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sets

Confirmed Cases Case report date (conCase). Date that the numbers were reported publicly for a confirmed case. Typically, reported data and this measure is most commonly reported and used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Number of confirmed cases (conf). Number of confirmed cases. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases episode date (confEp). Episode date is the earliest of onset, test or reported date for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed Cases Onset date (confOn). Earliest that symptoms were reported for a confirmed case. This data is often not known and reported. In lieu, episode is used. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Confirmed cases test date (conTest). Date that the covid-19 test was performed for a confirmed case. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Cumulative count (cuCo). Units for describing a population measure of a cumulative count of cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Daily count (daiCo). Units for describing a population measure of a daily count of new cases of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Deaths (death). Units for describing a population measure of patients who have died from a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

General Hospital Admissions (hosa). Hospital admissions or patients newly admitted to hospital. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Hospital Census (hosc). Hospital census or the number of people admitted with an ailment. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Intensive care unit patients (icu). Units for describing a population measure of patients who are in intensive care due to a given cause. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population of newly vaccinated persons (newVax). Population of newly vaccinated persons Status: active. First released: 2.0.0. Last updated: 2.0.0.

Percent positivity rate (pprt). Percent positivity rate of tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Number of positive tests (ptot). Number of positive tests conducted within a day in a given region. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Recovered patients (recov). Units for describing a population measure of patients who have recovered from a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Rate per 100,000 (rP100). Units for describing a population measure of a rate of case incidence per 100,000 people of a given disease. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Number of tests performed (test). Number of tests performed. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Population with 1 dose of vaccine (vax1). A measure of the population with a single dose of vaccine. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Population with 2 doses of vaccine (vax2). Population with 2 doses of vaccine Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Population with 3 doses of vaccine (vax3). Population with 3 doses of vaccine Status: active. First released: 2.0.0. Last updated: 2.0.0.

Precipitation unit set

(precipitationUnitSet)Unit set for percipitation measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Centimetres (cm). Unit part for the SI unit of centimetres. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres (mm). Unit part for the SI unit of millimetres. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Relative humidity unit set

(relHumidUnitSet)Unit set for relative humidity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Percent (perc). Percentage. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Specific humidity unit set

(specHumidUnitSet)Unit set for specific humidity. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Ratio (ratio). Ratio (unitless) Status: active. First released: 2.0.0. Last updated: 2.0.0.

Standard concentration unit set

(stdConcentrationUnitSet)Unit set for concentration measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Milligrams per litre (mgL). Milligrams per litre. Status: active. First released: 1.1.0. Last updated: 2.0.0.

parts per million (ppm). Parts per million. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Percent primary sludge (pps). Percentage of total solids, for primary sludge. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Temperature unit set

(temperatureUnitSet)Unit set for temperature measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Degrees Celcius (cel). Degrees Celsius. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Time unit set

(timeUnitSet)The unit set for measures of time. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Days (days). A unit for indicating a length of time in days. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Hours (hours). A unit for indicating a length of time in hours. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Minutes (minutes). A unit for indicating a length of time in minutes. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Sets

Months (months). A unit for indicating a length of time in months. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Years (years). A unit for indicating a length of time in years. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Turbidity unit set

(turbidityUnitSet)Unit set for turbidity measurements. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Nephelometric turbidity unit (ntu). Nephelometric Turbidity Units, a unit used in the measurement of turbidity (see “turbidity” under measureIDs). Status: active. First released: 2.0.0. Last updated: 2.0.0.

Unitless unit set

(unitlessUnitSet)Unit set for measurements that does not have units. Status: active. First Released: 2.0.0. Last updated: 2.0.0.

Unitless measure (unitless). A unit for unitless measures, like pH. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Volume unit set

(volumeUnitSet)Unit set of volume measurements. Status: active. First Released: 1.0.0. Last updated: 2.0.0.

Grams (grams). A unit of mass or weight. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Kilolitres (kl). Kilolitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Litres (l). Litres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millilitres (ml). Millilitres of volume Status: active. First released: 2.0.0. Last updated: 2.0.0.

Millimetres (mm). Unit part for the SI unit of millimetres. Status: active. First released: 1.1.0. Last updated: 2.0.0.

Percent positive (pp). Percent positive of sample measures. Status: active. First released: 1.0.0. Last updated: 2.0.0.

Sub-variant or lineage (subVar). A unit used to report a specific genetic lineage, or to specify that a reported variant is a sub-variant of another. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Frequency of variants detected (varFreq). A description of the frequency of a variants detected. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Wind speed unit set

(windSpeedUnitSet)Unit set for wind speed measurements. Status: deprecated. First Released: 1.0.0. Last updated: 2.0.0.

Metres per second (ms). metres per second. Status: active. First released: 2.0.0. Last updated: 2.0.0.

Tables

ODM version v2.0.0.

The ODM database model has 21 tables in its full relational database. The full model is commonly referred to as “long” tables as it stores data with one measurement per row. Each table has mandatory and optional headers.

i Data storage

Store all ODM data in the 21 recommended tables to ensure optimal interoperability and support by the [ODM validation toolkit](#) and other ODM libraries.

Users can, however, create custom tables or ‘views’ by merging fields from any table. A common format for data entry format are “wide” table format, with *one day per row* and multiple measures or attributes per row. It is important to note that this guide does not provide instructions on generating wide tables.

Table types

There are five table types:

1) **Regular report tables.** The **measures** and **samples** tables are used for daily reporting of new measurements and information on sample collection. There are three tables that support regular reporting:- **measureSets**

Tables

(optional) is used to link measures, `sampleRelationships` (optional) is used when samples are pooled or split, and `qualityReports` is used for detailed quality assurance and control measures.

2) **Contact information tables.** The `sites`, `organizations`, `contacts`, and `addresses` tables store contact information for the testing site and the person who performed the testing.

3) **Protocol or methods tables.** The `protocols` and `instruments` tables hold information on the methods used for sample collection or measurement. `protocolSteps` and `protocolRelationships` tables are used in conjunction with the `protocols` table.

4) **Other report tables.** The `polygons` table stores information on the sample area and the `datasets` table includes information on the data custodian and owner.

5) **Dictionary reference tables.** The `parts` table describes all elements of the ODM, including tables, table headers, measures, methods, categories, and units. The `sets` table stores collections of parts, such as units grouped together in a `unitSet`. The `languages` and `translations` tables support translations.

Header roles

Each table in the ODM has column headers, also known as table variables, fields, or entity relationship attributes. The header serves as the top row and contains the variable name.

Primary key (PK) - With the exception of the dictionary tables, all tables have one primary key. The primary key serves as a unique identifier for each row in the table, meaning that a primary key value cannot be repeated.

Mandatory, optional, mandatoryIf

Foreign key (FK) - Tables may have one or more foreign keys. A foreign key connects an entry to another table with a primary key. This allows for relationships to be established between tables, making it possible to link data between different tables.

Header (header) - This label is used in a table description when a header is not a primary or foreign key.

Mandatory, optional, mandatoryIf

Each header may be classified as either **mandatory** or **optional**. In the **measures** table, examples of mandatory headers are the **measureID** and **measurement date**. The ODM validation toolkit can be used to validate data, and will return an error if a table lacks a mandatory header or if there are missing row entries for mandatory fields.

Additionally, there are special headers called **mandatoryIf** headers. These headers are mandatory if certain conditions apply. For example, the **protocol steps** table has row entries that are either **measures** or **methods**. If the row entry is a **measure**, then a **unit** becomes **mandatory** since all measures must have a unit. Therefore, units are **mandatoryIf** in the protocol steps table.

Table list

Address table

Contact table

Dataset table

Instrument table

Language Look-up table

Tables

Measure report table

Measure set report table

Organization table

Parts Look-up table

Polygon table

Protocol relationships table

Protocols table

Protocol steps table

Quality reports table

Sample relationships table

Sample report table

Sets look-up table

Sites table

Translation look-up table

Address table

partID: addresses. The table that contains information about addresses. Addresses can be recorded for sites, organizations or contacts (individuals). The Sites, Organizations, and Contacts tables include a link to the addresses table through Address ID.

Address ID (addressID). A unique identifier for an address. Role: Primary Key. Requirement: mandatory. Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Address Line 1 (addL1). Line 1 (the street name, number and direction) for a given address.Role: Header.Requirement: mandatory.Data type: varchar.

Address Line 2 (addL2). Line 2 (the unit number) for a given address.Role: Header.Requirement: optional.Data type: varchar.

City (city). The city where a site or organization is located; part of the address.Role: Header.Requirement: mandatory.Data type: varchar.

State, Province, or Region (stateProvReg). The state, province, or region where a site or organization is located; part of the address.Role: Header.Requirement: mandatory.Data type: varchar.

Postal or Zip Code (pCode). The zip code or postal code for a given address, specifying a specific geographic area.Role: Header.Requirement: optional.Data type: varchar.

Country (country). The country where a site or organization is located; part of the address.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Tables

Contact table

partID: contacts. The table that contains information about a contact; a person who is the contact of a site or laboratory. adapt from version 1 documentation

Contact ID (contactID). A unique identifier for a given contact person.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

First name of contact (firstName). Specifies the first name of a given contact.Role: Header.Requirement: optional.Data type: varchar.

Last name of contact (lastName). Specifies the last name of a given contact.Role: Header.Requirement: optional.Data type: varchar.

Contact email (email). Contact e-mail address, for the lab.Role: Header.Requirement: mandatory.Data type: varchar.

Contact phone (phone). Contact phone number, for the lab.Role: Header.Requirement: optional.Data type: varchar.

Role of contact (role). Specifies the organizational role of a given contact.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Dataset table

partID: datasets. A report table for capturing details about data's parental data set and data custodians. Supplying attribution for data collectors.

Parent dataset ID (parDatasetID). The datasetID that is a parent to another datasetID.Role: Foreign Key.Requirement: optional.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset creation date (datasetDate). Specifies the date a given dataset was created.Role: Header.Requirement: recommended.Data type: date-time.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

License (license). The license of a dataset.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English.Role: Foreign Key.Requirement: optional.Data type: varchar.

Funder Contact ID (funderCont). A unique identifier for funders.Role: Foreign Key.Requirement: recommended.Data type: varchar.

Tables

Custody Contact ID (custodyCont). A unique identifier for data custodians.Role: Foreign Key.Requirement: recommended.Data type: varchar.

Funding agency ID (funderID). The funding agency of the dataset.Role: Foreign Key.Requirement: optional.Data type: varchar.

Data custodian ID (custodyID). The data custodian of a database.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Instrument table

partID: instruments. The table that contains information about instruments. adapt from version 1 documentation

Instrument ID (instrumentID). A unique identifier for an instrument.Role: Primary Key.Requirement: mandatory.Data type: categorical.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Model (model). Model number or version of the instrument.Role: Header.Requirement: mandatory.Data type: varchar.

Manufacturer (manufacturer). Manufacturer of an instrument.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Instrument Type (insType). Type of instrument used to perform the measurement.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Describe other instrument type, if applicable (insTypeOth). Description of the instrument in case it is not listed in instrumentType.Role: Header.Requirement: optional.Data type: varchar.

Index (index). Index number in case the measurement was taken multiple times.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Language Look-up table

partID: languages. Look up table for all languages, used to give structure to the translation table.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Language family (langFam). Specifies the language family of a given language for translation and language tracking purposes.Role: Header.Requirement: mandatory.Data type: varchar.

Language name (langName). Specifies the name of the language in roman alphabet characters for translation and language tracking purposes.Role: Header.Requirement: mandatory.Data type: varchar.

Native Name (natName). The native name of the language, i.e. what the language is called by its speakers.Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-1 (iso6391). The first part of the ISO 639 series of international standards for language codes. Part 1 covers the registration of two-letter codes. There are 183 two-letter codes registered as of June 2021. The registered codes cover the world's major languages.Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-2B (iso6392B). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding two-letter (ISO 639-1) codes where they exist. The 'B' specifies the bibliographic code (B code).Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-2T (iso6392T). A set of international standards that lists short codes for language names. These ISO639-2 are the three-letter codes defined in part two (ISO 639-2) of the standard, including the corresponding

two-letter (ISO 639-1) codes where they exist. The ‘T’ specifies the terminological code (T code).Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-3 (iso6393). A set of international standards that lists short codes for language names. ISO 639-3 extends the ISO 639-2 alpha-3 codes with an aim to cover all known natural languages.Role: Header.Requirement: mandatory.Data type: categorical.

ISO639-6 (iso6396). A set of international standards that lists short codes for language names. ISO 639-6 builds off ISO639-3 with the use of four-letter codes, and allowing users to differentiate between variants of languages and language families, such as historical vs. revived versions of languages.Role: Header.Requirement: mandatory.Data type: categorical.

First released version (firstReleased). The version in which a part was first released.Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: optional.Data type: varchar.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Measure report table

partID: measures. The table that contains information and details about a given measure adapt from version 1 documentation

Tables

Report ID (measureRepID). Unique identifier for a measurement.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Foreign Key.Requirement: optional.Data type: varchar.

Sample ID (sampleID). Unique identifier for a sample.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Purpose (purpose). The reason the measure or sample was taken.Role: Foreign Key.Requirement: optional.Data type: categorical.Measure, method or attribute set: purposeSet

Polygon ID (polygonID). Unique identifier for the polygon.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site ID (siteID). Unique identifier for the location where a sample was taken.Role: Foreign Key.Requirement: recommended.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Report set ID (meaureSetRepID). Unique identifier that links together a group of related measures.Role: Foreign Key.Requirement: optional.Data type: varchar.

Analysis date start (aDateStart). Date the measurement or analysis was started.Role: Header.Requirement: optional.Data type: datetime.

Analysis date end (aDateEnd). Date the measurement or analysis was completed.Role: Header.Requirement: mandatory.Data type: datetime.

Report date (reportDate). The date a measure was reported.Role: Header.Requirement: optional.Data type: datetime.

Specimen (specimen). The substance or thing upon which the observation was made.Role: Foreign Key.Requirement: mandatory.Data type: categorical.

Fraction analyzed (fraction). Fraction of the sample that is analyzed.Role: Foreign Key.Requirement: recommended.Data type: categorical.Measure, method or attribute set: fractionSet

Group (group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations.Role: Foreign Key.Requirement: optional.Data type: categorical.

Class (class). A unique identifier for a class, which is akin to a sub-group; it's a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class.Role: Foreign Key.Requirement: optional.Data type: categorical.

Measure (measure). A measurement or observation of any substance including a biological, physical or chemical substance.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Value (value). Value of a measure, observation or attribute.Role: Header.Requirement: mandatory.Data type: varchar.

Unit (unit). The units of a measurement.Role: Foreign Key.Requirement: mandatory.Data type: categorical.

Aggregation (aggregation). Statistical measures used to report a measure. Each aggregation has a corresponding value.Role: Foreign Key.Requirement: mandatory.Data type: categorical.

Nomenclature (nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases.Role: Foreign Key.Requirement: optional.Data type: categorical.

Index (index). Index number in case the measurement was taken multiple times.Role: Header.Requirement: optional.Data type: varchar.

Tables

Measure license (measureLic). Specifies the access and use licensing for a given single measurement.Role: Header.Requirement: optional.Data type: varchar.

Reportable (reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied.Role: Header.Requirement: optional.Data type: boolean.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Measure set report table

partID: measureSets. The table that identifies sets of measures. Examples of measure sets include a set of replicates, dilutions (used to generate a Ct curve) or variants that are identified in a single sample.

Report set ID (measureSetRepID). Unique identifier that links together a group of related measures.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Foreign Key.Requirement: optional.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Organization table

partID: organizations. The table that contains information about a laboratory. adapt from version 1 documentation

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Tables

Address ID (addressID). A unique identifier for an address.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization Type (orgType). Specifies the type or purpose of a given organization.Role: Header.Requirement: recommended.Data type: categorical.Measure, method or attribute set: orgTypeSet

Organization level (orgLevel). The geographic level of an organization.Role: Header.Requirement: recommended.Data type: categorical.Measure, method or attribute set: orgLevelSet

Organization sector (orgSector). The sector of an organization.Role: Header.Requirement: recommended.Data type: categorical.Measure, method or attribute set: orgSectorSet

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Parts Look-up table

partID: parts. Look up table containing all parts in of the data model. Contains all parts, including self-referential parts.

Part identifier (partID). The unique identify of any entity within the dictionary.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Label (partLabel). A human readable label of a part.Role: Header.Requirement: mandatory.Data type: varchar.

Part types (partType). Part types describe the purpose or use of the part.Role: Header.Requirement: mandatory.Data type: categorical.

Short names (shortName). Shortened names for tables and other important parts for use in wide names.Role: Header.Requirement: mandatory.Data type: varchar.

Part description (partDesc). The description of the part.Role: Header.Requirement: mandatory.Data type: varchar.

Part instruction (partInstr). Additional notes and instructions on how a part is used and/or defined.Role: Header.Requirement: mandatory.Data type: varchar.

Domain (domain). Domain is the highest level of describing of a measure.Role: Header.Requirement: mandatory.Data type: categorical.

Specimen sets (specimenSets). Sets of specimens.Role: Foreign Key.Requirement: mandatory.Data type: categorical.

Group (group). Unique identifier for a group of measures. Mostly applicable for measures, methods, units, and aggregations.Role: Header.Requirement: mandatory.Data type: categorical.

Class (class). A unique identifier for a class, which is akin to a subgroup; it's a way of grouping parts within a given group. A group can have one or more classes to describe different parts of the class.Role: Header.Requirement: mandatory.Data type: categorical.

Nomenclature (nomenclature). A classification system to report the measure class. Only applicable to variants, mutations, and diseases.Role: Header.Requirement: mandatory.Data type: categorical.

Ontology reference (ontologyRef). Ontology reference for a part.Role: Header.Requirement: mandatory.Data type: varchar.

LaTeX expression (latExp). LaTeX expression used to generate formulas, symbols, etc. Mainly relevant for units.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Measure, method, or attribute (mmaSet). The set for a measure, method, or attribute. Only applicable for categorical fields, and the categories that populate them.Role: Header.Requirement: mandatory.Data type: categorical.

Aggregation scale (aggragationScale). A scale used for an aggregation. Only applicable for measures and units.Role: Header.Requirement: mandatory.Data type: categorical.

Aggregation set (aggregationSet). The aggregation set for a unit.Role: Header.Requirement: mandatory.Data type: categorical.

Quality set (qualitySetID). The quality set that corresponds to a given part. Only applicable for samples and measures.Role: Header.Requirement: mandatory.Data type: varchar.

Missingness sets (missingnessSets). Missingness sets for measures, methods or attributes.Role: Header.Requirement: mandatory.Data type: categorical.

Status (status). Whether the part is still active and can be used in the most current ODM version. Values are ‘active’ or ‘inactive’.Role: Header.Requirement: mandatory.Data type: varchar.

First released version (firstReleased). The version in which a part was first released.Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: mandatory.Data type: varchar.

Protocol steps table (protocolSteps). The table for collecting metadata on individual steps in a protocol, methodological process, or assay.Role: Header.Requirement: mandatory.Data type: varchar.

Protocol Relationship table required headers (protocolRelationshipsRequired). Specifies the columns required in the Protocol Organization table.Role: Header.Requirement: mandatory.Data type: varchar.

Protocol steps table required headers (protocolStepsRequired). Specifies the columns required in a Protocol Steps table.Role: Header.Requirement: mandatory.Data type: varchar.

Protocol steps table column order (protocolStepsOrder). Specifies the order of the columns in a Protocol Steps table.Role: Header.Requirement: mandatory.Data type: integer.

Protocol relationships table (protocolRelationships). The table that contains the organizational information and details (such as order) about a given protocol, recorded as a series of protocol steps (protocolSteps).Role: Header.Requirement: mandatory.Data type: varchar.

Protocol relationships table column order (protocolRelationshipsOrder). Specifies the order of the columns in the Protocol Organization table.Role: Header.Requirement: mandatory.Data type: integer.

Measure report table (measures). The table that contains information and details about a given measureRole: Header.Requirement: mandatory.Data type: varchar.

Measures table required headers (measuresRequired). Specifies the columns required in a Measures table.Role: Header.Requirement: mandatory.Data type: varchar.

Measures table column order (measuresOrder). Specifies the order of the columns in a Measures table.Role: Header.Requirement: mandatory.Data type: integer.

Measure set report table (measureSets). The table that identifies sets of measures.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Measure sets table column order (measureSetsOrder). Specifies the order of the columns in a Measure Sets table.Role: Header.Requirement: mandatory.Data type: integer.

Measure Set table required headers (measureSetsRequired). Specifies the columns required in a Measure Sets table.Role: Header.Requirement: mandatory.Data type: varchar.

Dataset table (datasets). A report table for capturing details about data's parental data set and data custodians. Supplying attribution for data collectors.Role: Header.Requirement: mandatory.Data type: varchar.

Dataset table required headers (datasetsRequired). Specifies the columns required in a Datasets table.Role: Header.Requirement: mandatory.Data type: varchar.

Datasets table column order (datasetsOrder). Specifies the order of the columns in a Datasets table.Role: Header.Requirement: mandatory.Data type: integer.

Sites table (sites). The table that contains information about a site; the location where an environmental sample was taken.Role: Header.Requirement: mandatory.Data type: varchar.

Sites table required headers (sitesRequired). Specifies the columns required in a Sites table.Role: Header.Requirement: mandatory.Data type: varchar.

Sites table column order (sitesOrder). Specifies the order of the columns in a Sites table.Role: Header.Requirement: mandatory.Data type: integer.

Sample report table (samples). The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site.Role: Header.Requirement: mandatory.Data type: varchar.

Sample table required headers (samplesRequired). Specifies the columns required in a Samples table.Role: Header.Requirement: mandatory.Data type: varchar.

Samples table column order (samplesOrder). Specifies the order of the columns in a Samples table.Role: Header.Requirement: mandatory.Data type: integer.

Address table (addresses). The table that contains information about addresses. Addresses can be recorded for sites, organizations or contacts (individuals).Role: Header.Requirement: mandatory.Data type: varchar.

Address table required headers (addressesRequired). Specifies the columns required in a Addresses table.Role: Header.Requirement: mandatory.Data type: varchar.

Addresses table column order (addressesOrder). Specifies the order of the columns in the Addresses table.Role: Header.Requirement: mandatory.Data type: integer.

Contact table (contacts). The table that contains information about a contact; a person who is the contact of a site or laboratory.Role: Header.Requirement: mandatory.Data type: varchar.

Contact table required headers (contactsRequired). Specifies the columns required in a Contacts table.Role: Header.Requirement: mandatory.Data type: varchar.

Contacts table column order (contactsOrder). Specifies the order of the columns in a Contacts table.Role: Header.Requirement: mandatory.Data type: integer.

Organization table (organizations). The table that contains information about a laboratory.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Organiation table required headers (organizationsRequired). Specifies the columns required in a Organizations table.Role: Header.Requirement: mandatory.Data type: varchar.

Organizations table column order (organizationsOrder). Specifies the order of the columns in a Organizations table.Role: Header.Requirement: mandatory.Data type: integer.

Instrument table (instruments). The table that contains information about instruments.Role: Header.Requirement: mandatory.Data type: varchar.

Instrutment table required headers (instrumentsRequired). Specifies the columns required in the Instruments table.Role: Header.Requirement: mandatory.Data type: varchar.

Instruments table column order (instrumentsOrder). Specifies the order of the columns in the Instruments table.Role: Header.Requirement: mandatory.Data type: integer.

Polygon table (polygons). The table that contains information about the geometry of a geographic area.Role: Header.Requirement: mandatory.Data type: varchar.

Polygon table required headers (polygonsRequired). Specifies the columns required in a Polygons table.Role: Header.Requirement: mandatory.Data type: varchar.

Polygons table column order (polygonsOrder). Specifies the order of the columns in a Polygons table.Role: Header.Requirement: mandatory.Data type: integer.

Language Look-up table (languages). Look up table for all languages, used to give structure to the translation table.Role: Header.Requirement: mandatory.Data type: varchar.

Language table required headers (languagesRequired). Required headers in the Languages table.Role: Header.Requirement: mandatory.Data type: varchar.

Language table column order (languagesOrder). Specifies the order of the columns in the Languages table.Role: Header.Requirement: mandatory.Data type: integer.

Translation look-up table (translations). Look up table for translations of the description, label, and instruction for all parts.Role: Header.Requirement: mandatory.Data type: varchar.

Parts Look-up table (parts). Look up table containing all parts in of the data model.Role: Header.Requirement: mandatory.Data type: varchar.

Translation table required headers (translationsRequired). Required headers in the Translations table.Role: Header.Requirement: mandatory.Data type: varchar.

Parts table column order (partsOrder). Order of headers in the Parts table.Role: Header.Requirement: mandatory.Data type: varchar.

Translation table column order (translationsOrder). Specifies the order of the columns in the Translation table.Role: Header.Requirement: mandatory.Data type: integer.

Parts table required headers (partsRequired). Required headers in the Parts table.Role: Header.Requirement: mandatory.Data type: varchar.

Sets look-up table (sets). Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together.Role: Header.Requirement: mandatory.Data type: varchar.

Sets table column order (setsOrder). Specifies the order of the columns in the Sets table.Role: Header.Requirement: mandatory.Data type: integer.

Sets table required headers (setsRequired). Required headers in the Sets table.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

Quality reports table (qualityReports). The table for recording the various quality metrics and indicators for samples and measures.Role: Header.Requirement: mandatory.Data type: varchar.

Quality report table required headers (qualityReportsRequired). Specifies the columns required in the Quality Reports table.Role: Header.Requirement: mandatory.Data type: varchar.

Quality reports table column order (qualityReportsOrder). Specifies the order of the columns in the Quality Reports table.Role: Header.Requirement: mandatory.Data type: integer.

Sample relationships table (sampleRelationships). Table for recording the relationships between samples.Role: Header.Requirement: mandatory.Data type: varchar.

Sample relationships table required headers (sampleRelationshipsRequired). Specifies the columns required in a Sample Relationships table.Role: Header.Requirement: mandatory.Data type: varchar.

Sample relationships table column order (sampleRelationshipsOrder). Specifies the order of the columns in a Sample Relationships table.Role: Header.Requirement: mandatory.Data type: integer.

Protocols table (protocols). The table for protocols.Role: Header.Requirement: mandatory.Data type: varchar.

Protocols table required headers (protocolsRequired). Specifies the columns required in the Protocols table.Role: Header.Requirement: mandatory.Data type: varchar.

Protocols table column order (protocolsOrder). Specifies the order of the columns in the Protocols table.Role: Header.Requirement: mandatory.Data type: integer.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: mandatory.Data type: varchar.

Data types (dataTypes). The data type for a part.Role: Header.Requirement: mandatory.Data type: categorical.

Minimum value part support (minValue). The minimum value of part.Role: Header.Requirement: mandatory.Data type: varchar.

Maximum value part support (maxValue). The maximum value of a part.Role: Header.Requirement: mandatory.Data type: varchar.

Minimum length (minLength). The maximum value of measure.Role: Header.Requirement: mandatory.Data type: integer.

Maximum length (maxLength). The maximum length of the value of a part or measure.Role: Header.Requirement: mandatory.Data type: integer.

Polygon table

partID: polygons. The table that contains information about the geometry of a geographic area.

Polygon ID (polygonID). Unique identifier for the polygon.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Tables

Polygon Population (polyPop). An attribute of a polygon, which specifies the population of that polygon. A rough estimate of the number of human residents.Role: Header.Requirement: optional.Data type: integer.

Type of geography (geoType). Type of geography that is represented by the polygon.Role: Header.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: geoTypeSet

European Petroleum Survey Group Coordinates (geoEPSG). The unique EPSG code specifying a given geospatial area.Role: Header.Requirement: mandatory.Data type: float.

Well-known text (geoWKT). Well-known text of the polygonRole: Header.Requirement: mandatory.Data type: varchar.

File location of polygon (fileLocation). The location of the file containing the geometry of the polygon.Role: Header.Requirement: optional.Data type: blob.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Protocol relationships table

partID: protocolRelationships. The table that contains the organizational information and details (such as order) about a given protocol, recorded as a series of protocol steps (protocolSteps). Protocols are a group of related protocol steps. For example a protocol step can describe a sample concentration, extraction, inhibition or a measurement. Protocols can include both methodID (a procedure) and measureIDs (a measure). An example of a method is centrifugation; where the centrifugation rotation speed is described as a measure.

Protocol ID container (protocolIDContainer). Unique identifier for the protocol and the steps and other protocols that make it up.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Protocol ID object (protocolIDObj). The object of the relationship between one protocol and another protocol or step, within a given protocol ID container.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Step ID Object (stepIDObj). The object of the relationship between one protocol step and another step or protocol, within a given protocol ID container.Role: Foreign Key.Requirement: mandatoryIf.Data type: integer.

Sample relationship (sampleRelID). Attribute for specifying the relationship between a sample subject and object.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Protocol ID subject (protocolIDSub). The subject of the relationship between one protocol and another protocol or step, within a given protocol ID container.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Step ID Subject (stepIDSub). The subject of the relationship between one protocol step and another step or protocol, within a given protocol ID

Tables

container.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Protocols table

partID: protocols. The table for protocols.

Source Protocol ID (sourceProtocol). A protocol that served as a basis for another protocol. This may be due to version changes, or a referencing a protocol that was later adapted.Role: Foreign Key.Requirement: optional.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: optional.Data type: varchar.

Summary (summ). Short description of the assay and how it is different from the other assay methods.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Protocol version (protocolVersion). Specifies the version of a method set.Role: Header.Requirement: optional.Data type: integer.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Protocol steps table

partID: protocolSteps. The table for collecting metadata on individual steps in a protocol, methodological process, or assay.

Protocol Step ID (stepID). The unique identifier for a specific protocol step.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Method (method). A procedure for collecting a sample or performing a measure.Role: Foreign Key.Requirement: mandatoryIf.Data type: categorical.

Measure (measure). A measurement or observation of any substance including a biological, physical or chemical substance.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Summary (summ). Short description of the assay and how it is different from the other assay methods.Role: Header.Requirement: optional.Data type: varchar.

Tables

Protocol step source ID (sourceStep). Specifies the protocol step which serves as a basis for a given protocol step.Role: Foreign Key.Requirement: optional.Data type: varchar.

Protocol step version (stepVer). Specifies the version of a given protocol step.Role: Header.Requirement: optional.Data type: varchar.

Reference link (refLink). Link to the reference material for a part. May link to literature on a method, measure, etc.Role: Header.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Foreign Key.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Instrument ID (instrumentID). A unique identifier for an instrument.Role: Foreign Key.Requirement: optional.Data type: categorical.

Value (value). Value of a measure, observation or attribute.Role: Header.Requirement: mandatoryif.Data type: varchar.

Unit (unit). The units of a measurement.Role: Foreign Key.Requirement: mandatoryIf.Data type: categorical.

Aggregation (aggregation). Statistical measures used to report a measure. Each aggregation has a corresponding value.Role: Foreign Key.Requirement: mandatoryIf.Data type: categorical.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Quality reports table

partID: qualityReports. The table for recording the various quality metrics and indicators for samples and measures.

Quality report ID (qualityID). A unique identifier for a given quality report.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Report ID (measureRepID). Unique identifier for a measurement.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Sample ID (sampleID). Unique identifier for a sample.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Report set ID (measureSetRepID). Unique identifier that links together a group of related measures.Role: Foreign Key.Requirement: mandatoryIf.Data type: varchar.

Quality flag (qualityFlag). A field for reporting any quality concerns - of lack thereof - for a sample or measure.Role: Header.Requirement: mandatory.Data type: categorical.

Severity indicator (severity). An indicator of the severity or seriousness of a quality flag.Role: Header.Requirement: optional.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Tables

Sample relationships table

partID: sampleRelationships. Table for recording the relationships between samples. Samples can be pooled or split. The sample relationships table holds information on parent-child relationships between samples, and allow for tracking sample lineage for single and pooled samples.

Sample ID subject (sampleIDSubject). Populated by **sampleID** - this specifies the subject of a sample relationship.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Sample relationship (sampleRelID). Attribute for specifying the relationship between a sample subject and object.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Sample ID object (sampleIDObject). Populated by **sampleID** - this specifies the object of a sample relationship.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Sample report table

partID: samples. The table that contains information about a sample. A sample is defined as a representative volume of wastewater (or other forms of water or liquid), air, or surface area taken from a site. Samples can be combined, split, stored and reused. The **Sample relationships** (sampleRelationships) is used to describe the relationships between samples.

Sample ID (sampleID). Unique identifier for a sample.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Protocol ID (protocolID). A unique identifier for a given protocol.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site ID (siteID). Unique identifier for the location where a sample was taken.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Purpose (purpose). The reason the measure or sample was taken.Role: Foreign Key.Requirement: optional.Data type: categorical.Measure, method or attribute set: purposeSet

Sample material (saMaterial). Type of sample.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: sampleMatSet

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Sample origin (origin). An attribute of a sample specifying the origin.Role: Header.Requirement: optional.Data type: categorical.

Replicate Type (repType). Attribute of a sample, specifying whether the sample is unique, or a replicate. And if it is a replicate, what type of replicate is it?Role: Foreign Key.Requirement: optional.Data type: categorical.Measure, method or attribute set: replicateSet

Sample collection type (collType). The type of collection.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: collectSet

Tables

Collection period (collPer). Collection period. The time period over which the sample was collected, in hours. Alternatively, use collectionStart and collectionEnd.Role: Header.Requirement: mandatory.Data type: float.

Collection number (collNum). The number of subsamples that were combined to create the sample. Use NA for continuous, proportional or passive sampling.Role: Header.Requirement: mandatory.Data type: integer.

Collection number and period (collNumPer). Composite collection number.period.Role: Header.Requirement: mandatory.Data type: varchar.

Pooled (pooled). Is this a pooled sample, and therefore composed of multiple child samples obtained at different sites.Role: Header.Requirement: optional.Data type: boolean.

Collection date time (collDT). For grab samples this is the date, time and timezone the sample was taken.Role: Header.Requirement: mandatory.Data type: datetime.

Collection date time start (collDTStart). For integrated time averaged samples this is the date, time and timezone the sample was started being taken.Role: Header.Requirement: optional.Data type: datetime.

Collection date time end (collDTEnd). For integrated time average samples this is the date, time and timezone the sample was finished being taken.Role: Header.Requirement: optional.Data type: datetime.

Date sample was sent (sentDate). The date the sample was sent for analyses at a laboratory.Role: Header.Requirement: optional.Data type: datetime.

Date sample recieved (recDate). The date the sample was received at the laboratory for analysis.Role: Header.Requirement: optional.Data type: datetime.

Reportable (reportable). Flag for whether a measure is reportable or not, based on confidence in the measure and methods applied.Role: Header.Requirement: optional.Data type: boolean.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Sets look-up table

partID: sets. Look up table for all sets, managing how categorical inputs for various methods and attributes are grouped together. This table manages many-to-many relationships and category recycling within the ODM.

Set ID (setID). The unique identifier of a value set.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Set type (setType). The type of set. i.e. quality set, aggregation set, unit setRole: Header.Requirement: mandatory.

Set Value (setValue). The partID for the set value or category.Role: Foreign Key.Requirement: mandatory.

Label (partLabel). A human readable label of a part.Role: Header.Requirement: mandatory.Data type: varchar.

Part instruction (partInstr). Additional notes and instructions on how a part is used and/or defined.Role: Header.Requirement: mandatory.Data type: varchar.

Status (status). Whether the part is still active and can be used in the most current ODM version. Values are 'active' or 'inactive'.Role: Header.Requirement: mandatory.Data type: varchar.

Tables

First released version (firstReleased). The version in which a part was first released.Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: optional.Data type: varchar.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Sites table

partID: sites. The table that contains information about a site; the location where an environmental sample was taken. The site of an environmental sample. Information in the site table does not regularly change. Consider using the MeasureReport table if the information changes often.

Parent Site ID (parSiteID). The siteID that is a parent to another siteID; usually the sub-site will be contained in the parentSite.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site ID (siteID). Unique identifier for the location where a sample was taken.Role: Primary Key.Requirement: mandatory.Data type: varchar.

Dataset ID (datasetID). The name of the dataset that stores information for MeasureReport, SampleReport and other reporting tables.Role: Foreign Key.Requirement: optional.Data type: varchar.

Polygon ID (polygonID). Unique identifier for the polygon.Role: Foreign Key.Requirement: optional.Data type: varchar.

Site Type (siteType). Type of site or institution where sample was taken.Role: Foreign Key.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: siteTypeSet

Sample shed (sampleShed). A geographic area, physical space, or structure. A sample is taken from a sampleshed for a representative measurement of a substance(s).Role: Header.Requirement: mandatory.Data type: categorical.Measure, method or attribute set: shedSet

Address ID (addressID). A unique identifier for an address.Role: Foreign Key.Requirement: optional.Data type: varchar.

Organization ID (organizationID). A unique identifier for the organization to which the reporter is affiliated.Role: Header.Requirement: optional.Data type: varchar.

Contact ID (contactID). A unique identifier for a given contact person.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Name (name). Name of a domain, specimen, group, class, attribute, unit, nomenclature or other entity.Role: Header.Requirement: recommended.Data type: varchar.

Description (descr). A detailed description as an attribute for describing results or program elements.Role: Header.Requirement: optional.Data type: varchar.

Primary reporting authority (repOrg1). The primary or most responsible authority for routine surveillance reports or findings, or where the site is located.Role: Header.Requirement: optional.Data type: varchar.

Secondary reporting authority (repOrg2). The secondary, additional or alternative authority for routine surveillance reports or findings, or where the site is located.Role: Header.Requirement: optional.Data type: varchar.

Tables

Population Served (popServ). An attribute of a site, which specifies the population of/population served by a given site.Role: Header.Requirement: optional.Data type: categorical.

Latitude (geoLat). Geographical location, latitude in decimal coordinates, ie.: (45.424721)Role: Header.Requirement: mandatoryIf.Data type: float.

Longitude (geoLong). Geographical location, longitude in decimal coordinates, ie.: (-75.695000)Role: Header.Requirement: mandatoryIf.Data type: float.

European Petroleum Survey Group Coordinates (geoEPSG). The unique EPSG code specifying a given geospatial area.Role: Header.Requirement: mandatoryIf.Data type: float.

Last edited (lastEdited). The date the entry was last updated.Role: Header.Requirement: optional.Data type: datetime.

Notes (notes). A note used to describe details that are not captured in other attributesRole: Header.Requirement: optional.Data type: varchar.

Translation look-up table

partID: translations. Look up table for translations of the description, label, and instruction for all parts. The default language if a translation is not specified is English.

Language ID (lang). Language code for translation purposes. Specifies the language for each translation, other than the default English.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Part identifier (partID). The unique identify of any entity within the dictionary.Role: Foreign Key.Requirement: mandatory.Data type: varchar.

Label (partLabel). A human readable label of a part.Role: Header.Requirement: mandatory.Data type: varchar.

Part description (partDesc). The description of the part.Role: Header.Requirement: mandatory.Data type: varchar.

Part instruction (partInstr). Additional notes and instructions on how a part is used and/or defined.Role: Header.Requirement: optional.Data type: varchar.

First released version (firstReleased). The version in which a part was first released.Role: Header.Requirement: mandatory.Data type: varchar.

Last updated version (lastUpdated). The version in which the part was last updated.Role: Header.Requirement: mandatory.Data type: varchar.

Changes column (changes). A column for recording the changes from a previous version.Role: Header.Requirement: optional.Data type: varchar.

Notes (notes). A note used to describe details that are not captured in other attributes.Role: Header.Requirement: optional.Data type: varchar.

Page created on Fri Mar 17 00:00:00 2023.

